

Florida Department of Law Enforcement
Alcohol Testing Program

FEB 21 2016

FDLE
 Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: MARION COUNTY SO
 Time of Inspection: 15:56

Date of Inspection: 02/17/2016

Serial Number: 80-000813
 Software: 8100.27

| Check or Test | YES | NO |
|---|-----|----|
| Date and/or Time Adjusted | | No |
| Diagnostic Check (Pre-Inspection): OK | Yes | |
| Alcohol Free Subject Test: 0.000 | Yes | |
| Mouth Alcohol Test: Slope Not Met | Yes | |
| Interferent Detect Test: Interferent Detect | Yes | |
| Diagnostic Check (Post-Inspection): OK | Yes | |


| Alcohol Free Test (g/210L) | 0.05g/210L Test (g/210L) Lot#:201403F Exp: 03/25/2016 | 0.08g/210L Test (g/210L) Lot#:201502G Exp: 02/24/2017 | 0.20g/210L Test (g/210L) Lot#:201505A Exp: 05/12/2017 | 0.08 g/210L Dry Gas Std Test (g/210L) Lot#:83550 Exp: 09/26/2018 |
|----------------------------|---|---|---|--|
| 0.000 | 0.050 | 0.079 | 0.201 | 0.080 |
| 0.000 | 0.050 | 0.080 | 0.201 | 0.080 |
| 0.000 | 0.050 | 0.079 | 0.201 | 0.080 |

Number of Simulators Used: 5

Remarks:
 WITHIN COMPLIANCE

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

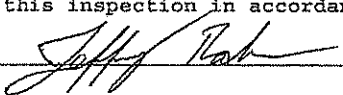

JEFFERY A RAKER

Signature and Printed Name

02/17/2016
 Date

99-80-000813

4:04 02/17/2016

MARION COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000813
02/17/2016
Software: 8100.27

| Test | g/210L | Time |
|--------------------|--------|-------|
| Air Blank | 0.000 | 15:00 |
| Control Test | 0.050 | 15:01 |
| Air Blank | 0.000 | 15:01 |
| Control Test | 0.049 | 15:02 |
| Air Blank | 0.000 | 15:02 |
| Control Test | 0.050 | 15:03 |
| Air Blank | 0.000 | 15:04 |
| Control Test Stats | | |
| Average | 0.0497 | |
| Std Dev | 0.0006 | |
| Rel Std Dev(%) | 1.1625 | |

[Signature]
Operator's Signature

MARION COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000813
02/17/2016
Software: 8100.27

| Test | g/210L | Time |
|--------------------|--------|-------|
| Air Blank | 0.000 | 15:05 |
| Control Test | 0.079 | 15:06 |
| Air Blank | 0.000 | 15:06 |
| Control Test | 0.079 | 15:07 |
| Air Blank | 0.000 | 15:07 |
| Control Test | 0.079 | 15:08 |
| Air Blank | 0.000 | 15:09 |
| Control Test Stats | | |
| Average | 0.0790 | |
| Std Dev | 0.0000 | |
| Rel Std Dev(%) | 0.0000 | |

[Signature]
Operator's Signature

MARION COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000813
02/17/2016
Software: 8100.27

| Test | g/210L | Time |
|--------------------|--------|-------|
| Air Blank | 0.000 | 15:11 |
| Control Test | 0.199 | 15:12 |
| Air Blank | 0.000 | 15:12 |
| Control Test | 0.200 | 15:13 |
| Air Blank | 0.000 | 15:13 |
| Control Test | 0.200 | 15:14 |
| Air Blank | 0.000 | 15:14 |
| Control Test Stats | | |
| Average | 0.1997 | |
| Std Dev | 0.0006 | |
| Rel Std Dev(%) | 0.2892 | |

[Signature]
Operator's Signature