



INSTRUMENT PROCESSING SHEET

Florida Department of
Law Enforcement

Agency Rockledge Police Department

S/N 80-001261

Date In 12/7/2020

DI Completion Date 12/7/2020

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake Performed By <u>DERR</u> ☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____	Quality Checks Performed By <u>DERR</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>268</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 104</u> 32 mm <u>0.105</u> (.139 - .169) 36 mm <u>0.128</u> (.156 - .190) 53 mm <u>0.203</u> (.228 - .278) 103 mm <u>0.484</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28663</u> ☒ Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> <tr> <td>0.050</td> <td>SD3967</td> <td>202010A 10/05/2022</td> </tr> <tr> <td>0.080</td> <td>SD3968</td> <td>202010B 10/05/2022</td> </tr> <tr> <td>0.200</td> <td>SD3969</td> <td>202010D 10/06/2022</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG003005 01/30/2022</td> </tr> </table>	Simulator	Serial #	Lot #/Exp	0.050	SD3967	202010A 10/05/2022	0.080	SD3968	202010B 10/05/2022	0.200	SD3969	202010D 10/06/2022	0.080 DGS	N/A	AG003005 01/30/2022	Flow Calibration Performed By <u>DERR</u> Flow Column # <u>ATP101</u> ☒ 5L/min - 17mm ☒ 15L/min - 53mm ☒ 30L/min - 103mm ☒ R-Value <u>224</u> ☒ Post Calibration Verification (L/s) Flow Column # <u>106</u> 32 mm <u>0.125/0.140</u> (.139 - .169) 36 mm <u>0.140/0.160</u> (.156 - .190) 53 mm <u>0.222/0.234</u> (.228 - .278) 103 mm <u>0.500/0.500</u> (.447 - .547)																																													
Simulator	Serial #	Lot #/Exp																																																												
0.050	SD3967	202010A 10/05/2022																																																												
0.080	SD3968	202010B 10/05/2022																																																												
0.200	SD3969	202010D 10/06/2022																																																												
0.080 DGS	N/A	AG003005 01/30/2022																																																												
Final Release Date FDLE Alcohol Testing Program	Digitally signed by FDLE Alcohol Testing Program Date: 2020.12.16 01:10:28 -05'00'																																																													
Calibration Adjustment Performed By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table>	Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection Performed By <u>DERR</u> Barometric Pressure ID# <u>28199</u> Gauge <u>1011</u> Instrument <u>1011</u> Mouth Alcohol Solution Lot # <u>2020A</u> Acetone Stock Solution Lot # <u>2019A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> <tr> <td>0.000</td> <td>SD3965</td> </tr> <tr> <td>Interferent</td> <td>SD3966</td> </tr> <tr> <td>0.050</td> <td>SD3967</td> </tr> <tr> <td>0.080</td> <td>SD3968</td> </tr> <tr> <td>0.200</td> <td>SD3969</td> </tr> </table>		Simulator	Serial Number	0.000	SD3965	Interferent	SD3966	0.050	SD3967	0.080	SD3968	0.200	SD3969
Simulator	Serial Number	Lot Number	Expiration																																																											
0.000		N/A	N/A																																																											
0.040																																																														
0.100																																																														
0.200																																																														
0.300																																																														
0.080 DGS	N/A																																																													
Simulator	Serial Number	Lot Number	Expiration																																																											
0.050																																																														
0.080																																																														
0.200																																																														
0.080 DGS	N/A																																																													
Simulator	Serial Number																																																													
0.000	SD3965																																																													
Interferent	SD3966																																																													
0.050	SD3967																																																													
0.080	SD3968																																																													
0.200	SD3969																																																													
Notes/Suggested Service: <u>Flow calibration twice to bring values closer to nominal. DERR</u> <u>Instrument would not upload, dial tone issue.</u> <u>Instrument was uploaded through direct connect andt was sent out for repair. DERR</u> _____ _____	Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☒ Form 40 ☐ Calibration Adjustment ☐ Other _____ ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☐ Return to/Place into Evidentiary Use ☒ Remain Out of Evidentiary Use ☐ Conduct an Agency Inspection Before Evidentiary Use																																																													
	<table style="width:100%;"> <tr> <td style="width:50%; text-align: center;"> 2020.12.15 15:52:57 -05'00' Michael S. Haughey Tech Review / Date </td> <td style="width:50%; text-align: center;"> 2020.12.16 01:09:41 [Signature] Admin Review / Date </td> </tr> </table>		2020.12.15 15:52:57 -05'00' Michael S. Haughey Tech Review / Date	2020.12.16 01:09:41 [Signature] Admin Review / Date																																																										
2020.12.15 15:52:57 -05'00' Michael S. Haughey Tech Review / Date	2020.12.16 01:09:41 [Signature] Admin Review / Date																																																													

Return Material Authorization

Ship to: ☐ CMI, Inc.

☒ Enforcement Electronics

Shipment to repair facility authorized by: Ivette Gomez on 12/09/2020

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-001261

Bill To Address:

Rockledge PD

Ivette Gomez

1776 JACK OATES BLVD

ROCKLEDGE FL 32955

Ship to Address:

FDLE ATP FMROC

Attn: Michael Haughey

4700 Terminal Dr. Suite 1

Fort Myers, FL 33907

Reason for Return:

Modem issue. Instrument responds "No Dial Tone" when attempting to upload. Tried at multiple locations/lines, all with the same response. Even responded with No Dial Tone on a line where another instrument had just uploaded.

Please choose one of the following options:

☐ 1. I _____, authorize all repairs.

☐ 2. I _____, authorize repairs up to \$ _____.

☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Ivette Gomez

Phone #: 321-508-762

Email: igomez@rockledgepolice.org

ATP Contact Name: Michael Haughey

ATP Email: michaelhaughey@fdle.state.fl.us

MX



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2729 Fort Knox Blvd.
Bldg. 2, Suite 1300
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001261, manufactured by CML, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number: 80-001261
Owning Agency: ROCKLEDGE P.D.
Calibration Date: 12/07/2020
Calibration Time: 20:13

UNCERTAINTY* ±	
0.050 g/ 210 L	0.004
0.080 g/ 210 L	0.005
0.200 g/ 210 L	0.007
0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration. *Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 April 2020
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1

Date 12/07/2020

DAVID REYES-RIVERA,
Department Inspector

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: ROCKLEDGE P.D.
Time of Inspection: 20:13

Date of Inspection: 12/07/2020

Serial Number: 80-001261
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202010A Exp: 10/05/2022	0.08g/210L Test (g/210L) Lot#:202010B Exp: 10/05/2022	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG003005 Exp: 01/30/2022
0.000	0.048	0.079	0.199	0.077
0.000	0.048	0.079	0.200	0.077
0.000	0.048	0.079	0.200	0.078
0.000	0.049	0.079	0.199	0.077
0.000	0.049	0.079	0.199	0.078
0.000	0.049	0.079	0.200	0.076
0.000	0.049	0.079	0.200	0.078
0.000	0.048	0.079	0.200	0.077
0.000	0.049	0.079	0.200	0.077
0.000	0.049	0.079	0.200	0.077

Standard Deviations	0.0005	0.0000	0.0004	0.0006
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

David E Reyes Rivera DAVID E REYES-RIVERA
Signature and Printed Name

12/07/2020
Date

MX
SR
2020.12.16
01:08:15 -05'00'

Type of Test	Serial Number	Agency	Date	Performed By
Stabilities	80-001261	Rockledge Police Department	12/7/2020	DERR <i>DELL</i>

0.05g/210L 0.047 to 0.053 ☒	0.08g/210L 0.077 to 0.083 ☒	0.20g/210L 0.194 to 0.206 ☒	DGS 0.08g/210L 0.077 to 0.083 ☒
ROCKLEDGE P.D. Intoxilyzer - Alcotest Analyzer Model 8000 SN 80-001261 12/07/2020 Software: 8100.27 Test g/210L Time Air Blank 0.000 17:42 Control Test 0.049 17:42 Air Blank 0.000 17:43 Control Test 0.048 17:44 Air Blank 0.000 17:44 Control Test 0.048 17:45 Air Blank 0.000 17:45 Control Test Stats Average 0.0483 Std Dev 0.0006 Rel Std Dev(%) 1.1945 Operator's Signature <i>DELL</i>	ROCKLEDGE P.D. Intoxilyzer - Alcotest Analyzer Model 8000 SN 80-001261 12/07/2020 Software: 8100.27 Test g/210L Time Air Blank 0.000 17:47 Control Test 0.079 17:47 Air Blank 0.000 17:48 Control Test 0.078 17:49 Air Blank 0.000 17:49 Control Test 0.078 17:50 Air Blank 0.000 17:50 Control Test Stats Average 0.0783 Std Dev 0.0006 Rel Std Dev(%) 0.7370 Operator's Signature <i>DELL</i>	ROCKLEDGE P.D. Intoxilyzer - Alcotest Analyzer Model 8000 SN 80-001261 12/07/2020 Software: 8100.27 Test g/210L Time Air Blank 0.000 17:51 Control Test 0.200 17:52 Air Blank 0.000 17:53 Control Test 0.198 17:53 Air Blank 0.000 17:54 Control Test 0.197 17:55 Air Blank 0.000 17:55 Control Test Stats Average 0.1983 Std Dev 0.0015 Rel Std Dev(%) 0.7702 Operator's Signature <i>DELL</i>	ROCKLEDGE P.D. Intoxilyzer - Alcotest Analyzer Model 8000 SN 80-001261 12/07/2020 Software: 8100.27 Test g/210L Time Air Blank 0.000 17:57 Control Test 0.078 17:57 Air Blank 0.000 17:58 Control Test 0.078 17:58 Air Blank 0.000 17:59 Control Test 0.078 17:59 Air Blank 0.000 17:59 Control Test Stats Average 0.0780 Std Dev 0.0000 Rel Std Dev(%) 0.0000 Operator's Signature <i>DELL</i>

ROCKLEDGE P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001261
12/07/2020
Software: 8100.27

Flow Rate Calibration*****
1: Rate (Liters/min) = 5
SQRT(Diff)) = 7.871
2: Rate (Liters/min) = 15
SQRT(Diff)) = 12.000
3: Rate (Liters/min) = 30
SQRT(Diff)) = 21.070
Dependent Data Scale Factor = 100000 L/min
Independent Data Scale Factor = 256
Rounded Slope = 724
Rounded Intercept = -864435
Correlation = 0.99521

Flow Cal 1
Printed on external
Printer
WLL

Flow Calibration	
SN:	80-001261
Agency:	Rockledge PD
Date:	12/7/2020
By:	DERR <i>WLL</i>

ROCKLEDGE P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001261
12/07/2020
Software: 8100.27

Flow Rate Calibration*****
1: Rate (Liters/min) = 5
SQRT(Diff)) = 7.141
2: Rate (Liters/min) = 15
SQRT(Diff)) = 11.957
3: Rate (Liters/min) = 30
SQRT(Diff)) = 21.141
Dependent Data Scale Factor = 100000 L/min
Independent Data Scale Factor = 256
Rounded Slope = 690
Rounded Intercept = -701465
Correlation = 0.99800

Flow Cal 2
WLL

2020.12.
16
01:07:12
-05'00"

WLL

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: ROCKLEDGE P.D.

Time of Inspection: 16:54

Date of Inspection: 12/07/2020

Serial Number: 80-001261

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

COMPLIANCE NOT DETERMINED, AI NOT CONDUCTED.

2020.12.
16
01:06:45
-05'00"

MX

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

David Reyes-Rivera

DAVID E REYES-RIVERA

Signature and Printed Name

12/07/2020
Date



INSTRUMENT PROCESSING SHEET

Florida Department of
Law Enforcement

Agency Rockledge PD


S/N 80-001261

Date In 3/3/2020

DI Completion Date 3/11/20

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake Performed By <u>DP</u> ☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ 	Quality Checks Performed By <u>SP</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>191</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP102</u> 32 mm <u>.144</u> (.139 - .169) 36 mm <u>.160</u> (.156 - .190) 53 mm <u>.230</u> (.228 - .278) 103 mm <u>.492</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>30793</u> ☒ Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td><u>SD1021</u></td> <td><u>201905A</u> <u>5-14-21</u></td> </tr> <tr> <td>0.080</td> <td><u>SD3962</u></td> <td><u>201905B</u> <u>5-14-21</u></td> </tr> <tr> <td>0.200</td> <td><u>SD3933</u></td> <td><u>201904D</u> <u>4-30-21</u></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td><u>AG931603</u> <u>11-12-21</u></td> </tr> </tbody> </table>	Simulator	Serial #	Lot #/Exp	0.050	<u>SD1021</u>	<u>201905A</u> <u>5-14-21</u>	0.080	<u>SD3962</u>	<u>201905B</u> <u>5-14-21</u>	0.200	<u>SD3933</u>	<u>201904D</u> <u>4-30-21</u>	0.080 DGS	N/A	<u>AG931603</u> <u>11-12-21</u>	Flow Calibration Performed By _____ Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance Performed By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ Temperature Checks Performed By <u>SP</u> ☒ Lab Temp °C <u>23.2</u> External Digital Therm. ID#: <u>300502</u> ☒ 34°C +/- .2 Serial #: <u>SD1021</u> ☒ 34°C +/- .2 Serial #: <u>SD3962</u> ☒ 34°C +/- .2 Serial #: <u>SD3933</u>																																	
Simulator	Serial #	Lot #/Exp																																																
0.050	<u>SD1021</u>	<u>201905A</u> <u>5-14-21</u>																																																
0.080	<u>SD3962</u>	<u>201905B</u> <u>5-14-21</u>																																																
0.200	<u>SD3933</u>	<u>201904D</u> <u>4-30-21</u>																																																
0.080 DGS	N/A	<u>AG931603</u> <u>11-12-21</u>																																																
Final Release Date FDLE MAR 16 2020 Alcohol Testing Program	Calibration Adjustment Performed By _____ Barometric Pressure Gauge ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>		Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A		
Simulator	Serial Number	Lot Number	Expiration																																															
0.000		N/A	N/A																																															
0.040																																																		
0.100																																																		
0.200																																																		
0.300																																																		
0.080 DGS	N/A																																																	
Simulator	Serial Number	Lot Number	Expiration																																															
0.050																																																		
0.080																																																		
0.200																																																		
0.080 DGS	N/A																																																	
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ 		Department Inspection Performed By <u>SP</u> Barometric Pressure ID# <u>30793</u> Gauge <u>1021</u> Instrument <u>1020</u> Mouth Alcohol Solution Lot # <u>2019-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td><u>G4444</u></td> </tr> <tr> <td>Interferent</td> <td><u>G16621</u></td> </tr> <tr> <td>0.050</td> <td><u>SD1021</u></td> </tr> <tr> <td>0.080</td> <td><u>SD3962</u></td> </tr> <tr> <td>0.200</td> <td><u>SD3933</u></td> </tr> </tbody> </table> Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☐ Other _____ ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use <u>ADPM 3/16/20</u> Tech Review / Date <u>Brett K. Kuhlman 3/30/2020</u> Admin Review / Date	Simulator	Serial Number	0.000	<u>G4444</u>	Interferent	<u>G16621</u>	0.050	<u>SD1021</u>	0.080	<u>SD3962</u>	0.200	<u>SD3933</u>																																				
Simulator	Serial Number																																																	
0.000	<u>G4444</u>																																																	
Interferent	<u>G16621</u>																																																	
0.050	<u>SD1021</u>																																																	
0.080	<u>SD3962</u>																																																	
0.200	<u>SD3933</u>																																																	

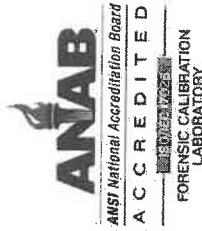

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2729 Fort Knox Blvd.
Bldg. 2, Suite 1300
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001261, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001261</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>ROCKLEDGE P.D.</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>03/11/2020</u>	0.080 g/ 210 L 0.005
Calibration Time:	<u>13:23</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

03/11/2020

Date

Shayla Platt

SHAYLA D PLATT,
Department Inspector

FDLE/ATP Form 69 January 2020

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1

CPM TSK 3/16/2020

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: ROCKLEDGE P.D.

Time of Inspection: 13:23

Date of Inspection: 03/11/2020

Serial Number: 80-001261

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG931603 Exp: 11/12/2021
0.000	0.047	0.078	0.200	0.077
0.000	0.048	0.079	0.201	0.077
0.000	0.048	0.078	0.202	0.077
0.000	0.048	0.079	0.201	0.077
0.000	0.048	0.079	0.202	0.077
0.000	0.048	0.078	0.201	0.077
0.000	0.048	0.078	0.202	0.077
0.000	0.048	0.078	0.202	0.077
0.000	0.049	0.078	0.201	0.078
0.000	0.049	0.078	0.202	0.077

Standard Deviations	0.0005	0.0004	0.0006	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

BK
3/16/2020

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT
Signature and Printed Name

03/11/2020
Date

Stability Checks #80-001261

ROCKLEDGE P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001261
03/11/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:21
Control Test	0.048	10:21
Air Blank	0.000	10:22
Control Test	0.047	10:23
Air Blank	0.000	10:23
Control Test	0.047	10:24
Air Blank	0.000	10:25
Control Test Stats		
Average	0.0473	
Std Dev	0.0006	
Rel Std Dev(%)	1.2198	

SP

Operator's Signature

ROCKLEDGE P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001261
03/11/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:16
Control Test	0.078	10:17
Air Blank	0.000	10:18
Control Test	0.078	10:18
Air Blank	0.000	10:19
Control Test	0.078	10:19
Air Blank	0.000	10:20
Control Test Stats		
Average	0.0780	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

SP

Operator's Signature

ROCKLEDGE P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001261
03/11/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:28
Control Test	0.200	10:29
Air Blank	0.000	10:29
Control Test	0.199	10:30
Air Blank	0.000	10:30
Control Test	0.200	10:31
Air Blank	0.000	10:32
Control Test Stats		
Average	0.1997	
Std Dev	0.0006	
Rel Std Dev(%)	0.2892	

SP

Operator's Signature

ROCKLEDGE P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001261
03/11/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:35
Control Test	0.078	10:35
Air Blank	0.000	10:36
Control Test	0.077	10:36
Air Blank	0.000	10:37
Control Test	0.077	10:37
Air Blank	0.000	10:37
Control Test Stats		
Average	0.0773	
Std Dev	0.0006	
Rel Std Dev(%)	0.7466	

SP

Operator's Signature

QAM YBK
3/16/2020