

# INSTRUMENT PROCESSING SHEET

Agency Alachua Police DepartmentS/N 80-001720

Florida Department of Law Enforcement

Date In 04/07/2020 DI Completion Date 4/15/20 Ship  P/U  H/D  CMI  EE

<b>Intake</b> Performed By <u>RAW</u> <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input checked="" type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input checked="" type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ _____ _____	<b>Quality Checks</b> Performed By <u>SP</u> <input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>187</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP-102</u> 32 mm <u>0.156</u> (.139 - .169) 36 mm <u>0.175</u> (.156 - .190) 53 mm <u>0.242</u> (.228 - .278) 103 mm <u>0.507</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>28427</u> <input checked="" type="checkbox"/> Stability Checks <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td><u>MP5088</u></td> <td><u>201905A</u> <u>05-14-2021</u></td> </tr> <tr> <td>0.080</td> <td><u>MP5089</u></td> <td><u>201905B</u> <u>05-19-2021</u></td> </tr> <tr> <td>0.200</td> <td><u>MP5090</u></td> <td><u>201904D</u> <u>04-30-2021</u></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td><u>AG931603</u> <u>11-12-2021</u></td> </tr> </tbody> </table>	Simulator	Serial #	Lot #/Exp	0.050	<u>MP5088</u>	<u>201905A</u> <u>05-14-2021</u>	0.080	<u>MP5089</u>	<u>201905B</u> <u>05-19-2021</u>	0.200	<u>MP5090</u>	<u>201904D</u> <u>04-30-2021</u>	0.080 DGS	N/A	<u>AG931603</u> <u>11-12-2021</u>	<b>Flow Calibration</b> Performed By _____ Flow Column # _____ <input type="checkbox"/> 5L/min - 17mm <input type="checkbox"/> 15L/min - 53mm <input type="checkbox"/> 30L/min - 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) <b>Maintenance</b> Performed By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____ <b>Temperature Checks</b> Performed By <u>SP</u> <input checked="" type="checkbox"/> Lab Temp °C <u>21.3</u> External Digital Therm. ID#: <u>300502</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5088</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5089</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5090</u>																																												
Simulator	Serial #	Lot #/Exp																																																											
0.050	<u>MP5088</u>	<u>201905A</u> <u>05-14-2021</u>																																																											
0.080	<u>MP5089</u>	<u>201905B</u> <u>05-19-2021</u>																																																											
0.200	<u>MP5090</u>	<u>201904D</u> <u>04-30-2021</u>																																																											
0.080 DGS	N/A	<u>AG931603</u> <u>11-12-2021</u>																																																											
<b>Calibration Adjustment</b> Performed By _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse; margin-bottom: 5px;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> <input type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>	Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			<b>Department Inspection</b> Performed By <u>SP</u> Barometric Pressure ID# <u>28421</u> Gauge <u>1014</u> Instrument <u>1012</u> Mouth Alcohol Solution Lot # <u>2019-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td><u>MP5086</u></td> </tr> <tr> <td>Interferent</td> <td><u>MP5087</u></td> </tr> <tr> <td>0.050</td> <td><u>MP5088</u></td> </tr> <tr> <td>0.080</td> <td><u>MP5089</u></td> </tr> <tr> <td>0.200</td> <td><u>MP5090</u></td> </tr> </tbody> </table> <b>Attachments</b> <input checked="" type="checkbox"/> Form 41 <input checked="" type="checkbox"/> Stability Checks <input checked="" type="checkbox"/> Calibration Certificate <input type="checkbox"/> Calibration Adjustment <input type="checkbox"/> Post-Stability Checks <input type="checkbox"/> Flow Calibration <input type="checkbox"/> Form 40 <input type="checkbox"/> Other _____	Simulator	Serial Number	0.000	<u>MP5086</u>	Interferent	<u>MP5087</u>	0.050	<u>MP5088</u>	0.080	<u>MP5089</u>	0.200	<u>MP5090</u>
Simulator	Serial Number	Lot Number	Expiration																																																										
0.000		N/A	N/A																																																										
0.040																																																													
0.100																																																													
0.200																																																													
0.300																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number	Lot Number	Expiration																																																										
0.050																																																													
0.080																																																													
0.200																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number																																																												
0.000	<u>MP5086</u>																																																												
Interferent	<u>MP5087</u>																																																												
0.050	<u>MP5088</u>																																																												
0.080	<u>MP5089</u>																																																												
0.200	<u>MP5090</u>																																																												
Notes/Suggested Service: _____ _____ _____ _____ _____	<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use <hr/> Tech Review / Date _____ Admin Review / Date _____																																																												



# Calibration Certificate

Florida Department of Law Enforcement  
Alcohol Testing Program  
2729 Fort Knox Blvd.  
Bldg. 2, Suite 1300  
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001720 , manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001720</u>	UNCERTAINTY * ±	
Owning Agency:	<u>ALACHUA PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>04/15/2020</u>	0.080 g/ 210 L	0.005
Calibration Time:	<u>12:24</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.  
Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.  
\*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

### TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

*Shayla Platt*

04/15/2020

Date

SHAYLA D PLATT,  
Department Inspector

FDLE/ATP Form 69 January 2020

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

# Florida Department of Law Enforcement Alcohol Testing Program

## DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: ALACHUA PD  
Time of Inspection: 12:24

Date of Inspection: 04/15/2020

Serial Number: 80-001720  
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG931603 Exp: 11/12/2021
0.000	0.049	0.079	0.201	0.078
0.000	0.049	0.079	0.201	0.078
0.000	0.049	0.079	0.201	0.078
0.000	0.049	0.079	0.201	0.079
0.000	0.049	0.080	0.200	0.078
0.000	0.049	0.080	0.200	0.078
0.000	0.049	0.080	0.200	0.078
0.000	0.049	0.080	0.200	0.078
0.000	0.049	0.080	0.200	0.078
0.000	0.049	0.080	0.200	0.078
0.000	0.049	0.080	0.200	0.078


Standard Deviations	0.0000	0.0005	0.0005	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (  ) does not comply (  ) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

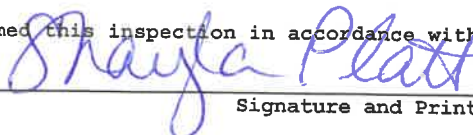

SHAYLA D PLATT

Signature and Printed Name

04/15/2020  
Date

# stability checks

ALACHUA PD  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-001720  
04/14/2020  
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:52
Control Test	0.049	07:53
Air Blank	0.000	07:54
Control Test	0.049	07:54
Air Blank	0.000	07:55
Control Test	0.049	07:56
Air Blank	0.000	07:56
Control Test Stats		
Average	0.0490	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

ALACHUA PD  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-001720  
04/14/2020  
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:58
Control Test	0.080	07:58
Air Blank	0.000	07:59
Control Test	0.080	08:00
Air Blank	0.000	08:00
Control Test	0.080	08:01
Air Blank	0.000	08:01
Control Test Stats		
Average	0.0800	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

ALACHUA PD  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-001720  
04/14/2020  
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:03
Control Test	0.202	08:03
Air Blank	0.000	08:04
Control Test	0.202	08:05
Air Blank	0.000	08:05
Control Test	0.202	08:06
Air Blank	0.000	08:06
Control Test Stats		
Average	0.2020	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

wet



Operator's Signature



Operator's Signature



Operator's Signature

ALACHUA PD  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-001720  
04/14/2020  
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:19
Control Test	0.078	08:19
Air Blank	0.000	08:19
Control Test	0.078	08:20
Air Blank	0.000	08:20
Control Test	0.078	08:21
Air Blank	0.000	08:21
Control Test Stats		
Average	0.0780	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Dry



Operator's Signature

**Return Material Authorization**

**Ship to:**  CMI, Inc.  
 Enforcement Electronics

Shipment to repair facility authorized by: J. Bennett on 10/21/19

**Items Returned:** Instrument  Supplies  Other  Describe: \_\_\_\_\_  
Instrument Model: Intoxilyzer 8000 Serial Number: 80-001720

<b>Bill To Address:</b> <u>Alachua PD</u> _____ _____ _____ _____	<b>Ship to Address:</b> <u>FDLE ATP Tallahassee</u> _____ _____ _____ _____
--	--

**Reason for Return:**  
No dial tone, unable to upload instrument via modem.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

**Please choose one of the following options:**

1. I \_\_\_\_\_, authorize all repairs.

2. I \_\_\_\_\_, authorize repairs up to \$ \_\_\_\_\_.

3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: J Bennett  
Phone #: (386) 462-1396 Email: jbennett@apd.cityofalachua.org  
ATP Contact Name: Shayla Platt ATP Email: shaylaplatt@fdle.state.fl.us





# INSTRUMENT PROCESSING SHEET

Agency Alachua PD S/N 80-001720

Florida Department of Law Enforcement

Date In 10/8/2019 DI Completion Date 10/17/19  Ship  P/U  H/D  CMI  EE

<b>Intake</b> Performed By <u>DP</u> <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input checked="" type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ _____ _____	<b>Quality Checks</b> Performed By <u>SP</u> <input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>180</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP103</u> 32 mm <u>.156</u> (.139 - .169) 36 mm <u>.167</u> (.156 - .190) 53 mm <u>.242</u> (.228 - .278) 103 mm <u>.503</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>26932</u> <input checked="" type="checkbox"/> Stability Checks	<b>Flow Calibration</b> Performed By _____ Flow Column # _____ <input type="checkbox"/> 5L/min - 17mm <input type="checkbox"/> 15L/min - 53mm <input type="checkbox"/> 30L/min - 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																												
<b>Final Release Date</b> _____ _____	<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>SD1012</td> <td>201905A 5-14-21</td> </tr> <tr> <td>0.080</td> <td>DR1279</td> <td>201905B 5-14-21</td> </tr> <tr> <td>0.200</td> <td>SD1013</td> <td>201904D 4-30-21</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG916501 6-14-21</td> </tr> </tbody> </table>	Simulator	Serial #	Lot #/Exp	0.050	SD1012	201905A 5-14-21	0.080	DR1279	201905B 5-14-21	0.200	SD1013	201904D 4-30-21	0.080 DGS	N/A	AG916501 6-14-21	<b>Maintenance</b> Performed By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____																																													
Simulator	Serial #	Lot #/Exp																																																												
0.050	SD1012	201905A 5-14-21																																																												
0.080	DR1279	201905B 5-14-21																																																												
0.200	SD1013	201904D 4-30-21																																																												
0.080 DGS	N/A	AG916501 6-14-21																																																												
<b>Calibration Adjustment</b> Performed By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> <input type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>	Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			<b>Department Inspection</b> Performed By <u>SP</u> Barometric Pressure ID# <u>26932</u> Gauge <u>1012</u> Instrument <u>1012</u> Mouth Alcohol Solution Lot # <u>2018-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>G2408</td> </tr> <tr> <td>Interferent</td> <td>G2882</td> </tr> <tr> <td>0.050</td> <td>SD1012</td> </tr> <tr> <td>0.080</td> <td>DR1279</td> </tr> <tr> <td>0.200</td> <td>SD1013</td> </tr> </tbody> </table>	Simulator	Serial Number	0.000	G2408	Interferent	G2882	0.050	SD1012	0.080	DR1279	0.200	SD1013	<b>Attachments</b> <input checked="" type="checkbox"/> Form 41 <input checked="" type="checkbox"/> Stability Checks <input checked="" type="checkbox"/> Calibration Certificate <input type="checkbox"/> Calibration Adjustment <input type="checkbox"/> Post-Stability Checks <input type="checkbox"/> Flow Calibration <input checked="" type="checkbox"/> Form 40 <input type="checkbox"/> Other _____
Simulator	Serial Number	Lot Number	Expiration																																																											
0.000		N/A	N/A																																																											
0.040																																																														
0.100																																																														
0.200																																																														
0.300																																																														
0.080 DGS	N/A																																																													
Simulator	Serial Number	Lot Number	Expiration																																																											
0.050																																																														
0.080																																																														
0.200																																																														
0.080 DGS	N/A																																																													
Simulator	Serial Number																																																													
0.000	G2408																																																													
Interferent	G2882																																																													
0.050	SD1012																																																													
0.080	DR1279																																																													
0.200	SD1013																																																													
Notes/Suggested Service: <u>Sending to repair for modern check. Unable to obtain dial tone for upload. Instrument uploaded via direct connect. SP</u> _____ _____	<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input checked="" type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use  <u>DP</u> 10/17/19 <u>Paul Kirkland</u> 10/23/19 Tech Review / Date      Admin Review / Date																																																													



# Calibration Certificate

Florida Department of Law Enforcement  
Alcohol Testing Program  
2729 Fort Knox Blvd.  
Bldg. 2, Suite 1300  
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001720, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001720</u>	UNCERTAINTY* ±	
Owning Agency:	<u>ALACHUA PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>10/17/2019</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>11:47</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

\*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

### TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

*Shayla Platt*

10/17/2019

Date

SHAYLA D PLATT,  
Department Inspector

FDLE/ATP Form 69 July 2018

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

*Good BK 10/23/19*

# Florida Department of Law Enforcement Alcohol Testing Program

## DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: ALACHUA PD

Time of Inspection: 11:47

Date of Inspection: 10/17/2019

Serial Number: 80-001720

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG916501 Exp: 06/14/2021
0.000	0.048	0.077	0.197	0.080
0.000	0.047	0.077	0.198	0.080
0.000	0.048	0.077	0.198	0.080
0.000	0.048	0.077	0.198	0.080
0.000	0.048	0.077	0.198	0.081
0.000	0.048	0.077	0.198	0.080
0.000	0.048	0.077	0.199	0.080
0.000	0.048	0.077	0.199	0.080
0.000	0.048	0.077	0.199	0.081
0.000	0.048	0.078	0.199	0.081

Standard Deviations	0.0003	0.0003	0.0006	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

PDM  
BK  
10/23/19

The above instrument complies (  ) does not comply (  ) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT

Signature and Printed Name

10/17/2019  
Date



ALACHUA PD  
 Intoxilyzer - Alcohol Analyzer  
 Model 8000 SN 80-001720  
 10/17/2019  
 Software: 8100.27

Stability  
 Checks  
 #80-001720

Test	g/210L	Time
Air Blank	0.000	09:26
Control Test	0.077	09:27
Air Blank	0.000	09:28
Control Test	0.077	09:28
Air Blank	0.000	09:29
Control Test	0.077	09:30
Air Blank	0.000	09:30
Control Test Stats		
Average	0.0770	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

SP

Operator's Signature

ALACHUA PD  
 Intoxilyzer - Alcohol Analyzer  
 Model 8000 SN 80-001720  
 10/17/2019  
 Software: 8100.27

ALACHUA PD  
 Intoxilyzer - Alcohol Analyzer  
 Model 8000 SN 80-001720  
 10/17/2019  
 Software: 8100.27

ALACHUA PD  
 Intoxilyzer - Alcohol Analyzer  
 Model 8000 SN 80-001720  
 10/17/2019  
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:31
Control Test	0.198	09:32
Air Blank	0.000	09:33
Control Test	0.199	09:33
Air Blank	0.000	09:34
Control Test	0.199	09:34
Air Blank	0.000	09:35
Control Test Stats		
Average	0.1987	
Std Dev	0.0006	
Rel Std Dev(%)	0.2906	

Test	g/210L	Time
Air Blank	0.000	09:43
Control Test	0.048	09:44
Air Blank	0.000	09:44
Control Test	0.048	09:45
Air Blank	0.000	09:45
Control Test	0.048	09:46
Air Blank	0.000	09:47
Control Test Stats		
Average	0.0480	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Test	g/210L	Time
Air Blank	0.000	09:36
Control Test	0.080	09:36
Air Blank	0.000	09:37
Control Test	0.080	09:37
Air Blank	0.000	09:37
Control Test	0.081	09:38
Air Blank	0.000	09:38
Control Test Stats		
Average	0.0803	
Std Dev	0.0006	
Rel Std Dev(%)	0.7187	

DGS

SP

Operator's Signature

SP

Operator's Signature

SP

Operator's Signature

PBam

BSK

10/23/19

# Florida Department of Law Enforcement Alcohol Testing Program

## AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: ALACHUA PD

Time of Inspection: 09:24

Date of Inspection: 10/17/2019

Serial Number: 80-001720

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: _____ Exp: _____	0.08g/210L Test (g/210L) Lot#: _____ Exp: _____	0.20g/210L Test (g/210L) Lot#: _____ Exp: _____	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: _____ Exp: _____

Number of Simulators Used: \_\_\_\_\_

Remarks:

BYPASSED AI TO OPERATE INSTRUMENT

Pg 1  
BK  
10/23/19

The above instrument complies (  ) does not comply (  ) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT

Signature and Printed Name

10/17/2019  
Date

# Alachua Police



# Department

---

ADAM B. BOUKARI  
CITY MANAGER

CHAD D. SCOTT  
CHIEF OF POLICE

## INTER-OFFICE COMMUNICATION

**DATE:** 10/07/2019

**TO:** BRETT KIRKLAND, PH.D. PROGRAM MANAGER

**FROM:** OFC. J. BENNETT 311

**RE:** Intox 8000 Repair and Annual Inspection

---

Greetings,

Please take a look at the phone modem. During the last inspection, the machine would not connect and send the inspection report. The building maintenance and the IT department inspected the phone outlets and were unable to find anything wrong.

The machine is also due for the annual inspection.

Thank you,

Ofc. J. Bennett 311

JB/jb

*efom*