




INSTRUMENT PROCESSING SHEET

Agency Osceola County Sheriff's Office S/N 80-003935

Florida Department of Law Enforcement

Date In 02/10/2020 DI Completion Date 2/21/20 Ship P/U H/D CMI EE

Intake Performed By <u>RAW</u> <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input checked="" type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ _____ _____	Quality Checks Performed By <u>SP</u> <input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>211</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP-103</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.226</u> (.228 - .278) 103 mm <u>0.406</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>26932</u> <input checked="" type="checkbox"/> Stability Checks	Flow Calibration Performed By <u>SP</u> Flow Column # <u>ATP105</u> <input checked="" type="checkbox"/> 5L/min - 17mm <input checked="" type="checkbox"/> 15L/min - 53mm <u>x3</u> <input checked="" type="checkbox"/> 30L/min - 103mm <input checked="" type="checkbox"/> R-Value <u>213</u> <input checked="" type="checkbox"/> Post Calibration Verification (L/s) Flow Column # <u>ATP103</u> 32 mm <u>.147</u> (.139 - .169) 36 mm <u>.187</u> (.156 - .190) 53 mm <u>.257</u> (.228 - .278) 103 mm <u>.472</u> (.447 - .547)																																																											
Final Release Date <p style="text-align: center; font-weight: bold;">FDLE</p> <p style="text-align: center; font-weight: bold;">FEB 25 2020</p> <p style="text-align: center; font-weight: bold;">Alcohol Testing Program</p>	<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td><u>SD1012</u></td> <td><u>201905A</u> <u>05-14-2021</u></td> </tr> <tr> <td>0.080</td> <td><u>DR1279</u></td> <td><u>201905B</u> <u>05-14-2021</u></td> </tr> <tr> <td>0.200</td> <td><u>SD1011</u></td> <td><u>201904D</u> <u>04-30-2021</u></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td><u>AG916501</u> <u>06-14-2021</u></td> </tr> </tbody> </table>	Simulator	Serial #	Lot #/Exp	0.050	<u>SD1012</u>	<u>201905A</u> <u>05-14-2021</u>	0.080	<u>DR1279</u>	<u>201905B</u> <u>05-14-2021</u>	0.200	<u>SD1011</u>	<u>201904D</u> <u>04-30-2021</u>	0.080 DGS	N/A	<u>AG916501</u> <u>06-14-2021</u>	Maintenance Performed By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____ Temperature Checks Performed By <u>SP</u> <input checked="" type="checkbox"/> Lab Temp °C <u>21.9</u> External Digital Therm. ID#: <u>300505</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5088</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5089</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5090</u>																																												
Simulator	Serial #	Lot #/Exp																																																											
0.050	<u>SD1012</u>	<u>201905A</u> <u>05-14-2021</u>																																																											
0.080	<u>DR1279</u>	<u>201905B</u> <u>05-14-2021</u>																																																											
0.200	<u>SD1011</u>	<u>201904D</u> <u>04-30-2021</u>																																																											
0.080 DGS	N/A	<u>AG916501</u> <u>06-14-2021</u>																																																											
Calibration Adjustment Performed By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> <input type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>	Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection Performed By <u>SP</u> Barometric Pressure ID# <u>30793</u> Gauge <u>1026</u> Instrument <u>1026</u> Mouth Alcohol Solution Lot # <u>2019-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td><u>MP5086</u></td> </tr> <tr> <td>Interferent</td> <td><u>MP5087</u></td> </tr> <tr> <td>0.050</td> <td><u>MP5088</u></td> </tr> <tr> <td>0.080</td> <td><u>MP5089</u></td> </tr> <tr> <td>0.200</td> <td><u>MP5090</u></td> </tr> </tbody> </table>	Simulator	Serial Number	0.000	<u>MP5086</u>	Interferent	<u>MP5087</u>	0.050	<u>MP5088</u>	0.080	<u>MP5089</u>	0.200	<u>MP5090</u>
Simulator	Serial Number	Lot Number	Expiration																																																										
0.000		N/A	N/A																																																										
0.040																																																													
0.100																																																													
0.200																																																													
0.300																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number	Lot Number	Expiration																																																										
0.050																																																													
0.080																																																													
0.200																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number																																																												
0.000	<u>MP5086</u>																																																												
Interferent	<u>MP5087</u>																																																												
0.050	<u>MP5088</u>																																																												
0.080	<u>MP5089</u>																																																												
0.200	<u>MP5090</u>																																																												
Notes/Suggested Service: _____ _____ _____ _____ _____	Attachments <input checked="" type="checkbox"/> Form 41 <input checked="" type="checkbox"/> Stability Checks <input checked="" type="checkbox"/> Calibration Certificate <input type="checkbox"/> Calibration Adjustment <input type="checkbox"/> Post-Stability Checks <input checked="" type="checkbox"/> Flow Calibration <u>x3</u> <input type="checkbox"/> Form 40 <input type="checkbox"/> Other _____																																																												
<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use <p style="text-align: center;"> <u>PPM 2/24/20</u> <u>Brett Kirkland 2/24/2020</u> Tech Review / Date Admin Review / Date </p>																																																													

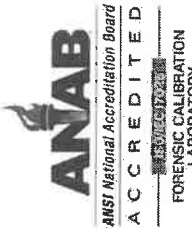

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2729 Fort Knox Blvd.
Bldg. 2, Suite 1300
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number **80-003935**, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-003935</u>	UNCERTAINTY* ±	
Owning Agency:	<u>OSCEOLA COUNTY S.O.</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>02/21/2020</u>	0.080 g/ 210 L	0.005
Calibration Time:	<u>11:03</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.
Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

02/21/2020

Date

Shayla Platt
SHAYLAD PLATT,
Department Inspector

FDLE/ATP Form 69 January 2020

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Handwritten: [Signature] 13K 2/24/2020

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: OSCEOLA COUNTY S.O.
Time of Inspection: 11:03

Date of Inspection: 02/21/2020

Serial Number: 80-003935
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG931603 Exp: 11/12/2021
0.000	0.049	0.081	0.201	0.080
0.000	0.049	0.081	0.202	0.080
0.000	0.050	0.080	0.202	0.080
0.000	0.050	0.080	0.201	0.079
0.000	0.049	0.081	0.201	0.080
0.000	0.049	0.081	0.201	0.080
0.000	0.049	0.081	0.201	0.080
0.000	0.049	0.081	0.201	0.080
0.000	0.049	0.082	0.202	0.079
0.000	0.049	0.082	0.202	0.080

Standard Deviations	0.0004	0.0006	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

RJM
TSK
2/24/2020

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT

Signature and Printed Name

02/21/2020
Date

Stability Checks

OSCEOLA COUNTY S.O.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003935
02/10/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:11
Control Test	0.048	14:12
Air Blank	0.000	14:12
Control Test	0.048	14:13
Air Blank	0.000	14:14
Control Test	0.048	14:14
Air Blank	0.000	14:15
Control Test Stats		
Average	0.0460	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator's Signature

OSCEOLA COUNTY S.O.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003935
02/10/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:16
Control Test	0.078	14:17
Air Blank	0.000	14:17
Control Test	0.078	14:18
Air Blank	0.000	14:18
Control Test	0.079	14:19
Air Blank	0.000	14:20
Control Test Stats		
Average	0.0783	
Std Dev	0.0006	
Rel Std Dev(%)	0.7370	

Wet


Operator's Signature

OSCEOLA COUNTY S.O.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003935
02/10/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:22
Control Test	0.199	14:23
Air Blank	0.000	14:23
Control Test	0.200	14:24
Air Blank	0.000	14:24
Control Test	0.198	14:25
Air Blank	0.000	14:26
Control Test Stats		
Average	0.1990	
Std Dev	0.0010	
Rel Std Dev(%)	0.5025	


Operator's Signature

OSCEOLA COUNTY S.O.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003935
02/10/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:27
Control Test	0.080	14:27
Air Blank	0.000	14:28
Control Test	0.080	14:28
Air Blank	0.000	14:29
Control Test	0.081	14:29
Air Blank	0.000	14:29
Control Test Stats		
Average	0.0803	
Std Dev	0.0006	
Rel Std Dev(%)	0.7187	

Dry


Operator's Signature

Open
TSK
2/24/2020

OSCEOLA COUNTY S.O.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003935
02/21/2020
Software: 8100.27

Flow Rate Calibration*****
1: Rate (Liters/min) = 5
SQRT(Diff) = 0.000
2: Rate (Liters/min) = 15
SQRT(Diff) = 0.000
3: Rate (Liters/min) = 30
SQRT(Diff) = 0.000
Dependent Data Scale Factor = 10000 L/min
Independent Data Scale Factor = 256
Rounded Slope = 0
Rounded Intercept = 0
Correlation = 0.00000

OSCEOLA COUNTY S.O.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003935
02/21/2020
Software: 8100.27

Flow Rate Calibration*****
1: Rate (Liters/min) = 5
SQRT(Diff) = 0.000
2: Rate (Liters/min) = 15
SQRT(Diff) = 1.000
3: Rate (Liters/min) = 30
SQRT(Diff) = 0.000
Dependent Data Scale Factor = 10000 L/min
Independent Data Scale Factor = 256
Rounded Slope = -976
Rounded Intercept = 1750000
Correlation = -0.11471

OSCEOLA COUNTY S.O.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003935
02/21/2020
Software: 8100.27

Flow Rate Calibration*****
1: Rate (Liters/min) = 5
SQRT(Diff) = 0.000
2: Rate (Liters/min) = 15
SQRT(Diff) = 9.590
3: Rate (Liters/min) = 30
SQRT(Diff) = 20.344
Dependent Data Scale Factor = 10000 L/min
Independent Data Scale Factor = 256
Rounded Slope = 481
Rounded Intercept = 437173
Correlation = 0.99665

FLOW
CAL
ADJUST
SP

Repeated flow cal adjustment
after noticing the flow rate
advanced prior to introducing
appropriate flow rate to
instrument during
adjustment. SP

PPGM
TK
2/24/2020

Return Material Authorization

Ship to: CMI, Inc.
 Enforcement Electronics

Shipment to repair facility authorized by: Dan Lyons on 10/22/19

Items Returned: Instrument Supplies Other Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: ~~80-001720~~ 80-003935 *SP*

Bill To Address:
Osceola County SO

Ship to Address:
FDLE ATP Tallahassee

Reason for Return:
Modem port pushed into instrument.

Please choose one of the following options:

- 1. I _____, authorize all repairs.
- 2. I _____, authorize repairs up to \$_____.
- 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.
Please contact: Name: Dan Lyons
Phone #: 407-269-1980 Email: dlyo@OSCEOLA.ORG
ATP Contact Name: Shayla Platt ATP Email: shaylplatt@fdle.state.fl.us

*ogm
FSK
2/24/2020*



INSTRUMENT PROCESSING SHEET

Agency Osceola CountyS/N 80-003935Florida Department of
Law EnforcementDate In 10/8/2019 DI Completion Date _____ Ship P/U H/D CMI EE

Intake Performed By <u>DP</u>		Quality Checks Performed By <u>SP</u>		Flow Calibration Performed By _____																
<input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ _____ _____		<input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>212</u> <input type="checkbox"/> Flow Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>210932</u> <input type="checkbox"/> Stability Checks		Flow Column # _____ <input type="checkbox"/> 5L/min - 17mm <input type="checkbox"/> 15L/min - 53mm <input type="checkbox"/> 30L/min - 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																
Final Release Date		<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr><td>0.050</td><td></td><td></td></tr> <tr><td>0.080</td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td></tr> </tbody> </table>		Simulator	Serial #	Lot #/Exp	0.050			0.080			0.200			0.080 DGS	N/A		Maintenance Performed By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____	
Simulator	Serial #	Lot #/Exp																		
0.050																				
0.080																				
0.200																				
0.080 DGS	N/A																			
				Temperature Checks Performed By _____ <input type="checkbox"/> Lab Temp °C _____ External Digital Therm. ID#: _____ <input type="checkbox"/> 34°C +/- .2 Serial #: _____ <input type="checkbox"/> 34°C +/- .2 Serial #: _____ <input type="checkbox"/> 34°C +/- .2 Serial #: _____																

Calibration Adjustment Performed By _____				Department Inspection Performed By _____																																			
Barometric Pressure Gauge _____ ID # _____				Barometric Pressure ID# _____																																			
<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr> <tr><td>0.040</td><td></td><td></td><td></td></tr> <tr><td>0.100</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.300</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table>				Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge _____ Instrument _____ Mouth Alcohol Solution Lot # _____ Acetone Stock Solution Lot # _____							
Simulator	Serial Number	Lot Number	Expiration																																				
0.000		N/A	N/A																																				
0.040																																							
0.100																																							
0.200																																							
0.300																																							
0.080 DGS	N/A																																						
<input type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr><td>0.050</td><td></td><td></td><td></td></tr> <tr><td>0.080</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table>				Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr><td>0.000</td><td></td></tr> <tr><td>Interferent</td><td></td></tr> <tr><td>0.050</td><td></td></tr> <tr><td>0.080</td><td></td></tr> <tr><td>0.200</td><td></td></tr> </tbody> </table>				Simulator	Serial Number	0.000		Interferent		0.050		0.080		0.200	
Simulator	Serial Number	Lot Number	Expiration																																				
0.050																																							
0.080																																							
0.200																																							
0.080 DGS	N/A																																						
Simulator	Serial Number																																						
0.000																																							
Interferent																																							
0.050																																							
0.080																																							
0.200																																							
Notes/Suggested Service: <u>Uploaded via direct connect</u> <u>Sending to repair due to broken</u> <u>modem port. sp</u>				<input type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use																																			
				Tech Review / Date _____ Admin Review / Date _____																																			

OSM
1/16
2/24/2020