



INSTRUMENT PROCESSING SHEET

Agency Highlands CSOS/N 80-007507Florida Department of
Law EnforcementDate In 06/09/2023 DI Completion Date 06/14/2023☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By TDG	Quality Checks	By TDG	Date <u>06/13/2023</u>	Flow Calibration	By	Date																											
<input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input checked="" type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		<input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>241</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.242</u> (.228 - .278) 103 mm <u>0.500</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>26932</u> <input checked="" type="checkbox"/> Stability Checks			Flow Column # _____ <input type="checkbox"/> 5L/min - 17mm <input type="checkbox"/> 15L/min - 53mm <input type="checkbox"/> 30L/min - 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																													
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5094</td><td>202201C 01/11/2024</td></tr><tr><td>0.080</td><td>MP5095</td><td>202201D 01/18/2024</td></tr><tr><td>0.200</td><td>MP5096</td><td>202201E 01/18/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG223802 08/26/2024</td></tr></tbody></table>		Simulator	Serial #	Lot #/Exp	0.050	MP5094	202201C 01/11/2024	0.080	MP5095	202201D 01/18/2024	0.200	MP5096	202201E 01/18/2024	0.080 DGS	N/A	AG223802 08/26/2024			Maintenance By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____ _____ _____ _____ _____ _____													
Simulator	Serial #	Lot #/Exp																																
0.050	MP5094	202201C 01/11/2024																																
0.080	MP5095	202201D 01/18/2024																																
0.200	MP5096	202201E 01/18/2024																																
0.080 DGS	N/A	AG223802 08/26/2024																																
Calibration Adjustment By _____		Department Inspection By TDG _____																																
Barometric Pressure Gauge ID # _____		Barometric Pressure ID# <u>26932</u>																																
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge <u>1014</u> Instrument <u>1011</u>				
Simulator	Serial #	Lot #	Expiration																															
0.000		N/A	N/A																															
0.040																																		
0.100																																		
0.200																																		
0.300																																		
0.080 DGS	N/A																																	
<input type="checkbox"/> Post Calibration Adjustment Stability Checks		Mouth Alcohol Solution Lot # <u>2021-D</u>																																
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Acetone Stock Solution Lot # <u>2022-B</u>												
Simulator	Serial #	Lot #	Expiration																															
0.050																																		
0.080																																		
0.200																																		
0.080 DGS	N/A																																	
		<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5092</td></tr><tr><td>Interferent</td><td>MP5093</td></tr><tr><td>0.050</td><td>MP5094</td></tr><tr><td>0.080</td><td>MP5095</td></tr><tr><td>0.200</td><td>MP5096</td></tr></tbody></table>		Simulator	Serial Number	0.000	MP5092	Interferent	MP5093	0.050	MP5094	0.080	MP5095	0.200	MP5096																			
Simulator	Serial Number																																	
0.000	MP5092																																	
Interferent	MP5093																																	
0.050	MP5094																																	
0.080	MP5095																																	
0.200	MP5096																																	
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		Attachments <input checked="" type="checkbox"/> Form 41 <input checked="" type="checkbox"/> Stability Checks <input checked="" type="checkbox"/> Calibration Certificate <input type="checkbox"/> Calibration Adjustment		<input type="checkbox"/> Post-Stability Checks <input type="checkbox"/> Flow Calibration <input type="checkbox"/> Form 40 <input type="checkbox"/> Other _____																														
		<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use																																
		Israel Soto <small>Digitally signed by Israel Soto Date: 2023.06.14 16:32:40 +04'00'</small>		Phil Nicodemo <small>Digitally signed by Phil Nicodemo Date: 2023.06.15 11:16:31 -04'00'</small>																														
		Tech Review / Date _____		Admin Review / Date _____																														

Type of Test	Serial Number	Agency	Date	Performed By
Stabilities	80-00 7507	Highlands CSU	06/13/2023	TDG MC

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
✓	✓	✓	✓
<p>HIGHLANDS COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-007507 06/13/2023 Software: 8100.27</p> <p>Test g/210L Time</p> <p>Air Blank 0.000 15:38</p> <p>Control Test 0.049 15:39</p> <p>Air Blank 0.000 15:39</p> <p>Control Test 0.049 15:40</p> <p>Air Blank 0.000 15:41</p> <p>Control Test 0.049 15:41</p> <p>Air Blank 0.000 15:42</p> <p>Control Test Stats</p> <p>Average 0.0490</p> <p>Std Dev 0.0000</p> <p>Rel Std Dev(%) 0.0000</p>	<p>HIGHLANDS COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-007507 06/13/2023 Software: 8100.27</p> <p>Test g/210L Time</p> <p>Air Blank 0.000 15:47</p> <p>Control Test 0.078 15:48</p> <p>Air Blank 0.000 15:48</p> <p>Control Test 0.078 15:49</p> <p>Air Blank 0.000 15:50</p> <p>Control Test 0.079 15:50</p> <p>Air Blank 0.000 15:51</p> <p>Control Test Stats</p> <p>Average 0.0783</p> <p>Std Dev 0.0006</p> <p>Rel Std Dev(%) 0.7370</p>	<p>HIGHLANDS COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-007507 06/13/2023 Software: 8100.27</p> <p>Test g/210L Time</p> <p>Air Blank 0.000 15:55</p> <p>Control Test 0.199 15:56</p> <p>Air Blank 0.000 15:56</p> <p>Control Test 0.199 15:57</p> <p>Air Blank 0.000 15:57</p> <p>Control Test 0.199 15:58</p> <p>Air Blank 0.000 15:59</p> <p>Control Test Stats</p> <p>Average 0.1990</p> <p>Std Dev 0.0000</p> <p>Rel Std Dev(%) 0.0000</p>	<p>HIGHLANDS COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-007507 06/13/2023 Software: 8100.27</p> <p>Test g/210L Time</p> <p>Air Blank 0.000 16:00</p> <p>Control Test 0.080 16:00</p> <p>Air Blank 0.000 16:01</p> <p>Control Test 0.080 16:01</p> <p>Air Blank 0.000 16:01</p> <p>Control Test 0.080 16:02</p> <p>Air Blank 0.000 16:02</p> <p>Control Test Stats</p> <p>Average 0.0800</p> <p>Std Dev 0.0000</p> <p>Rel Std Dev(%) 0.0000</p>
<p>MC</p> <p>Operator's Signature</p>	<p>MC</p> <p>Operator's Signature</p>	<p>MC</p> <p>Operator's Signature</p>	<p>MC</p> <p>Operator's Signature</p>

Comments:

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: HIGHLANDS COUNTY SO
Time of Inspection: 11:25

Date of Inspection: 06/14/2023

Serial Number: 80-007507
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG223802 Exp: 08/26/2024
0.000	0.048	0.078	0.199	0.081
0.000	0.049	0.078	0.199	0.080
0.000	0.051	0.078	0.200	0.079
0.000	0.052	0.078	0.200	0.081
0.000	0.053	0.078	0.200	0.080
0.000	0.053	0.078	0.200	0.080
0.000	0.053	0.078	0.200	0.081
0.000	0.054	0.078	0.200	0.080
0.000	0.055	0.078	0.200	0.080
0.000	0.055	0.078	0.200	0.080

Standard Deviations	0.0023	0.0000	0.0004	0.0006
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0008 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



TAYLOR D GUTSCHOW

Signature and Printed Name

06/14/2023
Date

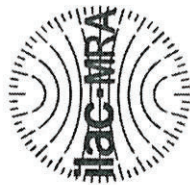

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-007507, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-007507</u>	UNCERTAINTY* \pm
Owning Agency:	<u>HIGHLANDS COUNTY SO</u>	0.050 g/210 L 0.004
Calibration Date:	<u>06/14/2023</u>	0.080 g/210 L 0.004
Calibration Time:	<u>11:25</u>	0.200 g/210 L 0.007
		0.080 g/210 L Dry Gas Control 0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

06/14/2023

Date


TAYLOR D GUTSCHOW,

Department Inspector

Service • Integrity • Respect • Quality