

Florida Department of Law Enforcement Alcohol Testing Program

RECEIVED
NOV 02 2018

AGENCY INSPECTION REPORT - INTOXILYZER 8000

FDLE
Alcohol Testing Program

Agency: CLAY COUNTY SO
Time of Inspection: 16:15

Date of Inspection: 12/23/2017

Serial Number: 80-801443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#:201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#:201702B Exp: 02/23/2019	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:22817080A4 Exp: 10/05/2019
0.000	0.048	0.077	0.198	0.000 / 0.080
0.000	0.048	0.077	0.198	0.000 / 0.080
0.000	0.048	0.078	0.198	0.000 / 0.080


Number of Simulators Used: 5

Remarks:

08:HOSE WAS NOT SEATED PROPERLY

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

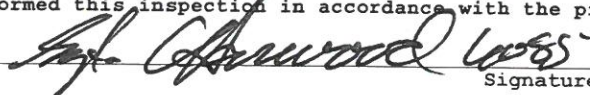

 Signature and Printed Name CHARLES M HARWOOD

12/23/2017
Date

(Handwritten initials)

AGENCY INSPECTION PROCEDURES – INTOXILYZER 8000

1. Prepare at least two simulators for use and allow them to warm up for at least thirty minutes prior to the first analysis. When changing solutions, allow the new solution to warm up for at least ten minutes after the heater light turns off for the first time. Ensure that each simulator maintains an air leak resistant seal and an operational temperature of 34 C (+/- 0.2 C).
2. Only distilled or deionized water must be used for the Alcohol Free Test and the Interferent Detect Test. Class A glassware must be used when measuring solutions.
3. Only approved and non-expired alcohol reference solution and non-expired dry gas standard from an approved source must be used during the applicable portions of the inspection.
4. The instrument must display READY MODE prior to beginning the inspection.
5. Press ESC twice to access the main menu. Enter Agency Inspector last name, first name and middle initial at USER prompts. Press the 2 key to access the Agency Inspector menu. Enter password. Scroll to Inspection Test and press ENTER. Enter Agency Inspector last name, first name and middle initial at INSPECTOR prompts.
6. All results are reported to three decimal places in g/210L. Results of all air blanks must display 0.000. The instrument will abort the inspection process if the air blank result is not 0.000. If any check or test is out of compliance, the instrument will prompt the Agency Inspector to REPEAT (Y/N) the check or test. Each check or test may only be repeated once. If a check or test must be repeated, the REASON must be entered when prompted and recorded in the Remarks section of FDLE/ATP Form 40 Agency Inspection Report – Intoxilyzer 8000.
 - ✓ Check DATE. Adjust if necessary. Press ENTER. Check TIME. Adjust if necessary. Press ENTER.
 - ✓ DIAGNOSTIC CHECK (Pre-Inspection). Press ENTER. All checks must display OK. If any check is not OK, the instrument will abort the inspection process.
 - ✓ NUMBER OF SIMULATORS USED. Enter the number of simulators used during the inspection.
 - ✓ ALCOHOL FREE SUBJECT/MOUTH ALCOHOL TEST. Press ENTER. When PROVIDE SAMPLE NOW is displayed, introduce an alcohol-free breath sample into the instrument. The instrument must display a result of 0.000. Rinse mouth with mouth alcohol solution. When PROVIDE SAMPLE NOW is again displayed introduce a breath sample into the instrument. The instrument must display SLOPE NOT MET.
 - ✓ ALCOHOL FREE TEST. Attach a simulator containing 500 mL distilled or deionized water to the instrument. Press ENTER. Conduct three (3) analyses. The results of all three (3) analyses must be displayed as 0.000.
 - ✓ INTERFERENT DETECT TEST. Attach a simulator containing 3 mL of acetone stock solution and 500 mL distilled or deionized water to the instrument. Press ENTER. Conduct three (3) analyses. The instrument must display INTERFERENT DETECT for each analysis.
 - ✓ 0.05 g/210L TEST. Attach a simulator containing 0.05 g/210L alcohol reference solution to the instrument. Press ENTER. Enter the lot number and expiration date of the alcohol reference solution used. Conduct three (3) analyses. The displayed results of all three (3) analyses must be within the acceptable range.
 - ✓ 0.08 g/210L TEST. Attach a simulator containing 0.08 g/210L alcohol reference solution to the instrument. Press ENTER. Enter the lot number and expiration date of the alcohol reference solution used. Conduct three (3) analyses. The displayed results of all three (3) analyses must be within the acceptable range.
 - ✓ 0.20 g/210L TEST. Attach a simulator containing 0.20 g/210L alcohol reference solution to the instrument. Press ENTER. Enter the lot number and expiration date of the alcohol reference solution used. Conduct three (3) analyses. The displayed results of all three (3) analyses must be within the acceptable range.
 - ✓ 0.08 g/210L DRY GAS STANDARD TEST. Attach a cylinder containing 0.08 g/210L dry gas standard to the instrument. Press ENTER. Enter the lot number and expiration date of the dry gas standard used. Conduct three (3) analyses. The displayed results of all three (3) analyses must be within the acceptable range.
 - ✓ DIAGNOSTIC CHECK (Post-Inspection). Press ENTER. All checks must display OK. If any check is not OK, the instrument will abort the inspection process.
 - ✓ REVIEW REMARKS. Enter Y/N.
 - ✓ IN COMPLIANCE. Enter Y/N to state whether the instrument complies or does not comply with the requirements of Chapter 11D-8, FAC. If the instrument does not comply with Chapter 11D-8, FAC, remove the instrument from service and notify the Department Inspector.
7. The results of the Agency Inspection must be recorded on FDLE/ATP Form 40 Agency Inspection Report – Intoxilyzer 8000.
8. For regulatory and administrative purposes only, the results of the agency inspection must be made electronically available to the Department within five (5) business days of completing the inspection.

Inspection Completed by Sgt. C. Howard 6085
 on 12-23-17
 - Sgt. Howard Coast

RECEIVED

Florida Department of Law Enforcement Alcohol Testing Program

NOV 27 2017

FDLE

Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: CLAY COUNTY SO
Time of Inspection: 12:07

Date of Inspection: 11/21/2017

Serial Number: 80-001443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#:201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#:201702B Exp: 02/23/2019	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:22817080A4 Exp: 10/05/2019
0.000	0.048	0.076	0.197	0.079
0.000	0.048	0.077	0.198	0.079
0.000	0.048	0.077	0.199	0.079

Number of Simulators Used: 5

Remarks:
MONTHLY INSPECTION

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

RICHARD L PATRONE

Signature and Printed Name

11/21/2017
Date

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000


Agency: CLAY COUNTY SO Date of Inspection: 10/21/2017 Serial Number: 80-001443
 Time of Inspection: 10:49 Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201702B Exp: 02/23/2019	0.08g/210L Test (g/210L) Lot#:201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#:201702B Exp: 02/23/2019	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:09117080A2 Exp: 06/05/2019
0.000	0.049	0.078	0.198	0.081
0.000	0.049	0.079	0.199	0.081
0.000	0.049	0.079	0.199	0.081

Number of Simulators Used: 5

Remarks:

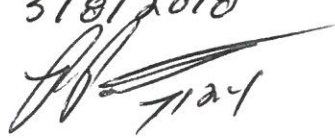MONTHLY INSPECTION
 .05 TEST LOT # ENTERED IN ERROR
 .05 TEST CORRECT # LOT 201603D
 EXP 3/8/2018


The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.


7124
Signature and Printed Name RICHARD L PATRONE

10/21/2017
 Date

Florida Department of Law Enforcement Alcohol Testing Program

RECEIVED

OCT 02 2017

FDLE

Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: CLAY COUNTY SO
Time of Inspection: 13:37

Date of Inspection: 09/24/2017

Serial Number: 80-001443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#:201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#:201702B Exp: 02/23/2019	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:09117080A2 Exp: 06/05/2019
0.000	0.048	0.078	0.200	0.080
0.000	0.048	0.079	0.199	0.080
0.000	0.048	0.079	0.200	0.080

Number of Simulators Used: 5

Remarks:
AGENCY INSPECTION

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

 7127
Signature and Printed Name

RICHARD L PATRONE

09/24/2017
Date

1443

RECEIVED

Florida Department of Law Enforcement Alcohol Testing Program

AUG 29 2017

FDLE

Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: CLAY COUNTY SO
Time of Inspection: 16:06

Date of Inspection: 08/22/2017

Serial Number: 80-001443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: 201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#: 201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#: 201702B Exp: 02/23/2019	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: 09117080A2 Exp: 06/05/2019
0.000	0.048	0.078	0.199	0.080
0.000	0.048	0.078	0.199	0.079
0.000	0.049	0.079	0.200	0.079

Number of Simulators Used: 5

Remarks:
MONTHLY INSPECTION

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

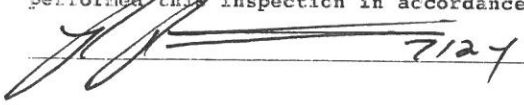

RICHARD L. PATRONE
Signature and Printed Name

08/22/2017
Date

1443

RECEIVED

Florida Department of Law Enforcement JUN 21 2017 Alcohol Testing Program

FDLE
Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: CLAY COUNTY SO
Time of Inspection: 12:26

Date of Inspection: 06/18/2017

Serial Number: 80-001443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#:201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#:201604C Exp: 04/05/2018	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:32316080A2 Exp: 01/05/2019
0.000	0.046	0.077	0.195	0.080
0.000	0.046	0.078	0.196	0.080
0.000	0.047	0.078	0.196	0.080

Number of Simulators Used: 5

Remarks:

MONTHLY INSPECTION RETURN FOR ANNUAL

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

RICHARD L PATRONE

Signature and Printed Name

06/18/2017

Date

1443

RECEIVED

Florida Department of Law Enforcement Alcohol Testing Program

MAY 22 2017

FDLE
Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: CLAY COUNTY SO
Time of Inspection: 11:32
Date of Inspection: 05/19/2017
Serial Number: 80-001443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#:201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#:201604C Exp: 04/05/2018	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:32316080A2 Exp: 01/05/2019
0.000	0.047	0.076	0.194	0.079
0.000	0.047	0.077	0.196	0.079
0.000	0.048	0.078	0.197	0.079

Number of Simulators Used: 5

Remarks: MONTHLY INSPECTION

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

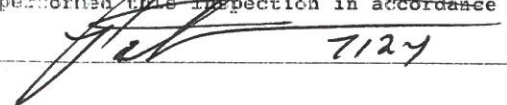 7127

RICHARD L. PATRONE
Signature and Printed Name

05/19/2017
Date

1443 (8)

Florida Department of Law Enforcement Alcohol Testing Program

RECEIVED

APR 27 2017

FDLE

AGENCY INSPECTION REPORT - INTOXILYZER 8000 Alcohol Testing Program

Agency: CLAY COUNTY SO
Time of Inspection: 08:14

Date of Inspection: 04/21/2017

Serial Number: 80-001443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#:201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#:201604C Exp: 04/05/2018	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:32316080A2 Exp: 01/05/2019
0.000	0.047	0.076	0.193	0.080
0.000	0.048	0.077	0.195	0.080
0.000	0.048	0.078	0.195	0.080

Number of Simulators Used: 5

Remarks:
MONTHLY INSPECTION

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

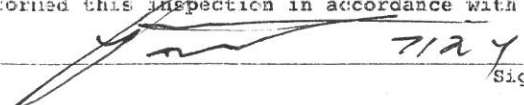

RICHARD L. PATRONE

Signature and Printed Name

04/21/2017
Date

Florida Department of Law Enforcement
Alcohol Testing Program

80-001443

AGENCY INSPECTION PROCEDURES – INTOXILYZER 8000

1. Prepare at least two simulators for use and allow them to warm up for at least thirty minutes prior to the first analysis. When changing solutions, allow the new solution to warm up for at least ten minutes after the heater light turns off for the first time. Ensure that each simulator maintains an air leak resistant seal and an operational temperature of 34 C (+/- 0.2 C).
2. Only distilled or deionized water must be used for the Alcohol Free Test and the Interferent Detect Test. Class A glassware must be used when measuring solutions.
3. Only approved and non-expired alcohol reference solution and non-expired dry gas standard from an approved source must be used during the applicable portions of the inspection.
4. The instrument must display READY MODE prior to beginning the inspection.
5. Press ESC twice to access the main menu. Enter Agency Inspector last name, first name and middle initial at USER prompts. Press the 2 key to access the Agency Inspector menu. Enter password. Scroll to Inspection Test and press ENTER. Enter Agency Inspector last name, first name and middle initial at INSPECTOR prompts.
6. All results are reported to three decimal places in g/210L. Results of all air blanks must display 0.000. The instrument will abort the inspection process if the air blank result is not 0.000. If any check or test is out of compliance, the instrument will prompt the Agency Inspector to REPEAT (Y/N) the check or test. Each check or test may only be repeated once. If a check or test must be repeated, the REASON must be entered when prompted and recorded in the Remarks section of FDLE/ATP Form 40 Agency Inspection Report – Intoxilyzer 8000.
 1. Check DATE. Adjust if necessary. Press ENTER. Check TIME. Adjust if necessary. Press ENTER.
 2. DIAGNOSTIC CHECK (Pre-Inspection). Press ENTER. All checks must display OK. If any check is not OK, the instrument will abort the inspection process.
 3. NUMBER OF SIMULATORS USED. Enter the number of simulators used during the inspection.
 4. ALCOHOL FREE SUBJECT/MOUTH ALCOHOL TEST. Press ENTER. When PROVIDE SAMPLE NOW is displayed, introduce an alcohol-free breath sample into the instrument. The instrument must display a result of 0.000. Rinse mouth with mouth alcohol solution. When PROVIDE SAMPLE NOW is again displayed introduce a breath sample into the instrument. The instrument must display SLOPE NOT MET.
 5. ALCOHOL FREE TEST. Attach a simulator containing 500 mL distilled or deionized water to the instrument. Press ENTER. Conduct three (3) analyses. The results of all three (3) analyses must be displayed as 0.000.
 6. INTERFERENT DETECT TEST. Attach a simulator containing 3 mL of acetone stock solution and 500 mL distilled or deionized water to the instrument. Press ENTER. Conduct three (3) analyses. The instrument must display INTERFERENT DETECT for each analysis.
 7. 0.05 g/210L TEST. Attach a simulator containing 0.05 g/210L alcohol reference solution to the instrument. Press ENTER. Enter the lot number and expiration date of the alcohol reference solution used. Conduct three (3) analyses. The displayed results of all three (3) analyses must be within the acceptable range.
 8. 0.08 g/210L TEST. Attach a simulator containing 0.08 g/210L alcohol reference solution to the instrument. Press ENTER. Enter the lot number and expiration date of the alcohol reference solution used. Conduct three (3) analyses. The displayed results of all three (3) analyses must be within the acceptable range.
 9. 0.20 g/210L TEST. Attach a simulator containing 0.20 g/210L alcohol reference solution to the instrument. Press ENTER. Enter the lot number and expiration date of the alcohol reference solution used. Conduct three (3) analyses. The displayed results of all three (3) analyses must be within the acceptable range.
 10. 0.08 g/210L DRY GAS STANDARD TEST. Attach a cylinder containing 0.08 g/210L dry gas standard to the instrument. Press ENTER. Enter the lot number and expiration date of the dry gas standard used. Conduct three (3) analyses. The displayed results of all three (3) analyses must be within the acceptable range.
7. DIAGNOSTIC CHECK (Post-Inspection). Press ENTER. All checks must display OK. If any check is not OK, the instrument will abort the inspection process.
8. REVIEW REMARKS. Enter Y/N.
9. IN COMPLIANCE. Enter Y/N to state whether the instrument complies or does not comply with the requirements of Chapter 11D-8, FAC. If the instrument does not comply with Chapter 11D-8, FAC, remove the instrument from service and notify the Department Inspector.
7. The results of the Agency Inspection must be recorded on FDLE/ATP Form 40 Agency Inspection Report – Intoxilyzer 8000.
8. For regulatory and administrative purposes only, the results of the agency inspection must be made electronically available to the Department within five (5) business days of completing the inspection.

Inspection Completed by Sgt. C Harwood
Sgt. C Harwood 3-18-17

RECEIVED

Florida Department of Law Enforcement APR 11 2017 Alcohol Testing Program

FDLE
Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: CLAY COUNTY SO
Time of Inspection: 11:31

Date of Inspection: 03/18/2017

Serial Number: 80-001443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#:201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#:201604C Exp: 04/05/2018	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:32316080A2 Exp: 01/05/2019
0.000	0.048	0.078	0.197	0.080
0.000	0.048	0.078	0.198	0.080
0.000	0.049	0.078	0.198	0.080

Number of Simulators Used: 5

Remarks:

Time-Date changed. A F / M A: Range Exceeded.

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Charles M. Harwood

Signature and Printed Name

CHARLES M. HARWOOD

03/18/2017
Date

Florida Department of Law Enforcement Alcohol Testing Program

1443
⑨

RECEIVED

FEB 15 2017

FDLE
Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: CLAY COUNTY SO
Time of Inspection: 14:30

Date of Inspection: 02/11/2017

Serial Number: 80-001443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: 201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#: 201611B Exp: 11/15/2018	0.20g/210L Test (g/210L) Lot#: 201604C Exp: 04/05/2018	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: 20016080A1 Exp: 09/05/2018
0.000	0.047	0.077	0.195	0.081
0.000	0.048	0.078	0.195	0.081
0.000	0.049	0.079	0.196	0.081

Number of Simulators Used: 5

Remarks:
MONTHLY INSPECTION

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.


RICHARD L. PATRONE
 Signature and Printed Name

02/11/2017
 Date

1443

RECEIVED

Florida Department of Law Enforcement JAN 26 2017

Alcohol Testing Program

FDLE
Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: CLAY COUNTY SO
Time of Inspection: 14:44

Date of Inspection: 01/21/2017

Serial Number: 80-001443
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK	Yes	
Alcohol Free Subject Test: 0.000	Yes	
Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes	
Diagnostic Check (Post-Inspection): OK	Yes	


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201603D Exp: 03/08/2018	0.08g/210L Test (g/210L) Lot#:201501F Exp: 01/26/2018	0.20g/210L Test (g/210L) Lot#:201604C Exp: 04/05/2018	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:20016080A1 Exp: 09/05/2018
0.000	0.048	0.078	0.195	0.080
0.000	0.048	0.078	0.197	0.080
0.000	0.048	0.078	0.196	0.080

Number of Simulators Used: 5

Remarks:
MONTHLY INSPECTION

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

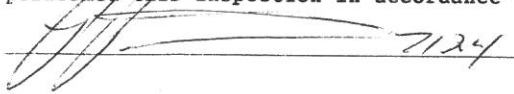 7124
Signature and Printed Name RICHARD L. PATRONE

01/21/2017
Date