



INSTRUMENT PROCESSING SHEET

Agency Lake Hamilton PDS/N 80-001866Florida Department of
Law EnforcementDate In 10/28/2019 DI Completion Date 11/5/19☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	Performed By <u>DP</u>	Quality Checks	Performed By <u>DP</u>	Flow Calibration	Performed By																																																										
☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>211</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP102</u> 32 mm <u>.142</u> (.139 - .169) 36 mm <u>.160</u> (.156 - .190) 53 mm <u>.238</u> (.228 - .278) 103 mm <u>.488</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28421</u> ☒ Stability Checks		Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																											
Final Release Date FDLE DEC 05 2019 Alcohol Testing Program		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td><u>SD1018</u></td><td><u>201905A</u> <u>5/14/21</u></td></tr><tr><td>0.080</td><td><u>SD3962</u></td><td><u>201905B</u> <u>5/14/21</u></td></tr><tr><td>0.200</td><td><u>G2078</u></td><td><u>201904D</u> <u>4/30/21</u></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td><u>A6916501</u> <u>6/14/21</u></td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	<u>SD1018</u>	<u>201905A</u> <u>5/14/21</u>	0.080	<u>SD3962</u>	<u>201905B</u> <u>5/14/21</u>	0.200	<u>G2078</u>	<u>201904D</u> <u>4/30/21</u>	0.080 DGS	N/A	<u>A6916501</u> <u>6/14/21</u>		Maintenance Performed By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ Temperature Checks Performed By <u>DP</u> ☒ Lab Temp °C <u>21.0</u> External Digital Therm. ID#: <u>300505</u> ☒ 34°C +/- .2 Serial #: <u>SD1018</u> ☒ 34°C +/- .2 Serial #: <u>SD3962</u> ☒ 34°C +/- .2 Serial #: <u>G2078</u>																																												
Simulator	Serial #	Lot #/Exp																																																													
0.050	<u>SD1018</u>	<u>201905A</u> <u>5/14/21</u>																																																													
0.080	<u>SD3962</u>	<u>201905B</u> <u>5/14/21</u>																																																													
0.200	<u>G2078</u>	<u>201904D</u> <u>4/30/21</u>																																																													
0.080 DGS	N/A	<u>A6916501</u> <u>6/14/21</u>																																																													
Calibration Adjustment Performed By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th><th>Lot Number</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th><th>Lot Number</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>	Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A				Department Inspection Performed By <u>DP</u> Barometric Pressure ID# <u>28421</u> Gauge <u>1020</u> Instrument <u>1016</u> Mouth Alcohol Solution Lot # <u>2018-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td><u>G4444</u></td></tr><tr><td>Interferent</td><td><u>G6621</u></td></tr><tr><td>0.050</td><td><u>SD1018</u></td></tr><tr><td>0.080</td><td><u>SD3962</u></td></tr><tr><td>0.200</td><td><u>G2078</u></td></tr></tbody></table> Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☐ Other _____ ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use <u>SP 11/7/19</u> <u>Brett Kirkland 12/4/19</u> Tech Review / Date Admin Review / Date	Simulator	Serial Number	0.000	<u>G4444</u>	Interferent	<u>G6621</u>	0.050	<u>SD1018</u>	0.080	<u>SD3962</u>	0.200	<u>G2078</u>	
Simulator	Serial Number	Lot Number	Expiration																																																												
0.000		N/A	N/A																																																												
0.040																																																															
0.100																																																															
0.200																																																															
0.300																																																															
0.080 DGS	N/A																																																														
Simulator	Serial Number	Lot Number	Expiration																																																												
0.050																																																															
0.080																																																															
0.200																																																															
0.080 DGS	N/A																																																														
Simulator	Serial Number																																																														
0.000	<u>G4444</u>																																																														
Interferent	<u>G6621</u>																																																														
0.050	<u>SD1018</u>																																																														
0.080	<u>SD3962</u>																																																														
0.200	<u>G2078</u>																																																														
Notes/Suggested Service: _____ _____ _____ _____ _____ _____																																																															

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: LAKE HAMILTON PD
~~US AIR FORCE MACDILL~~
 Time of Inspection: 12:21

Date of Inspection: 11/05/2019

Serial Number: 80-001866
 Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG916501 Exp: 06/14/2021
0.000	0.049	0.080	0.197	0.080
0.000	0.051	0.078	0.200	0.080
0.000	0.051	0.077	0.201	0.079
0.000	0.050	0.076	0.202	0.080
0.000	0.052	0.078	0.201	0.079
0.000	0.050	0.077	0.201	0.079
0.000	0.052	0.078	0.200	0.079
0.000	0.050	0.077	0.200	0.079
0.000	0.050	0.077	0.200	0.079
0.000	0.050	0.077	0.200	0.079

Standard Deviations	0.0009	0.0010	0.0013	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0009 Number of Simulators Used: 5

Remarks:

SP
BK
12/4/19

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Patrick J Murphy

PATRICK J MURPHY

Signature and Printed Name

11/05/2019
Date

US AIR FORCE MACDILL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001866
11/05/2019
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:53
Control Test	0.049	08:54
Air Blank	0.000	08:54
Control Test	0.049	08:55
Air Blank	0.000	08:55
Control Test	0.049	08:56
Air Blank	0.000	08:56
Control Test Stats		
Average	0.0490	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

P Murphy
Operator's Signature

US AIR FORCE MACDILL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001866
11/05/2019
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:58
Control Test	0.079	08:58
Air Blank	0.000	08:59
Control Test	0.079	09:00
Air Blank	0.000	09:00
Control Test	0.079	09:01
Air Blank	0.000	09:01
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

P Murphy
Operator's Signature

US AIR FORCE MACDILL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001866
11/05/2019
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:03
Control Test	0.196	09:04
Air Blank	0.000	09:05
Control Test	0.197	09:05
Air Blank	0.000	09:06
Control Test	0.200	09:06
Air Blank	0.000	09:07
Control Test Stats		
Average	0.1977	
Std Dev	0.0021	
Rel Std Dev(%)	1.0531	

P Murphy
Operator's Signature

US AIR FORCE MACDILL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001866
11/05/2019
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:15
Control Test	0.081	09:15
Air Blank	0.000	09:16
Control Test	0.080	09:16
Air Blank	0.000	09:16
Control Test	0.080	09:17
Air Blank	0.000	09:17
Control Test Stats		
Average	0.0803	
Std Dev	0.0006	
Rel Std Dev(%)	0.7187	

DGS

P Murphy
Operator's Signature

SP
B/K
12/4/19



ANAB
ACCREDITED
ISO/IEC 17025
FORENSIC CALIBRATION
LABORATORY

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2729 Fort Knox Blvd.
Bldg. 2, Suite 1300
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001866, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001866</u>	<u>LAKE HAMILTON P.D. QPM</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>US AIR FORCE MACDILL</u>		0.050 g/ 210 L 0.004
Calibration Date:	<u>11/05/2019</u>		0.080 g/ 210 L 0.004
Calibration Time:	<u>12:21</u>		0.200 g/ 210 L 0.007
			0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 July 2018

Issuing Authority: Alcohol Testing Program

11/05/2019

Date

PATRICK J MURPHY
PATRICK J MURPHY,
Department Inspector

Service • Integrity • Respect • Quality

Page 1 of 1

12/4/19

Return Material Authorization

Ship to: ☐ CMI, Inc.

☒ Enforcement Electronics

Shipment to repair facility authorized by: TJ Graham on 6/5/2019

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-001866

Bill To Address:
Lake Hamilton PD

Ship to Address:
FDLE Off-site Mail Facility

Alcohol Testing Program

813-B Lake Bradford Road

Tallahassee, FL 32304

Reason for Return:

Please check the instrument for analytical stability. Ethanol values vary when

I am running standards. A battery change will also be required. Please contact TJ Graham at ATP with a quote prior to performing repairs. This instrument was sent to ATP and is no longer registered.

Please choose one of the following options:

☐ 1. I _____, authorize all repairs.

☐ 2. I _____, authorize repairs up to \$ _____.

☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: TJ Graham

Phone #: 850-445-7938 Email: thomasgraham@fdle.state.fl.us

ATP Contact Name: TJ Graham ATP Email: thomasgraham@fdle.state.fl.us



INSTRUMENT PROCESSING SHEET

Agency US Air Force MacDill

S/N 80-001866

Florida Department of
Law Enforcement

Date In 01/25/2019 DI Completion Date _____

☐ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake Performed By <u>JO</u> ☒ Annual ☒ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____	Quality Checks Performed By <u>JO</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>209</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-105</u> 32 mm <u>.160</u> (.139 - .169) 36 mm <u>.175</u> (.156 - .190) 53 mm <u>.246</u> (.228 - .278) 103 mm <u>.523</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28421</u> ☒ Stability Checks <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> <tr> <td>0.050</td> <td>SD1012</td> <td>201707D 07/25/2019</td> </tr> <tr> <td>0.080</td> <td>DR1279</td> <td>201707E 07/25/2019</td> </tr> <tr> <td>0.200</td> <td>DR3856</td> <td>201707C 07/24/2019</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG805701 02/26/2020</td> </tr> </table>	Simulator	Serial #	Lot #/Exp	0.050	SD1012	201707D 07/25/2019	0.080	DR1279	201707E 07/25/2019	0.200	DR3856	201707C 07/24/2019	0.080 DGS	N/A	AG805701 02/26/2020	Flow Calibration Performed By _____ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance Performed By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ Temperature Checks Performed By _____ ☐ Lab Temp °C _____ External Digital Therm. ID#: _____ ☐ 34°C +/- .2 Serial #: _____ ☐ 34°C +/- .2 Serial #: _____ ☐ 34°C +/- .2 Serial #: _____																																	
Simulator	Serial #	Lot #/Exp																																																
0.050	SD1012	201707D 07/25/2019																																																
0.080	DR1279	201707E 07/25/2019																																																
0.200	DR3856	201707C 07/24/2019																																																
0.080 DGS	N/A	AG805701 02/26/2020																																																
Final Release Date _____ _____	Calibration Adjustment Performed By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table>		Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A		
Simulator	Serial Number	Lot Number	Expiration																																															
0.000		N/A	N/A																																															
0.040																																																		
0.100																																																		
0.200																																																		
0.300																																																		
0.080 DGS	N/A																																																	
Simulator	Serial Number	Lot Number	Expiration																																															
0.050																																																		
0.080																																																		
0.200																																																		
0.080 DGS	N/A																																																	
Notes/Suggested Service: <u>Instrument needs a battery change.</u> _____ _____ _____ _____ _____		Department Inspection Performed By _____ Barometric Pressure ID# _____ Gauge _____ Instrument _____ Mouth Alcohol Solution Lot # _____ Acetone Stock Solution Lot # _____ <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> <tr> <td>0.000</td> <td></td> </tr> <tr> <td>Interferent</td> <td></td> </tr> <tr> <td>0.050</td> <td></td> </tr> <tr> <td>0.080</td> <td></td> </tr> <tr> <td>0.200</td> <td></td> </tr> </table> Attachments ☐ Form 41 ☐ Stability Checks ☐ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____ ☐ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☐ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☐ Conduct an Agency Inspection Before Evidentiary Use Tech Review / Date _____ Admin Review / Date _____	Simulator	Serial Number	0.000		Interferent		0.050		0.080		0.200																																					
Simulator	Serial Number																																																	
0.000																																																		
Interferent																																																		
0.050																																																		
0.080																																																		
0.200																																																		

US AIR FORCE MACOILL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001866
01/28/2019
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:29
Control Test	0.046	09:30
Air Blank	0.000	09:30
Control Test	0.050	09:31
Air Blank	0.000	09:31
Control Test	0.050	09:32
Air Blank	0.000	09:32
Control Test Stats		
Average	0.0467	
Std Dev	0.0023	
Rel Std Dev(%)	4.7453	


Operator's Signature

US AIR FORCE MACOILL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 60-001866
01/28/2019
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:24
Control Test	0.081	09:25
Air Blank	0.000	09:25
Control Test	0.081	09:26
Air Blank	0.000	09:27
Control Test	0.082	09:27
Air Blank	0.000	09:28
Control Test Stats		
Average	0.0813	
Std Dev	0.0006	
Rel Std Dev(%)	0.7099	


Operator's Signature

US AIR FORCE MACOILL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001866
01/28/2019
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:38
Control Test	0.195	09:39
Air Blank	0.000	09:39
Control Test	0.196	09:40
Air Blank	0.000	09:40
Control Test	0.194	09:41
Air Blank	0.000	09:41
Control Test Stats		
Average	0.1950	
Std Dev	0.0010	
Rel Std Dev(%)	0.5128	


Operator's Signature

US AIR FORCE MACOILL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001866
01/28/2019
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:34
Control Test	0.080	09:34
Air Blank	0.000	09:35
Control Test	0.078	09:35
Air Blank	0.000	09:35
Control Test	0.080	09:36
Air Blank	0.000	09:36
Control Test Stats		
Average	0.0793	
Std Dev	0.0012	
Rel Std Dev(%)	1.4555	

DCS


Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: US AIR FORCE MACDILL

Time of Inspection: 08:52

Date of Inspection: 01/28/2019

Serial Number: 80-001866

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:


Number of Simulators Used: _____

Remarks:

BYPASSED AI TO OPERATE INSTRUMENT. NOT A COMPLIANCE CHECK.

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

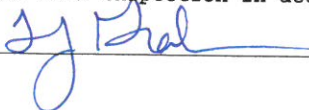

Signature and Printed Name

THOMAS J GRAHAM

01/28/2019
Date