




INSTRUMENT PROCESSING SHEET

Agency Pasco Co SOS/N 80-001074

Florida Department of Law Enforcement

Date In 08/03/2020DI Completion Date 8/4/2020Ship P/U H/D CMI EE

Intake Performed By <u>KAW</u> <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input checked="" type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ _____ _____	Quality Checks Performed By <u>[Signature]</u> <input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>107</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP-103</u> 32 mm <u>0.140</u> (.139 - .169) 36 mm <u>0.156</u> (.156 - .190) 53 mm <u>0.230</u> (.228 - .278) 103 mm <u>0.496</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>30793</u> <input checked="" type="checkbox"/> Stability Checks	Flow Calibration Performed By _____ Flow Column # _____ <input type="checkbox"/> 5L/min - 17mm <input type="checkbox"/> 15L/min - 53mm <input type="checkbox"/> 30L/min - 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																											
Final Release Date FDLE Alcohol Testing Program Digitally signed by FDLE Alcohol Testing Program Date: 2020.08.06 12:33:22 -04'00'	<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>MP5088</td> <td>201905A 05-17-2021</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> <td>201905B 05-17-2021</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> <td>201907D 07-30-2021</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG931603 11-12-2021</td> </tr> </tbody> </table>	Simulator	Serial #	Lot #/Exp	0.050	MP5088	201905A 05-17-2021	0.080	MP5089	201905B 05-17-2021	0.200	MP5090	201907D 07-30-2021	0.080 DGS	N/A	AG931603 11-12-2021	Maintenance Performed By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____ Temperature Checks Performed By <u>SP</u> <input checked="" type="checkbox"/> Lab Temp °C <u>21.6</u> External Digital Therm. ID#: <u>300505</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5088</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5089</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5090</u>																																												
Simulator	Serial #	Lot #/Exp																																																											
0.050	MP5088	201905A 05-17-2021																																																											
0.080	MP5089	201905B 05-17-2021																																																											
0.200	MP5090	201907D 07-30-2021																																																											
0.080 DGS	N/A	AG931603 11-12-2021																																																											
Calibration Adjustment Performed By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> <input type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>	Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection Performed By <u>SP</u> Barometric Pressure ID# <u>28421</u> Gauge <u>1014</u> Instrument <u>1011</u> Mouth Alcohol Solution Lot # <u>2019-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>MP5086</td> </tr> <tr> <td>Interferent</td> <td>MP287</td> </tr> <tr> <td>0.050</td> <td>MP5088</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> </tr> </tbody> </table>	Simulator	Serial Number	0.000	MP5086	Interferent	MP287	0.050	MP5088	0.080	MP5089	0.200	MP5090
Simulator	Serial Number	Lot Number	Expiration																																																										
0.000		N/A	N/A																																																										
0.040																																																													
0.100																																																													
0.200																																																													
0.300																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number	Lot Number	Expiration																																																										
0.050																																																													
0.080																																																													
0.200																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number																																																												
0.000	MP5086																																																												
Interferent	MP287																																																												
0.050	MP5088																																																												
0.080	MP5089																																																												
0.200	MP5090																																																												
Notes/Suggested Service: _____ _____ _____ _____ _____	<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use																																																												
	Tech Review / Date <u>Michael D. Haughey 2020.08.05 09:48:14 -04'00'</u>	Admin Review / Date <u>Brett Kirkland 2020.08.06 12:29:44 -04'00'</u>																																																											

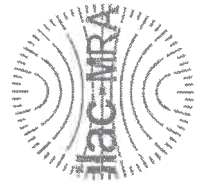

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2729 Fort Knox Blvd.
Bldg. 2, Suite 1300
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001074, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001074</u>	UNCERTAINTY* ±	
Owning Agency:	<u>PASCO COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>08/04/2020</u>	0.080 g/ 210 L	0.005
Calibration Time:	<u>14:25</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.
Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).
The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Shayla Platt

08/04/2020

Date

SHAYLA D PLATT,

Department Inspector

FDLE/ATP Form 69 April 2020
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

BK

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PASCO COUNTY SO
Time of Inspection: 14:25

Date of Inspection: 08/04/2020

Serial Number: 80-001074
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG931603 Exp: 11/12/2021
0.000	0.048	0.079	0.200	0.080
0.000	0.049	0.080	0.201	0.080
0.000	0.049	0.079	0.201	0.080
0.000	0.048	0.079	0.201	0.080
0.000	0.049	0.079	0.201	0.080
0.000	0.049	0.080	0.201	0.080
0.000	0.048	0.080	0.201	0.080
0.000	0.048	0.079	0.200	0.080
0.000	0.048	0.080	0.201	0.080
0.000	0.049	0.080	0.201	0.080
Standard Deviations	0.0005	0.0005	0.0004	0.0000

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

Int Det: Sequence Aborted *Instrument stalled on air blank after pressing green start test button the air blank restarted and completed as normal. All other air blanks during inspection continued without issue.-sp*

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT

Signature and Printed Name

08/04/2020
Date

MX
BK
2020.08.0
6 12:30:53
-04'00'

Stability Checks

PASCO COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-001074
 08/04/2020
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:50
Control Test	0.050	08:51
Air Blank	0.000	08:51
Control Test	0.049	08:52
Air Blank	0.000	08:52
Control Test	0.051	08:53
Air Blank	0.000	08:54
Control Test Stats		
Average	0.0500	
Std Dev	0.0010	
Rel Std Dev(%)	2.0000	

PASCO COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-001074
 08/04/2020
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:55
Control Test	0.080	08:55
Air Blank	0.000	08:56
Control Test	0.079	08:56
Air Blank	0.000	08:57
Control Test	0.079	08:58
Air Blank	0.000	08:58
Control Test Stats		
Average	0.0793	
Std Dev	0.0006	
Rel Std Dev(%)	0.7277	

PASCO COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-001074
 08/04/2020
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:59
Control Test	0.202	09:00
Air Blank	0.000	09:01
Control Test	0.201	09:01
Air Blank	0.000	09:02
Control Test	0.201	09:03
Air Blank	0.000	09:03
Control Test Stats		
Average	0.2013	
Std Dev	0.0006	
Rel Std Dev(%)	0.2868	

Wet

SS

Operator's Signature

SS

Operator's Signature

SS

Operator's Signature

PASCO COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-001074
 08/04/2020
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:04
Control Test	0.080	09:04
Air Blank	0.000	09:05
Control Test	0.081	09:05
Air Blank	0.000	09:05
Control Test	0.081	09:06
Air Blank	0.000	09:06
Control Test Stats		
Average	0.0807	
Std Dev	0.0006	
Rel Std Dev(%)	0.7157	

Dry

SS

Operator's Signature

MX