




INSTRUMENT PROCESSING SHEET

Agency Edgewater Police DepartmentS/N 80-001243Florida Department of
Law EnforcementDate In 04/07/2020DI Completion Date 4/20/2020☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake Performed By <u>RAW</u> ☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☐ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ 	Quality Checks Performed By <u>DA</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>205</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-103</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.226</u> (.228 - .278) 103 mm <u>0.492</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>30793</u> ☒ Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td><u>MP5088</u></td> <td><u>201905A</u> <u>05-14-2021</u></td> </tr> <tr> <td>0.080</td> <td><u>MP5089</u></td> <td><u>201905B</u> <u>05-14-2021</u></td> </tr> <tr> <td>0.200</td> <td><u>MP5090</u></td> <td><u>201904D</u> <u>04-30-2021</u></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td><u>AG931603</u> <u>11-12-2021</u></td> </tr> </tbody> </table>	Simulator	Serial #	Lot #/Exp	0.050	<u>MP5088</u>	<u>201905A</u> <u>05-14-2021</u>	0.080	<u>MP5089</u>	<u>201905B</u> <u>05-14-2021</u>	0.200	<u>MP5090</u>	<u>201904D</u> <u>04-30-2021</u>	0.080 DGS	N/A	<u>AG931603</u> <u>11-12-2021</u>	Flow Calibration Performed By <u>SP</u> Flow Column # <u>ATP102</u> ☒ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☒ R-Value <u>201</u> ☒ Post Calibration Verification (L/s) Flow Column # <u>ATP103</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.503</u> (.447 - .547) Maintenance Performed By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ Temperature Checks Performed By <u>SP</u> ☒ Lab Temp °C <u>21.5</u> External Digital Therm. ID#: <u>300502</u> ☒ 34°C +/-2 Serial #: <u>MP5088</u> ☒ 34°C +/-2 Serial #: <u>MP5089</u> ☒ 34°C +/-2 Serial #: <u>MP5090</u>																																													
Simulator	Serial #	Lot #/Exp																																																												
0.050	<u>MP5088</u>	<u>201905A</u> <u>05-14-2021</u>																																																												
0.080	<u>MP5089</u>	<u>201905B</u> <u>05-14-2021</u>																																																												
0.200	<u>MP5090</u>	<u>201904D</u> <u>04-30-2021</u>																																																												
0.080 DGS	N/A	<u>AG931603</u> <u>11-12-2021</u>																																																												
Final Release Date FDLE Alcohol Testing Program	Digitally signed by FDLE Alcohol Testing Program Date: 2020.04.20 16:12:57 -04'00'																																																													
Calibration Adjustment Performed By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>	Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection Performed By <u>SP</u> Barometric Pressure ID# <u>28421</u> Gauge <u>1003</u> Instrument <u>1001</u> Mouth Alcohol Solution Lot # <u>2019-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td><u>MP5086</u></td> </tr> <tr> <td>Interferent</td> <td><u>MP5087</u></td> </tr> <tr> <td>0.050</td> <td><u>MP5088</u></td> </tr> <tr> <td>0.080</td> <td><u>MP5089</u></td> </tr> <tr> <td>0.200</td> <td><u>MP5090</u></td> </tr> </tbody> </table> Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☒ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☐ Other _____ ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use <u>Michael D. Hagley</u> 2020.04.20 11:14:06 -04'00' <u>Brett Kirkland</u> 2020.04.20 16:04:30 -04'00' Tech Review / Date Admin Review / Date		Simulator	Serial Number	0.000	<u>MP5086</u>	Interferent	<u>MP5087</u>	0.050	<u>MP5088</u>	0.080	<u>MP5089</u>	0.200	<u>MP5090</u>
Simulator	Serial Number	Lot Number	Expiration																																																											
0.000		N/A	N/A																																																											
0.040																																																														
0.100																																																														
0.200																																																														
0.300																																																														
0.080 DGS	N/A																																																													
Simulator	Serial Number	Lot Number	Expiration																																																											
0.050																																																														
0.080																																																														
0.200																																																														
0.080 DGS	N/A																																																													
Simulator	Serial Number																																																													
0.000	<u>MP5086</u>																																																													
Interferent	<u>MP5087</u>																																																													
0.050	<u>MP5088</u>																																																													
0.080	<u>MP5089</u>																																																													
0.200	<u>MP5090</u>																																																													

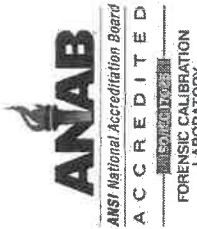

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2729 Fort Knox Blvd.
Bldg. 2, Suite 1300
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001243, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001243</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>EDGEWATER P.D.</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>04/20/2020</u>	0.080 g/ 210 L 0.005
Calibration Time:	<u>09:31</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

04/20/2020

Date

Shayla Platt

SHAYLA D PLATT,
Department Inspector

FDLE/ATP Form 69 January 2020

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1

2020.04.
20
16:05:21
-04'00'

114

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: EDGEWATER P.D.

Time of Inspection: 09:31

Date of Inspection: 04/20/2020

Serial Number: 80-001243

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG931603 Exp: 11/12/2021
0.000	0.047	0.077	0.195	0.078
0.000	0.047	0.077	0.195	0.078
0.000	0.047	0.077	0.194	0.078
0.000	0.048	0.078	0.195	0.078
0.000	0.048	0.077	0.195	0.078
0.000	0.047	0.078	0.195	0.078
0.000	0.047	0.077	0.195	0.078
0.000	0.047	0.078	0.194	0.078
0.000	0.047	0.077	0.193	0.078
0.000	0.047	0.077	0.195	0.078

Standard Deviations	0.0004	0.0004	0.0006	0.0000
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

M4

BK 2020.04.
20
16:05:58
-04'00'

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT

Signature and Printed Name

04/20/2020
Date

stability checks

EDGEWATER P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001243
04/14/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:48
Control Test	0.047	08:49
Air Blank	0.000	08:49
Control Test	0.048	08:50
Air Blank	0.000	08:50
Control Test	0.047	08:51
Air Blank	0.000	08:52
Control Test Stats		
Average	0.0473	
Std Dev	0.0006	
Rel Std Dev(%)	1.2198	

EDGEWATER P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001243
04/14/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:59
Control Test	0.079	09:00
Air Blank	0.000	09:00
Control Test	0.078	09:01
Air Blank	0.000	09:01
Control Test	0.079	09:02
Air Blank	0.000	09:03
Control Test Stats		
Average	0.0787	
Std Dev	0.0006	
Rel Std Dev(%)	0.7339	

EDGEWATER P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001243
04/14/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:05
Control Test	0.199	09:05
Air Blank	0.000	09:06
Control Test	0.198	09:07
Air Blank	0.000	09:07
Control Test	0.199	09:08
Air Blank	0.000	09:08
Control Test Stats		
Average	0.1987	
Std Dev	0.0006	
Rel Std Dev(%)	0.2906	

Wet



Operator's Signature



Operator's Signature



Operator's Signature

EDGEWATER P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001243
04/14/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:10
Control Test	0.077	09:10
Air Blank	0.000	09:11
Control Test	0.078	09:11
Air Blank	0.000	09:11
Control Test	0.078	09:12
Air Blank	0.000	09:12
Control Test Stats		
Average	0.0777	
Std Dev	0.0006	
Rel Std Dev(%)	0.7434	

Dry



Operator's Signature

MH

2020.04.
20
BK 16:06:33
-04'00'

EDGEWATER P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001243
04/20/2020
Software: 8100.27

Flow Rate Calibration*****
1: Rate (Liters/min) = 5
SQRT(Diff)) = 7.211
2: Rate (Liters/min) = 15
SQRT(Diff)) = 11.617
3: Rate (Liters/min) = 30
SQRT(Diff)) = 20.320
Dependent Data Scale Factor = 100000 L/min
Independent Data Scale Factor = 256
Rounded Slope = 735
Rounded Intercept = -788371
Correlation = 0.99740

MH

2020.04
BK .20
16:07:0
9-04'00'

SP

Return Material Authorization

Ship to: ☒ CMI, Inc.

☐ Enforcement Electronics

Shipment to repair facility authorized by: Stephen Binz on 1/17/2020

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-001243

Bill To Address:

Edgewater PD

Ship to Address:

FDLE ATP Tallahassee

Reason for Return:

Multiple INT Detect Exception Messages

Please choose one of the following options:

- ☐ 1. I _____, authorize all repairs.
- ☐ 2. I _____, authorize repairs up to \$ _____.
- ☐ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Stephen Binz

Phone #: _____ Email: sbinz@CITYOFEDGEWATER.ORG

ATP Contact Name: Shayla Platt ATP Email: shaylplatt@fdle.state.fl.us

MY
BK
2020.04.
20
16:07:51
-04'00'



INSTRUMENT PROCESSING SHEET

Florida Department of
Law Enforcement

Agency Edgewater PD

S/N 80-001243

Date In 12/27/2019

DI Completion Date 1/9/20

Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake Performed By SP ☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes:		Quality Checks Performed By SP ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value 201 ☒ Flow Verification (L/s) Flow Column # ATP105 32 mm .148 (.139 - .169) 36 mm .164 (.156 - .190) 53 mm .246 (.228 - .278) 103 mm .515 (.447 - .547) ☒ Barometric Pressure Check Gauge ID # 26932 ☒ Stability Checks		Flow Calibration Performed By Flow Column # ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value ☐ Post Calibration Verification (L/s) Flow Column # 32 mm (.139 - .169) 36 mm (.156 - .190) 53 mm (.228 - .278) 103 mm (.447 - .547)																																									
Final Release Date		<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>SD1012</td> <td>201905A 5/14/2021</td> </tr> <tr> <td>0.080</td> <td>DR1279</td> <td>201905B 5/14/2021</td> </tr> <tr> <td>0.200</td> <td>SD1013</td> <td>201904D 4/30/2021</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG916501 6-14-21</td> </tr> </tbody> </table>		Simulator	Serial #	Lot #/Exp	0.050	SD1012	201905A 5/14/2021	0.080	DR1279	201905B 5/14/2021	0.200	SD1013	201904D 4/30/2021	0.080 DGS	N/A	AG916501 6-14-21	Maintenance Performed By ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other																										
Simulator	Serial #	Lot #/Exp																																											
0.050	SD1012	201905A 5/14/2021																																											
0.080	DR1279	201905B 5/14/2021																																											
0.200	SD1013	201904D 4/30/2021																																											
0.080 DGS	N/A	AG916501 6-14-21																																											
Calibration Adjustment Performed By Barometric Pressure Gauge ID #		Department Inspection Performed By SP Barometric Pressure ID# 26932 Gauge 1030 Instrument 1026 Mouth Alcohol Solution Lot # 2018-B Acetone Stock Solution Lot # 2019-A		Temperature Checks Performed By SP ☒ Lab Temp °C 21.9 External Digital Therm. ID#: 300505 SP ☒ 34°C ±.2 Serial #: G2408 SD1012 ☒ 34°C ±.2 Serial #: DR1279 ☒ 34°C ±.2 Serial #: SD1013																																									
<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>		Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>G2408</td> </tr> <tr> <td>Interferent</td> <td>G2882</td> </tr> <tr> <td>0.050</td> <td>SD1012</td> </tr> <tr> <td>0.080</td> <td>DR1279</td> </tr> <tr> <td>0.200</td> <td>SD1013</td> </tr> </tbody> </table>		Simulator	Serial Number	0.000	G2408	Interferent	G2882	0.050	SD1012	0.080	DR1279	0.200	SD1013	Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other Additional Stabilities	
Simulator	Serial Number	Lot Number	Expiration																																										
0.000		N/A	N/A																																										
0.040																																													
0.100																																													
0.200																																													
0.300																																													
0.080 DGS	N/A																																												
Simulator	Serial Number																																												
0.000	G2408																																												
Interferent	G2882																																												
0.050	SD1012																																												
0.080	DR1279																																												
0.200	SD1013																																												
Notes/Suggested Service:		☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use		Tech Review / Date Admin Review / Date																																									



ANAB
ACCREDITED
FORENSIC CALIBRATION
LABORATORY

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2729 Fort Knox Blvd.
Bldg. 2, Suite 1300
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001243, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001243</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>EDGEWATER P.D.</u>	0.050 g/210 L	0.004
Calibration Date:	<u>01/09/2020</u>	0.080 g/210 L	0.005
Calibration Time:	<u>12:14</u>	0.200 g/210 L	0.007
		0.080 g/210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

2020.04.
20
16:09:24
-04'00'

Shayla Platt

01/09/2020

Date

SHAYLA D. PLATT,
Department Inspector

FDLE/ATP Form 69 July 2018

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: EDGEWATER P.D.

Time of Inspection: 12:14

Date of Inspection: 01/09/2020

Serial Number: 80-001243

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG916501 Exp: 06/14/2021
0.000	0.050	0.079 / 0.079	0.198	0.081
0.000	0.050	0.079 / 0.080	0.200	0.081
0.000	0.050	0.077 / 0.081	0.200	0.081
0.000	0.050	0.079 / 0.081	0.201	0.081
0.000	0.047	0.084 / 0.080	0.200	0.081
0.000	0.047	INT / 0.080	0.200	0.081
0.000	0.047	/ 0.080	0.202	0.083
0.000	0.048	/ 0.081	0.201	0.083
0.000	0.048	/ 0.081	0.201	0.083
0.000	0.048	/ 0.081	0.202	0.084
Standard Deviations	0.0013	/ 0.0006	0.0011	0.0011

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0010 Number of Simulators Used: 5

Remarks:

08: Interferent Detect. *cleared Room and blew out simulator tubing prior to retesting. will conduct additional stabilities to ensure accuracy and reliability. SP*

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

Signature and Printed Name

SHAYLA D PLATT

01/09/2020
Date

MY

Stability Checks #80-001243

EDGEWATER P.O.
Intoxilyzer - Alcohol Analyzer
Model 8000
01/09/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:57
Control Test	0.078	09:57
Air Blank	0.000	09:58
Control Test	0.080	09:59
Air Blank	0.000	09:59
Control Test	0.080	10:00
Air Blank	0.000	10:00
Control Test Stats		
Average	0.0793	
Std Dev	0.0012	
Rel Std Dev(%)	1.4555	

SP
Operator's Signature

EDGEWATER P.O.
Intoxilyzer - Alcohol Analyzer
Model 8000
01/09/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:50
Control Test	0.048	09:51
Air Blank	0.000	09:51
Control Test	0.048	09:52
Air Blank	0.000	09:52
Control Test	0.050	09:53
Air Blank	0.000	09:54
Control Test Stats		
Average	0.0487	
Std Dev	0.0012	
Rel Std Dev(%)	2.3727	

SP

EDGEWATER P.O.
Intoxilyzer - Alcohol Analyzer
Model 8000
01/09/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:01
Control Test	0.199	10:02
Air Blank	0.000	10:03
Control Test	0.199	10:03
Air Blank	0.000	10:04
Control Test	0.200	10:05
Air Blank	0.000	10:05
Control Test Stats		
Average	0.1993	
Std Dev	0.0006	
Rel Std Dev(%)	0.2896	

SP
Operator's Signature

EDGEWATER P.O.
Intoxilyzer - Alcohol Analyzer
Model 8000
01/09/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:06
Control Test	0.082	10:07
Air Blank	0.000	10:07
Control Test	0.081	10:07
Air Blank	0.000	10:08
Control Test	0.082	10:08
Air Blank	0.000	10:09
Control Test Stats		
Average	0.0817	
Std Dev	0.0006	
Rel Std Dev(%)	0.7070	

SP

SP
Operator's Signature

Stability
Checks
#80-001243

EDGEWATER P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001243
01/09/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	13:55
Control Test	0.080	13:56
Air Blank	0.000	13:56
Control Test	0.080	13:57
Air Blank	0.000	13:58
Control Test	0.081	13:58
Air Blank	0.000	13:59
Control Test	0.080	14:00
Air Blank	0.000	14:00
Control Test	0.081	14:01
Air Blank	0.000	14:01
Control Test	0.080	14:02
Air Blank	0.000	14:03
Control Test	0.081	14:03
Air Blank	0.000	14:04
Control Test	0.081	14:04
Air Blank	0.000	14:05
Control Test	0.081	14:06
Air Blank	0.000	14:06
Control Test	0.079	14:07
Air Blank	0.000	14:07
Control Test	0.081	14:08
Air Blank	0.000	14:09
Control Test	0.080	14:09
Air Blank	0.000	14:10
Control Test	0.081	14:11
Air Blank	0.000	14:11
Control Test	0.080	14:12
Air Blank	0.000	14:12
Control Test	0.081	14:13
Air Blank	0.000	14:14
Control Test	INT*	14:14
Air Blank	PUR**	14:15
Air Blank	PUR**	14:16

*Interferent Detect
**Purge Fail

MH

2020.04.
20
BK 16:11:35
-04'00'

SP

Operator's Signature

Return Material Authorization

Ship to:



CMI, Inc.



Enforcement Electronics

Shipment to repair facility authorized by: Stephen Binz on 1/17/2020

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-001243

Bill To Address:

Edgewater PD

Ship to Address:

FDLE ATP Tallahassee

Reason for Return:

Multiple INT Detect Exception Messages

Please choose one of the following options:

- ☐ 1. I _____, authorize all repairs.
- ☐ 2. I _____, authorize repairs up to \$_____.
- ☐ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Stephen Binz

Phone #: _____ Email: sbinz@CITYOFEDGEWATER.ORG

ATP Contact Name: Shayla Platt ATP Email: shaylplatt@fdle.state.fl.us

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: EDGEWATER P.D.

Time of Inspection: 12:14

Date of Inspection: 01/09/2020

Serial Number: 80-001243
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG916501 Exp: 06/14/2021
0.000	0.050	0.079 / 0.079	0.198	0.081
0.000	0.050	0.079 / 0.080	0.200	0.081
0.000	0.050	0.077 / 0.081	0.200	0.081
0.000	0.050	0.079 / 0.081	0.201	0.081
0.000	0.047	0.084 / 0.080	0.200	0.081
0.000	0.047	INT / 0.080	0.200	0.081
0.000	0.047	/ 0.080	0.202	0.083
0.000	0.048	/ 0.081	0.201	0.083
0.000	0.048	/ 0.081	0.201	0.083
0.000	0.048	/ 0.081	0.202	0.084
Standard Deviations	0.0013	/ 0.0006	0.0011	0.0011

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0010 Number of Simulators Used: 5

Remarks:

08: Interferent Detect. *Cleared Room and blew out simulator tubing prior to retesting. Will conduct additional stabilities to ensure accuracy and reliability.SP*

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT
Signature and Printed Name

01/09/2020
Date

Stability
Checks
#80-001243

EDGEWATER P.D.
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001243
01/09/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	13:55
Control Test	0.080	13:56
Air Blank	0.000	13:56
Control Test	0.080	13:57
Air Blank	0.000	13:58
Control Test	0.081	13:58
Air Blank	0.000	13:59
Control Test	0.080	14:00
Air Blank	0.000	14:00
Control Test	0.081	14:01
Air Blank	0.000	14:01
Control Test	0.080	14:02
Air Blank	0.000	14:03
Control Test	0.081	14:03
Air Blank	0.000	14:04
Control Test	0.081	14:04
Air Blank	0.000	14:05
Control Test	0.081	14:06
Air Blank	0.000	14:06
Control Test	0.079	14:07
Air Blank	0.000	14:07
Control Test	0.081	14:08
Air Blank	0.000	14:09
Control Test	0.080	14:09
Air Blank	0.000	14:10
Control Test	0.081	14:11
Air Blank	0.000	14:11
Control Test	0.080	14:12
Air Blank	0.000	14:12
Control Test	0.081	14:13
Air Blank	0.000	14:14
Control Test	INT*	14:14
Air Blank	PUR**	14:15
Air Blank	PUR**	14:16

*Interferent Detect
**Purge Fail



Operator's Signature