



INSTRUMENT PROCESSING SHEET

Agency Crestview PD

S/N 80-005058

Florida Department of
Law Enforcement

Date In 11-19-2020 DI Completion Date 11-19-2020

☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake Performed By <u>IS</u> ☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: _____		Quality Checks Performed By <u>IS</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>232</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-105</u> 32 mm <u>0.160</u> (.139 - .169) 36 mm <u>0.175</u> (.156 - .190) 53 mm <u>0.250</u> (.228 - .278) 103 mm <u>0.527</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>30793</u> ☒ Stability Checks		Flow Calibration Performed By _____ Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																							
Final Release Date FDLE Alcohol Testing Program Digitally signed by FDLE Alcohol Testing Program Date: 2020.11.20 12:25:39 -05'00'		<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>MP5088</td> <td>201905A 05-14-2021</td> </tr> <tr> <td>0.080</td> <td>MP5086</td> <td>201905B 05-14-2021</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> <td>201904D 04-30-2021</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG011102 04-20-2022</td> </tr> </tbody> </table>		Simulator	Serial #	Lot #/Exp	0.050	MP5088	201905A 05-14-2021	0.080	MP5086	201905B 05-14-2021	0.200	MP5090	201904D 04-30-2021	0.080 DGS	N/A	AG011102 04-20-2022	Maintenance Performed By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ Temperature Checks Performed By <u>IS</u> ☒ Lab Temp °C <u>21.16</u> External Digital Therm. ID#: <u>300503</u> ☒ 34°C +/- .2 Serial #: <u>MP5088</u> ☒ 34°C +/- .2 Serial #: <u>MP5086</u> ☒ 34°C +/- .2 Serial #: <u>MP5090</u>																								
Simulator	Serial #	Lot #/Exp																																									
0.050	MP5088	201905A 05-14-2021																																									
0.080	MP5086	201905B 05-14-2021																																									
0.200	MP5090	201904D 04-30-2021																																									
0.080 DGS	N/A	AG011102 04-20-2022																																									
Calibration Adjustment Performed By <u>IS</u> Barometric Pressure Gauge <u>1032</u> ID # <u>28421</u>		Department Inspection Performed By <u>IS</u> Barometric Pressure ID# <u>30793</u> Gauge <u>1030</u> Instrument <u>1031</u> Mouth Alcohol Solution Lot # <u>2019-B</u> Acetone Stock Solution Lot # <u>2019-A</u>																																									
<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>MP5091</td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td>MP5082</td> <td>20060</td> <td>02-10-2022</td> </tr> <tr> <td>0.100</td> <td>MP5083</td> <td>20420</td> <td>09-22-2022</td> </tr> <tr> <td>0.200</td> <td>MP5084</td> <td>20160</td> <td>03-18-2022</td> </tr> <tr> <td>0.300</td> <td>MP5085</td> <td>20030</td> <td>01-21-2022</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>08819080A1</td> <td>06-05-2021</td> </tr> </tbody> </table>		Simulator	Serial Number	Lot Number	Expiration	0.000	MP5091	N/A	N/A	0.040	MP5082	20060	02-10-2022	0.100	MP5083	20420	09-22-2022	0.200	MP5084	20160	03-18-2022	0.300	MP5085	20030	01-21-2022	0.080 DGS	N/A	08819080A1	06-05-2021	<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>G11621</td> </tr> <tr> <td>Interferent</td> <td>MP5087</td> </tr> <tr> <td>0.050</td> <td>MP5088</td> </tr> <tr> <td>0.080</td> <td>MP5086</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> </tr> </tbody> </table>		Simulator	Serial Number	0.000	G11621	Interferent	MP5087	0.050	MP5088	0.080	MP5086	0.200	MP5090
Simulator	Serial Number	Lot Number	Expiration																																								
0.000	MP5091	N/A	N/A																																								
0.040	MP5082	20060	02-10-2022																																								
0.100	MP5083	20420	09-22-2022																																								
0.200	MP5084	20160	03-18-2022																																								
0.300	MP5085	20030	01-21-2022																																								
0.080 DGS	N/A	08819080A1	06-05-2021																																								
Simulator	Serial Number																																										
0.000	G11621																																										
Interferent	MP5087																																										
0.050	MP5088																																										
0.080	MP5086																																										
0.200	MP5090																																										
☒ Post Calibration Adjustment Stability Checks		<table border="1"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>MP5088</td> <td>201905A</td> <td>05-14-2021</td> </tr> <tr> <td>0.080</td> <td>MP5086</td> <td>201905B</td> <td>05-14-2021</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> <td>201904D</td> <td>04-30-2021</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG011102</td> <td>04-20-2022</td> </tr> </tbody> </table>				Simulator	Serial Number	Lot Number	Expiration	0.050	MP5088	201905A	05-14-2021	0.080	MP5086	201905B	05-14-2021	0.200	MP5090	201904D	04-30-2021	0.080 DGS	N/A	AG011102	04-20-2022																		
Simulator	Serial Number	Lot Number	Expiration																																								
0.050	MP5088	201905A	05-14-2021																																								
0.080	MP5086	201905B	05-14-2021																																								
0.200	MP5090	201904D	04-30-2021																																								
0.080 DGS	N/A	AG011102	04-20-2022																																								
Notes/Suggested Service: <u>Pressure on instrument differed from barometric pressure gauge by more than 1%, performed optical bench cal adjust. Post stability 0.05 not nominal, tightened simulator lid and reran post stability, 0.05 nominal after tightening.</u>		☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Tech Review / Date <u>2020.11.19 14:09:58 -05'00'</u> Admin Review / Date <u>2020.11.20 12:24:41</u>																																									

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: CRESTVIEW PD
Time of Inspection: 12:05

Date of Inspection: 11/19/2020

Serial Number: 80-005058
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG011102 Exp: 04/20/2022
0.000	0.048	0.078	0.199	0.079
0.000	0.048	0.078	0.199	0.079
0.000	0.048	0.078	0.199	0.080
0.000	0.048	0.078	0.199	0.079
0.000	0.048	0.078	0.199	0.079
0.000	0.048	0.078	0.199	0.079
0.000	0.048	0.078	0.200	0.079
0.000	0.048	0.078	0.199	0.079
0.000	0.048	0.078	0.198	0.079
0.000	0.049	0.078	0.199	0.079
Standard Deviations	0.0003	0.0000	0.0004	0.0003

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL X SOTO

Signature and Printed Name

11/19/2020
Date

*No middle initial, required
to continue inspection*

MX

2020.11.20
12:23:49 -05'00'

Stability checks

CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:53
Control Test	0.050	08:54
Air Blank	0.000	08:55
Control Test	0.048	08:55
Air Blank	0.000	08:56
Control Test	0.048	08:56
Air Blank	0.000	08:57
Control Test Stats		
Average	0.0487	
Std Dev	0.0012	
Rel Std Dev(%)	2.3727	


CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:07
Control Test	0.079	09:08
Air Blank	0.000	09:09
Control Test	0.079	09:09
Air Blank	0.000	09:10
Control Test	0.079	09:10
Air Blank	0.000	09:11
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:12
Control Test	0.198	09:12
Air Blank	0.000	09:13
Control Test	0.198	09:14
Air Blank	0.000	09:14
Control Test	0.198	09:15
Air Blank	0.000	09:15
Control Test Stats		
Average	0.1980	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

wet

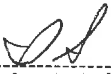

Operator's Signature



Operator's Signature




Operator's Signature

CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:25
Control Test	0.080	09:25
Air Blank	0.000	09:25
Control Test	0.080	09:26
Air Blank	0.000	09:26
Control Test	0.080	09:26
Air Blank	0.000	09:27
Control Test Stats		
Average	0.0800	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Dry

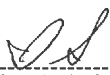

Operator's Signature



2020.11
.20
12:22:5
4 -05'00'

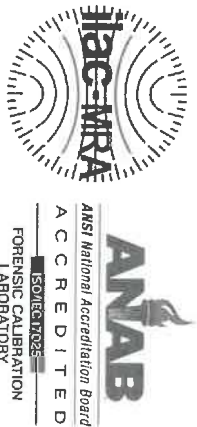

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2729 Fort Knox Blvd.
Bldg. 2, Suite 1300
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-005058, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-005058	UNCERTAINTY* $\pm$
Owning Agency:	CRESTVIEW PD	0.050 g/ 210 L
Calibration Date:	11/19/2020	0.080 g/ 210 L
Calibration Time:	12:05	0.200 g/ 210 L
		0.080 g/ 210 L Dry Gas Control
		0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration. *Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 April 2020
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Date

11/19/2020

ISRAEL SOTO,
Department Inspector

Page 1 of 1

Optical Bench cal. Adj.

Ed

CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020 09:28:51

Auto Calibration
Max Power Res Value = 36
Auto Range Res Value = 19

Sol Value = 0.000 g/210L ***
Fit value = 0.0000 mg/l %%%
Samples Taken = 4, Discarded = 1
Sum Io = 12607, Sum Io = 13344

Channel 1
Sample % Abs (% Abs Ref)
Sample #1 = 0.0650 (0.0010)
Sample #2 = 0.0820 (0.0360)
Sample #3 = 0.0860 (0.0640)
Sample #4 = 0.0750 (0.1050)
Avg % Abs = 0.0813 (0.0587)
STD DEV = 0.0050 (0.0352)
REL STD DEV = 6.188 (51.310)

Channel 2
Sample % Abs (% Abs Ref)
Sample #1 = 0.1280 (0.0080)
Sample #2 = 0.1540 (0.0120)
Sample #3 = 0.1680 (0.0190)
Sample #4 = 0.1550 (0.0330)
Avg % Abs = 0.1593 (0.0207)
STD DEV = 0.0076 (0.0116)
REL STD DEV = 4.752 (56.082)

Sol Value = 0.040 g/210L ***
Fit value = 0.1905 mg/l %%%
Samples Taken = 4, Discarded = 1
Sum Io = 12591, Sum Io = 13336

Channel 1
Sample % Abs (% Abs Ref)
Sample #1 = 0.8450 (-0.0150)
Sample #2 = 0.8410 (0.0000)
Sample #3 = 0.8530 (0.0050)
Sample #4 = 0.7910 (0.0390)
Avg % Abs = 0.8283 (0.0160)
STD DEV = 0.0329 (0.0204)
REL STD DEV = 3.970 (127.629)

Channel 2
Sample % Abs (% Abs Ref)
Sample #1 = 1.6050 (-0.0060)
Sample #2 = 1.5990 (-0.0100)
Sample #3 = 1.6110 (-0.0070)
Sample #4 = 1.5710 (0.0260)
Avg % Abs = 1.5937 (0.0030)
STD DEV = 0.0205 (0.0200)
REL STD DEV = 1.288 (665.833)

Sol Value = 0.100 g/210L ***
Fit value = 0.4762 mg/l %%%
Samples Taken = 4, Discarded = 1
Sum Io = 12585, Sum Io = 13333

Channel 1
Sample % Abs (% Abs Ref)
Sample #1 = 1.9380 (-0.0080)
Sample #2 = 1.9420 (0.0040)
Sample #3 = 1.8930 (0.0370)
Sample #4 = 1.9430 (0.0190)
Avg % Abs = 1.9260 (0.0200)
STD DEV = 0.0286 (0.0165)
REL STD DEV = 1.484 (82.614)

Channel 2
Sample % Abs (% Abs Ref)
Sample #1 = 3.7030 (-0.0060)
Sample #2 = 3.7010 (0.0040)
Sample #3 = 3.6690 (0.0180)
Sample #4 = 3.6950 (0.0070)
Avg % Abs = 3.6883 (0.0097)
STD DEV = 0.0170 (0.0074)
REL STD DEV = 0.461 (76.253)

Sol Value = 0.200 g/210L ***
Fit value = 0.9524 mg/l %%%
Samples Taken = 4, Discarded = 1
Sum Io = 12582, Sum Io = 13332

Channel 1
Sample % Abs (% Abs Ref)
Sample #1 = 3.7650 (-0.0060)
Sample #2 = 3.7580 (0.0170)
Sample #3 = 3.7150 (0.0440)
Sample #4 = 3.7470 (0.0280)
Avg % Abs = 3.7400 (0.0297)
STD DEV = 0.0223 (0.0135)
REL STD DEV = 0.597 (45.765)

Channel 2
Sample % Abs (% Abs Ref)
Sample #1 = 7.0820 (-0.0020)
Sample #2 = 7.0790 (0.0110)
Sample #3 = 7.0410 (0.0320)
Sample #4 = 7.0920 (-0.0040)
Avg % Abs = 7.0707 (0.0130)
STD DEV = 0.0265 (0.0181)
REL STD DEV = 0.375 (139.101)

Sol Value = 0.300 g/210L ***
Fit value = 1.4286 mg/l %%%
Samples Taken = 4, Discarded = 1
Sum Io = 12579, Sum Io = 13332

Channel 1
Sample % Abs (% Abs Ref)
Sample #1 = 5.4530 (-0.0060)
Sample #2 = 5.4060 (0.0310)
Sample #3 = 5.4240 (0.0130)
Sample #4 = 5.4450 (0.0080)
Avg % Abs = 5.4250 (0.0173)
STD DEV = 0.0195 (0.0121)
REL STD DEV = 0.360 (69.789)

Channel 2
Sample % Abs (% Abs Ref)
Sample #1 = 10.2010 (-0.0230)
Sample #2 = 10.1520 (0.0360)
Sample #3 = 10.1670 (0.0260)
Sample #4 = 10.1730 (0.0200)
Avg % Abs = 10.1640 (0.0273)
STD DEV = 0.0108 (0.0081)
REL STD DEV = 0.106 (29.572)

***** AUTO CAL DATA *****
Channel 1
Sol Val = 0.0000 mg/l or 0.000 g/210L
% Abs = 0.081
Std Dev = 0.01 Rel Std Dev = 6.19
Sol Val = 0.1905 mg/l or 0.040 g/210L
% Abs = 0.628
Std Dev = 0.03 Rel Std Dev = 3.97
Sol Val = 0.4762 mg/l or 0.100 g/210L
% Abs = 1.926
Std Dev = 0.03 Rel Std Dev = 1.48
Sol Val = 0.9524 mg/l or 0.200 g/210L
% Abs = 3.740
Std Dev = 0.02 Rel Std Dev = 0.60
Sol Val = 1.4286 mg/l or 0.300 g/210L
% Abs = 5.425
Std Dev = 0.02 Rel Std Dev = 0.36
Zero Order Coef = -183.51
First Order Coef = 2483.02
Second Order Coef = 33.37
Standard Deviation = 33.016754

Channel 2
Sol Val = 0.0000 mg/l or 0.000 g/210L
% Abs = 0.159
Std Dev = 0.01 Rel Std Dev = 4.75
Sol Val = 0.1905 mg/l or 0.040 g/210L
% Abs = 1.594
Std Dev = 0.02 Rel Std Dev = 1.29
Sol Val = 0.4762 mg/l or 0.100 g/210L
% Abs = 3.688
Std Dev = 0.02 Rel Std Dev = 0.46
Sol Val = 0.9524 mg/l or 0.200 g/210L
% Abs = 7.071
Std Dev = 0.03 Rel Std Dev = 0.37
Sol Val = 1.4286 mg/l or 0.300 g/210L
% Abs = 10.164
Std Dev = 0.01 Rel Std Dev = 0.11
Zero Order Coef = -195.09
First Order Coef = 1281.05
Second Order Coef = 13.88
Standard Deviation = 30.644089

2020.11.20
12:20:50 -05'00'

MX

Optical bench Cal Adj.

DA

Solution Stats Quadratic Fit Chan 1		
Act	Fit	Residual
g/210L	g/210L	g/210L
0.000	0.000	-0.0004
0.040	0.040	0.0002
0.100	0.099	0.0008
0.200	0.201	-0.0010
0.300	0.300	0.0004

Solution Stats Quadratic Fit Chan 2		
Act	Fit	Residual
g/210L	g/210L	g/210L
0.000	0.000	-0.0004
0.040	0.040	0.0003
0.100	0.099	0.0007
0.200	0.201	-0.0009
0.300	0.300	0.0003

Sol Value = 0.080 g/210L ***
 Fit value = 0.3810 mg/l %%%
 Samples Taken = 4, Discarded = 1
 ***** CHANNEL 1
 Sample #1 = 3348.00
 Sample #2 = 3367.00
 Sample #3 = 3486.00
 Sample #4 = 3316.00
 Average Result = 3389.6667
 STD DEV = 87.2372
 REL STD DEV = 2.574

 ***** CHANNEL 2
 Sample #1 = 3363.00
 Sample #2 = 3371.00
 Sample #3 = 3381.00
 Sample #4 = 3319.00
 Average Result = 3357.0000
 STD DEV = 33.2866
 REL STD DEV = 0.992

 Dry Gas H2O Adjust Results *****
 Barometric Pressure = 1031
 3 um H2O Adjust (mg/l*10,000) = 420
 9 um H2O Adjust (mg/l*10,000) = 452
 **** AUTO CAL PASS

MX


SR

2020.11.20
 12:20:08 -05'00'

Post stability checks

CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:16
Control Test	0.046	10:16
Air Blank	0.000	10:17
Control Test	0.048	10:17
Air Blank	0.000	10:18
Control Test	0.048	10:19
Air Blank	0.000	10:19
Control Test Stats		
Average	0.0473	
Std Dev	0.0012	
Rel Std Dev(%)	2.4395	

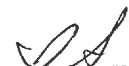

Operator's Signature

CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:20
Control Test	0.079	10:21
Air Blank	0.000	10:21
Control Test	0.079	10:22
Air Blank	0.000	10:23
Control Test	0.079	10:23
Air Blank	0.000	10:24
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

wet



Operator's Signature

CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:25
Control Test	0.048	10:25
Air Blank	0.000	10:26
Control Test	0.048	10:27
Air Blank	0.000	10:27
Control Test	0.048	10:28
Air Blank	0.000	10:28
Control Test Stats		
Average	0.0480	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

0.05 after
tightening lid



Operator's Signature

CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:29
Control Test	0.200	10:30
Air Blank	0.000	10:30
Control Test	0.199	10:31
Air Blank	0.000	10:32
Control Test	0.199	10:32
Air Blank	0.000	10:33
Control Test Stats		
Average	0.1993	
Std Dev	0.0006	
Rel Std Dev(%)	0.2896	



Operator's Signature

CRESTVIEW PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005058
11/19/2020
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:34
Control Test	0.079	10:34
Air Blank	0.000	10:35
Control Test	0.080	10:35
Air Blank	0.000	10:36
Control Test	0.080	10:36
Air Blank	0.000	10:36
Control Test Stats		
Average	0.0797	
Std Dev	0.0006	
Rel Std Dev(%)	0.7247	

Dry



Operator's Signature



2020.11.
20
12:17:20
-05'00'