





INSTRUMENT PROCESSING SHEET

Agency FHP Panama CityS/N 80-006760

Florida Department of Law Enforcement

Date In 05/05/2020 DI Completion Date 5/15/20 Ship P/U H/D CMI EE

| Intake Performed By <u>RAW</u> <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ _____ _____ | Quality Checks Performed By <u>DA</u> <input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>212</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP-103</u> 32 mm <u>0.160</u> (.139 - .169) 36 mm <u>0.175</u> (.156 - .190) 53 mm <u>0.246</u> (.228 - .278) 103 mm <u>0.511</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>28427</u> <input checked="" type="checkbox"/> Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td><u>MP5088</u></td> <td><u>201905A</u> <u>05-17-2021</u></td> </tr> <tr> <td>0.080</td> <td><u>MP5089</u></td> <td><u>201905B</u> <u>05-17-2021</u></td> </tr> <tr> <td>0.200</td> <td><u>MP5090</u></td> <td><u>201907D</u> <u>04-30-2021</u></td> </tr> <tr> <td>0.080 DGS</td> <td><u>N/A</u></td> <td><u>AG031603</u> <u>11-12-2021</u></td> </tr> </tbody> </table> | Simulator | Serial # | Lot #/Exp | 0.050 | <u>MP5088</u> | <u>201905A</u> <u>05-17-2021</u> | 0.080 | <u>MP5089</u> | <u>201905B</u> <u>05-17-2021</u> | 0.200 | <u>MP5090</u> | <u>201907D</u> <u>04-30-2021</u> | 0.080 DGS | <u>N/A</u> | <u>AG031603</u> <u>11-12-2021</u> | Flow Calibration Performed By _____ Flow Column # _____ <input type="checkbox"/> 5L/min - 17mm <input type="checkbox"/> 15L/min - 53mm <input type="checkbox"/> 30L/min - 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance Performed By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____ Temperature Checks Performed By <u>SP</u> <input checked="" type="checkbox"/> Lab Temp °C <u>21.24</u> External Digital Therm. ID#: <u>300502</u> <input type="checkbox"/> 34°C +/-2 Serial #: <u>MP5088</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5089</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP6090</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
|---|---|--------------------------------------|---------------|---------------|---------------|---------------|-------------------------------------|-------|---------------|-------------------------------------|---------------|---------------|-------------------------------------|--|------------|--------------------------------------|---|--|-------|--|--|--|-------|--|--|--|-----------|------------|--|--|-----------|---------------|------------|------------|-------|--|--|--|-------|--|--|--|-------|--|--|--|-----------|------------|--|--|
| Simulator | Serial # | Lot #/Exp | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.050 | <u>MP5088</u> | <u>201905A</u> <u>05-17-2021</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.080 | <u>MP5089</u> | <u>201905B</u> <u>05-17-2021</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.200 | <u>MP5090</u> | <u>201907D</u> <u>04-30-2021</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.080 DGS | <u>N/A</u> | <u>AG031603</u> <u>11-12-2021</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Final Release Date _____ _____ | Calibration Adjustment Performed By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td><u>N/A</u></td> <td></td> <td></td> </tr> </tbody> </table> <input type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td><u>N/A</u></td> <td></td> <td></td> </tr> </tbody> </table> Notes/Suggested Service: _____ _____ _____ _____ _____ | | Simulator | Serial Number | Lot Number | Expiration | 0.000 | | N/A | N/A | 0.040 | | | | 0.100 | | | | 0.200 | | | | 0.300 | | | | 0.080 DGS | <u>N/A</u> | | | Simulator | Serial Number | Lot Number | Expiration | 0.050 | | | | 0.080 | | | | 0.200 | | | | 0.080 DGS | <u>N/A</u> | | |
| Simulator | Serial Number | Lot Number | Expiration | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.000 | | N/A | N/A | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.040 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.100 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.200 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.300 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.080 DGS | <u>N/A</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Simulator | Serial Number | Lot Number | Expiration | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.050 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.080 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.200 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.080 DGS | <u>N/A</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Department Inspection Performed By <u>SP</u> Barometric Pressure ID# <u>26932</u> Gauge <u>1017</u> Instrument <u>1017</u> Mouth Alcohol Solution Lot # <u>2019-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td><u>MP5086</u></td> </tr> <tr> <td>Interferent</td> <td><u>MP5087</u></td> </tr> <tr> <td>0.050</td> <td><u>MP5088</u></td> </tr> <tr> <td>0.080</td> <td><u>MP5089</u></td> </tr> <tr> <td>0.200</td> <td><u>MP5090</u></td> </tr> </tbody> </table> Attachments <input checked="" type="checkbox"/> Form 41 <input checked="" type="checkbox"/> Stability Checks <input checked="" type="checkbox"/> Calibration Certificate <input type="checkbox"/> Calibration Adjustment <input type="checkbox"/> Post-Stability Checks <input type="checkbox"/> Flow Calibration <input checked="" type="checkbox"/> Form 40 <input type="checkbox"/> Other _____ | | Simulator | Serial Number | 0.000 | <u>MP5086</u> | Interferent | <u>MP5087</u> | 0.050 | <u>MP5088</u> | 0.080 | <u>MP5089</u> | 0.200 | <u>MP5090</u> | <input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use Tech Review / Date _____ Admin Review / Date _____ | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Simulator | Serial Number | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.000 | <u>MP5086</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Interferent | <u>MP5087</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.050 | <u>MP5088</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.080 | <u>MP5089</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 0.200 | <u>MP5090</u> | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |

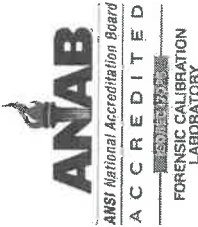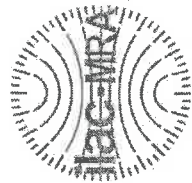

Florida Department of Law Enforcement
 Alcohol Testing Program
 2729 Fort Knox Blvd.
 Bldg. 2, Suite 1300
 Tallahassee, FL 32308

Calibration Certificate

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006760, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

| | | | |
|-------------------|-------------------|--------------------------------|-------|
| Serial Number: | <u>80-006760</u> | UNCERTAINTY* ± | |
| Owning Agency: | <u>FHP</u> | 0.050 g/ 210 L | 0.004 |
| Calibration Date: | <u>05/15/2020</u> | 0.080 g/ 210 L | 0.005 |
| Calibration Time: | <u>10:57</u> | 0.200 g/ 210 L | 0.007 |
| | | 0.080 g/ 210 L Dry Gas Control | 0.005 |

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
 *Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).
 The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
 This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

05/15/2020 Date
Shayla Platt
 SHAYLA D PLATT,
 Department Inspector

FDLE/ATP Form 69 April 2020
 Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: FHP

Time of Inspection: 10:57

Date of Inspection: 05/15/2020

Serial Number: 80-006760

Software: 8100.27

| Check or Test | YES | NO | Check or Test | YES | NO |
|---|-----|----|--|-----|----|
| Diagnostic Check (Pre-Inspection): OK | Yes | | Date and/or Time Adjusted | | No |
| Minimum Sample Volume Check: OK | Yes | | Barometric Pressure Sensor Check: OK | Yes | |
| Alcohol Free Subject Test: 0.000 | Yes | | Mouth Alcohol Test: Slope Not Met | Yes | |
| Interferent Detect Test: Interferent Detect | Yes | | Diagnostic Check (Post-Inspection): OK | Yes | |

| Alcohol Free Test (g/210L) | 0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021 | 0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021 | 0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021 | 0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG931603 Exp: 11/12/2021 |
|----------------------------|---|---|---|---|
| 0.000 | 0.048 | 0.079 | 0.197 | 0.078 |
| 0.000 | 0.048 | 0.079 | 0.198 | 0.078 |
| 0.000 | 0.048 | 0.078 | 0.198 | 0.078 |
| 0.000 | 0.048 | 0.078 | 0.198 | 0.078 |
| 0.000 | 0.048 | 0.079 | 0.198 | 0.078 |
| 0.000 | 0.048 | 0.078 | 0.198 | 0.078 |
| 0.000 | 0.048 | 0.079 | 0.198 | 0.078 |
| 0.000 | 0.049 | 0.079 | 0.198 | 0.078 |
| 0.000 | 0.049 | 0.079 | 0.197 | 0.078 |
| 0.000 | 0.048 | 0.079 | 0.197 | 0.078 |

| | | | | |
|---------------------|--------|--------|--------|--------|
| Standard Deviations | 0.0004 | 0.0004 | 0.0004 | 0.0000 |
|---------------------|--------|--------|--------|--------|

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT

Signature and Printed Name

05/15/2020
Date

stability checks

FHP
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006760
05/12/2020
Software: 8100.27

| Test | g/210L | Time |
|--------------------|--------|-------|
| Air Blank | 0.000 | 07:29 |
| Control Test | 0.049 | 07:29 |
| Air Blank | 0.000 | 07:30 |
| Control Test | 0.049 | 07:31 |
| Air Blank | 0.000 | 07:31 |
| Control Test | 0.048 | 07:32 |
| Air Blank | 0.000 | 07:32 |
| Control Test Stats | | |
| Average | 0.0487 | |
| Std Dev | 0.0006 | |
| Rel Std Dev(%) | 1.1863 | |



Operator's Signature

FHP
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006760
05/12/2020
Software: 8100.27

| Test | g/210L | Time |
|--------------------|--------|-------|
| Air Blank | 0.000 | 07:33 |
| Control Test | 0.079 | 07:34 |
| Air Blank | 0.000 | 07:35 |
| Control Test | 0.079 | 07:35 |
| Air Blank | 0.000 | 07:36 |
| Control Test | 0.079 | 07:37 |
| Air Blank | 0.000 | 07:37 |
| Control Test Stats | | |
| Average | 0.0790 | |
| Std Dev | 0.0000 | |
| Rel Std Dev(%) | 0.0000 | |

wet



Operator's Signature

FHP
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006760
05/12/2020
Software: 8100.27

| Test | g/210L | Time |
|--------------------|--------|-------|
| Air Blank | 0.000 | 07:39 |
| Control Test | 0.197 | 07:40 |
| Air Blank | 0.000 | 07:40 |
| Control Test | 0.198 | 07:41 |
| Air Blank | 0.000 | 07:41 |
| Control Test | 0.198 | 07:42 |
| Air Blank | 0.000 | 07:43 |
| Control Test Stats | | |
| Average | 0.1977 | |
| Std Dev | 0.0006 | |
| Rel Std Dev(%) | 0.2921 | |



Operator's Signature

FHP
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006760
05/12/2020
Software: 8100.27

| Test | g/210L | Time |
|--------------------|--------|-------|
| Air Blank | 0.000 | 07:45 |
| Control Test | 0.081 | 07:45 |
| Air Blank | 0.000 | 07:45 |
| Control Test | 0.081 | 07:46 |
| Air Blank | 0.000 | 07:46 |
| Control Test | 0.080 | 07:47 |
| Air Blank | 0.000 | 07:47 |
| Control Test Stats | | |
| Average | 0.0807 | |
| Std Dev | 0.0006 | |
| Rel Std Dev(%) | 0.7157 | |

Dry



Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: FHP

Time of Inspection: 06:47

Date of Inspection: 05/12/2020

Serial Number: 80-006760

Software: 8100.27

| Check or Test | YES | NO |
|---|-----|----|
| Date and/or Time Adjusted | | No |
| Diagnostic Check (Pre-Inspection): OK | | No |
| Alcohol Free Subject Test: 0.000 | | No |
| Mouth Alcohol Test: Slope Not Met | | No |
| Interferent Detect Test: Interferent Detect | | No |
| Diagnostic Check (Post-Inspection): OK | | No |

| Alcohol Free Test (g/210L) | 0.05g/210L Test (g/210L) Lot#: _____ Exp: _____ | 0.08g/210L Test (g/210L) Lot#: _____ Exp: _____ | 0.20g/210L Test (g/210L) Lot#: _____ Exp: _____ | 0.08 g/210L Dry Gas Std Test (g/210L) Lot#: _____ Exp: _____ |
|----------------------------|---|---|---|--|
| | | | | |
| | | | | |
| | | | | |

Number of Simulators Used: _____

Remarks:

BYPASSED AI TO OPERATE INSTRUMENT


 N/A compliancy not determined

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL SOTO

Signature and Printed Name

05/12/2020
Date