



# INSTRUMENT PROCESSING SHEET

Agency Fort Lauderdale PDS/N 80-007009

Florida Department of Law Enforcement

Date In 4/10/2020DI Completion Date 4/14/2020 Ship  P/U  H/D  CMI  EE

<b>Intake</b> Performed By <u>MH</u> <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input checked="" type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ _____ _____	<b>Quality Checks</b> Performed By <u>MH</u> <input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>211</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP 101</u> 32 mm <u>0.171</u> (.139 - .169) 36 mm <u>0.191</u> (.156 - .190) 53 mm <u>0.253</u> (.228 - .278) 103 mm <u>0.519</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>68639</u> <input checked="" type="checkbox"/> Stability Checks <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>MP4863</td> <td>201905A 05/14/2021</td> </tr> <tr> <td>0.080</td> <td>MP4864</td> <td>201905B 05/14/2021</td> </tr> <tr> <td>0.200</td> <td>MP5097</td> <td>201904D 04/30/2021</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG931603 11/12/2021</td> </tr> </tbody> </table>	Simulator	Serial #	Lot #/Exp	0.050	MP4863	201905A 05/14/2021	0.080	MP4864	201905B 05/14/2021	0.200	MP5097	201904D 04/30/2021	0.080 DGS	N/A	AG931603 11/12/2021	<b>Flow Calibration</b> Performed By <u>MH</u> Flow Column # <u>ATP 104</u> <input checked="" type="checkbox"/> 5L/min - 17mm <input checked="" type="checkbox"/> 15L/min - 53mm <input checked="" type="checkbox"/> 30L/min - 103mm <input checked="" type="checkbox"/> R-Value <u>211</u> <input checked="" type="checkbox"/> Post Calibration Verification (L/s) Flow Column # <u>ATP 101</u> 32 mm <u>0.152</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.238</u> (.228 - .278) 103 mm <u>0.507</u> (.447 - .547)																																												
Simulator	Serial #	Lot #/Exp																																																											
0.050	MP4863	201905A 05/14/2021																																																											
0.080	MP4864	201905B 05/14/2021																																																											
0.200	MP5097	201904D 04/30/2021																																																											
0.080 DGS	N/A	AG931603 11/12/2021																																																											
<b>Final Release Date</b> <div style="text-align: center; font-size: 1.2em; font-weight: bold;">FDLE</div> <div style="text-align: center; margin-top: 10px;">APR 30 2020</div> <div style="text-align: center; margin-top: 10px;">Alcohol Testing Program</div>	<b>Maintenance</b> Performed By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____	<b>Temperature Checks</b> Performed By <u>MH</u> <input checked="" type="checkbox"/> Lab Temp °C <u>22.84</u> External Digital Therm. ID#: <u>300504</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP4863</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP4864</u> <input checked="" type="checkbox"/> 34°C +/-2 Serial #: <u>MP5097</u>																																																											
<b>Calibration Adjustment</b> Performed By _____ Barometric Pressure Gauge ID # _____ <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> <input type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>	Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			<b>Department Inspection</b> Performed By <u>MH</u> Barometric Pressure ID# <u>28663</u> Gauge <u>1017</u> Instrument <u>1016</u> Mouth Alcohol Solution Lot # <u>2019-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1" style="width:100%; border-collapse: collapse; margin-top: 5px;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>SD1014</td> </tr> <tr> <td>Interferent</td> <td>SD1015</td> </tr> <tr> <td>0.050</td> <td>MP4863</td> </tr> <tr> <td>0.080</td> <td>MP4864</td> </tr> <tr> <td>0.200</td> <td>MP5097</td> </tr> </tbody> </table>	Simulator	Serial Number	0.000	SD1014	Interferent	SD1015	0.050	MP4863	0.080	MP4864	0.200	MP5097
Simulator	Serial Number	Lot Number	Expiration																																																										
0.000		N/A	N/A																																																										
0.040																																																													
0.100																																																													
0.200																																																													
0.300																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number	Lot Number	Expiration																																																										
0.050																																																													
0.080																																																													
0.200																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number																																																												
0.000	SD1014																																																												
Interferent	SD1015																																																												
0.050	MP4863																																																												
0.080	MP4864																																																												
0.200	MP5097																																																												
Notes/Suggested Service: <u>E-mailed</u> <div style="margin-top: 10px;"> <input checked="" type="checkbox"/> <span style="color: red; font-size: 1.2em; font-weight: bold;">APPROVED</span> <u>04/15/2020</u> </div>	<b>Attachments</b> <input checked="" type="checkbox"/> Form 41 <input type="checkbox"/> Post-Stability Checks <input checked="" type="checkbox"/> Stability Checks <input checked="" type="checkbox"/> Flow Calibration <input checked="" type="checkbox"/> Calibration Certificate <input checked="" type="checkbox"/> Form 40 <input type="checkbox"/> Calibration Adjustment <input type="checkbox"/> Other _____																																																												
_____ _____ _____	<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use <div style="margin-top: 10px;"> <u>SP 4/22/20</u> <u>Brett K...</u> <u>4/23/2020</u>            Tech Review / Date Admin Review / Date         </div>																																																												

# Florida Department of Law Enforcement Alcohol Testing Program

## DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: FORT LAUDERDALE PD  
Time of Inspection: 14:16

Date of Inspection: 04/14/2020

Serial Number: 80-007009  
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG931603 Exp: 11/12/2021
0.000	0.049	0.079	0.198	0.080
0.000	0.049	0.079	0.199	0.079
0.000	0.049	0.079	0.199	0.080
0.000	0.049	0.079	0.199	0.080
0.000	0.049	0.079	0.200	0.080
0.000	0.049	0.079	0.199	0.080
0.000	0.049	0.080	0.199	0.079
0.000	0.049	0.080	0.200	0.079
0.000	0.049	0.080	0.200	0.079
0.000	0.050	0.079	0.200	0.079

Standard Deviations	0.0003	0.0004	0.0006	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

SP  
BSK  
4/23/2020

The above instrument complies (  ) does not comply (  ) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Michael D Haughey MICHAEL D HAUGHEY  
Signature and Printed Name

04/14/2020  
Date

TYPE OF TEST	SERIAL NUMBER	AGENCY	DATE	PERFORMED BY
Stabilities	80-007009	Fort Lauderdale PD	04/14/2020	MXK

0.05g/210L	0.08g/210L	0.20g/210L	0.077 to 0.083	0.083g/210L	0.194 to 0.206	0.077 to 0.083
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>

g/210L	g/210L	g/210L	g/210L	g/210L	g/210L																																																																																																																																														
0.047 to 0.053	0.077 to 0.083	0.20g/210L	0.077 to 0.083	0.083g/210L	0.194 to 0.206																																																																																																																																														
<p>FORT LAUDERDALE PD Intoxilyzer - Alcohol Analyzer Model: 8000 04/14/2020 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>10:42</td></tr> <tr><td>Control Test</td><td>0.049</td><td>10:43</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:43</td></tr> <tr><td>Control Test</td><td>0.049</td><td>10:44</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:44</td></tr> <tr><td>Control Test</td><td>0.049</td><td>10:45</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:46</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0490</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	10:42	Control Test	0.049	10:43	Air Blank	0.000	10:43	Control Test	0.049	10:44	Air Blank	0.000	10:44	Control Test	0.049	10:45	Air Blank	0.000	10:46	Control Test Stats			Average	0.0490		Std Dev	0.0000		Rel Std Dev(%)	0.0000		<p>FORT LAUDERDALE PD Intoxilyzer - Alcohol Analyzer Model: 8000 04/14/2020 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>10:48</td></tr> <tr><td>Control Test</td><td>0.079</td><td>10:49</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:49</td></tr> <tr><td>Control Test</td><td>0.079</td><td>10:50</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:51</td></tr> <tr><td>Control Test</td><td>0.079</td><td>10:51</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:52</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0790</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	10:48	Control Test	0.079	10:49	Air Blank	0.000	10:49	Control Test	0.079	10:50	Air Blank	0.000	10:51	Control Test	0.079	10:51	Air Blank	0.000	10:52	Control Test Stats			Average	0.0790		Std Dev	0.0000		Rel Std Dev(%)	0.0000		<p>FORT LAUDERDALE PD Intoxilyzer - Alcohol Analyzer Model: 8000 04/14/2020 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>10:55</td></tr> <tr><td>Control Test</td><td>0.201</td><td>10:56</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:56</td></tr> <tr><td>Control Test</td><td>0.200</td><td>10:57</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:58</td></tr> <tr><td>Control Test</td><td>0.200</td><td>10:58</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:59</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.2003</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2882</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	10:55	Control Test	0.201	10:56	Air Blank	0.000	10:56	Control Test	0.200	10:57	Air Blank	0.000	10:58	Control Test	0.200	10:58	Air Blank	0.000	10:59	Control Test Stats			Average	0.2003		Std Dev	0.0006		Rel Std Dev(%)	0.2882		<p>FORT LAUDERDALE PD Intoxilyzer - Alcohol Analyzer Model: 8000 04/14/2020 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>11:01</td></tr> <tr><td>Control Test</td><td>0.080</td><td>11:02</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:02</td></tr> <tr><td>Control Test</td><td>0.080</td><td>11:03</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:03</td></tr> <tr><td>Control Test</td><td>0.080</td><td>11:03</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:04</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0800</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	11:01	Control Test	0.080	11:02	Air Blank	0.000	11:02	Control Test	0.080	11:03	Air Blank	0.000	11:03	Control Test	0.080	11:03	Air Blank	0.000	11:04	Control Test Stats			Average	0.0800		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:42																																																																																																																																																	
Control Test	0.049	10:43																																																																																																																																																	
Air Blank	0.000	10:43																																																																																																																																																	
Control Test	0.049	10:44																																																																																																																																																	
Air Blank	0.000	10:44																																																																																																																																																	
Control Test	0.049	10:45																																																																																																																																																	
Air Blank	0.000	10:46																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0490																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:48																																																																																																																																																	
Control Test	0.079	10:49																																																																																																																																																	
Air Blank	0.000	10:49																																																																																																																																																	
Control Test	0.079	10:50																																																																																																																																																	
Air Blank	0.000	10:51																																																																																																																																																	
Control Test	0.079	10:51																																																																																																																																																	
Air Blank	0.000	10:52																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0790																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:55																																																																																																																																																	
Control Test	0.201	10:56																																																																																																																																																	
Air Blank	0.000	10:56																																																																																																																																																	
Control Test	0.200	10:57																																																																																																																																																	
Air Blank	0.000	10:58																																																																																																																																																	
Control Test	0.200	10:58																																																																																																																																																	
Air Blank	0.000	10:59																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.2003																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.2882																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	11:01																																																																																																																																																	
Control Test	0.080	11:02																																																																																																																																																	
Air Blank	0.000	11:02																																																																																																																																																	
Control Test	0.080	11:03																																																																																																																																																	
Air Blank	0.000	11:03																																																																																																																																																	
Control Test	0.080	11:03																																																																																																																																																	
Air Blank	0.000	11:04																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0800																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		

DGS



SP  
TSK  
4/23/2020

MXK  
Operator's Signature

MXK  
Operator's Signature

MXK  
Operator's Signature

MXK  
Operator's Signature

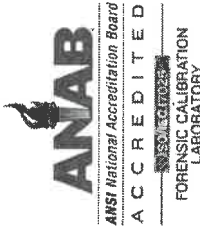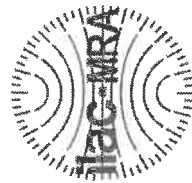

Florida Department of Law Enforcement  
 Alcohol Testing Program  
 4700 Terminal Drive, Suite 1  
 Ft. Myers, FL 33907

# Calibration Certificate

This is to certify the calibration of Intoxilyzer 8000 serial number 80-007009, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-007009</u>	UNCERTAINTY* ±
Owning Agency:	<u>FORT LAUDERDALE PD</u>	0.050 g/210 L      0.004
Calibration Date:	<u>04/14/2020</u>	0.080 g/210 L      0.005
Calibration Time:	<u>14:16</u>	0.200 g/210 L      0.007
		0.080 g/210 L Dry Gas Control    0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.  
 \*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

### TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 January 2020

Issuing Authority: Alcohol Testing Program

  
04/14/2020      Date  
**MICHAEL D HAUGHEY,**  
**Department Inspector**

SP  
 13K  
 4/23/2020

Service • Integrity • Respect • Quality

FORT LAUDERDALE PD  
Intoxilyzer - Alcohol Analyzer  
Model 8000  
04/14/2020  
SN 80-007009  
Software: 8100.27

Flow Rate Calibration\*\*\*\*\*  
1: Rate (Liters/min) = 5  
SQRT(DIFF) = 7.277  
2: Rate (Liters/min) = 15  
SQRT(DIFF) = 11.914  
3: Rate (Liters/min) = 30  
SQRT(DIFF) = 21.000  
Dependent Data Scale Factor = 100000 L/min  
Independent Data Scale Factor = 256  
Rounded Slope = 702  
Rounded Intercept = -742371  
Correlation = 0.99754

# Flow Calibration

80-607009  
04/14/2020

SP  
BK  
4/23/2020

MK

# Florida Department of Law Enforcement Alcohol Testing Program

## AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: FORT LAUDERDALE PD  
Time of Inspection: 09:45

Date of Inspection: 04/14/2020

Serial Number: 80-007009  
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: _____ Exp: _____	0.08g/210L Test (g/210L) Lot#: _____ Exp: _____	0.20g/210L Test (g/210L) Lot#: _____ Exp: _____	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: _____ Exp: _____

Number of Simulators Used: \_\_\_\_\_

Remarks:  
BY-PASSED AI. COMPLIANCE NOT DETERMINED.

SP  
12/14/20  
4/23/2020

The above instrument complies (  ) does not comply (  ) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Michael D Haughey

MICHAEL D HAUGHEY

Signature and Printed Name

04/14/2020  
Date