

Return Material Authorization

Ship to: CMI, Inc.
 Enforcement Electronics

Shipment to repair facility authorized by: James Daniels on 5/5/2020

Items Returned: Instrument Supplies Other Describe: _____
Instrument Model: Intoxilyzer 8000 Serial Number: 80-007159

Bill To Address: <u>Pensacola PD</u> _____ _____ _____	Ship to Address: <u>FDLE-ATP Tallahassee</u> _____ _____ _____
---	---

Reason for Return:
Instrument will not upload via phone line. Was able to upload via direct connect.

Please choose one of the following options:

1. I _____, authorize all repairs.

2. I _____, authorize repairs up to \$ _____.

3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: James Daniels
Phone #: 850-454-6018 Email: jdaniels@cityofpensacola.com

ATP Contact Name: Shayla Platt ATP Email: shaylplatt@fdle.state.fl.us

MH

BK
2020.05
.08
16:12:4
5-04'00





INSTRUMENT PROCESSING SHEET

Agency Pensacola Police DepartmentS/N 80-007159

Florida Department of Law Enforcement

Date In 04/23/2020 DI Completion Date 5/5/2020 Ship P/U H/D CMI EE

Intake Performed By <u>RAW</u> <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input checked="" type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ _____ _____	Quality Checks Performed By <u>RDW</u> <input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>197</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP-102</u> 32 mm <u>0.164</u> (.139 - .169) 36 mm <u>0.175</u> (.156 - .190) 53 mm <u>0.250</u> (.228 - .278) 103 mm <u>0.515</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>26932</u> <input checked="" type="checkbox"/> Stability Checks	Flow Calibration Performed By _____ Flow Column # _____ <input type="checkbox"/> 5L/min - 17mm <input type="checkbox"/> 15L/min - 53mm <input type="checkbox"/> 30L/min - 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																											
Final Release Date FDLE Alcohol Testing Program Digitally signed by FDLE Alcohol Testing Program Date: 2020.05.08 16:18:00 -04'00'	<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td><u>MP5088</u></td> <td><u>201905A</u> <u>05-14-2021</u></td> </tr> <tr> <td>0.080</td> <td><u>MP5089</u></td> <td><u>201905B</u> <u>05-14-2021</u></td> </tr> <tr> <td>0.200</td> <td><u>MP5090</u></td> <td><u>201904D</u> <u>04-30-2021</u></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td><u>AG931603</u> <u>11-12-2021</u></td> </tr> </tbody> </table>	Simulator	Serial #	Lot #/Exp	0.050	<u>MP5088</u>	<u>201905A</u> <u>05-14-2021</u>	0.080	<u>MP5089</u>	<u>201905B</u> <u>05-14-2021</u>	0.200	<u>MP5090</u>	<u>201904D</u> <u>04-30-2021</u>	0.080 DGS	N/A	<u>AG931603</u> <u>11-12-2021</u>	Maintenance Performed By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____ Temperature Checks Performed By <u>SP</u> <input checked="" type="checkbox"/> Lab Temp °C <u>21.1</u> External Digital Therm. ID#: <u>300502</u> <input checked="" type="checkbox"/> 34°C +-2 Serial #: <u>MP5088</u> <input checked="" type="checkbox"/> 34°C +-2 Serial #: <u>MP5089</u> <input checked="" type="checkbox"/> 34°C +-2 Serial #: <u>MP5090</u>																																												
Simulator	Serial #	Lot #/Exp																																																											
0.050	<u>MP5088</u>	<u>201905A</u> <u>05-14-2021</u>																																																											
0.080	<u>MP5089</u>	<u>201905B</u> <u>05-14-2021</u>																																																											
0.200	<u>MP5090</u>	<u>201904D</u> <u>04-30-2021</u>																																																											
0.080 DGS	N/A	<u>AG931603</u> <u>11-12-2021</u>																																																											
Calibration Adjustment Performed By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> <input type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> <th>Lot Number</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>	Simulator	Serial Number	Lot Number	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial Number	Lot Number	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection Performed By <u>SP</u> Barometric Pressure ID# <u>28427</u> Gauge <u>1013</u> Instrument <u>1012</u> Mouth Alcohol Solution Lot # <u>2019-B</u> Acetone Stock Solution Lot # <u>2019-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td><u>MP5086</u></td> </tr> <tr> <td>Interferent</td> <td><u>MP5087</u></td> </tr> <tr> <td>0.050</td> <td><u>MP5088</u></td> </tr> <tr> <td>0.080</td> <td><u>MP5089</u></td> </tr> <tr> <td>0.200</td> <td><u>MP5090</u></td> </tr> </tbody> </table>	Simulator	Serial Number	0.000	<u>MP5086</u>	Interferent	<u>MP5087</u>	0.050	<u>MP5088</u>	0.080	<u>MP5089</u>	0.200	<u>MP5090</u>
Simulator	Serial Number	Lot Number	Expiration																																																										
0.000		N/A	N/A																																																										
0.040																																																													
0.100																																																													
0.200																																																													
0.300																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number	Lot Number	Expiration																																																										
0.050																																																													
0.080																																																													
0.200																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number																																																												
0.000	<u>MP5086</u>																																																												
Interferent	<u>MP5087</u>																																																												
0.050	<u>MP5088</u>																																																												
0.080	<u>MP5089</u>																																																												
0.200	<u>MP5090</u>																																																												
Notes/Suggested Service: <u>Sending to Repair to check modern. Unable to upload instrument via phone line, downloaded via direct connect sp</u>	<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input type="checkbox"/> Return to/Place into Evidentiary Use <input checked="" type="checkbox"/> Remain Out of Evidentiary Use <input type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use <table style="width:100%;"> <tr> <td style="width:50%;"><u>Michael D. Houghay</u> 2020.05.07 14:55:32 -04'00'</td> <td style="width:50%;"><u>Brett Kirkland</u> 2020.05.08 16:13:21 -04'00'</td> </tr> </table>		<u>Michael D. Houghay</u> 2020.05.07 14:55:32 -04'00'	<u>Brett Kirkland</u> 2020.05.08 16:13:21 -04'00'																																																									
<u>Michael D. Houghay</u> 2020.05.07 14:55:32 -04'00'	<u>Brett Kirkland</u> 2020.05.08 16:13:21 -04'00'																																																												
	Tech Review / Date	Admin Review / Date																																																											

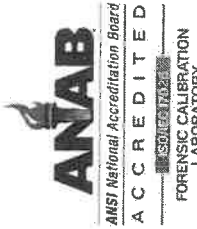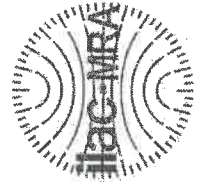

Florida Department of Law Enforcement
 Alcohol Testing Program
 2729 Fort Knox Blvd.
 Bldg. 2, Suite 1300
 Tallahassee, FL 32308

Calibration Certificate

This is to certify the calibration of Intoxilyzer 8000 serial number 80-007159, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-007159</u>	UNCERTAINTY* ±	
Owning Agency:	<u>PENSACOLA P.D.</u>	0.050 g/210 L	0.004
Calibration Date:	<u>05/05/2020</u>	0.080 g/210 L	0.005
Calibration Time:	<u>09:39</u>	0.200 g/210 L	0.007
		0.080 g/210 L Dry Gas Control	0.005

All results are reported in g/210 L.
 Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
 *Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).
 The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Shayla Platt
 SHAYLA D PLATT,
 Department Inspector

05/05/2020
 Date

FDLE/ATP Form 69 April 2020
 Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

MH
 BK
 2020.05.08
 16:14:01
 -04'00'

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PENSACOLA P.D.
Time of Inspection: 09:39

Date of Inspection: 05/05/2020

Serial Number: 80-007159
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:201905A Exp: 05/14/2021	0.08g/210L Test (g/210L) Lot#:201905B Exp: 05/14/2021	0.20g/210L Test (g/210L) Lot#:201904D Exp: 04/30/2021	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG931603 Exp: 11/12/2021
0.000	0.049	0.080	0.200	0.079
0.000	0.049	0.080	0.200	0.078
0.000	0.049	0.080	0.200	0.078
0.000	0.049	0.080	0.201	0.078
0.000	0.049	0.080	0.201	0.079
0.000	0.049	0.080	0.201	0.078
0.000	0.049	0.080	0.201	0.078
0.000	0.049	0.080	0.201	0.078
0.000	0.049	0.080	0.201	0.078
0.000	0.049	0.080	0.201	0.078
0.000	0.050	0.080	0.201	0.078

Standard Deviations	0.0003	0.0000	0.0004	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

Remarks:

MH

2020.05.
08
BK 16:14:43
-04'00'

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT

Signature and Printed Name

05/05/2020
Date

Stability Checks

PENSACOLA P.D.
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-007159
 05/05/2020
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	05:50
Control Test	0.050	05:51
Air Blank	0.000	05:51
Control Test	0.050	05:52
Air Blank	0.000	05:53
Control Test	0.050	05:53
Air Blank	0.000	05:54
Control Test Stats		
Average	0.0500	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

SA

Operator's Signature

PENSACOLA P.D.
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-007159
 05/05/2020
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	05:55
Control Test	0.080	05:56
Air Blank	0.000	05:56
Control Test	0.080	05:57
Air Blank	0.000	05:57
Control Test	0.080	05:58
Air Blank	0.000	05:59
Control Test Stats		
Average	0.0800	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

wet

SA

Operator's Signature

PENSACOLA P.D.
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-007159
 05/05/2020
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	06:00
Control Test	0.201	06:00
Air Blank	0.000	06:01
Control Test	0.201	06:02
Air Blank	0.000	06:02
Control Test	0.201	06:03
Air Blank	0.000	06:04
Control Test Stats		
Average	0.2010	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

SA

Operator's Signature

PENSACOLA P.D.
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-007159
 05/05/2020
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	06:05
Control Test	0.080	06:05
Air Blank	0.000	06:05
Control Test	0.078	06:06
Air Blank	0.000	06:06
Control Test	0.078	06:07
Air Blank	0.000	06:07
Control Test Stats		
Average	0.0787	
Std Dev	0.0012	
Rel Std Dev(%)	1.4678	

MA

BK 2020.05.0
 8 16:15:22
 -04'00"

Dry

SA

Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: PENSACOLA P.D.
Time of Inspection: 05:44

Date of Inspection: 05/05/2020

Serial Number: 80-007159
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: _____ Exp: _____	0.08g/210L Test (g/210L) Lot#: _____ Exp: _____	0.20g/210L Test (g/210L) Lot#: _____ Exp: _____	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: _____ Exp: _____

Number of Simulators Used: _____

Remarks:
BYPASSED AI TO OPERATE INSTRUMENT

MA
2020.05.
08
BK 16:16:08
-04'00

SA
N/A Compliancy not determined

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL SOTO

Signature and Printed Name

05/05/2020
Date



THE UPSIDE *of* FLORIDA

OFFICE OF THE CHIEF OF POLICE

Shayla,

This instrument will not transmit records and the software may need updated. I tried to transmit it multiple times over several days, to include going to Gulf Breeze Police Department and attempting transfer via their phone lines. Thus, the last Agency Inspection is still on the instrument.

Additionally, can you correct the phone number? It needs to be changed to: 9,18664253069.

Please let me know if you have any questions. Thank you.

A handwritten signature in black ink, appearing to be "MJH".

James Daniels
Pensacola Police Department