



INSTRUMENT PROCESSING SHEET

Agency Fort Myers PD

S/N 80-000914

Florida Department of
Law Enforcement

Date In 6/16/2021

DI Completion Date 6/17/2021

☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By MH	Quality Checks	By MH	Date 06/17/2021	Flow Calibration	By	Date																																																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>155</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 106</u> 32 mm <u>0.140</u> (.139 - .169) 36 mm <u>0.156</u> (.156 - .190) 53 mm <u>0.230</u> (.228 - .278) 103 mm <u>0.500</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28199</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																														
		<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>MP4863</td> <td>202010A 10/05/2022</td> </tr> <tr> <td>0.080</td> <td>MP4864</td> <td>202010B 10/05/2022</td> </tr> <tr> <td>0.200</td> <td>MP5097</td> <td>202010D 10/06/2022</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG026705 09/23/2022</td> </tr> </tbody> </table>		Simulator	Serial #	Lot #/Exp	0.050	MP4863	202010A 10/05/2022	0.080	MP4864	202010B 10/05/2022	0.200	MP5097	202010D 10/06/2022	0.080 DGS	N/A	AG026705 09/23/2022	Maintenance By MH ☒ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ DI Temp. Checks By MH ☒ Lab Temp °C <u>22.78</u> External Digital Therm. ID#: <u>300503</u> ☒ 34°C +/- .2 Serial #: <u>MP4863</u> ☒ 34°C +/- .2 Serial #: <u>MP4864</u> ☒ 34°C +/- .2 Serial #: <u>MP5097</u>																																																
Simulator	Serial #	Lot #/Exp																																																																	
0.050	MP4863	202010A 10/05/2022																																																																	
0.080	MP4864	202010B 10/05/2022																																																																	
0.200	MP5097	202010D 10/06/2022																																																																	
0.080 DGS	N/A	AG026705 09/23/2022																																																																	
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr> <tr><td>0.040</td><td></td><td></td><td></td></tr> <tr><td>0.100</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.300</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr><td>0.050</td><td></td><td></td><td></td></tr> <tr><td>0.080</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection By MH Barometric Pressure ID# <u>68639</u> Gauge <u>1016</u> Instrument <u>1016</u> Mouth Alcohol Solution Lot # <u>2021-A</u> Acetone Stock Solution Lot # <u>2020-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr><td>0.000</td><td>SD1014</td></tr> <tr><td>Interferent</td><td>SD1015</td></tr> <tr><td>0.050</td><td>MP4863</td></tr> <tr><td>0.080</td><td>MP4864</td></tr> <tr><td>0.200</td><td>MP5097</td></tr> </tbody> </table>						Simulator	Serial Number	0.000	SD1014	Interferent	SD1015	0.050	MP4863	0.080	MP4864	0.200	MP5097
Simulator	Serial #	Lot #	Expiration																																																																
0.000		N/A	N/A																																																																
0.040																																																																			
0.100																																																																			
0.200																																																																			
0.300																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial #	Lot #	Expiration																																																																
0.050																																																																			
0.080																																																																			
0.200																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial Number																																																																		
0.000	SD1014																																																																		
Interferent	SD1015																																																																		
0.050	MP4863																																																																		
0.080	MP4864																																																																		
0.200	MP5097																																																																		
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____ ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2021.06.17 12:17:10 -0400 17 15:16:25 04/00																																																																	
		Tech Review / Date		Admin Review / Date																																																															

TYPE OF TEST	SERIAL NUMBER	AGENCY	DATE	PERFORMED BY
Stabilities	80-000914	Fort Myers FD	06/07/21	MX

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																																
0.047 to 0.053 ☒	0.077 to 0.083 ☒	0.194 to 0.206 ☒	0.077 to 0.083 ☒																																																																																																																																																
FORT MYERS PD Intoxilyzer - Alcohol Analyzer Model 8000 06/17/2021 Software: 8100.27 SN 80-000914 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>10:00</td></tr><tr><td>Control Test</td><td>0.047</td><td>10:01</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:01</td></tr><tr><td>Control Test</td><td>0.049</td><td>10:02</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:02</td></tr><tr><td>Control Test</td><td>0.048</td><td>10:03</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:03</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0480</td><td></td></tr><tr><td>Std Dev</td><td>0.0010</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>2.0833</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	10:00	Control Test	0.047	10:01	Air Blank	0.000	10:01	Control Test	0.049	10:02	Air Blank	0.000	10:02	Control Test	0.048	10:03	Air Blank	0.000	10:03	Control Test Stats			Average	0.0480		Std Dev	0.0010		Rel Std Dev(%)	2.0833		FORT MYERS PD Intoxilyzer - Alcohol Analyzer Model 8000 06/17/2021 Software: 8100.27 SN 80-000914 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>10:05</td></tr><tr><td>Control Test</td><td>0.078</td><td>10:06</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:06</td></tr><tr><td>Control Test</td><td>0.079</td><td>10:07</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:08</td></tr><tr><td>Control Test</td><td>0.079</td><td>10:08</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:09</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0787</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7339</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	10:05	Control Test	0.078	10:06	Air Blank	0.000	10:06	Control Test	0.079	10:07	Air Blank	0.000	10:08	Control Test	0.079	10:08	Air Blank	0.000	10:09	Control Test Stats			Average	0.0787		Std Dev	0.0006		Rel Std Dev(%)	0.7339		FORT MYERS PD Intoxilyzer - Alcohol Analyzer Model 8000 06/17/2021 Software: 8100.27 SN 80-000914 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>10:10</td></tr><tr><td>Control Test</td><td>0.196</td><td>10:11</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:11</td></tr><tr><td>Control Test</td><td>0.197</td><td>10:12</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:13</td></tr><tr><td>Control Test</td><td>0.198</td><td>10:13</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:14</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.1970</td><td></td></tr><tr><td>Std Dev</td><td>0.0010</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.5076</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	10:10	Control Test	0.196	10:11	Air Blank	0.000	10:11	Control Test	0.197	10:12	Air Blank	0.000	10:13	Control Test	0.198	10:13	Air Blank	0.000	10:14	Control Test Stats			Average	0.1970		Std Dev	0.0010		Rel Std Dev(%)	0.5076		FORT MYERS PD Intoxilyzer - Alcohol Analyzer Model 8000 06/17/2021 Software: 8100.27 SN 80-000914 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>10:16</td></tr><tr><td>Control Test</td><td>0.078</td><td>10:16</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:16</td></tr><tr><td>Control Test</td><td>0.078</td><td>10:17</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:17</td></tr><tr><td>Control Test</td><td>0.077</td><td>10:17</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:18</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0777</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7434</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	10:16	Control Test	0.078	10:16	Air Blank	0.000	10:16	Control Test	0.078	10:17	Air Blank	0.000	10:17	Control Test	0.077	10:17	Air Blank	0.000	10:18	Control Test Stats			Average	0.0777		Std Dev	0.0006		Rel Std Dev(%)	0.7434	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:00																																																																																																																																																	
Control Test	0.047	10:01																																																																																																																																																	
Air Blank	0.000	10:01																																																																																																																																																	
Control Test	0.049	10:02																																																																																																																																																	
Air Blank	0.000	10:02																																																																																																																																																	
Control Test	0.048	10:03																																																																																																																																																	
Air Blank	0.000	10:03																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0480																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	2.0833																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:05																																																																																																																																																	
Control Test	0.078	10:06																																																																																																																																																	
Air Blank	0.000	10:06																																																																																																																																																	
Control Test	0.079	10:07																																																																																																																																																	
Air Blank	0.000	10:08																																																																																																																																																	
Control Test	0.079	10:08																																																																																																																																																	
Air Blank	0.000	10:09																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0787																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7339																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:10																																																																																																																																																	
Control Test	0.196	10:11																																																																																																																																																	
Air Blank	0.000	10:11																																																																																																																																																	
Control Test	0.197	10:12																																																																																																																																																	
Air Blank	0.000	10:13																																																																																																																																																	
Control Test	0.198	10:13																																																																																																																																																	
Air Blank	0.000	10:14																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1970																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	0.5076																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:16																																																																																																																																																	
Control Test	0.078	10:16																																																																																																																																																	
Air Blank	0.000	10:16																																																																																																																																																	
Control Test	0.078	10:17																																																																																																																																																	
Air Blank	0.000	10:17																																																																																																																																																	
Control Test	0.077	10:17																																																																																																																																																	
Air Blank	0.000	10:18																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0777																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7434																																																																																																																																																		
MX Operator's Signature	MX Operator's Signature	MX Operator's Signature	MX Operator's Signature																																																																																																																																																

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: FORT MYERS PD
Time of Inspection: 12:02

Date of Inspection: 06/17/2021

Serial Number: 80-000914
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202010A Exp: 10/05/2022	0.08g/210L Test (g/210L) Lot#:202010B Exp: 10/05/2022	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG026705 Exp: 09/23/2022
0.000	0.047	0.078	0.197	0.079
0.000	0.048	0.078	0.197	0.079
0.000	0.048	0.078	0.198	0.078
0.000	0.048	0.079	0.198	0.078
0.000	0.048	0.078	0.198	0.078
0.000	0.048	0.078	0.198	0.078
0.000	0.048	0.078	0.198	0.078
0.000	0.048	0.079	0.198	0.078
0.000	0.048	0.078	0.198	0.078
0.000	0.048	0.078	0.198	0.078

Standard Deviations	0.0003	0.0004	0.0004	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:



The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.


 Signature and Printed Name

MICHAEL D HAUGHEY

06/17/2021
 Date

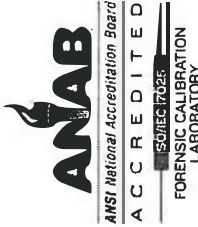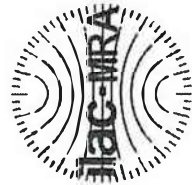

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000914, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000914</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>FORT MYERS PD</u>	0.050 g/ 210 L	0.005
Calibration Date:	<u>06/17/2021</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>12:02</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 January 2021

Issuing Authority: Alcohol Testing Program

06/17/2021

Date

MICHAEL D HAUGHEY,

Department Inspector

Service • Integrity • Respect • Quality