

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: CLEWISTON P.D.

Time of Inspection: 11:24

Date of Inspection: 06/15/2021

Serial Number: 80-000955

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202010A Exp: 10/05/2022	0.08g/210L Test (g/210L) Lot#:202010B Exp: 10/05/2022	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG026705 Exp: 09/23/2022
0.000	0.049	0.079	0.200	0.080
0.000	0.049	0.079	0.200	0.081
0.000	0.049	0.080	0.200	0.080
0.000	0.049	0.080	0.200	0.081
0.000	0.049	0.080	0.200	0.081
0.000	0.049	0.079	0.201	0.081
0.000	0.049	0.080	0.201	0.080
0.000	0.050	0.080	0.200	0.081
0.000	0.050	0.080	0.200	0.080
0.000	0.050	0.080	0.201	0.081


Standard Deviations	0.0004	0.0004	0.0004	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

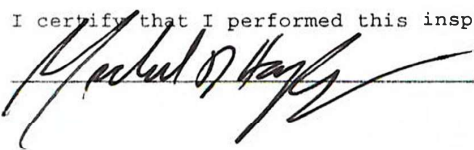

MICHAEL D HAUGHEY
Signature and Printed Name

06/15/2021
Date

TYPE OF TEST	SERIAL NUMBER	AGENCY	DATE	PERFORMED BY
Stabilities	80-000955	Clewiston PD	06/15/21	MX

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																																
0.047 to 0.053 ☒	0.077 to 0.083 ☒	0.194 to 0.206 ☒	0.077 to 0.083 ☒																																																																																																																																																
CLEWISTON P.D. Intoxilyzer - Alcohol Analyzer Model 8000 06/15/2021 Software: 8100.27 SN 80-000955 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>08:54</td></tr><tr><td>Control Test</td><td>0.050</td><td>08:55</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:56</td></tr><tr><td>Control Test</td><td>0.049</td><td>08:56</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:57</td></tr><tr><td>Control Test</td><td>0.050</td><td>08:58</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:59</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0497</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>1.1625</td><td></td></tr></table> Operator's Signature <u>MX</u>	Test	g/210L	Time	Air Blank	0.000	08:54	Control Test	0.050	08:55	Air Blank	0.000	08:56	Control Test	0.049	08:56	Air Blank	0.000	08:57	Control Test	0.050	08:58	Air Blank	0.000	08:59	Control Test Stats			Average	0.0497		Std Dev	0.0006		Rel Std Dev(%)	1.1625		CLEWISTON P.D. Intoxilyzer - Alcohol Analyzer Model 8000 06/15/2021 Software: 8100.27 SN 80-000955 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>09:00</td></tr><tr><td>Control Test</td><td>0.081</td><td>09:00</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:01</td></tr><tr><td>Control Test</td><td>0.079</td><td>09:02</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:02</td></tr><tr><td>Control Test</td><td>0.081</td><td>09:03</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:03</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0803</td><td></td></tr><tr><td>Std Dev</td><td>0.0012</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>1.4374</td><td></td></tr></table> Operator's Signature <u>MX</u>	Test	g/210L	Time	Air Blank	0.000	09:00	Control Test	0.081	09:00	Air Blank	0.000	09:01	Control Test	0.079	09:02	Air Blank	0.000	09:02	Control Test	0.081	09:03	Air Blank	0.000	09:03	Control Test Stats			Average	0.0803		Std Dev	0.0012		Rel Std Dev(%)	1.4374		CLEWISTON P.D. Intoxilyzer - Alcohol Analyzer Model 8000 06/15/2021 Software: 8100.27 SN 80-000955 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>09:05</td></tr><tr><td>Control Test</td><td>0.200</td><td>09:06</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:06</td></tr><tr><td>Control Test</td><td>0.199</td><td>09:07</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:08</td></tr><tr><td>Control Test</td><td>0.199</td><td>09:08</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:09</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.1993</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.2896</td><td></td></tr></table> Operator's Signature <u>MX</u>	Test	g/210L	Time	Air Blank	0.000	09:05	Control Test	0.200	09:06	Air Blank	0.000	09:06	Control Test	0.199	09:07	Air Blank	0.000	09:08	Control Test	0.199	09:08	Air Blank	0.000	09:09	Control Test Stats			Average	0.1993		Std Dev	0.0006		Rel Std Dev(%)	0.2896		CLEWISTON P.D. Intoxilyzer - Alcohol Analyzer Model 8000 06/15/2021 Software: 8100.27 SN 80-000955 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>09:10</td></tr><tr><td>Control Test</td><td>0.080</td><td>09:10</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:11</td></tr><tr><td>Control Test</td><td>0.081</td><td>09:11</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:12</td></tr><tr><td>Control Test</td><td>0.081</td><td>09:12</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:12</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0807</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7157</td><td></td></tr></table> Operator's Signature <u>MX</u>	Test	g/210L	Time	Air Blank	0.000	09:10	Control Test	0.080	09:10	Air Blank	0.000	09:11	Control Test	0.081	09:11	Air Blank	0.000	09:12	Control Test	0.081	09:12	Air Blank	0.000	09:12	Control Test Stats			Average	0.0807		Std Dev	0.0006		Rel Std Dev(%)	0.7157	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	08:54																																																																																																																																																	
Control Test	0.050	08:55																																																																																																																																																	
Air Blank	0.000	08:56																																																																																																																																																	
Control Test	0.049	08:56																																																																																																																																																	
Air Blank	0.000	08:57																																																																																																																																																	
Control Test	0.050	08:58																																																																																																																																																	
Air Blank	0.000	08:59																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0497																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	1.1625																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	09:00																																																																																																																																																	
Control Test	0.081	09:00																																																																																																																																																	
Air Blank	0.000	09:01																																																																																																																																																	
Control Test	0.079	09:02																																																																																																																																																	
Air Blank	0.000	09:02																																																																																																																																																	
Control Test	0.081	09:03																																																																																																																																																	
Air Blank	0.000	09:03																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0803																																																																																																																																																		
Std Dev	0.0012																																																																																																																																																		
Rel Std Dev(%)	1.4374																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	09:05																																																																																																																																																	
Control Test	0.200	09:06																																																																																																																																																	
Air Blank	0.000	09:06																																																																																																																																																	
Control Test	0.199	09:07																																																																																																																																																	
Air Blank	0.000	09:08																																																																																																																																																	
Control Test	0.199	09:08																																																																																																																																																	
Air Blank	0.000	09:09																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1993																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.2896																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	09:10																																																																																																																																																	
Control Test	0.080	09:10																																																																																																																																																	
Air Blank	0.000	09:11																																																																																																																																																	
Control Test	0.081	09:11																																																																																																																																																	
Air Blank	0.000	09:12																																																																																																																																																	
Control Test	0.081	09:12																																																																																																																																																	
Air Blank	0.000	09:12																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0807																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7157																																																																																																																																																		

Return Material Authorization

Ship to:

☒ CMI, Inc.

☐ Enforcement Electronics

Shipment to repair facility authorized by: Sgt. Rosario on 01/14/2021

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-000955

Bill To Address:

Clewiston PD

Ship to Address:

ATP FMROC

ATTN: Michael Haughey

4700 Terminal Dr.

Fort Myers FL 33907

Reason for Return:

Multiple DSP Fails. Flagged 0.08 ARS as interferent.

Please choose one of the following options:

☐ 1. I _____, authorize all repairs.

☐ 2. I _____, authorize repairs up to \$_____.

☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.


Please contact: Name: Sgt. Rosario

Phone #: (863) 599-4944

Email: Fausto.rosario@clewiston-fl.gov

ATP Contact Name: Michael Haughey

ATP Email: MichaelHaughey@fdle.state.fl

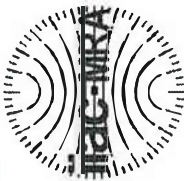

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000955, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000955</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>CLEWISTON P.D.</u>	0.050 g/ 210 L	0.005
Calibration Date:	<u>06/15/2021</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>11:24</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

06/15/2021

Date


MICHAEL D HAUGHEY,

Department Inspector

FDLE/ATP Form 69 January 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality