



INSTRUMENT PROCESSING SHEET

Agency Okeechobee PD

S/N 80-001322

Florida Department of
Law Enforcement

Date In 4/1/2021

DI Completion Date 4/8/2021

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By MH	Quality Checks	By MH	Date	04/08/21	Flow Calibration	By	Date																																																													
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>183</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 101</u> 32 mm <u>0.164</u> (.139 - .169) 36 mm <u>0.179</u> (.156 - .190) 53 mm <u>0.257</u> (.228 - .278) 103 mm <u>0.531</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28199</u> ☒ Stability Checks				Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																															
						Maintenance By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____																																																															
						DI Temp. Checks By MH _____ ☒ Lab Temp °C <u>21.85</u> External Digital Therm. ID#: <u>300504</u> ☒ 34°C +/-2 Serial #: <u>MP4863</u> ☒ 34°C +/-2 Serial #: <u>MP4864</u> ☒ 34°C +/-2 Serial #: <u>MP5097</u>																																																															
Calibration Adjustment					Department Inspection																																																																
Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr> <tr><td>0.040</td><td></td><td></td><td></td></tr> <tr><td>0.100</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.300</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr><td>0.050</td><td></td><td></td><td></td></tr> <tr><td>0.080</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table>					Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Barometric Pressure Gauge ID# <u>68639</u> Gauge <u>1018</u> Instrument <u>1017</u> Mouth Alcohol Solution Lot # <u>2020-A</u> Acetone Stock Solution Lot # <u>2020-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr><td>0.000</td><td>SD1014</td></tr> <tr><td>Interferent</td><td>SD1015</td></tr> <tr><td>0.050</td><td>MP4863</td></tr> <tr><td>0.080</td><td>MP4864</td></tr> <tr><td>0.200</td><td>MP5097</td></tr> </tbody> </table>					Simulator	Serial Number	0.000	SD1014	Interferent	SD1015	0.050	MP4863	0.080	MP4864	0.200	MP5097
Simulator	Serial #	Lot #	Expiration																																																																		
0.000		N/A	N/A																																																																		
0.040																																																																					
0.100																																																																					
0.200																																																																					
0.300																																																																					
0.080 DGS	N/A																																																																				
Simulator	Serial #	Lot #	Expiration																																																																		
0.050																																																																					
0.080																																																																					
0.200																																																																					
0.080 DGS	N/A																																																																				
Simulator	Serial Number																																																																				
0.000	SD1014																																																																				
Interferent	SD1015																																																																				
0.050	MP4863																																																																				
0.080	MP4864																																																																				
0.200	MP5097																																																																				
Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment					☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 ☐ Other _____																																																																
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____ _____					☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use 2021.04. Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2021.04.09 09:31:38 -0400 09 15:42:16 Tech Review / Date Admin Review / Date																																																																

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: OKEECHOBEE P.D.

Time of Inspection: 14:43

Date of Inspection: 04/08/2021

Serial Number: 80-001322

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202010A Exp: 10/05/2022	0.08g/210L Test (g/210L) Lot#:202010B Exp: 10/05/2022	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG026705 Exp: 09/23/2022
0.000	0.049	0.079	0.197	0.079
0.000	0.049	0.078	0.197	0.079
0.000	0.049	0.079	0.197	0.078
0.000	0.049	0.079	0.197	0.078
0.000	0.049	0.079	0.197	0.078
0.000	0.049	0.078	0.197	0.078
0.000	0.050	0.079	0.197	0.078
0.000	0.050	0.078	0.197	0.079
0.000	0.050	0.078	0.197	0.079
0.000	0.050	0.079	0.197	0.078


Standard Deviations	0.0005	0.0005	0.0000	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

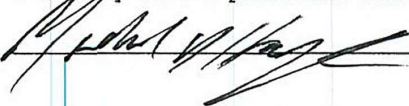

MICHAEL D HAUGHEY

Signature and Printed Name

04/08/2021
Date

TYPE OF TEST	SERIAL NUMBER	AGENCY	DATE	PERFORMED BY
Stabilities	80-001322	<i>C. Keechobee P.D.</i>	<i>4/8/21</i>	<i>MK</i>

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																				
0.047 to 0.053 ☒	0.077 to 0.083 ☒	0.194 to 0.206 ☒	0.077 to 0.083 ☒																																																																																																																																				
OKEECHOBEE P.D. Intoxilyzer - Alconol Analyzer Model 8000 04/08/2021 Software: 8100.27 SN 80-001322 Test g/210L Time <table><tr><td>Air Blank</td><td>0.000</td><td>11:03</td></tr><tr><td>Control Test</td><td>0.050</td><td>11:04</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:05</td></tr><tr><td>Control Test</td><td>0.049</td><td>11:05</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:06</td></tr><tr><td>Control Test</td><td>0.050</td><td>11:06</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:07</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0497</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>1.1625</td><td></td></tr></table>	Air Blank	0.000	11:03	Control Test	0.050	11:04	Air Blank	0.000	11:05	Control Test	0.049	11:05	Air Blank	0.000	11:06	Control Test	0.050	11:06	Air Blank	0.000	11:07	Control Test Stats			Average	0.0497		Std Dev	0.0006		Rel Std Dev(%)	1.1625		OKEECHOBEE P.D. Intoxilyzer - Alconol Analyzer Model 8000 04/08/2021 Software: 8100.27 SN 80-001322 Test g/210L Time <table><tr><td>Air Blank</td><td>0.000</td><td>11:09</td></tr><tr><td>Control Test</td><td>0.079</td><td>11:09</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:10</td></tr><tr><td>Control Test</td><td>0.078</td><td>11:11</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:11</td></tr><tr><td>Control Test</td><td>0.079</td><td>11:12</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:13</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0787</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7339</td><td></td></tr></table>	Air Blank	0.000	11:09	Control Test	0.079	11:09	Air Blank	0.000	11:10	Control Test	0.078	11:11	Air Blank	0.000	11:11	Control Test	0.079	11:12	Air Blank	0.000	11:13	Control Test Stats			Average	0.0787		Std Dev	0.0006		Rel Std Dev(%)	0.7339		OKEECHOBEE P.D. Intoxilyzer - Alconol Analyzer Model 8000 04/08/2021 Software: 8100.27 SN 80-001322 Test g/210L Time <table><tr><td>Air Blank</td><td>0.000</td><td>11:14</td></tr><tr><td>Control Test</td><td>0.198</td><td>11:15</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:15</td></tr><tr><td>Control Test</td><td>0.196</td><td>11:16</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:17</td></tr><tr><td>Control Test</td><td>0.196</td><td>11:17</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:18</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.1967</td><td></td></tr><tr><td>Std Dev</td><td>0.0012</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.5871</td><td></td></tr></table>	Air Blank	0.000	11:14	Control Test	0.198	11:15	Air Blank	0.000	11:15	Control Test	0.196	11:16	Air Blank	0.000	11:17	Control Test	0.196	11:17	Air Blank	0.000	11:18	Control Test Stats			Average	0.1967		Std Dev	0.0012		Rel Std Dev(%)	0.5871		OKEECHOBEE P.D. Intoxilyzer - Alconol Analyzer Model 8000 04/08/2021 Software: 8100.27 SN 80-001322 Test g/210L Time <table><tr><td>Air Blank</td><td>0.000</td><td>11:19</td></tr><tr><td>Control Test</td><td>0.079</td><td>11:20</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:20</td></tr><tr><td>Control Test</td><td>0.078</td><td>11:21</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:21</td></tr><tr><td>Control Test</td><td>0.079</td><td>11:21</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:22</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0787</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7339</td><td></td></tr></table>	Air Blank	0.000	11:19	Control Test	0.079	11:20	Air Blank	0.000	11:20	Control Test	0.078	11:21	Air Blank	0.000	11:21	Control Test	0.079	11:21	Air Blank	0.000	11:22	Control Test Stats			Average	0.0787		Std Dev	0.0006		Rel Std Dev(%)	0.7339	
Air Blank	0.000	11:03																																																																																																																																					
Control Test	0.050	11:04																																																																																																																																					
Air Blank	0.000	11:05																																																																																																																																					
Control Test	0.049	11:05																																																																																																																																					
Air Blank	0.000	11:06																																																																																																																																					
Control Test	0.050	11:06																																																																																																																																					
Air Blank	0.000	11:07																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0497																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	1.1625																																																																																																																																						
Air Blank	0.000	11:09																																																																																																																																					
Control Test	0.079	11:09																																																																																																																																					
Air Blank	0.000	11:10																																																																																																																																					
Control Test	0.078	11:11																																																																																																																																					
Air Blank	0.000	11:11																																																																																																																																					
Control Test	0.079	11:12																																																																																																																																					
Air Blank	0.000	11:13																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0787																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.7339																																																																																																																																						
Air Blank	0.000	11:14																																																																																																																																					
Control Test	0.198	11:15																																																																																																																																					
Air Blank	0.000	11:15																																																																																																																																					
Control Test	0.196	11:16																																																																																																																																					
Air Blank	0.000	11:17																																																																																																																																					
Control Test	0.196	11:17																																																																																																																																					
Air Blank	0.000	11:18																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.1967																																																																																																																																						
Std Dev	0.0012																																																																																																																																						
Rel Std Dev(%)	0.5871																																																																																																																																						
Air Blank	0.000	11:19																																																																																																																																					
Control Test	0.079	11:20																																																																																																																																					
Air Blank	0.000	11:20																																																																																																																																					
Control Test	0.078	11:21																																																																																																																																					
Air Blank	0.000	11:21																																																																																																																																					
Control Test	0.079	11:21																																																																																																																																					
Air Blank	0.000	11:22																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0787																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.7339																																																																																																																																						

MK
Operator's Signature

MK
Operator's Signature

MK
Operator's Signature

MK
Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: OKEECHOBEE P.D.

Time of Inspection: 09:54

Date of Inspection: 04/08/2021

Serial Number: 80-001322

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

COMPLIANCE NOT DETERMINED.AI BY-PASS.

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.


I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



MICHAEL D HAUGHEY

Signature and Printed Name

04/08/2021
Date

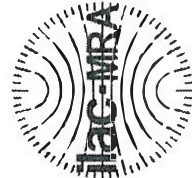

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001322 , manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001322</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>OKEECHOBEE P.D.</u>	0.050 g/ 210 L	0.005
Calibration Date:	<u>04/08/2021</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>14:43</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full,
without written approval of the Florida Department of
Law Enforcement Alcohol Testing Program.

04/08/2021

Date

MICHAEL D HAUGHEY,

Department Inspector

FDLE/ATP Form 69 January 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1