



INSTRUMENT PROCESSING SHEET

Agency Miami Gardens PD

S/N 80-002988

Florida Department of
Law Enforcement

Date In 5/24/2021 DI Completion Date 5/24/2021

☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By MH	Quality Checks	By MH	Date 05/24/2021	Flow Calibration	By	Date																																								
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>184</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 106</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.515</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28 199</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																										
					Maintenance By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____																																										
					DI Temp. Checks By MH _____ ☒ Lab Temp °C <u>21.83</u> External Digital Therm. ID#: <u>300503</u> ☒ 34°C +/- .2 Serial #: <u>MP4863</u> ☒ 34°C +/- .2 Serial #: <u>MP4864</u> ☒ 34°C +/- .2 Serial #: <u>MP5097</u>																																										
Calibration Adjustment By _____ Barometric Pressure Gauge ID # _____				Department Inspection By MH _____ Barometric Pressure ID# <u>28663</u> Gauge <u>1021</u> Instrument <u>1017</u> Mouth Alcohol Solution Lot # <u>2021-A</u> Acetone Stock Solution Lot # <u>2020-A</u>																																											
<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr> <tr><td>0.040</td><td></td><td></td><td></td></tr> <tr><td>0.100</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.300</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table>				Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr><td>0.000</td><td>SD1014</td></tr> <tr><td>Interferent</td><td>SD1015</td></tr> <tr><td>0.050</td><td>MP4863</td></tr> <tr><td>0.080</td><td>MP4864</td></tr> <tr><td>0.200</td><td>MP5097</td></tr> </tbody> </table>				Simulator	Serial Number	0.000	SD1014	Interferent	SD1015	0.050	MP4863	0.080	MP4864	0.200	MP5097
Simulator	Serial #	Lot #	Expiration																																												
0.000		N/A	N/A																																												
0.040																																															
0.100																																															
0.200																																															
0.300																																															
0.080 DGS	N/A																																														
Simulator	Serial Number																																														
0.000	SD1014																																														
Interferent	SD1015																																														
0.050	MP4863																																														
0.080	MP4864																																														
0.200	MP5097																																														
☐ Post Calibration Adjustment Stability Checks				Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____																																											
<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr><td>0.050</td><td></td><td></td><td></td></tr> <tr><td>0.080</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table>				Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																							
Simulator	Serial #	Lot #	Expiration																																												
0.050																																															
0.080																																															
0.200																																															
0.080 DGS	N/A																																														
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____				Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2021.05.24 12:58:11 -04'00' 2021.05.24 16:19:29 -04'00' 																																											
				Tech Review / Date _____ Admin Review / Date _____																																											

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI GARDENS PD
Time of Inspection: 11:57

Date of Inspection: 05/24/2021

Serial Number: 80-002988
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202010A Exp: 10/05/2022	0.08g/210L Test (g/210L) Lot#:202010B Exp: 10/05/2022	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG026705 Exp: 09/23/2022
0.000	0.047	0.078	0.195	0.080
0.000	0.048	0.078	0.195	0.081
0.000	0.048	0.077	0.195	0.080
0.000	0.048	0.078	0.196	0.081
0.000	0.048	0.078	0.196	0.080
0.000	0.048	0.078	0.196	0.080
0.000	0.048	0.078	0.195	0.080
0.000	0.048	0.078	0.195	0.080
0.000	0.048	0.078	0.196	0.079
0.000	0.048	0.078	0.195	0.079


Standard Deviations	0.0003	0.0003	0.0005	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

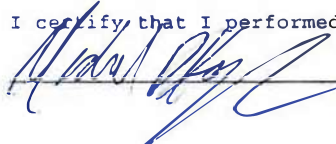



MICHAEL D HAUGHEY

Signature and Printed Name

05/24/2021
Date

TYPE OF TEST	SERIAL NUMBER	AGENCY	DATE	PERFORMED BY
Stabilities	80-002988	Miami Gardens	05/24/21	NJK

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																																
0.047 to 0.053 ☐	0.077 to 0.083 ☐	0.194 to 0.206 ☐	0.077 to 0.083 ☐																																																																																																																																																
MIAMI GARDENS PD Intoxilyzer - Alcotest Analyzer Model 8000 05/24/2021 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>09:56</td></tr><tr><td>Control Test</td><td>0.049</td><td>09:57</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:57</td></tr><tr><td>Control Test</td><td>0.048</td><td>09:58</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:59</td></tr><tr><td>Control Test</td><td>0.049</td><td>09:59</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:00</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0483</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>1.1945</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	09:56	Control Test	0.049	09:57	Air Blank	0.000	09:57	Control Test	0.048	09:58	Air Blank	0.000	09:59	Control Test	0.049	09:59	Air Blank	0.000	10:00	Control Test Stats			Average	0.0483		Std Dev	0.0006		Rel Std Dev(%)	1.1945		MIAMI GARDENS PD Intoxilyzer - Alcotest Analyzer Model 8000 05/24/2021 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>10:01</td></tr><tr><td>Control Test</td><td>0.078</td><td>10:01</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:02</td></tr><tr><td>Control Test</td><td>0.078</td><td>10:03</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:03</td></tr><tr><td>Control Test</td><td>0.078</td><td>10:04</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:04</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0780</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	10:01	Control Test	0.078	10:01	Air Blank	0.000	10:02	Control Test	0.078	10:03	Air Blank	0.000	10:03	Control Test	0.078	10:04	Air Blank	0.000	10:04	Control Test Stats			Average	0.0780		Std Dev	0.0000		Rel Std Dev(%)	0.0000		MIAMI GARDENS PD Intoxilyzer - Alcotest Analyzer Model 8000 05/24/2021 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>10:06</td></tr><tr><td>Control Test</td><td>0.195</td><td>10:06</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:07</td></tr><tr><td>Control Test</td><td>0.194</td><td>10:08</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:08</td></tr><tr><td>Control Test</td><td>0.195</td><td>10:09</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:09</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.1947</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.2966</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	10:06	Control Test	0.195	10:06	Air Blank	0.000	10:07	Control Test	0.194	10:08	Air Blank	0.000	10:08	Control Test	0.195	10:09	Air Blank	0.000	10:09	Control Test Stats			Average	0.1947		Std Dev	0.0006		Rel Std Dev(%)	0.2966		MIAMI GARDENS PD Intoxilyzer - Alcotest Analyzer Model 8000 05/24/2021 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>10:11</td></tr><tr><td>Control Test</td><td>0.080</td><td>10:11</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:11</td></tr><tr><td>Control Test</td><td>0.080</td><td>10:12</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:12</td></tr><tr><td>Control Test</td><td>0.080</td><td>10:13</td></tr><tr><td>Air Blank</td><td>0.000</td><td>10:13</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0800</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	10:11	Control Test	0.080	10:11	Air Blank	0.000	10:11	Control Test	0.080	10:12	Air Blank	0.000	10:12	Control Test	0.080	10:13	Air Blank	0.000	10:13	Control Test Stats			Average	0.0800		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	09:56																																																																																																																																																	
Control Test	0.049	09:57																																																																																																																																																	
Air Blank	0.000	09:57																																																																																																																																																	
Control Test	0.048	09:58																																																																																																																																																	
Air Blank	0.000	09:59																																																																																																																																																	
Control Test	0.049	09:59																																																																																																																																																	
Air Blank	0.000	10:00																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0483																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	1.1945																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:01																																																																																																																																																	
Control Test	0.078	10:01																																																																																																																																																	
Air Blank	0.000	10:02																																																																																																																																																	
Control Test	0.078	10:03																																																																																																																																																	
Air Blank	0.000	10:03																																																																																																																																																	
Control Test	0.078	10:04																																																																																																																																																	
Air Blank	0.000	10:04																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0780																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:06																																																																																																																																																	
Control Test	0.195	10:06																																																																																																																																																	
Air Blank	0.000	10:07																																																																																																																																																	
Control Test	0.194	10:08																																																																																																																																																	
Air Blank	0.000	10:08																																																																																																																																																	
Control Test	0.195	10:09																																																																																																																																																	
Air Blank	0.000	10:09																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1947																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.2966																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:11																																																																																																																																																	
Control Test	0.080	10:11																																																																																																																																																	
Air Blank	0.000	10:11																																																																																																																																																	
Control Test	0.080	10:12																																																																																																																																																	
Air Blank	0.000	10:12																																																																																																																																																	
Control Test	0.080	10:13																																																																																																																																																	
Air Blank	0.000	10:13																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0800																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Operator's Signature	Operator's Signature	Operator's Signature	Operator's Signature																																																																																																																																																

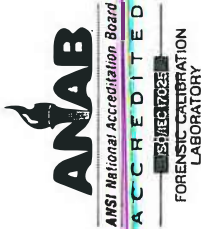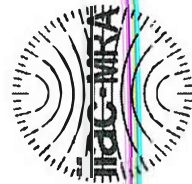

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-002988, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-002988</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>MIAMI GARDENS PD</u>	0.050 g/ 210 L	0.005
Calibration Date:	<u>05/24/2021</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>11:57</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Thermometer temperatures are checked with NIST traceable Eutechnics 4400 digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the uses of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

05/24/2021

Date


MICHAEL D HAUGHEY,
Department Inspector

FDLE/ATP Form 69 January 2021
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1