



INSTRUMENT PROCESSING SHEET

Agency Columbia County Sheriff's OfficeS/N 80-000775Florida Department of
Law EnforcementDate In 12/19/2022 DI Completion Date 12/22/2022☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By PN	Quality Checks	By PN	Date <u>12/22/2022</u>	Flow Calibration	By	Date																																							
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>250</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 105</u> 32 mm <u>.160</u> (.139 - .169) 36 mm <u>.167</u> (.156 - .190) 53 mm <u>.242</u> (.228 - .278) 103 mm <u>.519</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28662</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																									
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202201C 01/11/2024</td></tr><tr><td>0.080</td><td>MP5089</td><td>202201D 01/18/2024</td></tr><tr><td>0.200</td><td>MP5084</td><td>202201E 01/18/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG229803 10/25/2024</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP5088	202201C 01/11/2024	0.080	MP5089	202201D 01/18/2024	0.200	MP5084	202201E 01/18/2024	0.080 DGS	N/A	AG229803 10/25/2024																													
Simulator	Serial #	Lot #/Exp																																												
0.050	MP5088	202201C 01/11/2024																																												
0.080	MP5089	202201D 01/18/2024																																												
0.200	MP5084	202201E 01/18/2024																																												
0.080 DGS	N/A	AG229803 10/25/2024																																												
Calibration Adjustment		By _____		Department Inspection			By <u>PN</u>																																							
Barometric Pressure Gauge _____ ID # _____				Barometric Pressure ID# <u>28662</u> Gauge <u>1008</u> Instrument <u>1001</u> Mouth Alcohol Solution Lot # <u>2022-A</u> Acetone Stock Solution Lot # <u>2022-B</u>																																										
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A					<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5086</td></tr><tr><td>Interferent</td><td>MP5087</td></tr><tr><td>0.050</td><td>MP5088</td></tr><tr><td>0.080</td><td>MP5089</td></tr><tr><td>0.200</td><td>MP5084</td></tr></tbody></table>			Simulator	Serial Number	0.000	MP5086	Interferent	MP5087	0.050	MP5088	0.080	MP5089	0.200	MP5084
Simulator	Serial #	Lot #	Expiration																																											
0.000		N/A	N/A																																											
0.040																																														
0.100																																														
0.200																																														
0.300																																														
0.080 DGS	N/A																																													
Simulator	Serial Number																																													
0.000	MP5086																																													
Interferent	MP5087																																													
0.050	MP5088																																													
0.080	MP5089																																													
0.200	MP5084																																													
☐ Post Calibration Adjustment Stability Checks				Attachments																																										
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A					<table border="1"><tbody><tr><td>☒ Form 41</td><td>☐ Post-Stability Checks</td></tr><tr><td>☒ Stability Checks</td><td>☐ Flow Calibration</td></tr><tr><td>☒ Calibration Certificate</td><td>☐ Form 40</td></tr><tr><td>☐ Calibration Adjustment</td><td>☐ Other _____</td></tr></tbody></table>			☒ Form 41	☐ Post-Stability Checks	☒ Stability Checks	☐ Flow Calibration	☒ Calibration Certificate	☐ Form 40	☐ Calibration Adjustment	☐ Other _____												
Simulator	Serial #	Lot #	Expiration																																											
0.050																																														
0.080																																														
0.200																																														
0.080 DGS	N/A																																													
☒ Form 41	☐ Post-Stability Checks																																													
☒ Stability Checks	☐ Flow Calibration																																													
☒ Calibration Certificate	☐ Form 40																																													
☐ Calibration Adjustment	☐ Other _____																																													
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____ _____				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																										
		David Eliezer Reyes-Rivera		5 12:23:51 Tech Review / Date Admin Review / Date																																										

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: COLUMBIA COUNTY SO
Time of Inspection: 15:09

Date of Inspection: 12/22/2022

Serial Number: 80-000775
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG229803 Exp: 10/25/2024
0.000	0.050	0.079	0.200	0.080
0.000	0.050	0.079	0.201	0.080
0.000	0.050	0.079	0.200	0.079
0.000	0.049	0.079	0.201	0.079
0.000	0.050	0.078	0.201	0.080
0.000	0.049	0.079	0.200	0.080
0.000	0.050	0.079	0.200	0.079
0.000	0.050	0.079	0.200	0.079
0.000	0.049	0.079	0.200	0.080
0.000	0.050	0.079	0.200	0.080


Standard Deviations	0.0004	0.0003	0.0004	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

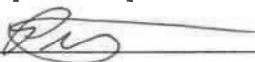

PHIL NICODEMO
Signature and Printed Name

12/22/2022
Date

STABILITY CHECKS 80-000775

COLUMBIA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000775
12/22/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:50
Control Test	0.050	12:50
Air Blank	0.000	12:51
Control Test	0.051	12:52
Air Blank	0.000	12:52
Control Test	0.050	12:53
Air Blank	0.000	12:53
Control Test Stats		
Average	0.0503	
Std Dev	0.0006	
Rel Std Dev(%)	1.1471	


Operator's Signature

COLUMBIA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000775
12/22/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:54
Control Test	0.079	12:55
Air Blank	0.000	12:56
Control Test	0.079	12:56
Air Blank	0.000	12:57
Control Test	0.079	12:58
Air Blank	0.000	12:58
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator's Signature

COLUMBIA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000775
12/22/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	13:00
Control Test	0.202	13:01
Air Blank	0.000	13:01
Control Test	0.200	13:02
Air Blank	0.000	13:02
Control Test	0.201	13:03
Air Blank	0.000	13:04
Control Test Stats		
Average	0.2010	
Std Dev	0.0010	
Rel Std Dev(%)	0.4975	


Operator's Signature

COLUMBIA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000775
12/22/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	13:05
Control Test	0.080	13:05
Air Blank	0.000	13:06
Control Test	0.081	13:06
Air Blank	0.000	13:07
Control Test	0.080	13:07
Air Blank	0.000	13:07
Control Test Stats		
Average	0.0803	
Std Dev	0.0006	
Rel Std Dev(%)	0.7187	

D6S


Operator's Signature



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000775, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000775</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>COLUMBIA COUNTY SO</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>12/22/2022</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>15:09</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Phil Nicodemo
Digitally signed by Phil
Nicodemo
Date: 2022.12.22 15:49:55 -05'00'

12/22/2022

Date

PHIL NICODEMO,
Department Inspector

FDLE/ATP Form 69 March 2022
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality