

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI-DADE PD
Time of Inspection: 08:43

Date of Inspection: 03/01/2022

Serial Number: 80-000880
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202010A Exp: 10/05/2022	0.08g/210L Test (g/210L) Lot#:202010B Exp: 10/05/2022	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG115904 Exp: 06/08/2023
0.000	0.050	0.080	0.200	0.077
0.000	0.050	0.081	0.200	0.077
0.000	0.050	0.081	0.199	0.077
0.000	0.050	0.080	0.200	0.077
0.000	0.050	0.080	0.200	0.076
0.000	0.050	0.081	0.200	0.076
0.000	0.050	0.080	0.200	0.077
0.000	0.050	0.080	0.200	0.076
0.000	0.049	0.080	0.200	0.076
0.000	0.049	0.080	0.200	0.076

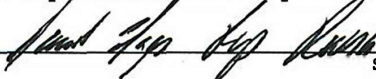
Standard Deviations	0.0004	0.0004	0.0003	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

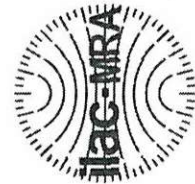

 Signature and Printed Name

DAVID E REYES-RIVERA

03/01/2022
 Date

Type of Test	Serial Number	Agency	Date	Performed By
Stabilities	80-000880	Miami-Dade Police Department	3/1/2022	DERR <i>MDL</i>

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																																			
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083																																																																																																																																																			
MIAMI-DADE PD Intoxilyzer - Alconol Analyzer Model: 8100 SN 80-000880 03/01/2022 Software: 8100.27 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>06:52</td></tr><tr><td>Control Test</td><td>0.051</td><td>06:53</td></tr><tr><td>Air Blank</td><td>0.000</td><td>06:53</td></tr><tr><td>Control Test</td><td>0.050</td><td>06:54</td></tr><tr><td>Air Blank</td><td>0.000</td><td>06:54</td></tr><tr><td>Control Test</td><td>0.050</td><td>06:55</td></tr><tr><td>Air Blank</td><td>0.000</td><td>06:55</td></tr><tr><td>Control Test</td><td>0.050</td><td>06:56</td></tr><tr><td>Control Test State</td><td></td><td></td></tr><tr><td>Average</td><td>0.0503</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel. Std Dev(%)</td><td>1.1471</td><td></td></tr></table>	Test	g/210L	Time	Air Blank	0.000	06:52	Control Test	0.051	06:53	Air Blank	0.000	06:53	Control Test	0.050	06:54	Air Blank	0.000	06:54	Control Test	0.050	06:55	Air Blank	0.000	06:55	Control Test	0.050	06:56	Control Test State			Average	0.0503		Std Dev	0.0006		Rel. Std Dev(%)	1.1471		MIAMI-DADE PD Intoxilyzer - Alconol Analyzer Model: 8100 SN 80-000880 03/01/2022 Software: 8100.27 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>06:57</td></tr><tr><td>Control Test</td><td>0.081</td><td>06:57</td></tr><tr><td>Air Blank</td><td>0.000</td><td>06:58</td></tr><tr><td>Control Test</td><td>0.080</td><td>06:59</td></tr><tr><td>Air Blank</td><td>0.000</td><td>06:59</td></tr><tr><td>Control Test</td><td>0.080</td><td>07:00</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:01</td></tr><tr><td>Control Test State</td><td></td><td></td></tr><tr><td>Average</td><td>0.0807</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel. Std Dev(%)</td><td>0.7157</td><td></td></tr></table>	Test	g/210L	Time	Air Blank	0.000	06:57	Control Test	0.081	06:57	Air Blank	0.000	06:58	Control Test	0.080	06:59	Air Blank	0.000	06:59	Control Test	0.080	07:00	Air Blank	0.000	07:01	Control Test State			Average	0.0807		Std Dev	0.0006		Rel. Std Dev(%)	0.7157		MIAMI-DADE PD Intoxilyzer - Alconol Analyzer Model: 8100 SN 80-000880 03/01/2022 Software: 8100.27 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>07:02</td></tr><tr><td>Control Test</td><td>0.201</td><td>07:03</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:03</td></tr><tr><td>Control Test</td><td>0.199</td><td>07:04</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:04</td></tr><tr><td>Control Test</td><td>0.199</td><td>07:05</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:06</td></tr><tr><td>Control Test State</td><td></td><td></td></tr><tr><td>Average</td><td>0.1997</td><td></td></tr><tr><td>Std Dev</td><td>0.0012</td><td></td></tr><tr><td>Rel. Std Dev(%)</td><td>0.5783</td><td></td></tr></table>	Test	g/210L	Time	Air Blank	0.000	07:02	Control Test	0.201	07:03	Air Blank	0.000	07:03	Control Test	0.199	07:04	Air Blank	0.000	07:04	Control Test	0.199	07:05	Air Blank	0.000	07:06	Control Test State			Average	0.1997		Std Dev	0.0012		Rel. Std Dev(%)	0.5783		MIAMI-DADE PD Intoxilyzer - Alconol Analyzer Model: 8100 SN 80-000880 03/01/2022 Software: 8100.27 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>07:07</td></tr><tr><td>Control Test</td><td>0.078</td><td>07:07</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:08</td></tr><tr><td>Control Test</td><td>0.078</td><td>07:08</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:09</td></tr><tr><td>Control Test</td><td>0.078</td><td>07:09</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:09</td></tr><tr><td>Control Test State</td><td></td><td></td></tr><tr><td>Average</td><td>0.0780</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel. Std Dev(%)</td><td>0.0000</td><td></td></tr></table>	Test	g/210L	Time	Air Blank	0.000	07:07	Control Test	0.078	07:07	Air Blank	0.000	07:08	Control Test	0.078	07:08	Air Blank	0.000	07:09	Control Test	0.078	07:09	Air Blank	0.000	07:09	Control Test State			Average	0.0780		Std Dev	0.0000		Rel. Std Dev(%)	0.0000	
Test	g/210L	Time																																																																																																																																																				
Air Blank	0.000	06:52																																																																																																																																																				
Control Test	0.051	06:53																																																																																																																																																				
Air Blank	0.000	06:53																																																																																																																																																				
Control Test	0.050	06:54																																																																																																																																																				
Air Blank	0.000	06:54																																																																																																																																																				
Control Test	0.050	06:55																																																																																																																																																				
Air Blank	0.000	06:55																																																																																																																																																				
Control Test	0.050	06:56																																																																																																																																																				
Control Test State																																																																																																																																																						
Average	0.0503																																																																																																																																																					
Std Dev	0.0006																																																																																																																																																					
Rel. Std Dev(%)	1.1471																																																																																																																																																					
Test	g/210L	Time																																																																																																																																																				
Air Blank	0.000	06:57																																																																																																																																																				
Control Test	0.081	06:57																																																																																																																																																				
Air Blank	0.000	06:58																																																																																																																																																				
Control Test	0.080	06:59																																																																																																																																																				
Air Blank	0.000	06:59																																																																																																																																																				
Control Test	0.080	07:00																																																																																																																																																				
Air Blank	0.000	07:01																																																																																																																																																				
Control Test State																																																																																																																																																						
Average	0.0807																																																																																																																																																					
Std Dev	0.0006																																																																																																																																																					
Rel. Std Dev(%)	0.7157																																																																																																																																																					
Test	g/210L	Time																																																																																																																																																				
Air Blank	0.000	07:02																																																																																																																																																				
Control Test	0.201	07:03																																																																																																																																																				
Air Blank	0.000	07:03																																																																																																																																																				
Control Test	0.199	07:04																																																																																																																																																				
Air Blank	0.000	07:04																																																																																																																																																				
Control Test	0.199	07:05																																																																																																																																																				
Air Blank	0.000	07:06																																																																																																																																																				
Control Test State																																																																																																																																																						
Average	0.1997																																																																																																																																																					
Std Dev	0.0012																																																																																																																																																					
Rel. Std Dev(%)	0.5783																																																																																																																																																					
Test	g/210L	Time																																																																																																																																																				
Air Blank	0.000	07:07																																																																																																																																																				
Control Test	0.078	07:07																																																																																																																																																				
Air Blank	0.000	07:08																																																																																																																																																				
Control Test	0.078	07:08																																																																																																																																																				
Air Blank	0.000	07:09																																																																																																																																																				
Control Test	0.078	07:09																																																																																																																																																				
Air Blank	0.000	07:09																																																																																																																																																				
Control Test State																																																																																																																																																						
Average	0.0780																																																																																																																																																					
Std Dev	0.0000																																																																																																																																																					
Rel. Std Dev(%)	0.0000																																																																																																																																																					
MDL Operator's Signature	MDL Operator's Signature	MDL Operator's Signature	MDL Operator's Signature																																																																																																																																																			



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000880, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000880</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>MIAMI-DADE PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>03/01/2022</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>08:43</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

03/01/2022

Date


DAVID E REYES-RIVERA,
Department Inspector

FDLE/ATP Form 69 December 2021
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1