



INSTRUMENT PROCESSING SHEET

Agency FFWCCS/N 80-000900Florida Department of
Law EnforcementDate In 03/07/2022DI Completion Date 04-04-2022☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By <u>RAW</u>	Quality Checks	By <u>IS</u>	Date <u>04-04-2022</u>	Flow Calibration	By _____	Date _____																														
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: Instrument came with dry gas canister in pelican case AI was notified and will be picking it up.		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>227</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-102</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.492</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>30793</u> ☒ Stability Checks			☐ Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202201C 01-11-2024</td></tr><tr><td>0.080</td><td>MP6292</td><td>202201D 01-18-2024</td></tr><tr><td>0.200</td><td>MP6293</td><td>202201E 01-18-2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG113403 05-14-2023</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP6291	202201C 01-11-2024	0.080	MP6292	202201D 01-18-2024	0.200	MP6293	202201E 01-18-2024	0.080 DGS	N/A	AG113403 05-14-2023	<table border="1"><thead><tr><th colspan="2">Maintenance</th><th>By _____</th></tr></thead><tbody><tr><td>☐ Battery Replacement</td><td colspan="2"></td></tr><tr><td>☐ Dry Gas Regulator Replacement</td><td colspan="2"></td></tr><tr><td>☐ Breath Tube Replacement</td><td colspan="2"></td></tr><tr><td>☐ Other _____</td><td colspan="2"></td></tr></tbody></table>			Maintenance		By _____	☐ Battery Replacement			☐ Dry Gas Regulator Replacement			☐ Breath Tube Replacement			☐ Other _____		
Simulator	Serial #	Lot #/Exp																																			
0.050	MP6291	202201C 01-11-2024																																			
0.080	MP6292	202201D 01-18-2024																																			
0.200	MP6293	202201E 01-18-2024																																			
0.080 DGS	N/A	AG113403 05-14-2023																																			
Maintenance		By _____																																			
☐ Battery Replacement																																					
☐ Dry Gas Regulator Replacement																																					
☐ Breath Tube Replacement																																					
☐ Other _____																																					
Calibration Adjustment		By _____		Department Inspection			By <u>IS</u>																														
Barometric Pressure Gauge _____ ID # _____				Barometric Pressure ID# <u>28421</u>																																	
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A					Gauge <u>1014</u> Instrument <u>1012</u>					
Simulator	Serial #	Lot #	Expiration																																		
0.000		N/A	N/A																																		
0.040																																					
0.100																																					
0.200																																					
0.300																																					
0.080 DGS	N/A																																				
☐ Post Calibration Adjustment Stability Checks				Mouth Alcohol Solution Lot # <u>2021-D</u>																																	
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A					Acetone Stock Solution Lot # <u>2021-C</u>													
Simulator	Serial #	Lot #	Expiration																																		
0.050																																					
0.080																																					
0.200																																					
0.080 DGS	N/A																																				
				<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6289</td></tr><tr><td>Interferent</td><td>MP6290</td></tr><tr><td>0.050</td><td>MP6291</td></tr><tr><td>0.080</td><td>MP6292</td></tr><tr><td>0.200</td><td>MP6293</td></tr></tbody></table>			Simulator	Serial Number	0.000	MP6289	Interferent	MP6290	0.050	MP6291	0.080	MP6292	0.200	MP6293																			
Simulator	Serial Number																																				
0.000	MP6289																																				
Interferent	MP6290																																				
0.050	MP6291																																				
0.080	MP6292																																				
0.200	MP6293																																				
Notes/Suggested Service: _____				Attachments																																	
				☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment			☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____																														
				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC																																	
				☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use																																	
				☒ Conduct an Agency Inspection Before Evidentiary Use																																	
				David Eliezer Reyes-Rivera Tech Review / Date			05 10:44:08 -0400																														

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: FFWCC

Time of Inspection: 12:03

Date of Inspection: 04/04/2022

Serial Number: 80-000900

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG113403 Exp: 05/14/2023
0.000	0.049	0.078	0.198	0.079
0.000	0.049	0.079	0.198	0.080
0.000	0.049	0.079	0.198	0.079
0.000	0.049	0.078	0.198	0.079
0.000	0.049	0.079	0.198	0.079
0.000	0.050	0.079	0.198	0.080
0.000	0.049	0.079	0.199	0.079
0.000	0.049	0.079	0.198	0.079
0.000	0.050	0.079	0.198	0.079
0.000	0.049	0.078	0.198	0.080

Standard Deviations	0.0004	0.0004	0.0003	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL SOTO

Signature and Printed Name

04/04/2022

Date

stability checks

FFWCC

Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000900
04/04/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:36
Control Test	0.049	09:36
Air Blank	0.000	09:37
Control Test	0.049	09:38
Air Blank	0.000	09:38
Control Test	0.049	09:39
Air Blank	0.000	09:39
Control Test Stats		
Average	0.0490	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator's Signature

FFWCC

Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000900
04/04/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:40
Control Test	0.079	09:41
Air Blank	0.000	09:42
Control Test	0.078	09:42
Air Blank	0.000	09:43
Control Test	0.078	09:44
Air Blank	0.000	09:44
Control Test Stats		
Average	0.0783	
Std Dev	0.0006	
Rel Std Dev(%)	0.7370	

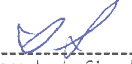
wet


Operator's Signature

FFWCC

Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000900
04/04/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:45
Control Test	0.199	09:46
Air Blank	0.000	09:46
Control Test	0.199	09:47
Air Blank	0.000	09:48
Control Test	0.199	09:48
Air Blank	0.000	09:49
Control Test Stats		
Average	0.1990	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

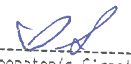

Operator's Signature

FFWCC

Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000900
04/04/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:51
Control Test	0.080	09:51
Air Blank	0.000	09:52
Control Test	0.080	09:52
Air Blank	0.000	09:52
Control Test	0.080	09:53
Air Blank	0.000	09:53
Control Test Stats		
Average	0.0800	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Dry


Operator's Signature



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000900, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000900</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>FFWCC</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>04/04/2022</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>12:03</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Israel Soto
Soto
Digitally signed by Israel Soto
Date: 2022.04.04 13:49:34 -04'00'

04/04/2022

Date

ISRAEL SOTO,

Department Inspector

FDLE/ATP Form 69 March 2022

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1