

INSTRUMENT PROCESSING SHEET

Agency Collier County Sheriff's Office

S/N 80-000942

Florida Department of
Law Enforcement


Date In 4/15/2022

DI Completion Date 4/18/2022

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By DERR	Quality Checks	By DERR	Date 4/18/2022	Flow Calibration	By	Date																																																												
☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>204</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP106</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.171</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.492</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28199</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																														
		<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td rowspan="2">MP6286</td> <td>202201C 01/11/2024</td> </tr> <tr> <td>0.080</td> <td>202201D 01/18/2024</td> </tr> <tr> <td>0.200</td> <td rowspan="2">MP6288</td> <td>202201E 01/18/2024</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A AG115904 06/08/2023</td> </tr> </tbody> </table>		Simulator	Serial #	Lot #/Exp	0.050	MP6286	202201C 01/11/2024	0.080	202201D 01/18/2024	0.200	MP6288	202201E 01/18/2024	0.080 DGS	N/A AG115904 06/08/2023	<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Maintenance</th> <th>By</th> </tr> </thead> <tbody> <tr> <td>☐ Battery Replacement</td> <td></td> </tr> <tr> <td>☐ Dry Gas Regulator Replacement</td> <td></td> </tr> <tr> <td>☐ Breath Tube Replacement</td> <td></td> </tr> <tr> <td>☐ Other _____</td> <td></td> </tr> <tr> <td>_____</td> <td></td> </tr> <tr> <td>_____</td> <td></td> </tr> <tr> <td>_____</td> <td></td> </tr> <tr> <td>_____</td> <td></td> </tr> </tbody> </table>		Maintenance	By	☐ Battery Replacement		☐ Dry Gas Regulator Replacement		☐ Breath Tube Replacement		☐ Other _____		_____		_____		_____		_____																																
Simulator	Serial #	Lot #/Exp																																																																	
0.050	MP6286	202201C 01/11/2024																																																																	
0.080		202201D 01/18/2024																																																																	
0.200	MP6288	202201E 01/18/2024																																																																	
0.080 DGS		N/A AG115904 06/08/2023																																																																	
Maintenance	By																																																																		
☐ Battery Replacement																																																																			
☐ Dry Gas Regulator Replacement																																																																			
☐ Breath Tube Replacement																																																																			
☐ Other _____																																																																			

Calibration Adjustment By _____				Department Inspection By DERR _____																																																															
Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>				Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Barometric Pressure ID# <u>68639</u> Gauge <u>1013</u> Instrument <u>1013</u> Mouth Alcohol Solution Lot # <u>2021-D</u> Acetone Stock Solution Lot # <u>2021-C</u> <table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>MP6284</td> </tr> <tr> <td>Interferent</td> <td>MP6285</td> </tr> <tr> <td>0.050</td> <td>MP6286</td> </tr> <tr> <td>0.080</td> <td>MP6287</td> </tr> <tr> <td>0.200</td> <td>MP6288</td> </tr> </tbody> </table>				Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288
Simulator	Serial #	Lot #	Expiration																																																																
0.000		N/A	N/A																																																																
0.040																																																																			
0.100																																																																			
0.200																																																																			
0.300																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial #	Lot #	Expiration																																																																
0.050																																																																			
0.080																																																																			
0.200																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial Number																																																																		
0.000	MP6284																																																																		
Interferent	MP6285																																																																		
0.050	MP6286																																																																		
0.080	MP6287																																																																		
0.200	MP6288																																																																		
Notes/Suggested Service: I had to go back to the level 2 menu after the Department Inspection to use the communications function to upload the data. DERR _____ _____ _____ _____ _____ _____ _____ _____ _____				Attachments <table border="1" style="width: 100%; border-collapse: collapse;"> <tbody> <tr> <td>☒ Form 41</td> <td>☐ Post-Stability Checks</td> </tr> <tr> <td>☒ Stability Checks</td> <td>☐ Flow Calibration</td> </tr> <tr> <td>☒ Calibration Certificate</td> <td>☐ Form 40</td> </tr> <tr> <td>☐ Calibration Adjustment</td> <td>☒ Other <u>Form 51</u></td> </tr> </tbody> </table> ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Israel Soto Digitally signed by Israel Soto Date: 2022.04.20 09:32:23 -0400				☒ Form 41	☐ Post-Stability Checks	☒ Stability Checks	☐ Flow Calibration	☒ Calibration Certificate	☐ Form 40	☐ Calibration Adjustment	☒ Other <u>Form 51</u>																																																				
☒ Form 41	☐ Post-Stability Checks																																																																		
☒ Stability Checks	☐ Flow Calibration																																																																		
☒ Calibration Certificate	☐ Form 40																																																																		
☐ Calibration Adjustment	☒ Other <u>Form 51</u>																																																																		
				Tech Review / Date _____ Admin Review / Date _____																																																															

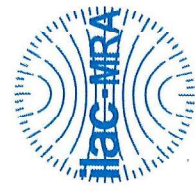

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000942 , manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000942</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>COLLIER COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>04/18/2022</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>09:29</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration. *Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

04/18/2022


DAVID E REYES-RIVERA,

Date

Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: COLLIER COUNTY SO
Time of Inspection: 09:29

Date of Inspection: 04/18/2022

Serial Number: 80-000942
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG115904 Exp: 06/08/2023
0.000	0.049	0.079	0.201	0.080
0.000	0.049	0.079	0.202	0.080
0.000	0.049	0.079	0.202	0.080
0.000	0.049	0.079	0.202	0.079
0.000	0.048	0.079	0.201	0.080
0.000	0.049	0.079	0.202	0.080
0.000	0.049	0.079	0.202	0.080
0.000	0.049	0.079	0.201	0.079
0.000	0.048	0.079	0.201	0.080
0.000	0.049	0.079	0.201	0.079


Standard Deviations	0.0004	0.0000	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

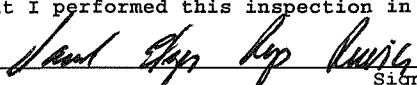

 Signature and Printed Name

DAVID E REYES-RIVERA

04/18/2022
 Date

Type of Test	Serial Number	Agency	Date	Performed By
Stabilities	80-000942	Collier County Sheriff's Office	4/18/2022	DERR <i>DM</i>

0.05g/210L 0.047 to 0.053	0.08g/210L 0.077 to 0.083	0.20g/210L 0.194 to 0.206	DGS 0.08g/210L 0.077 to 0.083
COLLIER COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 04/18/2022 Software: 8100.27 SN 80-000942 Test g/210L Time Air Blank 0.000 07:27 Control Test 0.049 07:28 Air Blank 0.000 07:28 Control Test 0.048 07:29 Air Blank 0.000 07:29 Control Test 0.048 07:30 Air Blank 0.000 07:31 Control Test Stats Average 0.0483 Std Dev 0.0006 Rel Std Dev(%) 1.1945	COLLIER COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 04/18/2022 Software: 8100.27 SN 80-000942 Test g/210L Time Air Blank 0.100 07:32 Control Test 0.178 07:33 Air Blank 0.100 07:33 Control Test 0.178 07:34 Air Blank 0.100 07:34 Control Test 0.178 07:35 Air Blank 0.100 07:36 Control Test Stats Average 0.0780 Std Dev 0.0000 Rel Std Dev(%) 0.0000	COLLIER COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 04/18/2022 Software: 8100.27 SN 80-000942 Test g/210L Time Air Blank 0.000 07:36 Control Test 0.201 07:37 Air Blank 0.000 07:38 Control Test 0.201 07:38 Air Blank 0.000 07:39 Control Test 0.200 07:40 Air Blank 0.000 07:40 Control Test Stats Average 0.2007 Std Dev 0.0005 Rel Std Dev(%) 0.2877	COLLIER COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 04/18/2022 Software: 8100.27 SN 80-000942 Test g/210L Time Air Blank 0.000 07:41 Control Test 0.079 07:42 Air Blank 0.000 07:42 Control Test 0.079 07:42 Air Blank 0.000 07:43 Control Test 0.079 07:43 Air Blank 0.000 07:44 Control Test Stats Average 0.0790 Std Dev 0.0000 Rel Std Dev(%) 0.0000
<i>DM</i> Operator's Signature	<i>DM</i> Operator's Signature	<i>DM</i> Operator's Signature	<i>DM</i> Operator's Signature

Return Material Authorization

Ship to:

☒ CMI, Inc.

☐ Enforcement Electronics

Shipment to repair facility authorized by: David Reyes on 2/28/2022

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000__ Serial Number: 80-000942_____

Bill To Address:

FDLE

4700 Terminal Drive

FT. Myers, FL 33907

Ship to Address:

FDLE

4700 Terminal Drive Suite 1

Fort Myers, FL 33907

Reason for Return:

Warranty repair, came back from CMI and according to the DVM the cell temp is 0

Please choose one of the following options:

☐ 1. I _____, authorize all repairs.

☐ 2. I _____, authorize repairs up to \$_____.

☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: _____

Phone #: _____ Email: _____

ATP Contact Name: David Reyes-Rivera Email: Davidreyes@fdle.state.fl.us

Return Material Authorization

Ship to:



CMI, Inc.



Enforcement Electronics

Shipment to repair facility authorized by: Thomas Graham _____ on 6/15/2021_

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000__ Serial Number: 80-000942_____

Bill To Address:

FDLE _____

4700 Terminal Drive _____

FT. Myers, FL 33907 _____

Apply any service credits towards this instrument from the FDLE account.

Ship to Address:

FDLE _____

4700 Terminal Drive Suite 1 _____

Fort Myers, FL 33907 _____

Reason for Return:

Instrument displays ERR 012 Contact CMI Repairs will be handled by FDLE, please apply to service credit.

Please choose one of the following options:

☐ 1. I _____, authorize all repairs.

☐ 2. I _____, authorize repairs up to \$ _____.

☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: David Whitt _____

Phone #: _____ Email: David.Whitt@colliersheriff.org _____

ATP Contact Name: David Reyes-Rivera Email: Davidreyes@fdle.state.fl.us