

INSTRUMENT PROCESSING SHEET

Agency Pinellas Park PD

S/N 80-001042

Florida Department of
Law Enforcement

Date In 1/21/2022 DI Completion Date 1/27/2022

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By TDG	Quality Checks	By TDG	Date 1/27/2022	Flow Calibration	By	Date																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value 192 ☒ Flow Verification (L/s) Flow Column # ATP104 32 mm 0.164 (.139 - .169) 36 mm 0.171 (.156 - .190) 53 mm 0.246 (.228 - .278) 103 mm 0.507 (.447 - .547) ☒ Barometric Pressure Check Gauge ID # 68639 ☒ Stability Checks			Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																														
					Maintenance By _____																														
					☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____																														
Calibration Adjustment By _____				Department Inspection By TDG _____																															
Barometric Pressure Gauge _____ ID # _____				Barometric Pressure ID# 28663																															
<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr> </thead> <tbody> <tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr> <tr><td>0.040</td><td></td><td></td><td></td></tr> <tr><td>0.100</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.300</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table>				Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge 1019 Instrument 1018 Mouth Alcohol Solution Lot # 2021-B Acetone Stock Solution Lot # 2020-A			
Simulator	Serial #	Lot #	Expiration																																
0.000		N/A	N/A																																
0.040																																			
0.100																																			
0.200																																			
0.300																																			
0.080 DGS	N/A																																		
☐ Post Calibration Adjustment Stability Checks				<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th><th>Serial Number</th></tr> </thead> <tbody> <tr><td>0.000</td><td>SD3963</td></tr> <tr><td>Interferent</td><td>SD1017</td></tr> <tr><td>0.050</td><td>MP5092</td></tr> <tr><td>0.080</td><td>MP4864</td></tr> <tr><td>0.200</td><td>MP5094</td></tr> </tbody> </table>				Simulator	Serial Number	0.000	SD3963	Interferent	SD1017	0.050	MP5092	0.080	MP4864	0.200	MP5094																
Simulator	Serial Number																																		
0.000	SD3963																																		
Interferent	SD1017																																		
0.050	MP5092																																		
0.080	MP4864																																		
0.200	MP5094																																		
<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr> </thead> <tbody> <tr><td>0.050</td><td></td><td></td><td></td></tr> <tr><td>0.080</td><td></td><td></td><td></td></tr> <tr><td>0.200</td><td></td><td></td><td></td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr> </tbody> </table>				Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other Extra Stabilities											
Simulator	Serial #	Lot #	Expiration																																
0.050																																			
0.080																																			
0.200																																			
0.080 DGS	N/A																																		
Notes/Suggested Service: Unknown Interferent detected during 0.08 DGS portion of DI. Let area air out for a few minutes and repeated successfully. Ran extra 20 Stabilities afterward using the same lot of DGS without incident (TDG)				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																															
				Israel Soto Digitally signed by Israel Soto Date: 2022.01.27 13:39:55 +05'00' Tech Review / Date Admin Review / Date																															

Type of Test	Serial Number	Agency	Date	Performed By
Stabilities	80-00 1042	Pinellas Park PD	01/27/2022	TDG <i>ML</i>

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
PINELLAS PARK PD Intoxilyzer - Alcohol Analyzer Model 8000 01/27/2022 Software: 8100.27 SN 80-001042 Test g/210L Time Air Blank 0.000 09:26 Control Test 0.048 09:27 Air Blank 0.000 09:27 Control Test 0.049 09:28 Air Blank 0.000 09:28 Control Test 0.049 09:29 Air Blank 0.000 09:30 Control Test Stats Average 0.0487 Std Dev 0.0006 Rel Std Dev(%) 1.1863	PINELLAS PARK PD Intoxilyzer - Alcohol Analyzer Model 8000 01/27/2022 Software: 8100.27 SN 80-001042 Test g/210L Time Air Blank 0.000 09:33 Control Test 0.079 09:34 Air Blank 0.000 09:34 Control Test 0.079 09:35 Air Blank 0.000 09:35 Control Test 0.079 09:36 Air Blank 0.000 09:36 Control Test Stats Average 0.0790 Std Dev 0.0000 Rel Std Dev(%) 0.0000	PINELLAS PARK PD Intoxilyzer - Alcohol Analyzer Model 8000 01/27/2022 Software: 8100.27 SN 80-001042 Test g/210L Time Air Blank 0.000 09:40 Control Test 0.198 09:41 Air Blank 0.000 09:41 Control Test 0.198 09:42 Air Blank 0.000 09:42 Control Test 0.198 09:43 Air Blank 0.000 09:43 Control Test Stats Average 0.1980 Std Dev 0.0000 Rel Std Dev(%) 0.0000	PINELLAS PARK PD Intoxilyzer - Alcohol Analyzer Model 8000 01/27/2022 Software: 8100.27 SN 80-001042 Test g/210L Time Air Blank 0.000 09:46 Control Test 0.078 09:46 Air Blank 0.000 09:46 Control Test 0.078 09:47 Air Blank 0.000 09:47 Control Test 0.077 09:47 Air Blank 0.000 09:48 Control Test Stats Average 0.0777 Std Dev 0.0006 Rel Std Dev(%) 0.7434
<i>ML</i> Operator's Signature	<i>ML</i> Operator's Signature	<i>ML</i> Operator's Signature	<i>ML</i> Operator's Signature

Comments:

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PINELLAS PARK PD
Time of Inspection: 12:20

Date of Inspection: 01/27/2022

Serial Number: 80-001042
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202010A Exp: 10/05/2022	0.08g/210L Test (g/210L) Lot#:202010B Exp: 10/05/2022	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG115904 Exp: 06/08/2023
0.000	0.048	0.077	0.197	0.078 / 0.078
0.000	0.048	0.078	0.197	0.077 / 0.077
0.000	0.048	0.078	0.197	0.076 / 0.077
0.000	0.048	0.079	0.197	0.077 / 0.077
0.000	0.048	0.078	0.197	0.076 / 0.077
0.000	0.048	0.078	0.197	0.076 / 0.077
0.000	0.047	0.079	0.197	0.076 / 0.076
0.000	0.047	0.079	0.197	0.076 / 0.076
0.000	0.048	0.079	0.197	INT / 0.077
0.000	0.048	0.079	0.197	/ 0.077

Standard Deviations	0.0004	0.0006	0.0000	/ 0.0005
---------------------	--------	--------	--------	----------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5


Remarks:

08: Interferent Detect. UNKNOWN. LET AREA AIR OUT

0.08 DGS. UNKNOWN. LET AREA AIR OUT.

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

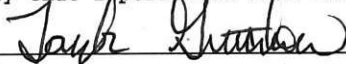

TAYLOR D GUTSCHOW

Signature and Printed Name

01/27/2022
Date

Extra Stabilities

MG 1/27/2022

PINELLAS PARK PD

Intoxilyzer - Alcohol Analyzer

Model 8000

SN 80-001042

01/27/2022


Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:22
Control Test	0.077	12:22
Air Blank	0.000	12:23
Control Test	0.078	12:23
Air Blank	0.000	12:24
Control Test	0.078	12:24
Air Blank	0.000	12:24
Control Test	0.078	12:25
Air Blank	0.000	12:25
Control Test	0.078	12:25
Air Blank	0.000	12:26
Control Test	0.078	12:26
Air Blank	0.000	12:27
Control Test	0.077	12:27
Air Blank	0.000	12:28
Control Test	0.078	12:28
Air Blank	0.000	12:28
Control Test	0.078	12:29
Air Blank	0.000	12:29
Control Test	0.077	12:30
Air Blank	0.000	12:30
Control Test	0.078	12:30
Air Blank	0.000	12:31
Control Test	0.078	12:31
Air Blank	0.000	12:32
Control Test	0.078	12:32
Air Blank	0.000	12:32
Control Test	0.078	12:33
Air Blank	0.000	12:33
Control Test	0.078	12:34
Air Blank	0.000	12:34
Control Test	0.078	12:34
Air Blank	0.000	12:35
Control Test	0.077	12:35
Air Blank	0.000	12:36
Control Test	0.078	12:36
Air Blank	0.000	12:36
Control Test	0.078	12:37
Air Blank	0.000	12:37
Control Test	0.078	12:38
Air Blank	0.000	12:38
Control Test Stats		
Average	0.0778	
Std Dev	0.0004	
Rel Std Dev(%)	0.5275	

Extra Stabilities
Page 1/2

MB

Operator's Signature

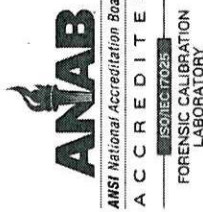

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001042, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001042</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>PINELLAS PARK PD</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>01/27/2022</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>12:20</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.


TAYLOR D GUTSCHOW,
Department Inspector

01/27/2022

Date

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality