



INSTRUMENT PROCESSING SHEET

Agency Clearwater PDS/N 80-001068Florida Department of
Law EnforcementDate In 3/22/2022 DI Completion Date 3/29/2022☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By TDG	Quality Checks	By TDG	Date <u>3/29/2022</u>	Flow Calibration	By	Date																																							
☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>187</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.156</u> (.139 - .169) 36 mm <u>0.171</u> (.156 - .190) 53 mm <u>0.242</u> (.228 - .278) 103 mm <u>0.503</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>68639</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																									
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202010A 10/05/2022</td></tr><tr><td>0.080</td><td>MP6287</td><td>202010B 10/05/2022</td></tr><tr><td>0.200</td><td>MP6288</td><td>202010D 10/06/2022</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG115904 06/08/2023</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP6286	202010A 10/05/2022	0.080	MP6287	202010B 10/05/2022	0.200	MP6288	202010D 10/06/2022	0.080 DGS	N/A	AG115904 06/08/2023	Maintenance By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____ _____																										
Simulator	Serial #	Lot #/Exp																																												
0.050	MP6286	202010A 10/05/2022																																												
0.080	MP6287	202010B 10/05/2022																																												
0.200	MP6288	202010D 10/06/2022																																												
0.080 DGS	N/A	AG115904 06/08/2023																																												
Calibration Adjustment By _____		Department Inspection By <u>TDG</u>																																												
Barometric Pressure Gauge _____ ID # _____		Barometric Pressure ID# <u>28663</u> Gauge <u>1022</u> Instrument <u>1023</u> Mouth Alcohol Solution Lot # <u>2021-B</u> Acetone Stock Solution Lot # <u>2021-C</u>																																												
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table>					Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288
Simulator	Serial #	Lot #	Expiration																																											
0.000		N/A	N/A																																											
0.040																																														
0.100																																														
0.200																																														
0.300																																														
0.080 DGS	N/A																																													
Simulator	Serial Number																																													
0.000	MP6284																																													
Interferent	MP6285																																													
0.050	MP6286																																													
0.080	MP6287																																													
0.200	MP6288																																													
☐ Post Calibration Adjustment Stability Checks		Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☒ Form 40 ☐ Calibration Adjustment ☒ Other <u>Form 51 and 2021 IPS</u>																																												
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																								
Simulator	Serial #	Lot #	Expiration																																											
0.050																																														
0.080																																														
0.200																																														
0.080 DGS	N/A																																													
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		David Eliezer Reyes-Rivera Digitally signed by David Eliezer Reyes-Rivera Date: 2022.03.30 10:18:36 -0400 13:15:3 6 Tech Review / Date _____ Admin Review / Date _____																																												

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: CLEARWATER PD
Time of Inspection: 09:49

Date of Inspection: 03/29/2022

Serial Number: 80-001068
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

BYPASSED AI TO OPERATE INSTRUMENT

Not determined 716
03/29/2022

The above instrument complies (☒) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Taylor D Gutschow

TAYLOR D GUTSCHOW

Signature and Printed Name

03/29/2022
Date

Type of Test	Serial Number	Agency	Date	Performed By
Stabilities	80-00 1068	Clearwater PD	03/29/2022	TDG MC

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
CLEARWATER PD Intoxilyzer - Alcohol Analyzer Model 8000 03/29/2022 Software: 8100.27 SN 80-001068 Test g/210L Time Air Blank 0.000 10:58 Control Test 0.051 10:59 Air Blank 0.000 11:00 Control Test 0.049 11:00 Air Blank 0.000 11:01 Control Test 0.050 11:01 Air Blank 0.000 11:02 Control Test Stats Average 0.0500 Std Dev 0.0010 Rel Std Dev(%) 2.0000 Operator's Signature <i>MC</i>	CLEARWATER PD Intoxilyzer - Alcohol Analyzer Model 8000 03/29/2022 Software: 8100.27 SN 80-001068 Test g/210L Time Air Blank 0.000 11:07 Control Test 0.080 11:07 Air Blank 0.000 11:08 Control Test 0.080 11:09 Air Blank 0.000 11:09 Control Test 0.080 11:10 Air Blank 0.000 11:10 Control Test Stats Average 0.0800 Std Dev 0.0000 Rel Std Dev(%) 0.0000 Operator's Signature <i>MC</i>	CLEARWATER PD Intoxilyzer - Alcohol Analyzer Model 8000 03/29/2022 Software: 8100.27 SN 80-001068 Test g/210L Time Air Blank 0.000 11:15 Control Test 0.198 11:16 Air Blank 0.000 11:16 Control Test 0.198 11:17 Air Blank 0.000 11:18 Control Test 0.198 11:18 Air Blank 0.000 11:19 Control Test Stats Average 0.1980 Std Dev 0.0000 Rel Std Dev(%) 0.0000 Operator's Signature <i>MC</i>	CLEARWATER PD Intoxilyzer - Alcohol Analyzer Model 8000 03/29/2022 Software: 8100.27 SN 80-001068 Test g/210L Time Air Blank 0.000 10:53 Control Test 0.077 10:53 Air Blank 0.000 10:54 Control Test 0.077 10:54 Air Blank 0.000 10:55 Control Test 0.079 10:55 Air Blank 0.000 10:55 Control Test Stats Average 0.0777 Std Dev 0.0012 Rel Std Dev(%) 1.4867 Operator's Signature <i>MC</i>

Comments:

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: CLEARWATER PD
Time of Inspection: 14:06

Date of Inspection: 03/29/2022

Serial Number: 80-001068
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202010A Exp: 10/05/2022	0.08g/210L Test (g/210L) Lot#:202010B Exp: 10/05/2022	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG115904 Exp: 06/08/2023
0.000	0.050	0.080	0.199	0.078
0.000	0.050	0.079	0.199	0.078
0.000	0.050	0.080	0.199	0.078
0.000	0.049	0.079	0.199	0.078
0.000	0.049	0.079	0.199	0.078
0.000	0.049	0.080	0.199	0.078
0.000	0.049	0.079	0.199	0.078
0.000	0.049	0.079	0.199	0.078
0.000	0.049	0.079	0.199	0.078
0.000	0.049	0.079	0.199	0.078
0.000	0.049	0.079	0.198	0.077

Standard Deviations	0.0004	0.0004	0.0003	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5


Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

 TAYLOR D GUTSCHOW
Signature and Printed Name

03/29/2022
Date

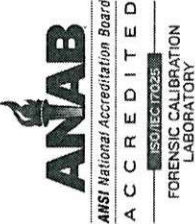

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001068, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001068</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>CLEARWATER PD</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>03/29/2022</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>14:06</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

03/29/2022

Date



TAYLOR D GUTSCHOW,

Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality



INSTRUMENT PROCESSING SHEET

Agency Clearwater PD

S/N 80-001068

Florida Department of
Law Enforcement

Date In 7/16/2021

DI Completion Date _____

☐ Ship

☐ P/U

☐ H/D

☐ CMI

☒ EE

Intake	By TDG	Quality Checks	By TDG	Date 7/19/2021	Flow Calibration	By	Date																																																										
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☐ Instrument Set Up Verified ☐ R-Value _____ ☐ Flow Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) ☐ Barometric Pressure Check Gauge ID # _____ ☐ Stability Checks			Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																												
		<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> </tr> </tbody> </table>		Simulator	Serial #	Lot #/Exp	0.050			0.080			0.200			0.080 DGS	N/A		Maintenance By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ DI Temp. Checks By _____ ☐ Lab Temp °C _____ External Digital Therm. ID#: _____ ☐ 34°C +/- .2 Serial #: _____ ☐ 34°C +/- .2 Serial #: _____ ☐ 34°C +/- .2 Serial #: _____																																														
Simulator	Serial #	Lot #/Exp																																																															
0.050																																																																	
0.080																																																																	
0.200																																																																	
0.080 DGS	N/A																																																																
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection By _____ Barometric Pressure ID# _____ Gauge _____ Instrument _____ Mouth Alcohol Solution Lot # _____ Acetone Stock Solution Lot # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> </tr> <tr> <td>Interferent</td> <td></td> </tr> <tr> <td>0.050</td> <td></td> </tr> <tr> <td>0.080</td> <td></td> </tr> <tr> <td>0.200</td> <td></td> </tr> </tbody> </table>				Simulator	Serial Number	0.000		Interferent		0.050		0.080		0.200	
Simulator	Serial #	Lot #	Expiration																																																														
0.000		N/A	N/A																																																														
0.040																																																																	
0.100																																																																	
0.200																																																																	
0.300																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial #	Lot #	Expiration																																																														
0.050																																																																	
0.080																																																																	
0.200																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial Number																																																																
0.000																																																																	
Interferent																																																																	
0.050																																																																	
0.080																																																																	
0.200																																																																	
Notes/Suggested Service: <u>Sending to repair. Instrument will not power on. (TDG)</u> _____ _____ _____ _____ _____ _____ _____ _____		Attachments ☐ Form 41 ☐ Stability Checks ☐ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____ ☐ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☐ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☐ Conduct an Agency Inspection Before Evidentiary Use Tech Review / Date _____ Admin Review / Date _____																																																															

Return Material Authorization

Ship to:

☐ CMI, Inc.

☒ Enforcement Electronics

Shipment to repair facility authorized by: Dave Nugent on 7/19/2021

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-001068

Bill To Address:

Clearwater Police Department

Attn: Dave Nugent

645 Pierce Street

Clearwater, FL 33756

Ship to Address:

Florida Department of Law Enforcement

Fort Myers Regional Operations Center

Attn: Taylor Gutschow

4700 Terminal Drive, Suite 1

Fort Myers, FL 33907

Reason for Return:

Instrument will not power on. After pressing the green Start Test button, there is a long pause
before the light turns from red to orange, then there is a clicking sound and the light blinks off.
After another clicking sound, the red light turns back on.

Please choose one of the following options:

☐ 1. I _____, authorize all repairs.

☐ 2. I _____, authorize repairs up to \$_____.

☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Dave Nugent

Phone #: 727-492-0058

Email: Dave.Nugent@MyClearwater.com

ATP Contact Name: Taylor Gutschow

ATP Email: TaylorGutschow@fdle.state.fl.us



INSTRUMENT PROCESSING SHEET

Agency Clearwater PDS/N 80-001068Florida Department of
Law EnforcementDate In 7/29/2022DI Completion Date 7/29/2022☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By TDG	Quality Checks	By TDG	Date <u>7/29/2022</u>	Flow Calibration	By	Date																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: AI reports the instrument failed the Temp Regulation part of the Diagnostics and won't go into Ready Mode.		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>175</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP106</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.507</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>68639</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																														
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5092</td><td>202201C 01/11/2024</td></tr><tr><td>0.080</td><td>MP5093</td><td>202201D 01/18/2024</td></tr><tr><td>0.200</td><td>MP5094</td><td>202201E 01/18/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>00521080A2 02/05/2023</td></tr></tbody></table>					Simulator	Serial #	Lot #/Exp	0.050	MP5092	202201C 01/11/2024	0.080	MP5093	202201D 01/18/2024	0.200	MP5094	202201E 01/18/2024	0.080 DGS	N/A	00521080A2 02/05/2023	Maintenance By <u>TDG</u> ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☒ Breath Tube Replacement ☐ Other _____															
Simulator	Serial #	Lot #/Exp																																	
0.050	MP5092	202201C 01/11/2024																																	
0.080	MP5093	202201D 01/18/2024																																	
0.200	MP5094	202201E 01/18/2024																																	
0.080 DGS	N/A	00521080A2 02/05/2023																																	
Calibration Adjustment By _____					Department Inspection By <u>TDG</u>																														
Barometric Pressure Gauge _____ ID # _____					Barometric Pressure ID# <u>28199</u>																														
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>					Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge <u>1019</u> Instrument <u>1021</u> Mouth Alcohol Solution Lot # <u>2021-D</u> Acetone Stock Solution Lot # <u>2021-C</u>		
Simulator	Serial #	Lot #	Expiration																																
0.000		N/A	N/A																																
0.040																																			
0.100																																			
0.200																																			
0.300																																			
0.080 DGS	N/A																																		
☐ Post Calibration Adjustment Stability Checks					<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5095</td></tr><tr><td>Interferent</td><td>SD3966</td></tr><tr><td>0.050</td><td>MP5092</td></tr><tr><td>0.080</td><td>MP5093</td></tr><tr><td>0.200</td><td>MP5094</td></tr></tbody></table>			Simulator	Serial Number	0.000	MP5095	Interferent	SD3966	0.050	MP5092	0.080	MP5093	0.200	MP5094																
Simulator	Serial Number																																		
0.000	MP5095																																		
Interferent	SD3966																																		
0.050	MP5092																																		
0.080	MP5093																																		
0.200	MP5094																																		
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>					Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☐ Other _____										
Simulator	Serial #	Lot #	Expiration																																
0.050																																			
0.080																																			
0.200																																			
0.080 DGS	N/A																																		
Notes/Suggested Service: <u>Breath tube was replaced prior to Quality Checks. The DVM screen indicated it wasn't heating up properly. (TDG)</u>					☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																														
David Eliezer Reyes-Rivera Digitally signed by David Eliezer Reyes-Rivera Date: 2022.07.29 14:34:34 -04'00'					Date: 2022.07.29 15:29:50 -04'00'																														
Tech Review / Date					Admin Review / Date																														

Type of Test	Serial Number	Agency	Date	Performed By
Stabilities	80-001068	Clearwater PD	07/29/2022	TDG MG

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
✓	✓	✓	✓
CLEARWATER PD Intoxilyzer - Alcohol Analyzer Model 8000 07/29/2022 Software: 8100.27 SN 80-001068 Test g/210L Time Air Blank 0.000 09:56 Control Test 0.050 09:57 Air Blank 0.000 09:58 Control Test 0.050 09:58 Air Blank 0.000 09:59 Control Test 0.050 10:00 Air Blank 0.000 10:00 Control Test Stats Average 0.0500 Std Dev 0.0000 Rel Std Dev(%) 0.0000	CLEARWATER PD Intoxilyzer - Alcohol Analyzer Model 8000 07/29/2022 Software: 8100.27 SN 80-001068 Test g/210L Time Air Blank 0.000 10:01 Control Test 0.079 10:02 Air Blank 0.000 10:02 Control Test 0.079 10:03 Air Blank 0.000 10:04 Control Test 0.079 10:04 Air Blank 0.000 10:05 Control Test Stats Average 0.0790 Std Dev 0.0000 Rel Std Dev(%) 0.0000	CLEARWATER PD Intoxilyzer - Alcohol Analyzer Model 8000 07/29/2022 Software: 8100.27 SN 80-001068 Test g/210L Time Air Blank 0.000 10:06 Control Test 0.199 10:06 Air Blank 0.000 10:07 Control Test 0.199 10:08 Air Blank 0.000 10:08 Control Test 0.199 10:09 Air Blank 0.000 10:09 Control Test Stats Average 0.1990 Std Dev 0.0000 Rel Std Dev(%) 0.0000	CLEARWATER PD Intoxilyzer - Alcohol Analyzer Model 8000 07/29/2022 Software: 8100.27 SN 80-001068 Test g/210L Time Air Blank 0.000 10:14 Control Test 0.079 10:14 Air Blank 0.000 10:14 Control Test 0.079 10:15 Air Blank 0.000 10:15 Control Test 0.078 10:16 Air Blank 0.000 10:16 Control Test Stats Average 0.0787 Std Dev 0.0006 Rel Std Dev(%) 0.7339
Operator's Signature MG	Operator's Signature MG	Operator's Signature MG	Operator's Signature MG

Comments:

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: CLEARWATER PD
Time of Inspection: 12:00

Date of Inspection: 07/29/2022

Serial Number: 80-001068
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:00521080A2 Exp: 02/05/2023
0.000	0.049	0.079	0.199	0.079
0.000	0.049	0.079	0.200	0.078
0.000	0.049	0.079	0.200	0.079
0.000	0.050	0.079	0.200	0.079
0.000	0.050	0.079	0.199	0.079
0.000	0.049	0.079	0.199	0.078
0.000	0.049	0.078	0.200	0.079
0.000	0.049	0.079	0.199	0.078
0.000	0.049	0.078	0.199	0.078
0.000	0.049	0.079	0.200	0.078


Standard Deviations	0.0004	0.0004	0.0005	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

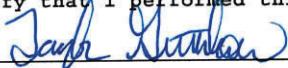

TAYLOR D GUTSCHOW

Signature and Printed Name

07/29/2022
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001068, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-001068	UNCERTAINTY* $\pm$
Owning Agency:	CLEARWATER PD	0.050 g/ 210 L 0.004
Calibration Date:	07/29/2022	0.080 g/ 210 L 0.004
Calibration Time:	12:00	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.



07/29/2022

Date

TAYLOR D GUTSCHOW,

Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality