



INSTRUMENT PROCESSING SHEET

Agency Ponce Inlet PD

S/N 80-001158

Florida Department of
Law Enforcement

Date In 02/18/22

DI Completion Date 02-23-2022

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>RAW</u> ☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>Agency indicated</u> <u>difficulties uploading instrument</u> <u>and that instrument may have</u> <u>a modem issue. IS</u>	Quality Checks By <u>IS</u> Date <u>02-23-2022</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>214</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-102</u> 32 mm <u>0.156</u> (.139 - .169) 36 mm <u>0.167</u> (.156 - .190) 53 mm <u>0.238</u> (.228 - .278) 103 mm <u>0.500</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks	Flow Calibration By _____ Date _____ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																													
	<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>MP5088</td> <td>202201C 01-11-2024</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> <td>202201D 01-18-2024</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> <td>202010D 10-06-2022</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG113403 05-14-2023</td> </tr> </tbody> </table>	Simulator	Serial #	Lot #/Exp	0.050	MP5088	202201C 01-11-2024	0.080	MP5089	202201D 01-18-2024	0.200	MP5090	202010D 10-06-2022	0.080 DGS	N/A	AG113403 05-14-2023	Maintenance By <u>IS</u> ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☒ Other <u>Diagnostic Check to check</u> <u>modem.</u>																																														
Simulator	Serial #	Lot #/Exp																																																													
0.050	MP5088	202201C 01-11-2024																																																													
0.080	MP5089	202201D 01-18-2024																																																													
0.200	MP5090	202010D 10-06-2022																																																													
0.080 DGS	N/A	AG113403 05-14-2023																																																													
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection By <u>IS</u> Barometric Pressure ID# <u>28421</u> Gauge <u>1025</u> Instrument <u>1024</u> Mouth Alcohol Solution Lot # <u>2021-D</u> Acetone Stock Solution Lot # <u>2021-C</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>MP5086</td> </tr> <tr> <td>Interferent</td> <td>MP5087</td> </tr> <tr> <td>0.050</td> <td>MP5088</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> </tr> </tbody> </table>		Simulator	Serial Number	0.000	MP5086	Interferent	MP5087	0.050	MP5088	0.080	MP5089	0.200	MP5090
Simulator	Serial #	Lot #	Expiration																																																												
0.000		N/A	N/A																																																												
0.040																																																															
0.100																																																															
0.200																																																															
0.300																																																															
0.080 DGS	N/A																																																														
Simulator	Serial #	Lot #	Expiration																																																												
0.050																																																															
0.080																																																															
0.200																																																															
0.080 DGS	N/A																																																														
Simulator	Serial Number																																																														
0.000	MP5086																																																														
Interferent	MP5087																																																														
0.050	MP5088																																																														
0.080	MP5089																																																														
0.200	MP5090																																																														
Notes/Suggested Service: <u>Instrument passed all modem</u> <u>checks during Diagnostic Checks. Uploaded instrument</u> <u>through phone line with no issue. IS</u>		Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☒ Other <u>Diagnostic Check</u>																																																													
		☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																																													
		David Eliezer Reyes-Rivera Digitally signed by David Eliezer Reyes-Rivera Date: 2022.02.23 14:29:47 -05'00' Tech Review / Date _____ Admin Review / Date _____																																																													

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PONCE INLET PD
Time of Inspection: 13:21

Date of Inspection: 02/23/2022

Serial Number: 80-001158
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG113403 Exp: 05/14/2023
0.000	0.049	0.078	0.199	0.079
0.000	0.049	0.078	0.199	0.080
0.000	0.049	0.078	0.198	0.080
0.000	0.049	0.078	0.198	0.079
0.000	0.049	0.078	0.199	0.079
0.000	0.049	0.078	0.199	0.079
0.000	0.049	0.078	0.199	0.080
0.000	0.049	0.078	0.198	0.079
0.000	0.049	0.078	0.199	0.079
0.000	0.049	0.078	0.198	0.079


Standard Deviations	0.0000	0.0000	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

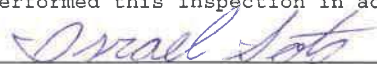

Signature and Printed Name

ISRAEL SOTO

02/23/2022
Date

Stability Checks

PONCE INLET PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001156
02/23/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:37
Control Test	0.048	10:37
Air Blank	0.000	10:38
Control Test	0.049	10:38
Air Blank	0.000	10:39
Control Test	0.049	10:40
Air Blank	0.000	10:40
Control Test Stats		
Average	0.0497	
Std Dev	0.0006	
Rel Std Dev(%)	1.1863	


Operator's Signature

PONCE INLET PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001158
02/23/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:54
Control Test	0.079	10:55
Air Blank	0.000	10:55
Control Test	0.079	10:56
Air Blank	0.000	10:56
Control Test	0.079	10:56
Air Blank	0.000	10:57
Control Test Stats		
Average	0.0798	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Dry


Operator's Signature

PONCE INLET PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001158
02/23/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:49
Control Test	0.199	10:50
Air Blank	0.000	10:50
Control Test	0.199	10:51
Air Blank	0.000	10:51
Control Test	0.199	10:52
Air Blank	0.000	10:53
Control Test Stats		
Average	0.1990	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator's Signature

PONCE INLET PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001158
02/23/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:44
Control Test	0.079	10:45
Air Blank	0.000	10:45
Control Test	0.078	10:46
Air Blank	0.000	10:46
Control Test	0.079	10:47
Air Blank	0.000	10:47
Control Test Stats		
Average	0.0783	
Std Dev	0.0006	
Rel Std Dev(%)	0.7370	

Dry
wet


Operator's Signature

Diagnostic Check to Check Modem



PONCE INLET PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001158
02/23/2022
Software: 8100.27

DIAGNOSTICS

Voltage/Current Test	OK
RAM Test	OK
EEPROM Checksum Test	OK
Real Time Clock Test	OK
DSP Test	OK
Analytical Stability Test	OK
Internal Printer Test	OK
Modem Test	OK
Temperature Regulation Test	OK



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2729 Fort Knox Blvd.
Bldg. 2, Suite 1300
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001158, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001158</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>PONCE INLET PD</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>02/23/2022</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>13:21</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Israel Soto
Soto
Date: 2022.02.23 14:02:58
-05'00'

Date

02/23/2022

ISRAEL SOTO,

Department Inspector

FDLE/ATP Form 69 January 2022

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1