



INSTRUMENT PROCESSING SHEET

Agency Santa Rosa Sheriff's OfficeS/N 80-001302Florida Department of
Law EnforcementDate In 10-17-2022 DI Completion Date 10-18-2022☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By IS	Quality Checks	By IS	Date <u>10-18-2022</u>	Flow Calibration	By	Date																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>133</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-105</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.503</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28427</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																														
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202201C 01-11-2024</td></tr><tr><td>0.080</td><td>MP6292</td><td>202201D 01-18-2024</td></tr><tr><td>0.200</td><td>MP6293</td><td>202201E 01-18-2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>08121080A1 05-05-2023</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP6291	202201C 01-11-2024	0.080	MP6292	202201D 01-18-2024	0.200	MP6293	202201E 01-18-2024	0.080 DGS	N/A	08121080A1 05-05-2023																		
Simulator	Serial #	Lot #/Exp																																	
0.050	MP6291	202201C 01-11-2024																																	
0.080	MP6292	202201D 01-18-2024																																	
0.200	MP6293	202201E 01-18-2024																																	
0.080 DGS	N/A	08121080A1 05-05-2023																																	
Calibration Adjustment		By IS		Department Inspection			By IS																												
Barometric Pressure Gauge <u>1010/1011</u> ID # <u>28421</u>				Barometric Pressure ID# <u>30793</u>																															
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td>MP5091</td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td>MP5082</td><td>21070</td><td>03-01-2023</td></tr><tr><td>0.100</td><td>MP5083</td><td>21380</td><td>09-13-2023</td></tr><tr><td>0.200</td><td>MP6297</td><td>20510</td><td>12-03-2022</td></tr><tr><td>0.300</td><td>MP5085</td><td>21420</td><td>10-20-2023</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG113403</td><td>05-14-2023</td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000	MP5091	N/A	N/A	0.040	MP5082	21070	03-01-2023	0.100	MP5083	21380	09-13-2023	0.200	MP6297	20510	12-03-2022	0.300	MP5085	21420	10-20-2023	0.080 DGS	N/A	AG113403	05-14-2023			Gauge <u>1010</u> Instrument <u>1010</u> Mouth Alcohol Solution Lot # <u>2022-A</u> Acetone Stock Solution Lot # <u>2022-B</u>			
Simulator	Serial #	Lot #	Expiration																																
0.000	MP5091	N/A	N/A																																
0.040	MP5082	21070	03-01-2023																																
0.100	MP5083	21380	09-13-2023																																
0.200	MP6297	20510	12-03-2022																																
0.300	MP5085	21420	10-20-2023																																
0.080 DGS	N/A	AG113403	05-14-2023																																
☒ Post Calibration Adjustment Stability Checks				<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6289</td></tr><tr><td>Interferent</td><td>MP6290</td></tr><tr><td>0.050</td><td>MP6291</td></tr><tr><td>0.080</td><td>MP6292</td></tr><tr><td>0.200</td><td>MP6293</td></tr></tbody></table>			Simulator	Serial Number	0.000	MP6289	Interferent	MP6290	0.050	MP6291	0.080	MP6292	0.200	MP6293			Attachments														
Simulator	Serial Number																																		
0.000	MP6289																																		
Interferent	MP6290																																		
0.050	MP6291																																		
0.080	MP6292																																		
0.200	MP6293																																		
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202201C</td><td>01-11-2024</td></tr><tr><td>0.080</td><td>MP6292</td><td>202201D</td><td>01-18-2024</td></tr><tr><td>0.200</td><td>MP6293</td><td>202201E</td><td>01-18-2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>08121080A1</td><td>05-05-2023</td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.050	MP6291	202201C	01-11-2024	0.080	MP6292	202201D	01-18-2024	0.200	MP6293	202201E	01-18-2024	0.080 DGS	N/A	08121080A1	05-05-2023			☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☒ Calibration Adjustment X2			☒ Post-Stability Checks X2 ☐ Flow Calibration ☐ Form 40 ☐ Other _____								
Simulator	Serial #	Lot #	Expiration																																
0.050	MP6291	202201C	01-11-2024																																
0.080	MP6292	202201D	01-18-2024																																
0.200	MP6293	202201E	01-18-2024																																
0.080 DGS	N/A	08121080A1	05-05-2023																																
Notes/Suggested Service: <u>Pressure readout from instrument more than 1% different than barometric pressure gauge. Performed Optical Bench Calibration Adjustment. Post Stability results not nominal. Repeated Optical Bench Calibration Adjust using same simulator serial #s and solutions. IS</u>				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																															
				Phil Nicodemo Tech Review / Date			24 14:26:14 -04'00" Admin Review / Date																												

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: SANTA ROSA COUNTY SO
Time of Inspection: 13:11

Date of Inspection: 10/18/2022

Serial Number: 80-001302
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:08121080A1 Exp: 05/05/2023
0.000	0.050	0.079	0.201	0.079
0.000	0.049	0.079	0.201	0.079
0.000	0.049	0.080	0.201	0.079
0.000	0.049	0.079	0.201	0.079
0.000	0.050	0.079	0.201	0.079
0.000	0.050	0.080	0.201	0.079
0.000	0.049	0.079	0.200	0.079
0.000	0.050	0.080	0.200	0.079
0.000	0.049	0.080	0.201	0.079
0.000	0.050	0.080	0.201	0.079

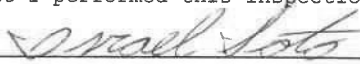
Standard Deviations	0.0005	0.0005	0.0004	0.0000
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

 ISRAEL SOTO
Signature and Printed Name

10/18/2022
Date

Stability Checks

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:10
Control Test	0.050	07:11
Air Blank	0.000	07:12
Control Test	0.050	07:12
Air Blank	0.000	07:13
Control Test	0.049	07:14
Air Blank	0.000	07:14
Control Test Stats		
Average	0.0497	
Std Dev	0.0006	
Rel Std Dev(%)	1.1625	

Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:15
Control Test	0.079	07:16
Air Blank	0.000	07:16
Control Test	0.079	07:17
Air Blank	0.000	07:18
Control Test	0.080	07:18
Air Blank	0.000	07:19
Control Test Stats		
Average	0.0793	
Std Dev	0.0006	
Rel Std Dev(%)	0.7277	

wet

Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:21
Control Test	0.201	07:22
Air Blank	0.000	07:23
Control Test	0.202	07:23
Air Blank	0.000	07:24
Control Test	0.200	07:24
Air Blank	0.000	07:25
Control Test Stats		
Average	0.2010	
Std Dev	0.0010	
Rel Std Dev(%)	0.4975	

Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:27
Control Test	0.080	07:28
Air Blank	0.000	07:28
Control Test	0.080	07:29
Air Blank	0.000	07:29
Control Test	0.080	07:29
Air Blank	0.000	07:30
Control Test Stats		
Average	0.0800	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

DF/

Operator's Signature



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001302, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001302</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>SANTA ROSA COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>10/18/2022</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>13:11</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Israel Soto
Soto
Date: 2022.10.18
14:19:24 -04'00'

10/18/2022

Date

ISRAEL SOTO,

Department Inspector

FDLE/ATP Form 69 March 2022
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022 08:01:46

Optical Bench
Calibration Adjustment #1


Auto Calibration

pg 1 of 2

	<<<<<	3um	>>>>>	<<<<<	9um	>>>>>

Solution = 0.000 g/210L or 0.0000 mg/l, Samples = 4, Discarded = 1						
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	0.1010	(-0.0390)		0.0920	(-0.0160)	
Sample #2	0.0990	(0.0000)		0.0960	(-0.0180)	
Sample #3	0.1020	(-0.0160)		0.1220	(-0.0300)	
Sample #4	0.0970	(0.0130)		0.0790	(-0.0030)	
Avg % Abs	0.0993	(-0.0010)		0.0990	(-0.0170)	
STD DEV	0.0025	(0.0145)		0.0217	(0.0135)	
REL STD DEV	2.534	(1452.584)		21.875	(79.575)	

Solution = 0.040 g/210L or 0.1905 mg/l, Samples = 4, Discarded = 1						
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	0.8780	(-0.0290)		1.5550	(-0.0120)	
Sample #2	0.8240	(-0.0070)		1.5390	(-0.0080)	
Sample #3	0.8640	(0.0000)		1.5490	(-0.0030)	
Sample #4	0.8540	(0.0020)		1.5580	(-0.0070)	
Avg % Abs	0.8473	(-0.0017)		1.5487	(-0.0060)	
STD DEV	0.0208	(0.0047)		0.0095	(0.0026)	
REL STD DEV	2.457	(283.549)		0.614	(44.096)	

Solution = 0.100 g/210L or 0.4762 mg/l, Samples = 4, Discarded = 1						
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	1.9270	(-0.0160)		3.6720	(-0.0110)	
Sample #2	1.9630	(-0.0190)		3.6850	(-0.0130)	
Sample #3	1.9320	(-0.0190)		3.6870	(-0.0370)	
Sample #4	1.9420	(0.0080)		3.6790	(-0.0030)	
Avg % Abs	1.9457	(-0.0100)		3.6837	(-0.0177)	
STD DEV	0.0158	(0.0156)		0.0042	(0.0175)	
REL STD DEV	0.813	(155.885)		0.113	(98.908)	

Solution = 0.200 g/210L or 0.9524 mg/l, Samples = 4, Discarded = 1						
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	3.7070	(-0.0020)		6.9900	(-0.0150)	
Sample #2	3.7240	(-0.0110)		7.0020	(-0.0120)	
Sample #3	3.7080	(0.0170)		6.9700	(-0.0030)	
Sample #4	3.6990	(-0.0020)		7.0030	(-0.0240)	
Avg % Abs	3.7103	(0.0013)		6.9917	(-0.0130)	
STD DEV	0.0127	(0.0143)		0.0188	(0.0105)	
REL STD DEV	0.341	(1072.089)		0.268	(81.043)	

Solution = 0.300 g/210L or 1.4286 mg/l, Samples = 4, Discarded = 1						
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	5.4040	(-0.0150)		10.1420	(-0.0050)	
Sample #2	5.3640	(-0.0070)		10.1140	(0.0020)	
Sample #3	5.3870	(0.0020)		10.1040	(0.0120)	
Sample #4	5.3570	(0.0000)		10.1060	(0.0000)	
Avg % Abs	5.3693	(-0.0017)		10.1080	(0.0047)	
STD DEV	0.0157	(0.0047)		0.0053	(0.0064)	
REL STD DEV	0.292	(283.549)		0.052	(137.766)	

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022 08:01:46

Optical Bench
calibration Adjust
#1

Auto Calibration

pg 2 of 2

<<<< 3um >>>>

Zero Order Coef -237.46
First Order Coef 2483.45
Second Order Coef 41.00

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.000	-0.0002
0.040	0.040	0.0002
0.100	0.100	0.0003
0.200	0.200	-0.0004
0.300	0.300	0.0001

<<<< 9um >>>>

-119.63
1273.62
15.00

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.000	-0.0001
0.040	0.040	0.0003
0.100	0.100	-0.0003
0.200	0.200	0.0001
0.300	0.300	-0.0000

<<<< 3um >>>> <<<< 9um >>>>

Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1
Sample
Sample #1 3275.00 3552.00
Sample #2 3296.00 3538.00
Sample #3 3268.00 3541.00
Sample #4 3210.00 3514.00
Avg 3258.0000 3531.0000
STD DEV 43.8634 14.7986
REL STD DEV 1.346 0.419
H2O adjust (mg/l*10k) 551 278

Barometric Pressure = 1010

*****CALIBRATION SUCCESSFUL*****

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022 10:11:16

Optical Bench
Calibration Adjustment #2

JS

Auto Calibration

pg 1 of 2

<<<<< 3um >>>>>			<<<<< 9um >>>>>		

Solution = 0.000 g/210L or 0.0000 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	0.0840	(-0.0140)	0.0640	(-0.0110)	
Sample #2	0.1190	(0.0060)	0.0970	(-0.0130)	
Sample #3	0.1060	(0.0040)	0.0940	(-0.0250)	
Sample #4	0.0930	(0.0400)	0.0830	(-0.0070)	
Avg % Abs	0.1060	(0.0167)	0.0913	(-0.0150)	
STD DEV	0.0130	(0.0202)	0.0074	(0.0092)	
REL STD DEV	12.264	(121.392)	8.071	(61.101)	

Solution = 0.040 g/210L or 0.1905 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	0.8160	(0.0090)	1.5310	(0.0000)	
Sample #2	0.8040	(0.0290)	1.5180	(0.0050)	
Sample #3	0.8350	(0.0100)	1.5360	(0.0030)	
Sample #4	0.8060	(0.0350)	1.5250	(-0.0020)	
Avg % Abs	0.8150	(0.0247)	1.5263	(0.0020)	
STD DEV	0.0173	(0.0131)	0.0091	(0.0036)	
REL STD DEV	2.129	(52.910)	0.594	(180.278)	

Solution = 0.100 g/210L or 0.4762 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	1.9320	(-0.0060)	3.6910	(-0.0050)	
Sample #2	1.9330	(0.0260)	3.6790	(0.0070)	
Sample #3	1.9270	(0.0220)	3.6640	(0.0240)	
Sample #4	1.9050	(0.0490)	3.6280	(0.0470)	
Avg % Abs	1.9217	(0.0323)	3.6570	(0.0260)	
STD DEV	0.0147	(0.0146)	0.0262	(0.0201)	
REL STD DEV	0.767	(45.067)	0.717	(77.211)	

Solution = 0.200 g/210L or 0.9524 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	3.7270	(-0.0070)	7.0430	(0.0000)	
Sample #2	3.7100	(0.0410)	7.0030	(0.0480)	
Sample #3	3.6900	(0.0220)	6.9860	(0.0490)	
Sample #4	3.6500	(0.0620)	6.9250	(0.0970)	
Avg % Abs	3.6833	(0.0417)	6.9713	(0.0647)	
STD DEV	0.0306	(0.0200)	0.0410	(0.0280)	
REL STD DEV	0.829	(48.020)	0.588	(43.308)	

Solution = 0.300 g/210L or 1.4286 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	5.4230	(-0.0160)	10.2260	(-0.0040)	
Sample #2	5.3860	(0.0110)	10.1480	(0.0690)	
Sample #3	5.3590	(0.0400)	10.1490	(0.0810)	
Sample #4	5.3890	(0.0240)	10.1770	(0.0480)	
Avg % Abs	5.3780	(0.0250)	10.1580	(0.0660)	
STD DEV	0.0165	(0.0145)	0.0165	(0.0167)	
REL STD DEV	0.307	(58.103)	0.162	(25.308)	

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022 10:11:16

Optical Bench
Calibration Adjustment #2
SF

Auto Calibration

pg 2 of 2

<<<< 3um >>>>

Zero Order Coef -243.16
First Order Coef 2558.71
Second Order Coef 26.31

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.001	-0.0006
0.040	0.039	0.0009
0.100	0.100	-0.0002
0.200	0.200	-0.0003
0.300	0.300	0.0002

<<<< 9um >>>>

-116.68
1297.00
11.94

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.000	-0.0000
0.040	0.040	0.0003
0.100	0.101	-0.0005
0.200	0.200	0.0004
0.300	0.300	-0.0001

<<<< 3um >>>> <<<< 9um >>>>

Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1
Sample
Sample #1 3287.00 3548.00
Sample #2 3374.00 3566.00
Sample #3 3348.00 3588.00
Sample #4 3284.00 3561.00
Avg 3335.3333 3571.6667
STD DEV 46.3177 14.3643
REL STD DEV 1.389 0.402
H2O adjust (mg/l*10k) 474 238

Barometric Pressure = 1011

*****CALIBRATION SUCCESSFUL*****

Post stability Checks #1

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:56
Control Test	0.047	08:57
Air Blank	0.000	08:58
Control Test	0.047	08:58
Air Blank	0.000	08:59
Control Test	0.047	09:00
Air Blank	0.000	09:00
Control Test Stats		
Average	0.0470	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:01
Control Test	0.077	09:01
Air Blank	0.000	09:02
Control Test	0.076	09:03
Air Blank	0.000	09:03
Control Test	0.077	09:04
Air Blank	0.000	09:04
Control Test Stats		
Average	0.0767	
Std Dev	0.0006	
Rel Std Dev(%)	0.7531	

Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:05
Control Test	0.197	09:06
Air Blank	0.000	09:07
Control Test	0.196	09:07
Air Blank	0.000	09:08
Control Test	0.197	09:09
Air Blank	0.000	09:09
Control Test Stats		
Average	0.1973	
Std Dev	0.0006	
Rel Std Dev(%)	0.2926	

Operator's Signature

wet

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:15
Control Test	0.079	09:16
Air Blank	0.000	09:16
Control Test	0.079	09:17
Air Blank	0.000	09:17
Control Test	0.079	09:17
Air Blank	0.000	09:18
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

Dry

Post stability checks #2

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:00
Control Test	0.049	11:01
Air Blank	0.000	11:01
Control Test	0.049	11:02
Air Blank	0.000	11:03
Control Test	0.049	11:03
Air Blank	0.000	11:04
Control Test Stats		
Average	0.0490	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:07
Control Test	0.079	11:07
Air Blank	0.000	11:08
Control Test	0.079	11:09
Air Blank	0.000	11:09
Control Test	0.079	11:10
Air Blank	0.000	11:10
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:14
Control Test	0.201	11:14
Air Blank	0.000	11:15
Control Test	0.199	11:16
Air Blank	0.000	11:16
Control Test	0.199	11:17
Air Blank	0.000	11:17
Control Test Stats		
Average	0.1997	
Std Dev	0.0012	
Rel Std Dev(%)	0.5783	

Operator's Signature

wet

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
10/18/2022
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:25
Control Test	0.079	11:26
Air Blank	0.000	11:26
Control Test	0.079	11:26
Air Blank	0.000	11:27
Control Test	0.080	11:27
Air Blank	0.000	11:28
Control Test Stats		
Average	0.0793	
Std Dev	0.0006	
Rel Std Dev(%)	0.7277	

Dry

Operator's Signature