




INSTRUMENT PROCESSING SHEET

Agency Miami Police DepartmentS/N 80-006457Florida Department of
Law EnforcementDate In 08/30/2022 DI Completion Date 09/14/2022☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By DERR	Quality Checks	By DERR	Date	Flow Calibration	By	Date																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>200</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP101</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.230</u> (.228 - .278) 103 mm <u>0.472</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28199</u> ☒ Stability Checks		<u>09/01/2022</u>	☐ Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																														
					Maintenance By _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____ _____																														
Calibration Adjustment By DERR _____				Department Inspection By DERR _____																															
Barometric Pressure Gauge <u>1016/1013</u> ID # <u>68639</u>				Barometric Pressure ID# <u>26932</u>																															
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td>MP5099</td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td>MP5096</td><td>21070</td><td>03/01/2023</td></tr><tr><td>0.100</td><td>MP5098</td><td>21080</td><td>03/08/2023</td></tr><tr><td>0.200</td><td>MP5100</td><td>20510</td><td>12/03/2022</td></tr><tr><td>0.300</td><td>MP5101</td><td>21030</td><td>02/02/2023</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG115904</td><td>06/08/2023</td></tr></tbody></table>				Simulator	Serial #	Lot #	Expiration	0.000	MP5099	N/A	N/A	0.040	MP5096	21070	03/01/2023	0.100	MP5098	21080	03/08/2023	0.200	MP5100	20510	12/03/2022	0.300	MP5101	21030	02/02/2023	0.080 DGS	N/A	AG115904	06/08/2023	Gauge <u>1013</u> Instrument <u>1013</u> Mouth Alcohol Solution Lot # <u>2021-D</u> Acetone Stock Solution Lot # <u>2021-C</u>			
Simulator	Serial #	Lot #	Expiration																																
0.000	MP5099	N/A	N/A																																
0.040	MP5096	21070	03/01/2023																																
0.100	MP5098	21080	03/08/2023																																
0.200	MP5100	20510	12/03/2022																																
0.300	MP5101	21030	02/02/2023																																
0.080 DGS	N/A	AG115904	06/08/2023																																
☒ Post Calibration Adjustment Stability Checks				<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5095</td></tr><tr><td>Interferent</td><td>SD3966</td></tr><tr><td>0.050</td><td>MP5092</td></tr><tr><td>0.080</td><td>MP5093</td></tr><tr><td>0.200</td><td>MP5094</td></tr></tbody></table>				Simulator	Serial Number	0.000	MP5095	Interferent	SD3966	0.050	MP5092	0.080	MP5093	0.200	MP5094																
Simulator	Serial Number																																		
0.000	MP5095																																		
Interferent	SD3966																																		
0.050	MP5092																																		
0.080	MP5093																																		
0.200	MP5094																																		
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td>MP5092</td><td>202201C</td><td>01/11/2024</td></tr><tr><td>0.080</td><td>MP5093</td><td>202201D</td><td>01/18/2024</td></tr><tr><td>0.200</td><td>MP5094</td><td>202201E</td><td>01/18/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>00521080A2</td><td>02/05/2023</td></tr></tbody></table>				Simulator	Serial #	Lot #	Expiration	0.050	MP5092	202201C	01/11/2024	0.080	MP5093	202201D	01/18/2024	0.200	MP5094	202201E	01/18/2024	0.080 DGS	N/A	00521080A2	02/05/2023	Attachments ☒ Form 41 ☒ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☒ Form 40 ☒ Calibration Adjustment ☐ Other _____											
Simulator	Serial #	Lot #	Expiration																																
0.050	MP5092	202201C	01/11/2024																																
0.080	MP5093	202201D	01/18/2024																																
0.200	MP5094	202201E	01/18/2024																																
0.080 DGS	N/A	00521080A2	02/05/2023																																
Notes/Suggested Service: _____ Instrument was calibrated to bring values closer to nominal. DERR _____ _____ _____ _____ _____ _____ _____ _____				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use 2022.09. Israel Soto Digitally signed by Israel Soto Date: 2022.09.14 12:39:46 +0400 14:22:46 -0400 Tech Review / Date Admin Review / Date																															

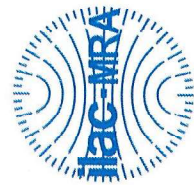

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006457, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-006457</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>MIAMI PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>09/14/2022</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>11:19</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

09/14/2022

Date

DAVID E R YES-RIVERA,

Department Inspector

Service • Integrity • Respect • Quality

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI PD
Time of Inspection: 11:19

Date of Inspection: 09/14/2022

Serial Number: 80-006457
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:00521080A2 Exp: 02/05/2023
0.000	0.048	0.077	0.196	0.080
0.000	0.048	0.077	0.197	0.079
0.000	0.049	0.076	0.197	0.079
0.000	0.048	0.076	0.197	0.080
0.000	0.048	0.076	0.197	0.079
0.000	0.048	0.077	0.198	0.079
0.000	0.049	0.077	0.197	0.079
0.000	0.048	0.077	0.197	0.079
0.000	0.048	0.077	0.198	0.079
0.000	0.048	0.076	0.198	0.079


Standard Deviations	0.0004	0.0005	0.0006	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

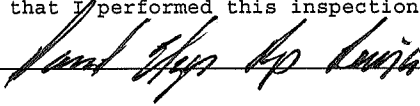

 Signature and Printed Name

DAVID E REYES-RIVERA

09/14/2022
 Date

Type of Test	Serial Number	Agency	Date	Performed By
Post Stabilities	80-006457	Miami Police Department	9/14/2022	DERR <i>me</i>

0.05g/210L 0.047 to 0.053	0.08g/210L 0.077 to 0.083	0.20g/210L 0.194 to 0.206	DGS 0.08g/210L 0.077 to 0.083																																																																																																																																																
MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/14/2022 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>07:55</td></tr><tr><td>Control Test</td><td>0.048</td><td>07:56</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:57</td></tr><tr><td>Control Test</td><td>0.049</td><td>07:57</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:58</td></tr><tr><td>Control Test</td><td>0.049</td><td>07:59</td></tr><tr><td>Air Blank</td><td>0.000</td><td>07:59</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0487</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>1.1863</td><td></td></tr></tbody></table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	07:55	Control Test	0.048	07:56	Air Blank	0.000	07:57	Control Test	0.049	07:57	Air Blank	0.000	07:58	Control Test	0.049	07:59	Air Blank	0.000	07:59	Control Test Stats			Average	0.0487		Std Dev	0.0006		Rel Std Dev(%)	1.1863		MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/14/2022 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>08:00</td></tr><tr><td>Control Test</td><td>0.077</td><td>08:01</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:02</td></tr><tr><td>Control Test</td><td>0.078</td><td>08:02</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:03</td></tr><tr><td>Control Test</td><td>0.077</td><td>08:03</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:04</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0773</td><td></td></tr><tr><td>Std Dev</td><td>0.0016</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7466</td><td></td></tr></tbody></table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	08:00	Control Test	0.077	08:01	Air Blank	0.000	08:02	Control Test	0.078	08:02	Air Blank	0.000	08:03	Control Test	0.077	08:03	Air Blank	0.000	08:04	Control Test Stats			Average	0.0773		Std Dev	0.0016		Rel Std Dev(%)	0.7466		MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/14/2022 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>08:06</td></tr><tr><td>Control Test</td><td>0.197</td><td>08:07</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:07</td></tr><tr><td>Control Test</td><td>0.197</td><td>08:08</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:08</td></tr><tr><td>Control Test</td><td>0.197</td><td>08:09</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:10</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.1970</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></tbody></table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	08:06	Control Test	0.197	08:07	Air Blank	0.000	08:07	Control Test	0.197	08:08	Air Blank	0.000	08:08	Control Test	0.197	08:09	Air Blank	0.000	08:10	Control Test Stats			Average	0.1970		Std Dev	0.0000		Rel Std Dev(%)	0.0000		MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/14/2022 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>08:11</td></tr><tr><td>Control Test</td><td>0.080</td><td>08:12</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:12</td></tr><tr><td>Control Test</td><td>0.080</td><td>08:13</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:13</td></tr><tr><td>Control Test</td><td>0.080</td><td>08:13</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:14</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0800</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></tbody></table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	08:11	Control Test	0.080	08:12	Air Blank	0.000	08:12	Control Test	0.080	08:13	Air Blank	0.000	08:13	Control Test	0.080	08:13	Air Blank	0.000	08:14	Control Test Stats			Average	0.0800		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	07:55																																																																																																																																																	
Control Test	0.048	07:56																																																																																																																																																	
Air Blank	0.000	07:57																																																																																																																																																	
Control Test	0.049	07:57																																																																																																																																																	
Air Blank	0.000	07:58																																																																																																																																																	
Control Test	0.049	07:59																																																																																																																																																	
Air Blank	0.000	07:59																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0487																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	1.1863																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	08:00																																																																																																																																																	
Control Test	0.077	08:01																																																																																																																																																	
Air Blank	0.000	08:02																																																																																																																																																	
Control Test	0.078	08:02																																																																																																																																																	
Air Blank	0.000	08:03																																																																																																																																																	
Control Test	0.077	08:03																																																																																																																																																	
Air Blank	0.000	08:04																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0773																																																																																																																																																		
Std Dev	0.0016																																																																																																																																																		
Rel Std Dev(%)	0.7466																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	08:06																																																																																																																																																	
Control Test	0.197	08:07																																																																																																																																																	
Air Blank	0.000	08:07																																																																																																																																																	
Control Test	0.197	08:08																																																																																																																																																	
Air Blank	0.000	08:08																																																																																																																																																	
Control Test	0.197	08:09																																																																																																																																																	
Air Blank	0.000	08:10																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1970																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	08:11																																																																																																																																																	
Control Test	0.080	08:12																																																																																																																																																	
Air Blank	0.000	08:12																																																																																																																																																	
Control Test	0.080	08:13																																																																																																																																																	
Air Blank	0.000	08:13																																																																																																																																																	
Control Test	0.080	08:13																																																																																																																																																	
Air Blank	0.000	08:14																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0800																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		

MIAMI PD
Intoxilyzer - Alcohol Analyzer
Model 8000
09/14/2022
SN 80-006457
07:03:55

Auto Calibration
Max Power Res Value = 93
Auto Range Res Value = 77

Channel 2
Sample % Abs (% Abs Ref)
Sample #1 = 1.5270 (0.0010)
Sample #2 = 1.5570 (-0.0080)
Sample #3 = 1.5750 (0.0000)
Sample #4 = 1.5520 (0.0230)
Avg % Abs = 1.5613 (0.0050)
STD DEV = 0.0121 (0.0161)
REL STD DEV = 0.775 (321.870)

Sol Value = 0.000 g/210L ***
Fit value = 0.0000 mg/l ****
Samples Taken = 4, Discarded = 1
Sum Io = 12612, Sum Io = 12635

Channel 1
Sample % Abs (% Abs Ref)
Sample #1 = 0.1540 (-0.0120)
Sample #2 = 0.1340 (0.0400)
Sample #3 = 0.1320 (0.0500)
Sample #4 = 0.1520 (0.0700)
Avg % Abs = 0.1393 (0.0533)
STD DEV = 0.0110 (0.0153)
REL STD DEV = 7.906 (28.641)

Channel 2
Sample % Abs (% Abs Ref)
Sample #1 = 0.1890 (-0.0050)
Sample #2 = 0.1630 (0.0260)
Sample #3 = 0.1650 (0.0320)
Sample #4 = 0.1810 (0.0320)
Avg % Abs = 0.1697 (0.0300)
STD DEV = 0.0099 (0.0035)
REL STD DEV = 5.815 (11.547)

Sol Value = 0.040 g/210L ***
Fit value = 0.1905 mg/l ****
Samples Taken = 4, Discarded = 1
Sum Io = 12595, Sum Io = 12627

Channel 1
Sample % Abs (% Abs Ref)
Sample #1 = 0.8410 (0.0000)
Sample #2 = 0.8570 (0.0000)
Sample #3 = 0.8780 (0.0270)
Sample #4 = 0.8420 (0.0550)
Avg % Abs = 0.8590 (0.0273)
STD DEV = 0.0181 (0.0275)
REL STD DEV = 2.105 (100.615)

Channel 1
Sol Val = 0.0000 mg/l or 0.000 g/210L
% Abs = 0.139
Std Dev = 0.01 Rel Std Dev = 7.91
Sol Val = 0.1905 mg/l or 0.040 g/210L
% Abs = 0.859
Std Dev = 0.02 Rel Std Dev = 2.11
Sol Val = 0.4762 mg/l or 0.100 g/210L
% Abs = 1.938
Std Dev = 0.00 Rel Std Dev = 0.24
Sol Val = 0.9524 mg/l or 0.200 g/210L
% Abs = 3.757
Std Dev = 0.02 Rel Std Dev = 0.41
Sol Val = 1.4286 mg/l or 0.300 g/210L
% Abs = 5.457
Std Dev = 0.03 Rel Std Dev = 0.46
Zero Order Coef = -328.13
First Order Coef = 2560.47
Second Order Coef = 20.82
Standard Deviation = 44.210552

Channel 2
Sol Val = 0.0000 mg/l or 0.000 g/210L
% Abs = 0.170
Std Dev = 0.01 Rel Std Dev = 5.81
Sol Val = 0.1905 mg/l or 0.040 g/210L
% Abs = 1.561
Std Dev = 0.01 Rel Std Dev = 0.77
Sol Val = 0.4762 mg/l or 0.100 g/210L
% Abs = 3.539
Std Dev = 0.01 Rel Std Dev = 0.30
Sol Val = 0.9524 mg/l or 0.200 g/210L
% Abs = 6.844
Std Dev = 0.16 Rel Std Dev = 2.32
Sol Val = 1.4286 mg/l or 0.300 g/210L
% Abs = 9.751
Std Dev = 0.02 Rel Std Dev = 0.17
Zero Order Coef = -197.16
First Order Coef = 1319.22
Second Order Coef = 16.68
Standard Deviation = 63.863907

Optical Calibration
SN: 80-006457
Agency: Miami PD
Date: 9/14/2022
Quadratic Fit: +/- 0.002g/210L
By: DERR

Solution Stats Quadratic Fit Chan 1
Act Fit Residual
g/210L g/210L g/210L
0.000 0.001 -0.0006
0.040 0.040 0.0004
0.100 0.099 0.0010
0.200 0.201 -0.0013
0.300 0.300 0.0005

Solution Stats Quadratic Fit Chan 2
Act Fit Residual
g/210L g/210L g/210L
0.000 0.001 -0.0006
0.040 0.040 0.0000
0.100 0.098 0.0017
0.200 0.202 -0.0019
0.300 0.299 0.0007

Sol Value = 0.080 g/210L ***
Fit value = 0.3810 mg/l ****
Samples Taken = 4, Discarded = 1

Channel 1
Sample #1 = 3027.00
Sample #2 = 3086.00
Sample #3 = 3033.00
Sample #4 = 3051.00
Average Result = 3056.6667
STD DEV = 26.9506
REL STD DEV = 0.882

Channel 2
Sample #1 = 3204.00
Sample #2 = 3221.00
Sample #3 = 3216.00
Sample #4 = 3208.00
Average Result = 3215.0000
STD DEV = 6.5574
REL STD DEV = 0.204

Dry Gas H2O Adjust Results *****
Barometric Pressure = 1013
3 um H2O Adjust (mg/l*10,000) = 753
9 um H2O Adjust (mg/l*10,000) = 594
**** AUTO CAL PASS

Type of Test	Serial Number	Agency	Date	Performed By
Post Stabilities	80-006457	Miami Police Department	9/1/2022	DERR <i>[Signature]</i>

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																				
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083																																																																																																																																				
MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/01/2022 Software: 8100.27 Test: g/210LTime <table><tr><td>Air Blank</td><td>0.000</td><td>13:11</td></tr><tr><td>Control Test</td><td>0.048</td><td>13:11</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:12</td></tr><tr><td>Control Test</td><td>0.048</td><td>13:12</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:13</td></tr><tr><td>Control Test</td><td>0.048</td><td>13:14</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:14</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0480</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></table> Operator's Signature	Air Blank	0.000	13:11	Control Test	0.048	13:11	Air Blank	0.000	13:12	Control Test	0.048	13:12	Air Blank	0.000	13:13	Control Test	0.048	13:14	Air Blank	0.000	13:14	Control Test Stats			Average	0.0480		Std Dev	0.0000		Rel Std Dev(%)	0.0000		MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/01/2022 Software: 8100.27 Test: g/210LTime <table><tr><td>Air Blank</td><td>0.000</td><td>13:16</td></tr><tr><td>Control Test</td><td>0.078</td><td>13:17</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:17</td></tr><tr><td>Control Test</td><td>0.077</td><td>13:18</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:19</td></tr><tr><td>Control Test</td><td>0.078</td><td>13:19</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:20</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0777</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7434</td><td></td></tr></table> Operator's Signature	Air Blank	0.000	13:16	Control Test	0.078	13:17	Air Blank	0.000	13:17	Control Test	0.077	13:18	Air Blank	0.000	13:19	Control Test	0.078	13:19	Air Blank	0.000	13:20	Control Test Stats			Average	0.0777		Std Dev	0.0006		Rel Std Dev(%)	0.7434		MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/01/2022 Software: 8100.27 Test: g/210LTime <table><tr><td>Air Blank</td><td>0.000</td><td>13:21</td></tr><tr><td>Control Test</td><td>0.194</td><td>13:22</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:22</td></tr><tr><td>Control Test</td><td>0.194</td><td>13:23</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:23</td></tr><tr><td>Control Test</td><td>0.193</td><td>13:24</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:24</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.1937</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.2981</td><td></td></tr></table> Operator's Signature	Air Blank	0.000	13:21	Control Test	0.194	13:22	Air Blank	0.000	13:22	Control Test	0.194	13:23	Air Blank	0.000	13:23	Control Test	0.193	13:24	Air Blank	0.000	13:24	Control Test Stats			Average	0.1937		Std Dev	0.0006		Rel Std Dev(%)	0.2981		MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/01/2022 Software: 8100.27 Test: g/210LTime <table><tr><td>Air Blank</td><td>0.000</td><td>13:27</td></tr><tr><td>Control Test</td><td>0.080</td><td>13:28</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:28</td></tr><tr><td>Control Test</td><td>0.080</td><td>13:29</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:29</td></tr><tr><td>Control Test</td><td>0.080</td><td>13:29</td></tr><tr><td>Air Blank</td><td>0.000</td><td>13:30</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0800</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></table> Operator's Signature	Air Blank	0.000	13:27	Control Test	0.080	13:28	Air Blank	0.000	13:28	Control Test	0.080	13:29	Air Blank	0.000	13:29	Control Test	0.080	13:29	Air Blank	0.000	13:30	Control Test Stats			Average	0.0800		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Air Blank	0.000	13:11																																																																																																																																					
Control Test	0.048	13:11																																																																																																																																					
Air Blank	0.000	13:12																																																																																																																																					
Control Test	0.048	13:12																																																																																																																																					
Air Blank	0.000	13:13																																																																																																																																					
Control Test	0.048	13:14																																																																																																																																					
Air Blank	0.000	13:14																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0480																																																																																																																																						
Std Dev	0.0000																																																																																																																																						
Rel Std Dev(%)	0.0000																																																																																																																																						
Air Blank	0.000	13:16																																																																																																																																					
Control Test	0.078	13:17																																																																																																																																					
Air Blank	0.000	13:17																																																																																																																																					
Control Test	0.077	13:18																																																																																																																																					
Air Blank	0.000	13:19																																																																																																																																					
Control Test	0.078	13:19																																																																																																																																					
Air Blank	0.000	13:20																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0777																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.7434																																																																																																																																						
Air Blank	0.000	13:21																																																																																																																																					
Control Test	0.194	13:22																																																																																																																																					
Air Blank	0.000	13:22																																																																																																																																					
Control Test	0.194	13:23																																																																																																																																					
Air Blank	0.000	13:23																																																																																																																																					
Control Test	0.193	13:24																																																																																																																																					
Air Blank	0.000	13:24																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.1937																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.2981																																																																																																																																						
Air Blank	0.000	13:27																																																																																																																																					
Control Test	0.080	13:28																																																																																																																																					
Air Blank	0.000	13:28																																																																																																																																					
Control Test	0.080	13:29																																																																																																																																					
Air Blank	0.000	13:29																																																																																																																																					
Control Test	0.080	13:29																																																																																																																																					
Air Blank	0.000	13:30																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0800																																																																																																																																						
Std Dev	0.0000																																																																																																																																						
Rel Std Dev(%)	0.0000																																																																																																																																						

Type of Test	Serial Number	Agency	Date	Performed By
Stabilities	80-006457	Miami Police Department	9/1/2022	DERR <i>[Signature]</i>

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																																
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083																																																																																																																																																
MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/01/2022 Software: 8100.27 SN 80-006457 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>08:43</td></tr><tr><td>Control Test</td><td>INT*</td><td>08:44</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:44</td></tr><tr><td>Control Test</td><td>1.055</td><td>08:45</td></tr><tr><td>Air Blank</td><td>1.000</td><td>08:45</td></tr><tr><td>Control Test</td><td>1.046</td><td>08:46</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:46</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0470</td><td></td></tr><tr><td>Std Dev</td><td>0.0075</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>16.0635</td><td></td></tr></table> *Interferent Detect  Operator's Signature	Test	g/210L	Time	Air Blank	0.000	08:43	Control Test	INT*	08:44	Air Blank	0.000	08:44	Control Test	1.055	08:45	Air Blank	1.000	08:45	Control Test	1.046	08:46	Air Blank	0.000	08:46	Control Test Stats			Average	0.0470		Std Dev	0.0075		Rel Std Dev(%)	16.0635		MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/01/2022 Software: 8100.27 SN 80-006457 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>08:52</td></tr><tr><td>Control Test</td><td>0.080</td><td>08:53</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:54</td></tr><tr><td>Control Test</td><td>0.079</td><td>08:54</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:55</td></tr><tr><td>Control Test</td><td>0.079</td><td>08:56</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:56</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0753</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7277</td><td></td></tr></table>  Operator's Signature	Test	g/210L	Time	Air Blank	0.000	08:52	Control Test	0.080	08:53	Air Blank	0.000	08:54	Control Test	0.079	08:54	Air Blank	0.000	08:55	Control Test	0.079	08:56	Air Blank	0.000	08:56	Control Test Stats			Average	0.0753		Std Dev	0.0006		Rel Std Dev(%)	0.7277		MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/01/2022 Software: 8100.27 SN 80-006457 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>08:58</td></tr><tr><td>Control Test</td><td>0.193</td><td>08:59</td></tr><tr><td>Air Blank</td><td>0.000</td><td>08:59</td></tr><tr><td>Control Test</td><td>0.193</td><td>09:00</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:00</td></tr><tr><td>Control Test</td><td>0.193</td><td>09:01</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:02</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.1930</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></table>  Operator's Signature	Test	g/210L	Time	Air Blank	0.000	08:58	Control Test	0.193	08:59	Air Blank	0.000	08:59	Control Test	0.193	09:00	Air Blank	0.000	09:00	Control Test	0.193	09:01	Air Blank	0.000	09:02	Control Test Stats			Average	0.1930		Std Dev	0.0000		Rel Std Dev(%)	0.0000		MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 09/01/2022 Software: 8100.27 SN 80-006457 <table><tr><th>Test</th><th>g/210L</th><th>Time</th></tr><tr><td>Air Blank</td><td>0.000</td><td>09:03</td></tr><tr><td>Control Test</td><td>0.000</td><td>09:04</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:04</td></tr><tr><td>Control Test</td><td>0.000</td><td>09:05</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:05</td></tr><tr><td>Control Test</td><td>0.000</td><td>09:06</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:06</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0000</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></table>  Operator's Signature	Test	g/210L	Time	Air Blank	0.000	09:03	Control Test	0.000	09:04	Air Blank	0.000	09:04	Control Test	0.000	09:05	Air Blank	0.000	09:05	Control Test	0.000	09:06	Air Blank	0.000	09:06	Control Test Stats			Average	0.0000		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	08:43																																																																																																																																																	
Control Test	INT*	08:44																																																																																																																																																	
Air Blank	0.000	08:44																																																																																																																																																	
Control Test	1.055	08:45																																																																																																																																																	
Air Blank	1.000	08:45																																																																																																																																																	
Control Test	1.046	08:46																																																																																																																																																	
Air Blank	0.000	08:46																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0470																																																																																																																																																		
Std Dev	0.0075																																																																																																																																																		
Rel Std Dev(%)	16.0635																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	08:52																																																																																																																																																	
Control Test	0.080	08:53																																																																																																																																																	
Air Blank	0.000	08:54																																																																																																																																																	
Control Test	0.079	08:54																																																																																																																																																	
Air Blank	0.000	08:55																																																																																																																																																	
Control Test	0.079	08:56																																																																																																																																																	
Air Blank	0.000	08:56																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0753																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7277																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	08:58																																																																																																																																																	
Control Test	0.193	08:59																																																																																																																																																	
Air Blank	0.000	08:59																																																																																																																																																	
Control Test	0.193	09:00																																																																																																																																																	
Air Blank	0.000	09:00																																																																																																																																																	
Control Test	0.193	09:01																																																																																																																																																	
Air Blank	0.000	09:02																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1930																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	09:03																																																																																																																																																	
Control Test	0.000	09:04																																																																																																																																																	
Air Blank	0.000	09:04																																																																																																																																																	
Control Test	0.000	09:05																																																																																																																																																	
Air Blank	0.000	09:05																																																																																																																																																	
Control Test	0.000	09:06																																																																																																																																																	
Air Blank	0.000	09:06																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0000																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI PD

Time of Inspection: 07:12

Date of Inspection: 09/01/2022

Serial Number: 80-006457

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:_____ Exp:_____	0.08g/210L Test (g/210L) Lot#:_____ Exp:_____	0.20g/210L Test (g/210L) Lot#:_____ Exp:_____	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:_____ Exp:_____


Number of Simulators Used: _____

Remarks:

COMPLIANCE NOT DETERMINED, AI NOT CONDUCTED.

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

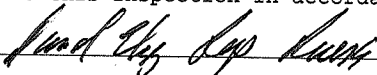 DAVID E REYES-RIVERA
Signature and Printed Name

09/01/2022
Date