



INSTRUMENT PROCESSING SHEET

Agency Tampa FHPS/N 80-006633Florida Department of
Law EnforcementDate In 1/18/2022 DI Completion Date 1/21/2022☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By TDG	Quality Checks	By TDG	Date <u>1/21/2022</u>	Flow Calibration	By	Date																																																																				
☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>Missing feet on DGS shelf</u> 		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>212</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.167</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.488</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>68639</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																																						
					<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:20%;">Simulator</th> <th style="width:20%;">Serial #</th> <th style="width:20%;">Lot #/Exp</th> </tr> <tr> <td>0.050</td> <td>MP5092</td> <td>202010A 10/05/2022</td> </tr> <tr> <td>0.080</td> <td>MP4864</td> <td>202010B 10/05/2022</td> </tr> <tr> <td>0.200</td> <td>MP5094</td> <td>202010D 10/06/2022</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG115904 06/08/2023</td> </tr> </table>			Simulator	Serial #	Lot #/Exp	0.050	MP5092	202010A 10/05/2022	0.080	MP4864	202010B 10/05/2022	0.200	MP5094	202010D 10/06/2022	0.080 DGS	N/A	AG115904 06/08/2023																																																					
Simulator	Serial #	Lot #/Exp																																																																									
0.050	MP5092	202010A 10/05/2022																																																																									
0.080	MP4864	202010B 10/05/2022																																																																									
0.200	MP5094	202010D 10/06/2022																																																																									
0.080 DGS	N/A	AG115904 06/08/2023																																																																									
					<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:50%;">Maintenance</th> <th style="width:50%;">By</th> </tr> <tr> <td> ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ </td> <td></td> </tr> </table>			Maintenance	By	☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ 																																																																	
Maintenance	By																																																																										
☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ 																																																																											
<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:40%;">Calibration Adjustment</th> <th style="width:60%;">By</th> </tr> <tr> <td> Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:20%;">Simulator</th> <th style="width:20%;">Serial #</th> <th style="width:20%;">Lot #</th> <th style="width:20%;">Expiration</th> </tr> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:20%;">Simulator</th> <th style="width:20%;">Serial #</th> <th style="width:20%;">Lot #</th> <th style="width:20%;">Expiration</th> </tr> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table> </td> <td></td> </tr> </table>				Calibration Adjustment	By	Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:20%;">Simulator</th> <th style="width:20%;">Serial #</th> <th style="width:20%;">Lot #</th> <th style="width:20%;">Expiration</th> </tr> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:20%;">Simulator</th> <th style="width:20%;">Serial #</th> <th style="width:20%;">Lot #</th> <th style="width:20%;">Expiration</th> </tr> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table>	Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A				<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:50%;">Department Inspection</th> <th style="width:50%;">By TDG</th> </tr> <tr> <td> Barometric Pressure ID# <u>28663</u> Gauge <u>1019</u> Instrument <u>1018</u> Mouth Alcohol Solution Lot # <u>2021-B</u> Acetone Stock Solution Lot # <u>2020-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:50%;">Simulator</th> <th style="width:50%;">Serial Number</th> </tr> <tr> <td>0.000</td> <td>SD3963</td> </tr> <tr> <td>Interferent</td> <td>SD1017</td> </tr> <tr> <td>0.050</td> <td>MP5092</td> </tr> <tr> <td>0.080</td> <td>MP4864</td> </tr> <tr> <td>0.200</td> <td>MP5094</td> </tr> </table> </td> <td></td> </tr> </table>				Department Inspection	By TDG	Barometric Pressure ID# <u>28663</u> Gauge <u>1019</u> Instrument <u>1018</u> Mouth Alcohol Solution Lot # <u>2021-B</u> Acetone Stock Solution Lot # <u>2020-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:50%;">Simulator</th> <th style="width:50%;">Serial Number</th> </tr> <tr> <td>0.000</td> <td>SD3963</td> </tr> <tr> <td>Interferent</td> <td>SD1017</td> </tr> <tr> <td>0.050</td> <td>MP5092</td> </tr> <tr> <td>0.080</td> <td>MP4864</td> </tr> <tr> <td>0.200</td> <td>MP5094</td> </tr> </table>	Simulator	Serial Number	0.000	SD3963	Interferent	SD1017	0.050	MP5092	0.080	MP4864	0.200	MP5094	
Calibration Adjustment	By																																																																										
Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:20%;">Simulator</th> <th style="width:20%;">Serial #</th> <th style="width:20%;">Lot #</th> <th style="width:20%;">Expiration</th> </tr> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:20%;">Simulator</th> <th style="width:20%;">Serial #</th> <th style="width:20%;">Lot #</th> <th style="width:20%;">Expiration</th> </tr> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table>	Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A																													
Simulator	Serial #	Lot #	Expiration																																																																								
0.000		N/A	N/A																																																																								
0.040																																																																											
0.100																																																																											
0.200																																																																											
0.300																																																																											
0.080 DGS	N/A																																																																										
Simulator	Serial #	Lot #	Expiration																																																																								
0.050																																																																											
0.080																																																																											
0.200																																																																											
0.080 DGS	N/A																																																																										
Department Inspection	By TDG																																																																										
Barometric Pressure ID# <u>28663</u> Gauge <u>1019</u> Instrument <u>1018</u> Mouth Alcohol Solution Lot # <u>2021-B</u> Acetone Stock Solution Lot # <u>2020-A</u> <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:50%;">Simulator</th> <th style="width:50%;">Serial Number</th> </tr> <tr> <td>0.000</td> <td>SD3963</td> </tr> <tr> <td>Interferent</td> <td>SD1017</td> </tr> <tr> <td>0.050</td> <td>MP5092</td> </tr> <tr> <td>0.080</td> <td>MP4864</td> </tr> <tr> <td>0.200</td> <td>MP5094</td> </tr> </table>	Simulator	Serial Number	0.000	SD3963	Interferent	SD1017	0.050	MP5092	0.080	MP4864	0.200	MP5094																																																															
Simulator	Serial Number																																																																										
0.000	SD3963																																																																										
Interferent	SD1017																																																																										
0.050	MP5092																																																																										
0.080	MP4864																																																																										
0.200	MP5094																																																																										
<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th style="width:50%;">Attachments</th> <th style="width:50%;">By</th> </tr> <tr> <td> ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment </td> <td> ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other <u>2021 IPS and Form 51</u> </td> </tr> </table>				Attachments	By	☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other <u>2021 IPS and Form 51</u>	☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																																																			
Attachments	By																																																																										
☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other <u>2021 IPS and Form 51</u>																																																																										
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____				Israel Soto Digitally signed by Israel Soto Date: 2022.01.21 16:19:46 +05'00' 25 14:45:09 05/00/																																																																							
Tech Review / Date _____				Admin Review / Date _____																																																																							

Type of Test	Serial Number	Agency	Date	Performed By
Stabilities	80-006633	FHP	01/21/2022	TDG

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
✓	✓	✓	✓
FHP Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006633 01/21/2022 Software: 8100.27 Test g/210L Time Air Blank 0.000 09:57 Control Test 0.048 09:57 Air Blank 0.000 09:58 Control Test 0.048 09:59 Air Blank 0.000 09:59 Control Test 0.048 10:00 Air Blank 0.000 10:00 Control Test Stats Average 0.0480 Std Dev 0.0000 Rel Std Dev(%) 0.0000	FHP Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006633 01/21/2022 Software: 8100.27 Test g/210L Time Air Blank 0.000 10:03 Control Test 0.079 10:04 Air Blank 0.000 10:04 Control Test 0.079 10:05 Air Blank 0.000 10:06 Control Test 0.078 10:06 Air Blank 0.000 10:07 Control Test Stats Average 0.0787 Std Dev 0.0006 Rel Std Dev(%) 0.7339	FHP Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006633 01/21/2022 Software: 8100.27 Test g/210L Time Air Blank 0.000 10:14 Control Test 0.198 10:14 Air Blank 0.000 10:15 Control Test 0.198 10:15 Air Blank 0.000 10:16 Control Test 0.198 10:17 Air Blank 0.000 10:17 Control Test Stats Average 0.1980 Std Dev 0.0000 Rel Std Dev(%) 0.0000	065 FHP Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006633 01/21/2022 Software: 8100.27 Test g/210L Time Air Blank 0.000 09:50 Control Test 0.080 09:50 Air Blank 0.000 09:51 Control Test 0.079 09:51 Air Blank 0.000 09:52 Control Test 0.079 09:52 Air Blank 0.000 09:52 Control Test Stats Average 0.0793 Std Dev 0.0006 Rel Std Dev(%) 0.7277
Operator's Signature _____	Operator's Signature _____	Operator's Signature _____	Operator's Signature _____

Comments:

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PHP

Time of Inspection: 12:37

Date of Inspection: 01/21/2022

Serial Number: 80-006633

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202010A Exp: 10/05/2022	0.08g/210L Test (g/210L) Lot#:202010B Exp: 10/05/2022	0.20g/210L Test (g/210L) Lot#:202010D Exp: 10/06/2022	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG115904 Exp: 06/08/2023
0.000	0.048	0.078	0.198	0.079
0.000	0.048	0.078	0.198	0.079
0.000	0.048	0.078	0.197	0.079
0.000	0.048	0.078	0.198	0.079
0.000	0.048	0.078	0.198	0.079
0.000	0.048	0.078	0.197	0.079
0.000	0.048	0.079	0.198	0.079
0.000	0.048	0.078	0.197	0.080
0.000	0.048	0.078	0.198	0.079
0.000	0.048	0.079	0.198	0.080


Standard Deviations	0.0000	0.0004	0.0004	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

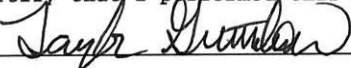

TAYLOR D GUTSCHOW

Signature and Printed Name

01/21/2022
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006633, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-006633</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>FHP</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>01/21/2022</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>12:37</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.



01/21/2022

Date

TAYLOR D GUTSCHOW,

Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality



INSTRUMENT PROCESSING SHEET

Agency Tampa FHPS/N 80-006633Florida Department of
Law EnforcementDate In 1/29/2021 DI Completion Date _____ ☒ Ship ☐ P/U ☐ H/D ☒ CMI ☐ EE

Intake	By <u>TDG</u>	Quality Checks	By _____	Date _____	Flow Calibration	By _____	Date _____																																																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>No feet on DGS shelf</u> _____ _____ _____ _____ _____ _____ _____ _____		☐ Breath Tube Screen ☐ Replace External O-Rings ☐ Instrument Set Up Verified ☐ R-Value _____ ☐ Flow Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) ☐ Barometric Pressure Check Gauge ID # _____ ☐ Stability Checks			Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																														
		<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> </tr> </tbody> </table>			Simulator	Serial #	Lot #/Exp	0.050			0.080			0.200			0.080 DGS	N/A		Maintenance By <u>TDG</u> ☒ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ DI Temp. Checks By _____ ☐ Lab Temp °C _____ External Digital Therm. ID#: _____ ☐ 34°C +/- .2 Serial #: _____ ☐ 34°C +/- .2 Serial #: _____ ☐ 34°C +/- .2 Serial #: _____																																															
Simulator	Serial #	Lot #/Exp																																																																	
0.050																																																																			
0.080																																																																			
0.200																																																																			
0.080 DGS	N/A																																																																		
Calibration Adjustment By _____		Department Inspection By _____																																																																	
Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Barometric Pressure ID# _____ Gauge _____ Instrument _____ Mouth Alcohol Solution Lot # _____ Acetone Stock Solution Lot # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> </tr> <tr> <td>Interferent</td> <td></td> </tr> <tr> <td>0.050</td> <td></td> </tr> <tr> <td>0.080</td> <td></td> </tr> <tr> <td>0.200</td> <td></td> </tr> </tbody> </table>						Simulator	Serial Number	0.000		Interferent		0.050		0.080		0.200	
Simulator	Serial #	Lot #	Expiration																																																																
0.000		N/A	N/A																																																																
0.040																																																																			
0.100																																																																			
0.200																																																																			
0.300																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial #	Lot #	Expiration																																																																
0.050																																																																			
0.080																																																																			
0.200																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial Number																																																																		
0.000																																																																			
Interferent																																																																			
0.050																																																																			
0.080																																																																			
0.200																																																																			
Attachments ☐ Form 41 ☐ Stability Checks ☐ Calibration Certificate ☐ Calibration Adjustment		☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____																																																																	
Notes/Suggested Service: <u>Battery replaced on 02/10/2021 but could not set Date/Time and check status of new battery. See Form 51. Will be sent to CMI for repair. (TDG 02/10/2021)</u> _____ _____ _____ _____ _____		☐ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☐ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☐ Conduct an Agency Inspection Before Evidentiary Use																																																																	
Israel Soto Digitally signed by Israel Soto Date: 2022.01.21 16:19:17 -0500		Tech Review / Date _____ Admin Review / Date _____																																																																	

Return Material Authorization

Ship to: ☒ CMI, Inc.

☐ Enforcement Electronics

Shipment to repair facility authorized by: Alain Hernandez on 02/10/2021

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-006633

Bill To Address:

Florida Highway Patrol

11305 N McKinley Drive

Tampa, FL 33612

Ship to Address:

Florida Department of Law Enforcement

Fort Myers Regional Operations Center

4700 Terminal Drive, Suite 1

Fort Myers, FL. 33907

Reason for Return:

Battery was replaced on 2/10/2021. After battery replacement, instrument began to give these 2 error messages upon startup: "EEPROM fault. Life data invalid. Baro cal invalid" and "Error 17: Contact CMI". Double-checked internal components were attached/secured. Found a tiny tan square piece freely moving inside instrument but could not find its source. Could not set the date/time to check status of new battery.

Please choose one of the following options:

- ☐ 1. I _____, authorize all repairs.
- ☐ 2. I _____, authorize repairs up to \$ _____.
- ☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Alain Hernandez

Phone #: 813-735-2304 Email: AlainHernandez-Ginarte@flhsmv.gov

ATP Contact Name: Taylor Gutschow ATP Email: TaylorGutschow@fdle.state.fl.us