



INSTRUMENT PROCESSING SHEET

Agency Pinellas County Sheriff's Office

S/N 80-000888

Florida Department of Law Enforcement

Date In 07/26/23

DI Completion Date 07-31-2023

Ship P/U H/D CMI EE

Intake, Quality Checks, Flow Calibration, Maintenance, Notes, and a table of simulator data.

Calibration Adjustment, Department Inspection, and Attachments sections.

Notes/Suggested Service, Tech Review, and signature blocks for Taylor Gutschow and Phil Nicodemo.

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PINELLAS COUNTY SO
Time of Inspection: 11:29

Date of Inspection: 07/31/2023

Serial Number: 80-000888
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:06723080A5 Exp: 04/05/2025
0.000	0.049	0.080	0.201	0.078
0.000	0.049	0.080	0.201	0.078
0.000	0.049	0.080	0.201	0.078
0.000	0.049	0.080	0.201	0.078
0.000	0.049	0.081	0.201	0.079
0.000	0.049	0.080	0.200	0.078
0.000	0.049	0.081	0.200	0.078
0.000	0.049	0.080	0.199	0.079
0.000	0.050	0.081	0.198	0.078
0.000	0.049	0.081	0.198	0.078

Standard Deviations	0.0003	0.0005	0.0012	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0006 Number of Simulators Used: 5

Remarks:

Int Det:REMADE SIM WITH NEW ACETONE BOTTLE.

First acetone test did not have Int Det for each measurement. Acetone solution was poured out, and a new acetone solution was made using a new acetone stock bottle of the same lot. Simulator allowed to heat 30 minutes prior to re-test, re-test was successful. ✓✓

The above instrument complies () does not comply () with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL SOTO

Signature and Printed Name

07/31/2023
Date

Stability checks

WELLS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000888
07/31/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:52
Control Test	0.048	07:52
Air Blank	0.000	07:53
Control Test	0.047	07:54
Air Blank	0.000	07:54
Control Test	0.048	07:55
Air Blank	0.000	07:55
Control Test Stats		
Average	0.0477	
Std Dev	0.0006	
Rel Std Dev(%)	1.2112	


Operator's Signature

WELLS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000888
07/31/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:57
Control Test	0.079	07:57
Air Blank	0.000	07:58
Control Test	0.080	07:59
Air Blank	0.000	07:59
Control Test	0.079	08:00
Air Blank	0.000	08:01
Control Test Stats		
Average	0.0793	
Std Dev	0.0006	
Rel Std Dev(%)	0.7277	

wet


Operator's Signature

WELLS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000888
07/31/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:09
Control Test	0.201	08:10
Air Blank	0.000	08:10
Control Test	0.198	08:11
Air Blank	0.000	08:12
Control Test	0.198	08:12
Air Blank	0.000	08:13
Control Test Stats		
Average	0.1990	
Std Dev	0.0017	
Rel Std Dev(%)	0.8704	


Operator's Signature

WELLS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000888
07/31/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	08:14
Control Test	0.079	08:14
Air Blank	0.000	08:15
Control Test	0.078	08:15
Air Blank	0.000	08:16
Control Test	0.078	08:16
Air Blank	0.000	08:17
Control Test Stats		
Average	0.0783	
Std Dev	0.0006	
Rel Std Dev(%)	0.7370	

dry


Operator's Signature



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000888, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000888</u>	UNCERTAINTY* ±	
Owning Agency:	<u>PINELLAS COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>07/31/2023</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>11:29</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Israel Soto
Soto
Digitally signed by Israel Soto
Date: 2023.07.31 11:53:34 -04'00'

07/31/2023

Date

ISRAEL SOTO,
Department Inspector

FDLE/ATP Form 69 March 2022

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality



INSTRUMENT PROCESSING SHEET

Agency Pinellas County Sheriff's OfficeS/N 80-000888

Florida Department of Law Enforcement

Date In 1-30-2023

DI Completion Date _____

 Ship P/U H/D CMI EE

Intake By PN _____ <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: _____ <u>Missing screw on front panel of instrument.</u> _____ _____ _____ _____ _____	Quality Checks By IS _____ Date <u>02-06-2023</u> <input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>204</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP-105</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.519</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>30793</u> <input checked="" type="checkbox"/> Stability Checks	Flow Calibration By _____ Date _____ Flow Column # _____ <input type="checkbox"/> 5L/min – 17mm <input type="checkbox"/> 15L/min – 53mm <input type="checkbox"/> 30L/min – 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																																								
<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td rowspan="2">MP6291</td> <td>202201C</td> </tr> <tr> <td></td> <td>01-11-2024</td> </tr> <tr> <td>0.080</td> <td rowspan="2">MP6292</td> <td>202201D</td> </tr> <tr> <td></td> <td>01-18-2024</td> </tr> <tr> <td>0.200</td> <td rowspan="2">MP6293</td> <td>202201E</td> </tr> <tr> <td></td> <td>01-18-2024</td> </tr> <tr> <td>0.080 DGS</td> <td rowspan="2">N/A</td> <td>AG113403</td> </tr> <tr> <td></td> <td>05-14-2023</td> </tr> </tbody> </table>			Simulator	Serial #	Lot #/Exp	0.050	MP6291	202201C		01-11-2024	0.080	MP6292	202201D		01-18-2024	0.200	MP6293	202201E		01-18-2024	0.080 DGS	N/A	AG113403		05-14-2023																																																	
Simulator	Serial #	Lot #/Exp																																																																								
0.050	MP6291	202201C																																																																								
		01-11-2024																																																																								
0.080	MP6292	202201D																																																																								
		01-18-2024																																																																								
0.200	MP6293	202201E																																																																								
		01-18-2024																																																																								
0.080 DGS	N/A	AG113403																																																																								
		05-14-2023																																																																								
<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th colspan="2">Calibration Adjustment</th> <th colspan="2">By IS _____</th> </tr> </thead> <tbody> <tr> <td colspan="4">Barometric Pressure Gauge <u>1019</u> ID # <u>28427</u></td> </tr> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> <tr> <td>0.000</td> <td>MP6294</td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td>MP6295</td> <td>21410</td> <td>09-30-2023</td> </tr> <tr> <td>0.100</td> <td>MP6296</td> <td>22310</td> <td>08-11-2024</td> </tr> <tr> <td>0.200</td> <td>MP6297</td> <td>22050</td> <td>02-07-2024</td> </tr> <tr> <td>0.300</td> <td>MP6298</td> <td>22220</td> <td>06-15-2024</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>08121080A1</td> <td>05-05-2023</td> </tr> <tr> <td colspan="4"><input checked="" type="checkbox"/> Post Calibration Adjustment Stability Checks</td> </tr> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> <tr> <td>0.050</td> <td>MP6291</td> <td>202201C</td> <td>01-11-2024</td> </tr> <tr> <td>0.080</td> <td>MP6292</td> <td>202201D</td> <td>01-18-2024</td> </tr> <tr> <td>0.200</td> <td>MP6293</td> <td>202201E</td> <td>01-18-2024</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG113403</td> <td>05-14-2023</td> </tr> </tbody> </table>			Calibration Adjustment		By IS _____		Barometric Pressure Gauge <u>1019</u> ID # <u>28427</u>				Simulator	Serial #	Lot #	Expiration	0.000	MP6294	N/A	N/A	0.040	MP6295	21410	09-30-2023	0.100	MP6296	22310	08-11-2024	0.200	MP6297	22050	02-07-2024	0.300	MP6298	22220	06-15-2024	0.080 DGS	N/A	08121080A1	05-05-2023	<input checked="" type="checkbox"/> Post Calibration Adjustment Stability Checks				Simulator	Serial #	Lot #	Expiration	0.050	MP6291	202201C	01-11-2024	0.080	MP6292	202201D	01-18-2024	0.200	MP6293	202201E	01-18-2024	0.080 DGS	N/A	AG113403	05-14-2023												
Calibration Adjustment		By IS _____																																																																								
Barometric Pressure Gauge <u>1019</u> ID # <u>28427</u>																																																																										
Simulator	Serial #	Lot #	Expiration																																																																							
0.000	MP6294	N/A	N/A																																																																							
0.040	MP6295	21410	09-30-2023																																																																							
0.100	MP6296	22310	08-11-2024																																																																							
0.200	MP6297	22050	02-07-2024																																																																							
0.300	MP6298	22220	06-15-2024																																																																							
0.080 DGS	N/A	08121080A1	05-05-2023																																																																							
<input checked="" type="checkbox"/> Post Calibration Adjustment Stability Checks																																																																										
Simulator	Serial #	Lot #	Expiration																																																																							
0.050	MP6291	202201C	01-11-2024																																																																							
0.080	MP6292	202201D	01-18-2024																																																																							
0.200	MP6293	202201E	01-18-2024																																																																							
0.080 DGS	N/A	AG113403	05-14-2023																																																																							
<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th colspan="2">Department Inspection</th> <th colspan="2">By _____</th> </tr> </thead> <tbody> <tr> <td colspan="4">Barometric Pressure ID# _____</td> </tr> <tr> <td colspan="2">Gauge _____</td> <td colspan="2">Instrument _____</td> </tr> <tr> <td colspan="4">Mouth Alcohol Solution Lot # _____</td> </tr> <tr> <td colspan="4">Acetone Stock Solution Lot # _____</td> </tr> <tr> <th>Simulator</th> <th colspan="3">Serial Number</th> </tr> <tr> <td>0.000</td> <td colspan="3"></td> </tr> <tr> <td>Interferent</td> <td colspan="3"></td> </tr> <tr> <td>0.050</td> <td colspan="3"></td> </tr> <tr> <td>0.080</td> <td colspan="3"></td> </tr> <tr> <td>0.200</td> <td colspan="3"></td> </tr> <tr> <th colspan="4">Attachments</th> </tr> <tr> <td><input type="checkbox"/> Form 41</td> <td colspan="3"><input checked="" type="checkbox"/> Post-Stability Checks</td> </tr> <tr> <td><input checked="" type="checkbox"/> Stability Checks</td> <td colspan="3"><input type="checkbox"/> Flow Calibration</td> </tr> <tr> <td><input type="checkbox"/> Calibration Certificate</td> <td colspan="3"><input type="checkbox"/> Form 40</td> </tr> <tr> <td><input checked="" type="checkbox"/> Calibration Adjustment</td> <td colspan="3"><input checked="" type="checkbox"/> Other <u>Form 51</u></td> </tr> <tr> <td colspan="4"> <input type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use </td> </tr> <tr> <td colspan="2">Tech Review / Date _____</td> <td colspan="2">Admin Review / Date _____</td> </tr> </tbody> </table>			Department Inspection		By _____		Barometric Pressure ID# _____				Gauge _____		Instrument _____		Mouth Alcohol Solution Lot # _____				Acetone Stock Solution Lot # _____				Simulator	Serial Number			0.000				Interferent				0.050				0.080				0.200				Attachments				<input type="checkbox"/> Form 41	<input checked="" type="checkbox"/> Post-Stability Checks			<input checked="" type="checkbox"/> Stability Checks	<input type="checkbox"/> Flow Calibration			<input type="checkbox"/> Calibration Certificate	<input type="checkbox"/> Form 40			<input checked="" type="checkbox"/> Calibration Adjustment	<input checked="" type="checkbox"/> Other <u>Form 51</u>			<input type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use				Tech Review / Date _____		Admin Review / Date _____	
Department Inspection		By _____																																																																								
Barometric Pressure ID# _____																																																																										
Gauge _____		Instrument _____																																																																								
Mouth Alcohol Solution Lot # _____																																																																										
Acetone Stock Solution Lot # _____																																																																										
Simulator	Serial Number																																																																									
0.000																																																																										
Interferent																																																																										
0.050																																																																										
0.080																																																																										
0.200																																																																										
Attachments																																																																										
<input type="checkbox"/> Form 41	<input checked="" type="checkbox"/> Post-Stability Checks																																																																									
<input checked="" type="checkbox"/> Stability Checks	<input type="checkbox"/> Flow Calibration																																																																									
<input type="checkbox"/> Calibration Certificate	<input type="checkbox"/> Form 40																																																																									
<input checked="" type="checkbox"/> Calibration Adjustment	<input checked="" type="checkbox"/> Other <u>Form 51</u>																																																																									
<input type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use																																																																										
Tech Review / Date _____		Admin Review / Date _____																																																																								
Notes/Suggested Service: <u>Instrument had Ambient and Purge Fail during 0.050 g/210L Stability Check. Ensured area round Intoxilyzer was clear and repeated 0.050 g/210L Stability Check with nominal results. 0.200 g/210L Stability Check results not nominal, instrument had Interferent Detect on 0.080 g/210 Dry Gas Stability Check. Performed Optical Bench Calibration Adjustment. Post Stability Checks still had Interferent Detect.</u> <u>Sending instrument to CMI. IS</u>																																																																										

Stability Checks

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023
 Software: 8100.27

Test	g/210L	Time
Air Blank	AMB*	09:32
Air Blank	PUR**	09:33
*Ambient Fail		
**Purge Fail		

First
 0.080 g/210L
 045

[Signature]
 Operator's Signature

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:34
Control Test	0.049	09:35
Air Blank	0.000	09:35
Control Test	0.048	09:36
Air Blank	0.000	09:37
Control Test	0.048	09:37
Air Blank	0.000	09:38
Control Test Stats		
Average	0.0483	
Std Dev	0.0006	
Rel Std Dev(%)	1.1945	

Second
 0.050 g/210L

[Signature]
 Operator's Signature

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:39
Control Test	0.078	09:39
Air Blank	0.000	09:40
Control Test	0.077	09:41
Air Blank	0.000	09:41
Control Test	0.077	09:42
Air Blank	0.000	09:42
Control Test Stats		
Average	0.0773	
Std Dev	0.0006	
Rel Std Dev(%)	0.7466	

wet

[Signature]
 Operator's Signature

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:43
Control Test	0.197	09:44
Air Blank	0.000	09:45
Control Test	0.194	09:45
Air Blank	0.000	09:46
Control Test	0.193	09:47
Air Blank	0.000	09:47
Control Test Stats		
Average	0.1947	
Std Dev	0.0021	
Rel Std Dev(%)	1.0693	

[Signature]
 Operator's Signature

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:49
Control Test	INT*	09:49
Air Blank	0.000	09:50
Control Test	INT*	09:50
Air Blank	0.000	09:50
Control Test	INT*	09:51
Air Blank	0.000	09:51
Control Test Stats		
Average	0.0800	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

*Interferent Detect

Dry
 0.080 g/210L

[Signature]
 Operator's Signature

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023 10:47:49

Optical Bench
 Calibration Adjustment

Auto Calibration

pg 2 of 2

OK

<<<<< 3um >>>>>

 Zero Order Coef -59.97
 First Order Coef 2585.68
 Second Order Coef 43.56

<<<<< 9um >>>>>

 -165.53
 1330.83
 13.99

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.000	-0.0004
0.040	0.039	0.0006
0.100	0.100	-0.0001
0.200	0.200	-0.0002
0.300	0.300	0.0001

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.000	-0.0004
0.040	0.040	0.0003
0.100	0.099	0.0007
0.200	0.201	-0.0009
0.300	0.300	0.0004

<<<<< 3um >>>>> <<<<< 9um >>>>>

 Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1
 Sample
 Sample #1 2607.00 3387.00
 Sample #2 2949.00 3362.00
 Sample #3 2953.00 3378.00
 Sample #4 2750.00 3464.00
 Avg 2884.0000 3401.3333
 STD DEV 116.0646 54.8574
 REL STD DEV 4.024 1.613
 H2O adjust (mg/l*10k) 925 408

Barometric Pressure = 1018


*****CALIBRATION SUCCESSFUL*****

Post stability checks

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:19
Control Test	INT*	14:20
Air Blank	0.000	14:20
Control Test	0.048	14:21
Air Blank	0.000	14:21
Control Test	0.049	14:22
Air Blank	0.000	14:23
Control Test Stats		
Average	0.0487	
Std Dev	0.0006	
Rel Std Dev(%)	1.1863	

*Interferent Detect

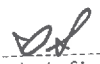

 Operator's Signature

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:24
Control Test	INT*	14:24
Air Blank	0.000	14:25
Control Test	0.078	14:26
Air Blank	0.000	14:26
Control Test	0.079	14:27
Air Blank	0.000	14:28
Control Test Stats		
Average	0.0787	
Std Dev	0.0006	
Rel Std Dev(%)	0.7339	

*Interferent Detect

wet



 Operator's Signature

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:29
Control Test	0.202	14:30
Air Blank	0.000	14:30
Control Test	0.199	14:31
Air Blank	0.000	14:32
Control Test	0.198	14:32
Air Blank	0.000	14:33
Control Test Stats		
Average	0.1997	
Std Dev	0.0021	
Rel Std Dev(%)	1.0426	



 Operator's Signature

PINELLAS COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-000888
 02/06/2023
 Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:35
Control Test	0.090	14:35
Air Blank	0.000	14:36
Control Test	0.081	14:36
Air Blank	0.000	14:37
Control Test	0.081	14:37
Air Blank	0.000	14:37
Control Test Stats		
Average	0.0807	
Std Dev	0.0006	
Rel Std Dev(%)	0.7157	

Dry



 Operator's Signature

Return Material Authorization

Ship to: CMI, Inc.
 Enforcement Electronics

Shipment to repair facility authorized by: Malcolm Deane on 02-15-2023

Items Returned: Instrument Supplies Other Describe: _____
Instrument Model: Intoxilyzer 8000 Serial Number: 80-000888

Bill To Address: <u>Pinellas County Sheriff's Office</u> <u>Florida</u> _____ _____ _____	Ship to Address: <u>Florida Department of Law Enforcement</u> <u>Alcohol Testing Program</u> <u>813 B Lake Bradford Road</u> <u>Tallahassee, FL 32304</u> _____ _____
--	---

Reason for Return:
Instrument is producing Ambient Fail and Purge Fail intermittently, and produces
Interferent Detect on alcohol reference solutions.

Please choose one of the following options:

1. I _____, authorize all repairs.

2. I _____, authorize repairs up to \$_____.

3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Malcolm Dean
Phone #: 727-464-6330 Email: mdeane@pcsonet.com
ATP Contact Name: Israel Soto ATP Email: israelsoto@fdle.state.fl.us