



Agency Florida Fish and Wildlife Conservation Commission S/N 80-000893

Date In 08/03/23 DI Completion Date 08-08-2023 ☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Stability Checks

FFWCC
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000893
08/07/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:03
Control Test	0.048	11:03
Air Blank	0.000	11:04
Control Test	0.049	11:05
Air Blank	0.000	11:05
Control Test	0.049	11:06
Air Blank	0.000	11:06
Control Test Stats		
Average	0.0487	
Std Dev	0.0006	
Rel Std Dev(%)	1.1863	

Benjamin Sidesey
Operator's Signature

FFWCC
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000893
08/07/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:08
Control Test	0.078	11:09
Air Blank	0.000	11:09
Control Test	0.077	11:10
Air Blank	0.000	11:11
Control Test	0.078	11:11
Air Blank	0.000	11:12
Control Test Stats		
Average	0.0777	
Std Dev	0.0006	
Rel Std Dev(%)	0.7434	

Benjamin Sidesey
Operator's Signature

FFWCC
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000893
08/07/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:14
Control Test	0.199	11:15
Air Blank	0.000	11:15
Control Test	0.196	11:16
Air Blank	0.000	11:17
Control Test	0.196	11:17
Air Blank	0.000	11:18
Control Test Stats		
Average	0.1970	
Std Dev	0.0017	
Rel Std Dev(%)	0.8792	

Benjamin Sidesey
Operator's Signature

FFWCC
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000893
08/07/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:20
Control Test	0.078	11:20
Air Blank	0.000	11:21
Control Test	0.078	11:21
Air Blank	0.000	11:21
Control Test	0.077	11:22
Air Blank	0.000	11:22
Control Test Stats		
Average	0.0777	
Std Dev	0.0006	
Rel Std Dev(%)	0.7434	

DGS

Benjamin Sidesey
Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: FFWCC

Time of Inspection: 09:51

Date of Inspection: 08/08/2023

Serial Number: 80-000893

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:06723080A5 Exp: 04/05/2025
0.000	0.049	0.079	0.200	0.078
0.000	0.049	0.079	0.200	0.078
0.000	0.049	0.079	0.200	0.077
0.000	0.049	0.079	0.199	0.077
0.000	0.049	0.079	0.198	0.077
0.000	0.049	0.079	0.198	0.077
0.000	0.049	0.079	0.197	0.077
0.000	0.049	0.079	0.197	0.077
0.000	0.049	0.079	0.197	0.077
0.000	0.049	0.079	0.197	0.077
0.000	0.050	0.079	0.197	0.077

Standard Deviations	0.0003	0.0000	0.0013	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0005 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL SOTO

Signature and Printed Name

08/08/2023
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000893, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000893</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>FFWCC</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>08/08/2023</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>09:51</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Israel Soto
Soto
Date: 2023.08.08 12:23:45
-04'00'

Date

ISRAEL SOTO,
Department Inspector

FDLE/ATP Form 69 March 2022

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality



Agency FFWCC

S/N 80-000893

Date In 01-05-2024 DI Completion Date 1/11/2024

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

FDLE/ATP Form 48 January 2022
Issuing Authority: Alcohol Testing Program

Page 1 of 1

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: FFWCC
Time of Inspection: 16:28

Date of Inspection: 01/11/2024

Serial Number: 80-000893
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG310901 Exp: 04/19/2025
0.000	0.050	0.079	0.197	0.077
0.000	0.049	0.078	0.197	0.076
0.000	0.049	0.079	0.197	0.077
0.000	0.050	0.079	0.197	0.077
0.000	0.050	0.078	0.197	0.077
0.000	0.050	0.079	0.197	0.077
0.000	0.050	0.079	0.197	0.077
0.000	0.049	0.079	0.197	0.076
0.000	0.049	0.078	0.197	0.076
0.000	0.049	0.079	0.196	0.076

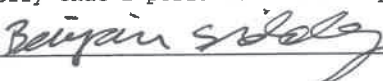
Standard Deviations	0.0005	0.0004	0.0003	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

 BENJAMIN W SIDDOWAY
Signature and Printed Name

01/11/2024
Date

Stability checks

FFWCC
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000893
01/08/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:31
Control Test	0.050	10:31
Air Blank	0.000	10:32
Control Test	0.050	10:33
Air Blank	0.000	10:33
Control Test	0.049	10:34
Air Blank	0.000	10:35
Control Test Stats		
Average	0.0497	
Std Dev	0.0006	
Rel Std Dev(%)	1.1625	

Operator's Signature

FFWCC
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000893
01/08/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:36
Control Test	0.080	10:37
Air Blank	0.000	10:37
Control Test	0.078	10:38
Air Blank	0.000	10:39
Control Test	0.079	10:39
Air Blank	0.000	10:40
Control Test Stats		
Average	0.0790	
Std Dev	0.0010	
Rel Std Dev(%)	1.2658	

Operator's Signature

FFWCC
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000893
01/08/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:41
Control Test	0.199	10:42
Air Blank	0.000	10:43
Control Test	0.197	10:43
Air Blank	0.000	10:44
Control Test	0.197	10:45
Air Blank	0.000	10:45
Control Test Stats		
Average	0.1977	
Std Dev	0.0012	
Rel Std Dev(%)	0.5842	

Operator's Signature

wet

FFWCC
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000893
01/08/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:48
Control Test	0.079	10:48
Air Blank	0.000	10:49
Control Test	0.079	10:49
Air Blank	0.000	10:50
Control Test	0.079	10:50
Air Blank	0.000	10:50
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

Dry



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000893, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000893</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>FFWCC</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>01/11/2024</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>16:28</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Benjamin
Siddoway
Digitally signed by
Benjamin Siddoway
Date: 2024.01.12 08:47:17
-05'00'

Date
01/11/2024
BENJAMIN W SIDDOWAY,
Department Inspector

FDLE/ATP Form 69 March 2022
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality