



Agency Pinellas Park PD

S/N 80-001042

Date In 11/6/2023

DI Completion Date 11-13-2023

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By <u>ALL</u>	Quality Checks	By <u>IS</u>	Date <u>11-9-2023</u>	Flow Calibration	By _____	Date _____																																
☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>199</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-105</u> 32 mm <u>0.164</u> (.139 - .169) 36 mm <u>0.175</u> (.156 - .190) 53 mm <u>0.257</u> (.228 - .278) 103 mm <u>0.527</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>30793</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																		
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202303K 03-29-2025</td></tr><tr><td>0.080</td><td>MP6292</td><td>202303L 03-29-2025</td></tr><tr><td>0.200</td><td>MP6293</td><td>202304C 04-05-2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG229803 10-25-2024</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP6291	202303K 03-29-2025	0.080	MP6292	202303L 03-29-2025	0.200	MP6293	202304C 04-05-2025	0.080 DGS	N/A	AG229803 10-25-2024	<table border="1"><thead><tr><th>Maintenance</th><th>By _____</th></tr></thead><tbody><tr><td>☐ Battery Replacement</td><td></td></tr><tr><td>☐ Dry Gas Regulator Replacement</td><td></td></tr><tr><td>☐ Breath Tube Replacement</td><td></td></tr><tr><td>☐ Other _____</td><td></td></tr><tr><td>_____</td><td></td></tr><tr><td>_____</td><td></td></tr><tr><td>_____</td><td></td></tr><tr><td>_____</td><td></td></tr></tbody></table>		Maintenance	By _____	☐ Battery Replacement		☐ Dry Gas Regulator Replacement		☐ Breath Tube Replacement		☐ Other _____		_____		_____		_____		_____	
Simulator	Serial #	Lot #/Exp																																					
0.050	MP6291	202303K 03-29-2025																																					
0.080	MP6292	202303L 03-29-2025																																					
0.200	MP6293	202304C 04-05-2025																																					
0.080 DGS	N/A	AG229803 10-25-2024																																					
Maintenance	By _____																																						
☐ Battery Replacement																																							
☐ Dry Gas Regulator Replacement																																							
☐ Breath Tube Replacement																																							
☐ Other _____																																							

Calibration Adjustment				By _____	
Barometric Pressure Gauge _____				ID # _____	
Simulator	Serial #	Lot #	Expiration		
0.000		N/A	N/A		
0.040					
0.100					
0.200					
0.300					
0.080 DGS	N/A				
☐ Post Calibration Adjustment Stability Checks					
Simulator	Serial #	Lot #	Expiration		
0.050					
0.080					
0.200					
0.080 DGS	N/A				
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____					

Department Inspection		By IS _____	
Barometric Pressure ID# 30793			
Gauge 1023	Instrument 1022		
Mouth Alcohol Solution Lot # 2023-A			
Acetone Stock Solution Lot # 2022-B			
Simulator	Serial Number		
0.000	MP6289		
Interferent	MP6290		
0.050	MP6291		
0.080	MP6292		
0.200	MP6293		
Attachments			
☒ Form 41		☐ Post-Stability Checks	
☒ Stability Checks		☐ Flow Calibration	
☒ Calibration Certificate		☐ Form 40	
☐ Calibration Adjustment		☐ Other _____	
☒ Instrument Complies with Chapter 11D-8, FAC			
☐ Instrument Does Not Comply with Chapter 11D-8, FAC			
☒ Return to/Place into Evidentiary Use			
☐ Remain Out of Evidentiary Use			
☒ Conduct an Agency Inspection Before Evidentiary Use			
Benjamin Siddoway Digitally signed by Benjamin Siddoway Date: 2023.11.13 14:54:09 -05'00'		Phil Nicodemo Digitally signed by Phil Nicodemo Date: 2023.11.15 10:27:02 -05'00'	
Tech Review / Date _____		Admin Review / Date _____	

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PINELLAS PARK PD
Time of Inspection: 10:35

Date of Inspection: 11/13/2023

Serial Number: 80-001042
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG229803 Exp: 10/25/2024
0.000	0.049	0.079	0.198	0.082
0.000	0.049	0.079	0.200	0.082
0.000	0.049	0.080	0.200	0.081
0.000	0.049	0.079	0.200	0.082
0.000	0.049	0.080	0.200	0.082
0.000	0.050	0.080	0.200	0.082
0.000	0.049	0.080	0.199	0.082
0.000	0.049	0.080	0.199	0.082
0.000	0.050	0.080	0.199	0.082
0.000	0.050	0.079	0.200	0.081

Standard Deviations	0.0004	0.0005	0.0007	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0005 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto ISRAEL SOTO
Signature and Printed Name

11/13/2023
Date

Stability Checks

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
11/09/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:15
Control Test	0.050	14:16
Air Blank	0.000	14:17
Control Test	0.049	14:17
Air Blank	0.000	14:18
Control Test	0.048	14:19
Air Blank	0.000	14:19
Control Test Stats		
Average	0.0490	
Std Dev	0.0010	
Rel Std Dev(%)	2.9408	

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
11/09/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:21
Control Test	0.078	14:21
Air Blank	0.000	14:22
Control Test	0.079	14:23
Air Blank	0.000	14:23
Control Test	0.078	14:24
Air Blank	0.000	14:24
Control Test Stats		
Average	0.0783	
Std Dev	0.0006	
Rel Std Dev(%)	0.7370	

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
11/09/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:25
Control Test	0.198	14:26
Air Blank	0.000	14:26
Control Test	0.196	14:27
Air Blank	0.000	14:28
Control Test	0.197	14:28
Air Blank	0.000	14:29
Control Test Stats		
Average	0.1970	
Std Dev	0.0010	
Rel Std Dev(%)	0.5076	

Wet


Operator's Signature


Operator's Signature


Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
11/09/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:30
Control Test	0.081	14:30
Air Blank	0.000	14:31
Control Test	0.081	14:31
Air Blank	0.000	14:32
Control Test	0.080	14:32
Air Blank	0.000	14:33
Control Test Stats		
Average	0.0807	
Std Dev	0.0006	
Rel Std Dev(%)	0.7157	

Dry


Operator's Signature



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001042, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001042</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>PINELLAS PARK PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>11/13/2023</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>10:35</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Israel Soto
Soto
Date: 2023.11.13
13:14:16 -05'00'

11/13/2023

Date

ISRAEL SOTO,
Department Inspector

FDLE/ATP Form 69 March 2022
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality



Agency Pinellas Park Police Department

S/N 80-001042

Date In 09-19-2023

DI Completion Date

☐ Ship

☐ P/U

☐ H/D

Intake	By <u>ALL</u>	Quality Checks	By <u>IS</u>	Date <u>09-21-2023</u>	Flow Calibration	By	Date																																																												
☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>196</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-105</u> 32 mm <u>0.164</u> (.139 - .169) 36 mm <u>0.179</u> (.156 - .190) 53 mm <u>0.246</u> (.228 - .278) 103 mm <u>0.511</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28427</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																														
		<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 15%;">Simulator</th> <th style="width: 20%;">Serial #</th> <th style="width: 25%;">Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td rowspan="2">0.050</td> <td rowspan="2">MP6291</td> <td>202303K</td> </tr> <tr> <td>03-29-2025</td> </tr> <tr> <td rowspan="2">0.080</td> <td rowspan="2">MP6292</td> <td>202303L</td> </tr> <tr> <td>03-29-2025</td> </tr> <tr> <td rowspan="2">0.200</td> <td rowspan="2">MP6293</td> <td>202304C</td> </tr> <tr> <td>04-05-2025</td> </tr> <tr> <td rowspan="2">0.080 DGS</td> <td rowspan="2">N/A</td> <td>06723080A5</td> </tr> <tr> <td>04-05-2025</td> </tr> </tbody> </table>			Simulator	Serial #	Lot #/Exp	0.050	MP6291	202303K	03-29-2025	0.080	MP6292	202303L	03-29-2025	0.200	MP6293	202304C	04-05-2025	0.080 DGS	N/A	06723080A5	04-05-2025	<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 60%;">Maintenance</th> <th style="width: 40%;">By</th> </tr> </thead> <tbody> <tr> <td>☐ Battery Replacement</td> <td></td> </tr> <tr> <td>☐ Dry Gas Regulator Replacement</td> <td></td> </tr> <tr> <td>☐ Breath Tube Replacement</td> <td></td> </tr> <tr> <td>☐ Other _____</td> <td></td> </tr> <tr><td> </td><td></td></tr> <tr><td> </td><td></td></tr> <tr><td> </td><td></td></tr> <tr><td> </td><td></td></tr> </tbody> </table>			Maintenance	By	☐ Battery Replacement		☐ Dry Gas Regulator Replacement		☐ Breath Tube Replacement		☐ Other _____																																
Simulator	Serial #	Lot #/Exp																																																																	
0.050	MP6291	202303K																																																																	
		03-29-2025																																																																	
0.080	MP6292	202303L																																																																	
		03-29-2025																																																																	
0.200	MP6293	202304C																																																																	
		04-05-2025																																																																	
0.080 DGS	N/A	06723080A5																																																																	
		04-05-2025																																																																	
Maintenance	By																																																																		
☐ Battery Replacement																																																																			
☐ Dry Gas Regulator Replacement																																																																			
☐ Breath Tube Replacement																																																																			
☐ Other _____																																																																			
Calibration Adjustment By <u>IS</u>				Department Inspection By <u>IS</u>																																																															
Barometric Pressure Gauge <u>1013</u> ID # <u>28662</u> <table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 15%;">Simulator</th> <th style="width: 20%;">Serial #</th> <th style="width: 25%;">Lot #</th> <th style="width: 40%;">Expiration</th> </tr> </thead> <tbody> <tr><td>0.000</td><td>MP6294</td><td>N/A</td><td>N/A</td></tr> <tr><td>0.040</td><td>MP6295</td><td>22460</td><td>12-28-2024</td></tr> <tr><td>0.100</td><td>MP6296</td><td>22430</td><td>11-30-2024</td></tr> <tr><td>0.200</td><td>MP6297</td><td>22400</td><td>10-12-2024</td></tr> <tr><td>0.300</td><td>MP6298</td><td>23070</td><td>03-06-2025</td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td>AG229803</td><td>10-25-2024</td></tr> </tbody> </table> ☒ Post Calibration Adjustment Stability Checks <table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 15%;">Simulator</th> <th style="width: 20%;">Serial #</th> <th style="width: 25%;">Lot #</th> <th style="width: 40%;">Expiration</th> </tr> </thead> <tbody> <tr><td>0.050</td><td>MP6291</td><td>202303K</td><td>03-29-2025</td></tr> <tr><td>0.080</td><td>MP6292</td><td>202303L</td><td>03-29-2025</td></tr> <tr><td>0.200</td><td>MP6293</td><td>202304C</td><td>04-05-2025</td></tr> <tr><td>0.080 DGS</td><td>N/A</td><td>06723080A5</td><td>04-05-2025</td></tr> </tbody> </table>				Simulator	Serial #	Lot #	Expiration	0.000	MP6294	N/A	N/A	0.040	MP6295	22460	12-28-2024	0.100	MP6296	22430	11-30-2024	0.200	MP6297	22400	10-12-2024	0.300	MP6298	23070	03-06-2025	0.080 DGS	N/A	AG229803	10-25-2024	Simulator	Serial #	Lot #	Expiration	0.050	MP6291	202303K	03-29-2025	0.080	MP6292	202303L	03-29-2025	0.200	MP6293	202304C	04-05-2025	0.080 DGS	N/A	06723080A5	04-05-2025	Barometric Pressure ID# <u>28427</u> Gauge <u>1013/1013</u> Instrument <u>1013/1013</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2022-B</u> <table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 50%;">Simulator</th> <th style="width: 50%;">Serial Number</th> </tr> </thead> <tbody> <tr><td>0.000</td><td>MP6289</td></tr> <tr><td>Interferent</td><td>MP6289</td></tr> <tr><td>0.050</td><td>MP6291</td></tr> <tr><td>0.080</td><td>MP6292</td></tr> <tr><td>0.200</td><td>MP6293</td></tr> </tbody> </table>				Simulator	Serial Number	0.000	MP6289	Interferent	MP6289	0.050	MP6291	0.080	MP6292	0.200	MP6293
Simulator	Serial #	Lot #	Expiration																																																																
0.000	MP6294	N/A	N/A																																																																
0.040	MP6295	22460	12-28-2024																																																																
0.100	MP6296	22430	11-30-2024																																																																
0.200	MP6297	22400	10-12-2024																																																																
0.300	MP6298	23070	03-06-2025																																																																
0.080 DGS	N/A	AG229803	10-25-2024																																																																
Simulator	Serial #	Lot #	Expiration																																																																
0.050	MP6291	202303K	03-29-2025																																																																
0.080	MP6292	202303L	03-29-2025																																																																
0.200	MP6293	202304C	04-05-2025																																																																
0.080 DGS	N/A	06723080A5	04-05-2025																																																																
Simulator	Serial Number																																																																		
0.000	MP6289																																																																		
Interferent	MP6289																																																																		
0.050	MP6291																																																																		
0.080	MP6292																																																																		
0.200	MP6293																																																																		
Notes/Suggested Service: <u>Instrument had Interferent Detect on dry gas test during Department Inspection. Optical Bench Calibration Adjustment performed. Second Department Inspection used same barometric pressure gauge, simulators, and mouth/acetone stocks. Instrument had Interferent Detect on dry gas test on second Department Inspection. Returning instrument to repair. IS</u>				<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 100%;">Attachments</th> </tr> </thead> <tbody> <tr> <td> ☒ Form 41 x2 ☒ Stability Checks ☐ Calibration Certificate ☒ Calibration Adjustment </td> <td> ☒ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other <u>Form 51</u> </td> </tr> </tbody> </table>				Attachments	☒ Form 41 x2 ☒ Stability Checks ☐ Calibration Certificate ☒ Calibration Adjustment	☒ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other <u>Form 51</u>																																																									
Attachments																																																																			
☒ Form 41 x2 ☒ Stability Checks ☐ Calibration Certificate ☒ Calibration Adjustment	☒ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other <u>Form 51</u>																																																																		
☐ Instrument Complies with Chapter 11D-8, FAC ☒ Instrument Does Not Comply with Chapter 11D-8, FAC ☐ Return to/Place into Evidentiary Use ☒ Remain Out of Evidentiary Use ☐ Conduct an Agency Inspection Before Evidentiary Use				<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 50%;">Tech Review / Date</th> <th style="width: 50%;">Admin Review / Date</th> </tr> </thead> <tbody> <tr> <td> </td> <td> </td> </tr> </tbody> </table>				Tech Review / Date	Admin Review / Date																																																										
Tech Review / Date	Admin Review / Date																																																																		

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PINELLAS PARK PD
Time of Inspection: 12:35

Date of Inspection: 09/21/2023

Serial Number: 80-001042
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:06723080A5 Exp: 04/05/2025
0.000	0.049	0.079	0.198	0.080 / 0.080
0.000	0.049	0.079	0.199	0.079 / 0.079
0.000	0.049	0.079	0.198	0.078 / 0.079
0.000	0.050	0.079	0.198	INT / INT
0.000	0.050	0.079	0.199	
0.000	0.050	0.079	0.199	
0.000	0.049	0.079	0.198	
0.000	0.049	0.079	0.199	
0.000	0.049	0.079	0.198	
0.000	0.049	0.079	0.198	

Standard Deviations	0.0005	0.0000	0.0005	
---------------------	--------	--------	--------	--

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: _____ Number of Simulators Used: 5

Remarks:

08: Interferent Detect ENSURED INLET SECURE, SWITCHED TANKS. Non-compliance: WILL PERFORM CAL ADJUST

Switched to a new tank of same lot #

The above instrument complies () does not comply (X) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL SOTO

Signature and Printed Name

09/21/2023
Date

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PINELLAS PARK PD
Time of Inspection: 14:20

Date of Inspection: 09/22/2023

Serial Number: 80-001042
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:06723080A5 Exp: 04/05/2025
0.000	0.049	0.079	0.198	0.081
0.000	0.050	0.080	0.198	INT
0.000	0.051	0.080	0.198	
0.000	0.050	0.081	0.198	
0.000	0.050	0.080	0.199	
0.000	0.050	0.081	0.199	
0.000	0.050	0.081	0.199	
0.000	0.050	0.081	0.199	
0.000	0.050	0.081	0.199	
0.000	0.050	0.081	0.199	

Standard Deviations	0.0004	0.0007	0.0005	
---------------------	--------	--------	--------	--

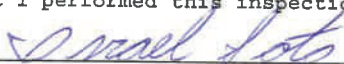
Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: _____ Number of Simulators Used: 5

Remarks:

08: Interferent Detect WILL SEND TO REPAIR. Non-compliance: INT DETECT ON DRY GAS.

The above instrument complies () does not comply (X) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

 ISRAEL SOTO
Signature and Printed Name

09/22/2023
Date

stability checks

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/21/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:52
Control Test	0.050	10:53
Air Blank	0.000	10:53
Control Test	0.049	10:54
Air Blank	0.000	10:55
Control Test	0.048	10:55
Air Blank	0.000	10:56
Control Test Stats		
Average	0.0490	
Std Dev	0.0010	
Rel Std Dev(%)	2.0408	

Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/21/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:57
Control Test	0.078	10:57
Air Blank	0.000	10:58
Control Test	0.078	10:59
Air Blank	0.000	10:59
Control Test	0.078	11:00
Air Blank	0.000	11:00
Control Test Stats		
Average	0.0780	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/21/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:01
Control Test	0.198	11:02
Air Blank	0.000	11:02
Control Test	0.197	11:03
Air Blank	0.000	11:03
Control Test	0.197	11:04
Air Blank	0.000	11:05
Control Test Stats		
Average	0.1973	
Std Dev	0.0006	
Rel Std Dev(%)	0.2926	

Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/21/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:06
Control Test	0.080	11:06
Air Blank	0.000	11:06
Control Test	0.079	11:07
Air Blank	0.000	11:07
Control Test	0.080	11:07
Air Blank	0.000	11:08
Control Test Stats		
Average	0.0797	
Std Dev	0.0006	
Rel Std Dev(%)	0.7247	

Operator's Signature

wet

Dry

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/22/2023 09:22:48

Optical Bench
Calibration Adjustment

Auto Calibration

pg 1 of 2

<<<<< 3um >>>>> <<<<< 9um >>>>>

Solution = 0.000 g/210L or 0.0000 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	0.2700	(-0.0440)	0.1780	(-0.0260)
Sample #2	0.2330	(0.0800)	0.1550	(0.0300)
Sample #3	0.2230	(0.1030)	0.1720	(0.0310)
Sample #4	0.2180	(0.1530)	0.1590	(0.0460)
Avg % Abs	0.2247	(0.1120)	0.1620	(0.0357)
STD DEV	0.0076	(0.0373)	0.0089	(0.0090)
REL STD DEV	3.400	(33.324)	5.487	(25.130)

Solution = 0.040 g/210L or 0.1905 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	0.9610	(-0.0280)	1.6470	(-0.0080)
Sample #2	0.9080	(0.0410)	1.6040	(0.0370)
Sample #3	0.8800	(0.0900)	1.6000	(0.0470)
Sample #4	0.8920	(0.1070)	1.6110	(0.0430)
Avg % Abs	0.8933	(0.0793)	1.6050	(0.0423)
STD DEV	0.0140	(0.0343)	0.0056	(0.0050)
REL STD DEV	1.572	(43.196)	0.347	(11.890)

Solution = 0.100 g/210L or 0.4762 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	1.9960	(-0.0120)	3.7550	(-0.0110)
Sample #2	1.9750	(0.0450)	3.7540	(0.0080)
Sample #3	1.9730	(0.0930)	3.7340	(0.0400)
Sample #4	1.9580	(0.1080)	3.7330	(0.0440)
Avg % Abs	1.9687	(0.0820)	3.7403	(0.0307)
STD DEV	0.0093	(0.0329)	0.0118	(0.0197)
REL STD DEV	0.472	(40.133)	0.317	(64.342)

Solution = 0.200 g/210L or 0.9524 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	3.6840	(-0.0180)	7.1010	(-0.0150)
Sample #2	3.6510	(0.0760)	7.0650	(0.0470)
Sample #3	3.6130	(0.1030)	7.0490	(0.0600)
Sample #4	3.6320	(0.1180)	7.0460	(0.0730)
Avg % Abs	3.6320	(0.0990)	7.0533	(0.0600)
STD DEV	0.0190	(0.0213)	0.0102	(0.0130)
REL STD DEV	0.523	(21.499)	0.145	(21.667)

Solution = 0.300 g/210L or 1.4286 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	5.3350	(-0.0150)	10.2650	(0.0000)
Sample #2	5.2800	(0.0610)	10.2120	(0.0680)
Sample #3	5.2540	(0.1220)	10.2030	(0.0970)
Sample #4	5.2480	(0.1170)	10.1770	(0.1030)
Avg % Abs	5.2607	(0.1000)	10.1973	(0.0893)
STD DEV	0.0170	(0.0339)	0.0182	(0.0187)
REL STD DEV	0.323	(33.867)	0.178	(20.952)

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/22/2023 09:22:48

Optical Bench
Calibration Adjustment

Auto Calibration

pg 2 of 2

<<<< 3um >>>>

Zero Order Coef -569.76
First Order Coef 2672.04
Second Order Coef 28.86

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.001	-0.0007
0.040	0.039	0.0014
0.100	0.101	-0.0009
0.200	0.200	0.0002
0.300	0.300	0.0000

<<<< 9um >>>>

-202.02
1280.20
13.80

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.000	-0.0001
0.040	0.040	0.0003
0.100	0.100	-0.0004
0.200	0.200	0.0002
0.300	0.300	-0.0000

<<<< 3um >>>>

Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1

Sample

Sample #1	2473.00	3210.00
Sample #2	2512.00	3190.00
Sample #3	2553.00	3187.00
Sample #4	2554.00	3216.00
Avg	2539.6667	3197.6667
STD DEV	23.9653	15.9478
REL STD DEV	0.944	0.499
H2O adjust (mg/l*10k)	1270	612

Barometric Pressure = 1013

*****CALIBRATION SUCCESSFUL*****

Post Stability checks

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/22/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:33
Control Test	0.050	10:34
Air Blank	0.000	10:34
Control Test	0.049	10:35
Air Blank	0.000	10:35
Control Test	0.049	10:36
Air Blank	0.000	10:37
Control Test Stats		
Average	0.0493	
Std Dev	0.0006	
Rel Std Dev(%)	1.1703	

Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/22/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:38
Control Test	0.079	10:39
Air Blank	0.000	10:39
Control Test	0.079	10:40
Air Blank	0.000	10:40
Control Test	0.079	10:41
Air Blank	0.000	10:42
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

wet

Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/22/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:44
Control Test	0.198	10:45
Air Blank	0.000	10:46
Control Test	0.196	10:46
Air Blank	0.000	10:47
Control Test	0.197	10:47
Air Blank	0.000	10:48
Control Test Stats		
Average	0.1970	
Std Dev	0.0010	
Rel Std Dev(%)	0.5076	

Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
09/22/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:49
Control Test	0.080	10:50
Air Blank	0.000	10:50
Control Test	0.080	10:50
Air Blank	0.000	10:51
Control Test	0.081	10:51
Air Blank	0.000	10:52
Control Test Stats		
Average	0.0803	
Std Dev	0.0006	
Rel Std Dev(%)	0.7187	

Dry

Operator's Signature

Return Material Authorization

Ship to: ☒ CMI, Inc.

☐ Enforcement Electronics

Shipment to repair facility authorized by: Israel Soto on 09-25-2023

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-001042

Bill To Address:

Returning to repair

Ship to Address:

FDLE Off-Site Mail Facility

c/o Florida Dept of Law Enforcement

Alcohol Testing Program

813 B Lake Bradford Road

Tallahassee, FL 32304

Reason for Return:

Instrument was shipped to us from repair. During calibration instrument produced an Interferent Detect twice during the dry 0.080 g/210L test. An optical bench calibration adjustment was successfully performed. On second calibration, after adjustment, instrument again produced Interferent Detect during the dry 0.080 g/210L test. See included paperwork.

Please choose one of the following options:

- ☐ 1. I _____, authorize all repairs.
- ☐ 2. I _____, authorize repairs up to \$_____.
- ☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Jeffrey DeTrano

Phone #: 727-369-7864 Email: jdetrano@pinellas-park.com

ATP Contact Name: Israel Soto ATP Email: israelsoto@fdle.state.fl.us



INSTRUMENT PROCESSING SHEET

Agency Pinellas Park Police DepartmentS/N 80-001042Florida Department of
Law EnforcementDate In 03-28-2023 DI Completion Date 05-01-2023☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By IS	Quality Checks	By IS	Date <u>04-03-2023</u>	Flow Calibration	By	Date																										
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☐ Handle ☒ Keyboard ☐ Dry Gas Shelf ☒ Feet ☐ Breath Tube ☒ Ports ☐ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>Exhaust port cover found broken and loose in box.</u>		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>187</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-102</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.167</u> (.156 - .190) 53 mm <u>0.238</u> (.228 - .278) 103 mm <u>0.496</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>30793</u> ☒ Stability Checks			☐ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																												
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202201C 01-11-2024</td></tr><tr><td>0.080</td><td>MP6292</td><td>202201D 01-18-2024</td></tr><tr><td>0.200</td><td>MP6293</td><td>202201E 01-18-2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG229803 10-25-2024</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP6291	202201C 01-11-2024	0.080	MP6292	202201D 01-18-2024	0.200	MP6293	202201E 01-18-2024	0.080 DGS	N/A	AG229803 10-25-2024																
Simulator	Serial #	Lot #/Exp																															
0.050	MP6291	202201C 01-11-2024																															
0.080	MP6292	202201D 01-18-2024																															
0.200	MP6293	202201E 01-18-2024																															
0.080 DGS	N/A	AG229803 10-25-2024																															
Calibration Adjustment		By _____		Department Inspection																													
Barometric Pressure Gauge _____ ID # _____				Barometric Pressure ID# <u>28662</u>																													
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A					Gauge <u>1006</u> Instrument <u>1005</u>	
Simulator	Serial #	Lot #	Expiration																														
0.000		N/A	N/A																														
0.040																																	
0.100																																	
0.200																																	
0.300																																	
0.080 DGS	N/A																																
☐ Post Calibration Adjustment Stability Checks				Mouth Alcohol Solution Lot # <u>2022-A</u>																													
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A					Acetone Stock Solution Lot # <u>2022-B</u>									
Simulator	Serial #	Lot #	Expiration																														
0.050																																	
0.080																																	
0.200																																	
0.080 DGS	N/A																																
Notes/Suggested Service: <u>Exhaust port cover replaced prior to Department Inspection. IS</u>				<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5086</td></tr><tr><td>Interferent</td><td>MP5087</td></tr><tr><td>0.050</td><td>MP5088</td></tr><tr><td>0.080</td><td>MP5089</td></tr><tr><td>0.200</td><td>MP5090</td></tr></tbody></table>		Simulator	Serial Number	0.000	MP5086	Interferent	MP5087	0.050	MP5088	0.080	MP5089	0.200	MP5090																
Simulator	Serial Number																																
0.000	MP5086																																
Interferent	MP5087																																
0.050	MP5088																																
0.080	MP5089																																
0.200	MP5090																																
				Attachments																													
				☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment																													
				☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 ☐ Other _____																													
				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC																													
				☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use																													
				☒ Conduct an Agency Inspection Before Evidentiary Use																													
		Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2023.05.02 13:34:21 -0400		Phil Nicodemo Digitally signed by Phil Nicodemo Date: 2023.05.04 15:39:00 -0400																													
		Tech Review / Date		Admin Review / Date																													

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PINELLAS PARK PD
Time of Inspection: 15:01

Date of Inspection: 05/01/2023

Serial Number: 80-001042
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG229803 Exp: 10/25/2024
0.000	0.047	0.078	0.198	0.076
0.000	0.048	0.078	0.198	0.076
0.000	0.048	0.078	0.198	0.076
0.000	0.048	0.077	0.198	0.076
0.000	0.048	0.077	0.198	0.076
0.000	0.049	0.080	0.198	0.076
0.000	0.049	0.080	0.199	0.076
0.000	0.049	0.081	0.199	0.076
0.000	0.049	0.082	0.199	0.076
0.000	0.049	0.082	0.199	0.077

Standard Deviations	0.0006	0.0019	0.0005	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0008 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



ISRAEL SOTO

Signature and Printed Name

05/01/2023
Date

stability checks

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
04/03/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:43
Control Test	0.048	14:43
Air Blank	0.000	14:44
Control Test	0.047	14:45
Air Blank	0.000	14:45
Control Test	0.047	14:46
Air Blank	0.000	14:46
Control Test Stats		
Average	0.0473	
Std Dev	0.0006	
Rel Std Dev(%)	1.2198	

Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
04/03/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:47
Control Test	0.077	14:48
Air Blank	0.000	14:49
Control Test	0.078	14:49
Air Blank	0.000	14:50
Control Test	0.078	14:50
Air Blank	0.000	14:51
Control Test Stats		
Average	0.0777	
Std Dev	0.0006	
Rel Std Dev(%)	0.7434	

wet

Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
04/03/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:52
Control Test	0.198	14:53
Air Blank	0.000	14:53
Control Test	0.198	14:54
Air Blank	0.000	14:55
Control Test	0.198	14:55
Air Blank	0.000	14:56
Control Test Stats		
Average	0.1980	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

PINELLAS PARK PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001042
04/03/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:57
Control Test	0.078	14:57
Air Blank	0.000	14:58
Control Test	0.077	14:58
Air Blank	0.000	14:58
Control Test	0.077	14:59
Air Blank	0.000	14:59
Control Test Stats		
Average	0.0773	
Std Dev	0.0006	
Rel Std Dev(%)	0.7466	

Dry

Operator's Signature



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001042 , manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001042</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>PINELLAS PARK PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>05/01/2023</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>15:01</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Israel Soto
Soto
Digitally signed by Israel Soto
Date: 2023.05.01 15:09:55 -04'00'

05/01/2023

Date

ISRAEL SOTO,
Department Inspector

FDLE/ATP Form 69 March 2022
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: PINELLAS PARK PD
Time of Inspection: 13:21

Date of Inspection: 05/01/2023

Serial Number: 80-001042
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

BYPASSED AI TO OPERATE INSTRUMENT

N/A Compliance not determined

The above instrument complies (☒) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL SOTO

Signature and Printed Name

05/01/2023
Date