



Agency St Johns County SO

S/N 80-001133

Date In 11/27/2023 DI Completion Date 11/29/2023

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

FDLE/ATP Form 48 January 2022
Issuing Authority: Alcohol Testing Program

Page 1 of 1

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: ST JOHNS COUNTY SO

Time of Inspection: 09:17

Date of Inspection: 11/28/2023

Serial Number: 80-001133

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

COMPLIANCE NOT DETERMINED, BYPASS AI TO OPERATE

N/A, compliance not determined

The above instrument ~~complies (X)~~ does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Benjamin W Siddoway

BENJAMIN W SIDDOWAY

Signature and Printed Name

11/28/2023
Date

Stability Checks and Flow Calibration Adjustment

ST JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
11/28/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:52
Control Test	0.050	10:53
Air Blank	0.000	10:54
Control Test	0.051	10:54
Air Blank	0.000	10:55
Control Test	0.050	10:55
Air Blank	0.000	10:56
Control Test Stats		
Average	0.0503	
Std Dev	0.0006	
Rel Std Dev(%)	1.1471	

Benjamin Saldy
Operator's Signature

ST JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
11/28/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:09
Control Test	0.079	11:10
Air Blank	0.000	11:10
Control Test	0.079	11:11
Air Blank	0.000	11:12
Control Test	0.079	11:12
Air Blank	0.000	11:13
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Benjamin Saldy
Operator's Signature

ST JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
11/28/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:17
Control Test	0.200	11:18
Air Blank	0.000	11:18
Control Test	0.200	11:19
Air Blank	0.000	11:20
Control Test	0.198	11:20
Air Blank	0.000	11:21
Control Test Stats		
Average	0.1993	
Std Dev	0.0012	
Rel Std Dev(%)	0.5793	

Benjamin Saldy
Operator's Signature

ST JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
11/28/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:22
Control Test	0.081	11:23
Air Blank	0.000	11:23
Control Test	0.081	11:23
Air Blank	0.000	11:24
Control Test	0.081	11:24
Air Blank	0.000	11:25
Control Test Stats		
Average	0.0810	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

DGS

Benjamin Saldy
Operator's Signature

ST JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
11/28/2023
Software: 8100.27

Flow Rate Calibration*****
1: Rate (Liters/min) = 5
SQRT(Diff) = 6.633
2: Rate (Liters/min) = 15
SQRT(Diff) = 11.445
3: Rate (Liters/min) = 30
SQRT(Diff) = 20.395
Dependent Data Scale Factor = 100000 L/min
Independent Data Scale Factor = 256
Rounded Slope = 703
Rounded Intercept = -640227
Correlation = 0.99838

Benjamin Saldy

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: ST JOHNS COUNTY SO
Time of Inspection: 10:42

Date of Inspection: 11/29/2023

Serial Number: 80-001133
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG325603 Exp: 09/13/2025
0.000	0.050	0.079	0.201	0.081
0.000	0.050	0.079	0.201	0.080
0.000	0.050	0.079	0.201	0.081
0.000	0.050	0.079	0.200	0.080
0.000	0.050	0.079	0.201	0.080
0.000	0.050	0.080	0.201	0.080
0.000	0.051	0.079	0.201	0.080
0.000	0.050	0.079	0.200	0.080
0.000	0.050	0.080	0.200	0.081
0.000	0.050	0.080	0.201	0.080


Standard Deviations	0.0003	0.0004	0.0004	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

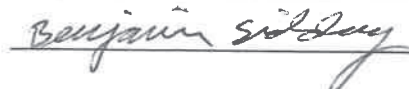

BENJAMIN W SIDDOWAY
Signature and Printed Name

11/29/2023
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001133 , manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001133</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>ST JOHNS COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>11/29/2023</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>10:42</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Benjamin
Siddoway
Digitally signed by Benjamin
Siddoway
Date: 2023.11.29 11:12:24
-05'00'

Date

11/29/2023

BENJAMIN W SIDDOWAY,
Department Inspector

FDLE/ATP Form 69 March 2022

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality



INSTRUMENT PROCESSING SHEET

Agency St John's County Sheriff's OfficeS/N 80-001133Florida Department of
Law EnforcementDate In 3/7/2023DI Completion Date 03-24-2023☐ Ship ☐ P/U ☐ H/D ☒ CMI ☐ EE

Intake	By PN	Quality Checks	By PN	Date 3/23/2023	Flow Calibration	By PN	Date 3/23/2023																																																																				
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: Prior inspection (2022) indicated instrument needed flow sensor and battery change. Instrument arrived with same items not addressed. PN 		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>92*</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 105</u> 32 mm <u>.140</u> (.139 - .169) 36 mm <u>.152*</u> (.156 - .190) 53 mm <u>.230</u> (.228 - .278) 103 mm <u>.515</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28662</u> ☒ Stability Checks			☒ Flow Column # <u>ATP 102</u> ☒ 5L/min - 17mm ☒ 15L/min - 53mm ☒ 30L/min - 103mm ☒ R-Value <u>92</u> ☒ Post Calibration Verification (L/s) Flow Column # <u>ATP 105</u> 32 mm <u>.140</u> (.139 - .169) 36 mm <u>.156</u> (.156 - .190) 53 mm <u>.238</u> (.228 - .278) 103 mm <u>.523</u> (.447 - .547)																																																																						
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202201C 01/11/2024</td></tr><tr><td>0.080</td><td>MP5089</td><td>202201D 01/18/2024</td></tr><tr><td>0.200</td><td>MP5090</td><td>202201E 01/18/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG229803 10-25-2024</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP5088	202201C 01/11/2024	0.080	MP5089	202201D 01/18/2024	0.200	MP5090	202201E 01/18/2024	0.080 DGS	N/A	AG229803 10-25-2024	<table border="1"><thead><tr><th colspan="2">Maintenance</th><th>By</th></tr></thead><tbody><tr><td>☐ Battery Replacement</td><td></td><td></td></tr><tr><td>☐ Dry Gas Regulator Replacement</td><td></td><td></td></tr><tr><td>☐ Breath Tube Replacement</td><td></td><td></td></tr><tr><td>☐ Other</td><td></td><td></td></tr><tr><td></td><td></td><td></td></tr><tr><td></td><td></td><td></td></tr><tr><td></td><td></td><td></td></tr><tr><td></td><td></td><td></td></tr></tbody></table>			Maintenance		By	☐ Battery Replacement			☐ Dry Gas Regulator Replacement			☐ Breath Tube Replacement			☐ Other																																								
Simulator	Serial #	Lot #/Exp																																																																									
0.050	MP5088	202201C 01/11/2024																																																																									
0.080	MP5089	202201D 01/18/2024																																																																									
0.200	MP5090	202201E 01/18/2024																																																																									
0.080 DGS	N/A	AG229803 10-25-2024																																																																									
Maintenance		By																																																																									
☐ Battery Replacement																																																																											
☐ Dry Gas Regulator Replacement																																																																											
☐ Breath Tube Replacement																																																																											
☐ Other																																																																											
Calibration Adjustment By _____				Department Inspection By <u>IS</u>																																																																							
Barometric Pressure Gauge _____ ID # _____ <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>				Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Barometric Pressure ID# <u>30793</u> Gauge <u>1016</u> Instrument <u>1008</u> Mouth Alcohol Solution Lot # <u>2022-A</u> Acetone Stock Solution Lot # <u>2022-B</u> <table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6289</td></tr><tr><td>Interferent</td><td>MP6290</td></tr><tr><td>0.050</td><td>MP6291</td></tr><tr><td>0.080</td><td>MP6292</td></tr><tr><td>0.200</td><td>MP6293</td></tr></tbody></table> Attachments <table border="1"><tbody><tr><td>☒ Form 41</td><td>☐ Post-Stability Checks</td></tr><tr><td>☒ Stability Checks</td><td>☒ Flow Calibration</td></tr><tr><td>☒ Calibration Certificate</td><td>☒ Form 40</td></tr><tr><td>☐ Calibration Adjustment</td><td>☒ Other <u>Form 51</u></td></tr></tbody></table>				Simulator	Serial Number	0.000	MP6289	Interferent	MP6290	0.050	MP6291	0.080	MP6292	0.200	MP6293	☒ Form 41	☐ Post-Stability Checks	☒ Stability Checks	☒ Flow Calibration	☒ Calibration Certificate	☒ Form 40	☐ Calibration Adjustment	☒ Other <u>Form 51</u>
Simulator	Serial #	Lot #	Expiration																																																																								
0.000		N/A	N/A																																																																								
0.040																																																																											
0.100																																																																											
0.200																																																																											
0.300																																																																											
0.080 DGS	N/A																																																																										
Simulator	Serial #	Lot #	Expiration																																																																								
0.050																																																																											
0.080																																																																											
0.200																																																																											
0.080 DGS	N/A																																																																										
Simulator	Serial Number																																																																										
0.000	MP6289																																																																										
Interferent	MP6290																																																																										
0.050	MP6291																																																																										
0.080	MP6292																																																																										
0.200	MP6293																																																																										
☒ Form 41	☐ Post-Stability Checks																																																																										
☒ Stability Checks	☒ Flow Calibration																																																																										
☒ Calibration Certificate	☒ Form 40																																																																										
☐ Calibration Adjustment	☒ Other <u>Form 51</u>																																																																										
Notes/Suggested Service: R-value below 100, sending instrument to repair. IS Tech Review: Instrument was meant to be sent to repair by Agency after 2022 Department Inspection. Instrument accidentally shipped back to ATP. Instrument being sent to repair directly by ATP. Instrument was not used between 2022 and 2023 Department Inspection. IS 04-13-2023				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☐ Return to/Place into Evidentiary Use ☒ Remain Out of Evidentiary Use ☐ Conduct an Agency Inspection Before Evidentiary Use Taylor Gutschow Tech Review / Date 09:10:08 -04'00' Admin Review / Date																																																																							

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: ST JOHNS COUNTY SO
Time of Inspection: 10:11

Date of Inspection: 03/24/2023

Serial Number: 80-001133
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG229803 Exp: 10/25/2024
0.000	0.049	0.077	0.200	0.077
0.000	0.049	0.078	0.200	0.078
0.000	0.049	0.078	0.200	0.078
0.000	0.049	0.078	0.200	0.077
0.000	0.049	0.078	0.200	0.077
0.000	0.049	0.078	0.200	0.078
0.000	0.049	0.078	0.200	0.078
0.000	0.049	0.077	0.200	0.078
0.000	0.049	0.076	0.200	0.077
0.000	0.049	0.075	0.200	0.078

Standard Deviations	0.0000	0.0010	0.0000	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto ISRAEL SOTO
Signature and Printed Name

03/24/2023
Date

ST JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
03/23/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	15:52
Control Test	0.048	15:53
Air Blank	0.000	15:53
Control Test	0.049	15:54
Air Blank	0.000	15:55
Control Test	0.048	15:55
Air Blank	0.000	15:56
Control Test Stats		
Average	0.0483	
Std Dev	0.0006	
Rel Std Dev(%)	1.1945	



Operator's Signature

ST JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
03/23/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	15:57
Control Test	0.078	15:57
Air Blank	0.000	15:58
Control Test	0.078	15:59
Air Blank	0.000	15:59
Control Test	0.078	16:00
Air Blank	0.000	16:01
Control Test Stats		
Average	0.0780	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	



Operator's Signature

ST JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
03/23/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	16:02
Control Test	0.200	16:03
Air Blank	0.000	16:03
Control Test	0.199	16:04
Air Blank	0.000	16:04
Control Test	0.198	16:05
Air Blank	0.000	16:06
Control Test Stats		
Average	0.1990	
Std Dev	0.0010	
Rel Std Dev(%)	0.5025	



Operator's Signature

ST JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
03/23/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	16:08
Control Test	0.078	16:08
Air Blank	0.000	16:09
Control Test	0.079	16:09
Air Blank	0.000	16:10
Control Test	0.079	16:10
Air Blank	0.000	16:11
Control Test Stats		
Average	0.0787	
Std Dev	0.0006	
Rel Std Dev(%)	0.7335	

DGS



Operator's Signature

STABILITY CHECKS



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001133, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001133</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>ST. JOHNS COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>03/24/2023</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>10:11</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Israel Soto
Soto
Date: 2023.03.24
10:47:10 -04'00'

03/24/2023

Date

ISRAEL SOTO,
Department Inspector

FDLE/ATP Form 69 March 2022
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

FLOW CALIBRATION ADJUSTMENT

ST. JOHNS COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001133
03/23/2023
Software: 8100.27

Flow Rate Calibration*****

1: Rate (Liters/min) = 5
SQRT(Diff) = 4.469

2: Rate (Liters/min) = 15
SQRT(Diff) = 9.379

3: Rate (Liters/min) = 30
SQRT(Diff) = 19.207

Dependent Data Scale Factor = 100000 L/min

Independent Data Scale Factor = 256

Rounded Slope = 653

Rounded Intercept = -175517

Correlation = 0.99716



Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: ST JOHNS COUNTY SO
Time of Inspection: 13:43

Date of Inspection: 03/23/2023

Serial Number: 80-001133
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:


Number of Simulators Used: _____

Remarks:

BYPASSED AI TO OPERATE INSTRUMENT, COMPLIANCE NOT DETERMINED

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

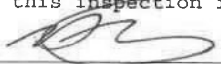

PHIL NICODEMO

Signature and Printed Name

03/23/2023
Date

Return Material Authorization

Ship to: ☒ CMI, Inc.

☐ Enforcement Electronics

Shipment to repair facility authorized by: Alex Ellis on 04-10-2023

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-001133

Bill To Address:

St. John's County Sheriff's Office

Ship to Address:

FDLE Off-Site Mail Facility

Alcohol Testing Program

813 B Lake Bradford Road

Tallahassee, FL 32304

Reason for Return:

Instrument R-Value is below 100.

Please choose one of the following options:

☐ 1. I _____, authorize all repairs.

☐ 2. I _____, authorize repairs up to \$_____.

☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Alex Ellis

Phone #: 904-755-5071

Email: aellis@sjso.org

ATP Contact Name: Phil Nicodemo ATP Email: philipnicodemo@fdle.state.fl.us