



Agency Santa Rosa County Sheriff's Office

S/N 80-001302

Date In 06/09/2023 DI Completion Date 06-13-2023

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Calibration Adjustment

By_

Barometric Pressure Gauge _____ ID # _____

❑ Post Calibration Adjustment Stability Checks

Notes/Suggested Service: Instrument displaying "Disabled Mode: Memory Full". Instrument was uploaded and cleared, no subject breath tests were present. Instrument continued displaying "Disabled Mode: Memory Full". Instrument entered Ready Mode after contacting CMI (see maintenance). IS

Admin Review: Removed extra "2" from Barometric Pressure ID#
under Department Inspection. IS 06-14-2023

Department Inspection

By ISBarometric Pressure ID# 284221 IS 06-14-2023Gauge 1011 Instrument 1010Mouth Alcohol Solution Lot # 2022-A

Acetone Stock Solution Lot # 2022-B

Attachments

- | | |
|---|--|
| ☒ Form 41 | ☐ Post-Stability Checks |
| ☒ Stability Checks | ☐ Flow Calibration |
| ☒ Calibration Certificate | ☐ Form 40 |
| ☐ Calibration Adjustment | ☐ Other _____ |

- ☒ Instrument Complies with Chapter 11D-8, FAC
 - ☐ Instrument Does Not Comply with Chapter 11D-8, FAC
 - ☒ Return to/Place into Evidentiary Use
 - ☐ Remain Out of Evidentiary Use
 - ☒ Conduct an Agency Inspection Before Evidentiary Use

Taylor
Gutschow

Digitally signed by Taylor
Gutschow
Date: 2023.06.14 10:59:45
-04'00'

Phil Nicodemo

Digitally signed by Phil Nicodemo
Date: 2023.06.14 12:07:18 -04'00'

Tech Review / Date

Admin Review / Date

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: SANTA ROSA COUNTY SO
Time of Inspection: 14:29

Date of Inspection: 06/13/2023

Serial Number: 80-001302
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG229803 Exp: 10/25/2024
0.000	0.049	0.078 / 0.080	0.196	0.079
0.000	0.047	0.078 / 0.080	0.196	0.079
0.000	0.048	0.079 / 0.078	0.195	0.078
0.000	0.048	0.078 / 0.078	0.196	0.078
0.000	0.048	0.079 / 0.078	0.195	0.079
0.000	0.049	0.078 / 0.077	0.196	0.079
0.000	0.048	0.078 / 0.078	0.196	0.078
0.000	0.047	0.078 / 0.078	0.195	0.078
0.000	0.047	/ 0.078	0.196	0.079
0.000	0.048	/ 0.078	0.195	0.079

Standard Deviations	0.0007	/ 0.0009	0.0005	0.0005
---------------------	--------	----------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0006 Number of Simulators Used: 5

Remarks:

08: RFI Detect COWORKER WALKED BY WITH CELLPHONE.

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL SOTO

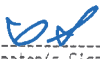
Signature and Printed Name

06/13/2023
Date

Stability checks

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
06/13/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:27
Control Test	0.048	10:28
Air Blank	0.000	10:29
Control Test	0.048	10:29
Air Blank	0.000	10:30
Control Test	0.048	10:31
Air Blank	0.000	10:31
Control Test Stats		
Average	0.0480	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
06/13/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:33
Control Test	0.080	10:33
Air Blank	0.000	10:34
Control Test	0.079	10:34
Air Blank	0.000	10:35
Control Test	0.079	10:36
Air Blank	0.000	10:36
Control Test Stats		
Average	0.0793	
Std Dev	0.0006	
Rel Std Dev(%)	0.7277	

wet


Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
06/13/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:37
Control Test	0.196	10:38
Air Blank	0.000	10:39
Control Test	0.196	10:39
Air Blank	0.000	10:40
Control Test	0.196	10:40
Air Blank	0.000	10:41
Control Test Stats		
Average	0.1960	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator's Signature

SANTA ROSA COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001302
06/13/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:42
Control Test	0.079	10:42
Air Blank	0.000	10:43
Control Test	0.079	10:43
Air Blank	0.000	10:44
Control Test	0.079	10:44
Air Blank	0.000	10:44
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Dry


Operator's Signature



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001302, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001302</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>SANTA ROSA COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>06/13/2023</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>14:29</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Israel Soto
Soto
Date: 2023.06.13 15:40:24
-04'00'

Date

ISRAEL SOTO,
Department Inspector

FDLE/ATP Form 69 March 2022

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality