



INSTRUMENT PROCESSING SHEET

Agency Florida Highway PatrolS/N 80-003407Florida Department of
Law EnforcementDate In 6/1/2023DI Completion Date 6/6/2023☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By PN	Quality Checks	By PN	Date <u>6/5/2023</u>	Flow Calibration	By	Date																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: <u>Instrument hand delivered.</u> <u>Missing two feet that sit directly</u> <u>under the dry gas shelf. PN</u> 		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>264</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 102</u> 32 mm <u>.148</u> (.139 - .169) 36 mm <u>.164</u> (.156 - .190) 53 mm <u>.238</u> (.228 - .278) 103 mm <u>.500</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>30793</u> ☒ Stability Checks			☐ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																														
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202201C 01/11/2024</td></tr><tr><td>0.080</td><td>MP6292</td><td>202201D 01/18/2024</td></tr><tr><td>0.200</td><td>MP6293</td><td>202201E 01/18/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>06723080A5 4/5/2025</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP6291	202201C 01/11/2024	0.080	MP6292	202201D 01/18/2024	0.200	MP6293	202201E 01/18/2024	0.080 DGS	N/A	06723080A5 4/5/2025																		
Simulator	Serial #	Lot #/Exp																																	
0.050	MP6291	202201C 01/11/2024																																	
0.080	MP6292	202201D 01/18/2024																																	
0.200	MP6293	202201E 01/18/2024																																	
0.080 DGS	N/A	06723080A5 4/5/2025																																	
Calibration Adjustment		By PN		Department Inspection			By PN																												
Barometric Pressure Gauge <u>1012/1011</u> ID # <u>28421/28421</u>				Barometric Pressure ID# <u>30793</u>																															
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td>MP5091</td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td>MP5082</td><td>22460</td><td>12/28/2024</td></tr><tr><td>0.100</td><td>MP6296</td><td>22430</td><td>11/30/2024</td></tr><tr><td>0.200</td><td>MP5084</td><td>22400</td><td>10/12/2024</td></tr><tr><td>0.300</td><td>MP5085</td><td>23070</td><td>3/6/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG229803</td><td>10/25/2024</td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000	MP5091	N/A	N/A	0.040	MP5082	22460	12/28/2024	0.100	MP6296	22430	11/30/2024	0.200	MP5084	22400	10/12/2024	0.300	MP5085	23070	3/6/2025	0.080 DGS	N/A	AG229803	10/25/2024			Gauge <u>1010</u> Instrument <u>1010</u> Mouth Alcohol Solution Lot # <u>2022-A</u> Acetone Stock Solution Lot # <u>2022-B</u>			
Simulator	Serial #	Lot #	Expiration																																
0.000	MP5091	N/A	N/A																																
0.040	MP5082	22460	12/28/2024																																
0.100	MP6296	22430	11/30/2024																																
0.200	MP5084	22400	10/12/2024																																
0.300	MP5085	23070	3/6/2025																																
0.080 DGS	N/A	AG229803	10/25/2024																																
☒ Post Calibration Adjustment Stability Checks				<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6289</td></tr><tr><td>Interferent</td><td>MP6290</td></tr><tr><td>0.050</td><td>MP6291</td></tr><tr><td>0.080</td><td>MP6292</td></tr><tr><td>0.200</td><td>MP6293</td></tr></tbody></table>			Simulator	Serial Number	0.000	MP6289	Interferent	MP6290	0.050	MP6291	0.080	MP6292	0.200	MP6293			Attachments														
Simulator	Serial Number																																		
0.000	MP6289																																		
Interferent	MP6290																																		
0.050	MP6291																																		
0.080	MP6292																																		
0.200	MP6293																																		
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202201C</td><td>01/11/2024</td></tr><tr><td>0.080</td><td>MP6292</td><td>202201D</td><td>01/18/2024</td></tr><tr><td>0.200</td><td>MP6293</td><td>202201E</td><td>01/18/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>06723080A5</td><td>4/5/2025</td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.050	MP6291	202201C	01/11/2024	0.080	MP6292	202201D	01/18/2024	0.200	MP6293	202201E	01/18/2024	0.080 DGS	N/A	06723080A5	4/5/2025			☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☒ Calibration Adjustment X2			☒ Post-Stability Checks X2 ☐ Flow Calibration ☐ Form 40 ☐ Other _____								
Simulator	Serial #	Lot #	Expiration																																
0.050	MP6291	202201C	01/11/2024																																
0.080	MP6292	202201D	01/18/2024																																
0.200	MP6293	202201E	01/18/2024																																
0.080 DGS	N/A	06723080A5	4/5/2025																																
Notes/Suggested Service: <u>Calibration adjustment was performed on 6/6.</u> <u>Post-stability checks resulted in the 0.08 g/210L</u> <u>wet bath and DGS averages exceeding 0.003 g/210L.</u> <u>Repeated calibration adjustment and post stability</u> <u>checks on 6/6 and all measurements were nominal.</u> <u>PN 6/6/2023</u>				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																															
				Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2023.06.07 10:20:52 -0400			Shayla Platt Digitally signed by Shayla Platt Date: 2023.06.12 10:12:36 -0400																												
				Tech Review / Date			Admin Review / Date																												

STABILITY CHECKS

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/05/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:14
Control Test	0.048	11:15
Air Blank	0.000	11:15
Control Test	0.047	11:16
Air Blank	0.000	11:16
Control Test	0.047	11:17
Air Blank	0.000	11:18
Control Test Stats		
Average	0.0473	
Std Dev	0.0006	
Rel Std Dev(%)	1.2198	

Operator's Signature

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/05/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:19
Control Test	0.077	11:20
Air Blank	0.000	11:20
Control Test	0.077	11:21
Air Blank	0.000	11:21
Control Test	0.077	11:22
Air Blank	0.000	11:23
Control Test Stats		
Average	0.0770	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/05/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:24
Control Test	0.198	11:25
Air Blank	0.000	11:25
Control Test	0.197	11:26
Air Blank	0.000	11:26
Control Test	0.197	11:27
Air Blank	0.000	11:27
Control Test Stats		
Average	0.1973	
Std Dev	0.0006	
Rel Std Dev(%)	0.2926	

Operator's Signature

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/05/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:29
Control Test	0.083	11:29
Air Blank	0.000	11:29
Control Test	0.083	11:30
Air Blank	0.000	11:30
Control Test	0.083	11:31
Air Blank	0.000	11:31
Control Test Stats		
Average	0.0830	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

DGS

Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: FL HIGHWAY PATROL
Time of Inspection: 15:27

Date of Inspection: 06/06/2023

Serial Number: 80-003407
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:06723080A5 Exp: 04/05/2025
0.000	0.049	0.078	0.197	0.080
0.000	0.049	0.078	0.198	0.080
0.000	0.049	0.078	0.197	0.079
0.000	0.049	0.078	0.197	0.080
0.000	0.049	0.078	0.197	0.080
0.000	0.049	0.078	0.197	0.080
0.000	0.049	0.078	0.197	0.080
0.000	0.049	0.078	0.197	0.080
0.000	0.049	0.078	0.197	0.080
0.000	0.049	0.078	0.197	0.080
Standard Deviations	0.0000	0.0000	0.0003	0.0003

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0001 Number of Simulators Used: 5


Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

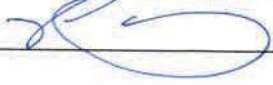 PHIL NICODEMO
Signature and Printed Name

06/06/2023
Date

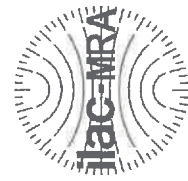

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-003407, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-003407</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>FL HIGHWAY PATROL</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>06/06/2023</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>15:27</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 March 2022

Issuing Authority: Alcohol Testing Program

06/06/2023

Date

PHIL NICODEMO,

Department Inspector

Service • Integrity • Respect • Quality

CALIBRATION ADJUSTMENT

FL HIGHWAY PATROL
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-003407
 06/06/2023 10:56:39

1

PN

Auto Calibration

pg 1 of 2

	<<<<<	3um	>>>>>	<<<<<	9um	>>>>>

Solution =	0.000 g/210L or 0.0000 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	0.0930	(-0.0140)		0.0890	(0.0110)	
Sample #2	0.0720	(0.0700)		0.1290	(0.0290)	
Sample #3	0.0700	(0.1280)		0.1200	(0.0570)	
Sample #4	0.0690	(0.1690)		0.1150	(0.0690)	
Avg % Abs	0.0703	(0.1223)		0.1213	(0.0517)	
STD DEV	0.0015	(0.0497)		0.0071	(0.0205)	
REL STD DEV	2.172	(40.662)		5.847	(39.729)	

Solution =	0.040 g/210L or 0.1905 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	0.7540	(-0.0050)		1.4540	(-0.0090)	
Sample #2	0.7180	(0.0310)		1.4450	(0.0020)	
Sample #3	0.7450	(0.0440)		1.4640	(-0.0070)	
Sample #4	0.7320	(0.0700)		1.4290	(0.0180)	
Avg % Abs	0.7317	(0.0483)		1.4460	(0.0043)	
STD DEV	0.0135	(0.0199)		0.0175	(0.0127)	
REL STD DEV	1.846	(41.085)		1.212	(292.206)	

Solution =	0.100 g/210L or 0.4762 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	1.7470	(0.0000)		3.4280	(-0.0070)	
Sample #2	1.7450	(0.0180)		3.4340	(0.0060)	
Sample #3	1.7620	(0.0250)		3.4510	(-0.0230)	
Sample #4	1.7550	(0.0480)		3.4420	(0.0150)	
Avg % Abs	1.7540	(0.0303)		3.4423	(-0.0007)	
STD DEV	0.0085	(0.0157)		0.0085	(0.0199)	
REL STD DEV	0.487	(51.742)		0.247	(2978.674)	

Solution =	0.200 g/210L or 0.9524 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	3.3560	(0.0000)		6.5750	(-0.0230)	
Sample #2	3.3730	(0.0050)		6.5530	(-0.0050)	
Sample #3	3.3330	(0.0420)		6.5370	(0.0130)	
Sample #4	3.3520	(0.0420)		6.5370	(0.0110)	
Avg % Abs	3.3527	(0.0297)		6.5423	(0.0063)	
STD DEV	0.0200	(0.0214)		0.0092	(0.0099)	
REL STD DEV	0.597	(72.007)		0.141	(155.775)	

Solution =	0.300 g/210L or 1.4286 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)		% Abs	(% Abs Ref)	
Sample #1	4.9560	(-0.0020)		9.5470	(0.0070)	
Sample #2	4.9550	(0.0230)		9.5470	(0.0240)	
Sample #3	4.9320	(0.0460)		9.5410	(0.0390)	
Sample #4	4.9640	(0.0380)		9.5570	(0.0190)	
Avg % Abs	4.9503	(0.0357)		9.5483	(0.0273)	
STD DEV	0.0165	(0.0117)		0.0081	(0.0104)	
REL STD DEV	0.333	(32.737)		0.085	(38.079)	

CALIBRATION ADJUSTMENT

#1

PN

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/06/2023 10:56:39

Auto Calibration

pg 2 of 2

```
<<<<<      3um      >>>>>
-----
Zero Order Coef   -192.23
First Order Coef  2808.95
Second Order Coef  23.75
```

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.000	-0.0001
0.040	0.039	0.0006
0.100	0.101	-0.0010
0.200	0.199	0.0007
0.300	0.300	-0.0002

```
<<<<<      9um      >>>>>
-----
Zero Order Coef   -164.81
First Order Coef  1398.71
Second Order Coef  12.07
```

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	0.000	-0.0001
0.040	0.040	0.0005
0.100	0.101	-0.0007
0.200	0.200	0.0004
0.300	0.300	-0.0001

```
<<<<<      3um      >>>>>
-----
Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1
Sample
Sample #1          3346.00          3266.00
Sample #2          3321.00          3318.00
Sample #3          3353.00          3339.00
Sample #4          3333.00          3310.00
Avg                3335.6667        3322.3333
STD DEV            16.1658          14.9778
REL STD DEV         0.485           0.451
H2O adjust (mg/l*10k) 474           487
```

Barometric Pressure = 1012

*****CALIBRATION SUCCESSFUL*****

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000
06/06/2023

SN 80-003407
12:47:59

CALIBRATION ADJUSTMENT

2

PN

Auto Calibration

pg 1 of 2

<<<< 3um >>>> <<<< 9um >>>>

Solution = 0.000 g/210L or 0.0000 mg/l, Samples = 4, Discarded = 1
Sample % Abs (% Abs Ref) % Abs (% Abs Ref)
Sample #1 0.0760 (-0.0060) 0.0860 (0.0140)
Sample #2 0.0740 (0.0120) 0.0990 (0.0030)
Sample #3 0.0500 (0.0510) 0.0890 (0.0250)
Sample #4 0.0410 (0.0600) 0.0860 (0.0280)
Avg % Abs 0.0550 (0.0410) 0.0913 (0.0187)
STD DEV 0.0171 (0.0255) 0.0068 (0.0137)
REL STD DEV 31.016 (62.231) 7.453 (73.127)

Solution = 0.040 g/210L or 0.1905 mg/l, Samples = 4, Discarded = 1
Sample % Abs (% Abs Ref) % Abs (% Abs Ref)
Sample #1 0.7600 (-0.0240) 1.4550 (0.0070)
Sample #2 0.7280 (0.0080) 1.4500 (-0.0020)
Sample #3 0.7440 (0.0080) 1.4480 (0.0210)
Sample #4 0.7260 (0.0260) 1.4380 (0.0430)
Avg % Abs 0.7327 (0.0140) 1.4453 (0.0207)
STD DEV 0.0099 (0.0104) 0.0064 (0.0225)
REL STD DEV 1.347 (74.231) 0.445 (108.880)

Solution = 0.100 g/210L or 0.4762 mg/l, Samples = 4, Discarded = 1
Sample % Abs (% Abs Ref) % Abs (% Abs Ref)
Sample #1 1.7760 (-0.0040) 3.4380 (-0.0040)
Sample #2 1.7260 (0.0240) 3.4330 (0.0250)
Sample #3 1.7500 (0.0220) 3.4570 (0.0200)
Sample #4 1.7170 (0.0390) 3.4220 (0.0270)
Avg % Abs 1.7310 (0.0283) 3.4373 (0.0240)
STD DEV 0.0171 (0.0093) 0.0179 (0.0036)
REL STD DEV 0.985 (32.794) 0.521 (15.023)

Solution = 0.200 g/210L or 0.9524 mg/l, Samples = 4, Discarded = 1
Sample % Abs (% Abs Ref) % Abs (% Abs Ref)
Sample #1 3.3690 (-0.0050) 6.5840 (0.0000)
Sample #2 3.3690 (0.0110) 6.5670 (0.0170)
Sample #3 3.3660 (0.0220) 6.5470 (0.0230)
Sample #4 3.3330 (0.0260) 6.5490 (0.0230)
Avg % Abs 3.3560 (0.0197) 6.5543 (0.0210)
STD DEV 0.0200 (0.0078) 0.0110 (0.0035)
REL STD DEV 0.595 (39.496) 0.168 (16.496)

Solution = 0.300 g/210L or 1.4286 mg/l, Samples = 4, Discarded = 1
Sample % Abs (% Abs Ref) % Abs (% Abs Ref)
Sample #1 4.9650 (0.0000) 9.5400 (0.0090)
Sample #2 4.9320 (0.0270) 9.5350 (0.0350)
Sample #3 4.9610 (0.0200) 9.5760 (0.0210)
Sample #4 4.9600 (0.0260) 9.5640 (0.0230)
Avg % Abs 4.9510 (0.0243) 9.5583 (0.0263)
STD DEV 0.0165 (0.0038) 0.0211 (0.0076)
REL STD DEV 0.332 (15.559) 0.221 (28.754)

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/06/2023 12:47:59

CALIBRATION ADJUSTMENT

#2 PN

Auto Calibration

pg 2 of 2

<<<< 3um >>>>

Zero Order Coef -159.77
First Order Coef 2808.18
Second Order Coef 22.24

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	-0.000	0.0001
0.040	0.040	-0.0001
0.100	0.100	-0.0001
0.200	0.200	0.0002
0.300	0.300	-0.0001

<<<< 9um >>>>

-129.46
1384.69
12.98

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	-0.000	0.0001
0.040	0.040	0.0001
0.100	0.100	-0.0005
0.200	0.200	0.0004
0.300	0.300	-0.0001

<<<< 3um >>>>

Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1
Sample

Sample #1	3357.00
Sample #2	3354.00
Sample #3	3382.00
Sample #4	3400.00
Avg	3378.6667
STD DEV	23.1805
REL STD DEV	0.686
H2O adjust (mg/l*10k)	431

<<<< 9um >>>>

3320.00
3325.00
3308.00
3314.00
3315.6667
8.6217
0.260
494

Barometric Pressure = 1010

*****CALIBRATION SUCCESSFUL*****

Post stability, check (#1)

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/06/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:17
Control Test	0.047	12:17
Air Blank	0.000	12:18
Control Test	0.048	12:18
Air Blank	0.000	12:19
Control Test	0.048	12:20
Air Blank	0.000	12:20
Control Test Stats		
Average	0.0477	
Std Dev	0.0006	
Rel Std Dev(%)	1.2112	

Operator's Signature

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/06/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:23
Control Test	0.077	12:24
Air Blank	0.000	12:25
Control Test	0.077	12:25
Air Blank	0.000	12:26
Control Test	0.077	12:27
Air Blank	0.000	12:27
Control Test Stats		
Average	0.0770	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/06/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:28
Control Test	0.197	12:29
Air Blank	0.000	12:29
Control Test	0.197	12:30
Air Blank	0.000	12:30
Control Test	0.197	12:31
Air Blank	0.000	12:32
Control Test Stats		
Average	0.1970	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

FL HIGHWAY PATROL
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-003407
06/06/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:33
Control Test	0.081	12:34
Air Blank	0.000	12:34
Control Test	0.080	12:34
Air Blank	0.000	12:35
Control Test	0.081	12:35
Air Blank	0.000	12:36
Control Test Stats		
Average	0.0807	
Std Dev	0.0006	
Rel Std Dev(%)	0.7157	


Operator's Signature

DGS

Post stability checks (#2)

FL HIGHWAY PATROL
Intoxilyzer - Alcotest Analyzer
Model 8000 SN 80-003407
06/06/2023
Software: 8100.27


Test	g/210L	Time
Air Blank	0.000	13:41
Control Test	0.048	13:42
Air Blank	0.000	13:42
Control Test	0.049	13:43
Air Blank	0.000	13:44
Control Test	0.049	13:44
Air Blank	0.000	13:45
Control Test Stats		
Average	0.0487	
Std Dev	0.0006	
Rel Std Dev(%)	1.1863	

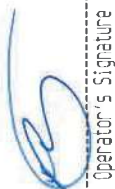

Operator's Signature

FL HIGHWAY PATROL
Intoxilyzer - Alcotest Analyzer
Model 8000 SN 80-003407
06/06/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	13:46
Control Test	0.078	13:47
Air Blank	0.000	13:47
Control Test	0.078	13:48
Air Blank	0.000	13:48
Control Test	0.078	13:49
Air Blank	0.000	13:50
Control Test Stats		
Average	0.0780	
Std Dev	0.0003	
Rel Std Dev(%)	0.0000	

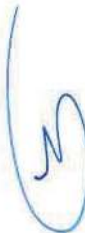

Operator's Signature

FL HIGHWAY PATROL
Intoxilyzer - Alcotest Analyzer
Model 8000 SN 80-003407
06/06/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	13:51
Control Test	0.198	13:51
Air Blank	0.000	13:52
Control Test	0.197	13:53
Air Blank	0.000	13:53
Control Test	0.197	13:54
Air Blank	0.000	13:54
Control Test Stats		
Average	0.1973	
Std Dev	0.0006	
Rel Std Dev(%)	0.2926	



Operator's Signature

FL HIGHWAY PATROL
Intoxilyzer - Alcotest Analyzer
Model 8000 SN 80-003407
06/06/2023
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	13:56
Control Test	0.080	13:56
Air Blank	0.000	13:57
Control Test	0.081	13:57
Air Blank	0.000	13:57
Control Test	0.080	13:58
Air Blank	0.000	13:58
Control Test Stats		
Average	0.0803	
Std Dev	0.0006	
Rel Std Dev(%)	0.7187	

DGS



Operator's Signature