







INSTRUMENT PROCESSING SHEET

Agency Hendry CSOS/N 80-006386Florida Department of
Law EnforcementDate In 10/11/2023 DI Completion Date 10/13/2023☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By TDG	Quality Checks	By TDG	Date <u>10/12/2023</u>	Flow Calibration	By	Date																																																											
<input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: Missing back two feet on DGS shelf.		<input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>192</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.167</u> (.139 - .169) 36 mm <u>0.183</u> (.156 - .190) 53 mm <u>0.257</u> (.228 - .278) 103 mm <u>0.523</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>26932</u> <input checked="" type="checkbox"/> Stability Checks			Flow Column # _____ <input type="checkbox"/> 5L/min - 17mm <input type="checkbox"/> 15L/min - 53mm <input type="checkbox"/> 30L/min - 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																													
				<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5094</td><td>202201C 01/11/2024</td></tr><tr><td>0.080</td><td>MP5095</td><td>202201D 01/18/2024</td></tr><tr><td>0.200</td><td>MP5096</td><td>202201E 01/18/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG223802 08/26/2024</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP5094	202201C 01/11/2024	0.080	MP5095	202201D 01/18/2024	0.200	MP5096	202201E 01/18/2024	0.080 DGS	N/A	AG223802 08/26/2024	Maintenance By _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____																																														
Simulator	Serial #	Lot #/Exp																																																																
0.050	MP5094	202201C 01/11/2024																																																																
0.080	MP5095	202201D 01/18/2024																																																																
0.200	MP5096	202201E 01/18/2024																																																																
0.080 DGS	N/A	AG223802 08/26/2024																																																																
Calibration Adjustment By TDG _____ Barometric Pressure Gauge <u>1010</u> ID # <u>28199</u> <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td>MP5097</td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td>MP5098</td><td>22460</td><td>12/28/2024</td></tr><tr><td>0.100</td><td>MP5099</td><td>22310</td><td>08/11/2024</td></tr><tr><td>0.200</td><td>MP5100</td><td>22050</td><td>02/07/2024</td></tr><tr><td>0.300</td><td>MP5101</td><td>22220</td><td>06/15/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG22220</td><td>08/10/2024</td></tr></tbody></table> <input checked="" type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td>MP5094</td><td>202201C</td><td>01/11/2024</td></tr><tr><td>0.080</td><td>MP5095</td><td>202201D</td><td>01/18/2024</td></tr><tr><td>0.200</td><td>MP5096</td><td>202201E</td><td>01/18/2024</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG223802</td><td>08/26/2024</td></tr></tbody></table>				Simulator	Serial #	Lot #	Expiration	0.000	MP5097	N/A	N/A	0.040	MP5098	22460	12/28/2024	0.100	MP5099	22310	08/11/2024	0.200	MP5100	22050	02/07/2024	0.300	MP5101	22220	06/15/2024	0.080 DGS	N/A	AG22220	08/10/2024	Simulator	Serial #	Lot #	Expiration	0.050	MP5094	202201C	01/11/2024	0.080	MP5095	202201D	01/18/2024	0.200	MP5096	202201E	01/18/2024	0.080 DGS	N/A	AG223802	08/26/2024	Department Inspection By TDG _____ Barometric Pressure ID# <u>26932</u> Gauge <u>1010</u> Instrument <u>1011</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2022-B</u> <table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP4863</td></tr><tr><td>Interferent</td><td>MP5093</td></tr><tr><td>0.050</td><td>MP5094</td></tr><tr><td>0.080</td><td>MP5095</td></tr><tr><td>0.200</td><td>MP5096</td></tr></tbody></table> Attachments <input checked="" type="checkbox"/> Form 41 <input checked="" type="checkbox"/> Post-Stability Checks <input checked="" type="checkbox"/> Stability Checks <input type="checkbox"/> Flow Calibration <input checked="" type="checkbox"/> Calibration Certificate <input type="checkbox"/> Form 40 <input checked="" type="checkbox"/> Calibration Adjustment <input type="checkbox"/> Other _____ <input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use Israel Soto <small>Digitally signed by Israel Soto Date: 2023.10.13 14:58:16 +04'00'</small> Phil Nicodemo <small>Digitally signed by Phil Nicodemo Date: 2023.10.13 15:14:20 -04'00'</small> Tech Review / Date Admin Review / Date			Simulator	Serial Number	0.000	MP4863	Interferent	MP5093	0.050	MP5094	0.080	MP5095	0.200	MP5096
Simulator	Serial #	Lot #	Expiration																																																															
0.000	MP5097	N/A	N/A																																																															
0.040	MP5098	22460	12/28/2024																																																															
0.100	MP5099	22310	08/11/2024																																																															
0.200	MP5100	22050	02/07/2024																																																															
0.300	MP5101	22220	06/15/2024																																																															
0.080 DGS	N/A	AG22220	08/10/2024																																																															
Simulator	Serial #	Lot #	Expiration																																																															
0.050	MP5094	202201C	01/11/2024																																																															
0.080	MP5095	202201D	01/18/2024																																																															
0.200	MP5096	202201E	01/18/2024																																																															
0.080 DGS	N/A	AG223802	08/26/2024																																																															
Simulator	Serial Number																																																																	
0.000	MP4863																																																																	
Interferent	MP5093																																																																	
0.050	MP5094																																																																	
0.080	MP5095																																																																	
0.200	MP5096																																																																	
Notes/Suggested Service: _____ _____ _____ _____ _____ _____																																																																		

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
<p>0.047 to 0.053 ✓</p> <p>HENDRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 10/12/2023 Software: 8100.27</p> <p>Test g/210L Time</p> <p>Air Blank 0.000 11:38 Control Test 0.049 11:38 Air Blank 0.000 11:39 Control Test 0.048 11:40 Air Blank 0.000 11:40 Control Test 0.048 11:41 Air Blank 0.000 11:42</p> <p>Control Test Stats Average 0.0483 Std Dev 0.0006 Rel Std Dev(%) 1.1945</p> <p>Operator's Signature </p>	<p>0.077 to 0.083 ✗</p> <p>HENDRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 10/12/2023 Software: 8100.27</p> <p>Test g/210L Time</p> <p>Air Blank 0.000 11:46 Control Test 0.078 11:46 Air Blank 0.000 11:47 Control Test 0.078 11:48 Air Blank 0.000 11:48 Control Test 0.076 11:49 Air Blank 0.000 11:49</p> <p>Control Test Stats Average 0.0773 Std Dev 0.0012 Rel Std Dev(%) 1.4932</p> <p>Operator's Signature </p>	<p>0.194 to 0.206 ✓</p> <p>HENDRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 10/12/2023 Software: 8100.27</p> <p>Test g/210L Time</p> <p>Air Blank 0.000 11:52 Control Test 0.198 11:52 Air Blank 0.000 11:53 Control Test 0.199 11:54 Air Blank 0.000 11:54 Control Test 0.199 11:55 Air Blank 0.000 11:56</p> <p>Control Test Stats Average 0.1987 Std Dev 0.0006 Rel Std Dev(%) 0.2906</p> <p>Operator's Signature </p>	<p>0.077 to 0.083 ✓</p> <p>0.003 of Wet ✓</p> <p>HENDRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 10/12/2023 Software: 8100.27</p> <p>Test g/210L Time</p> <p>Air Blank 0.000 11:58 Control Test 0.080 11:58 Air Blank 0.000 11:59 Control Test 0.080 11:59 Air Blank 0.000 12:00 Control Test 0.080 12:00 Air Blank 0.000 12:01</p> <p>Control Test Stats Average 0.0800 Std Dev 0.0000 Rel Std Dev(%) 0.0000</p> <p>Operator's Signature </p>

Q65

HENRY COUNTY SO
Intoxilizer - Alconol Analyzer
Model 8000
10/13/2023
SN 80-006386
09:09:48

<<<< CHANNEL 2 >>>>
Sample % Abs (% Abs Ref)
Sample #1 = 1.5850 (-0.0130)
Sample #2 = 1.5190 (-0.0300)
Sample #3 = 1.5970 (-0.0100)
Sample #4 = 1.5930 (-0.0060)
Avg % Abs = 1.6030 (-0.0153)
STD DEV = 0.0140 (0.0129)
REL STD DEV = 0.873 (83.658)

***** AUTO CAL DATA *****
<<<< CHANNEL 1 >>>>
Sol Val = 0.0000 mg/l or 0.000 g/210L
% Abs = 0.063
Std Dev = 0.03 Rel Std Dev = 48.93
Sol Val = 0.1905 mg/l or 0.040 g/210L
% Abs = 0.764
Std Dev = 0.02 Rel Std Dev = 3.07
Sol Val = 0.4762 mg/l or 0.100 g/210L
% Abs = 1.807
Std Dev = 0.03 Rel Std Dev = 1.70
Sol Val = 0.9524 mg/l or 0.200 g/210L
% Abs = 3.512
Std Dev = 0.03 Rel Std Dev = 0.89
Sol Val = 1.4286 mg/l or 0.300 g/210L
% Abs = 5.168
Std Dev = 0.02 Rel Std Dev = 0.43
Zero Order Coef = -164.34
First Order Coef = 2685.54
Second Order Coef = 21.29
Standard Deviation = 4.778313

Solution Stats Quadratic Fit Chan 2
Act Fit Residual
g/210L g/210L g/210L
0.000 -0.000 0.0002
0.040 -0.040 -0.0005
0.100 0.100 0.0005
0.200 0.200 -0.0002
0.300 0.300 0.0001

Auto Calibration
Max Power Res Value = 45
Auto Range Res Value = 28

Sol Value = 0.000 g/210L ***
Fit value = 0.0000 mg/l %%%
Samples Taken = 4, Discarded = 1
Sum Io = 12653, Sum Io = 13285

<<<< CHANNEL 1 >>>>
Sample % Abs (% Abs Ref)
Sample #1 = 1.7650 (-0.0010)
Sample #2 = 1.8190 (-0.0180)
Sample #3 = 1.7720 (0.0140)
Sample #4 = 1.8300 (0.0110)
Avg % Abs = 1.8070 (0.0023)
STD DEV = 0.0308 (0.0177)
REL STD DEV = 1.705 (757.412)

Sol Value = 0.300 g/210L ***
Fit value = 1.4286 mg/l %%%
Samples Taken = 4, Discarded = 1
Sum Io = 12629, Sum Io = 13273

<<<< CHANNEL 1 >>>>
Sample % Abs (% Abs Ref)
Sample #1 = 5.1540 (-0.0140)
Sample #2 = 5.1460 (0.0000)
Sample #3 = 5.1900 (-0.0040)
Sample #4 = 5.1680 (0.0260)
Avg % Abs = 5.1680 (0.0073)
STD DEV = 0.0220 (0.0163)
REL STD DEV = 0.426 (222.124)

Sol Value = 0.080 g/210L ***
Fit value = 0.3810 mg/l %%%
Samples Taken = 4, Discarded = 1
***** CHANNEL 1 *****
Sample #1 = 3469.00
Sample #2 = 3467.00
Sample #3 = 3557.00
Sample #4 = 3527.00
Average Result = 3517.0000
STD DEV = 45.8258
REL STD DEV = 1.303

***** CHANNEL 2 *****
Sample #1 = 3371.00
Sample #2 = 3309.00
Sample #3 = 3356.00
Sample #4 = 3356.00
Average Result = 3340.3333
STD DEV = 27.1355
REL STD DEV = 0.812

Dry Gas H2O Adjust Results *****
Barometric Pressure = 1011
3 um H2O Adjust (mg/l*10.000) = 292
9 um H2O Adjust (mg/l*10.000) = 469
**** AUTO CAL PASS

Solution Stats Quadratic Fit Chan 1
Act Fit Residual
g/210L g/210L g/210L
0.000 -0.000 -0.0001
0.040 0.040 0.0001
0.100 0.100 0.0001
0.200 0.200 -0.0001
0.300 0.300 0.0001

<<<< CHANNEL 2 >>>>
Sample % Abs (% Abs Ref)
Sample #1 = 3.6260 (-0.0090)
Sample #2 = 3.6390 (0.0010)
Sample #3 = 3.6090 (0.0060)
Sample #4 = 3.6480 (0.0140)
Avg % Abs = 3.6320 (0.0070)
STD DEV = 0.0204 (0.0066)
REL STD DEV = 0.562 (93.678)

Sol Value = 0.200 g/210L ***
Fit value = 0.9524 mg/l %%%
Samples Taken = 4, Discarded = 1
Sum Io = 12641, Sum Io = 13279

<<<< CHANNEL 1 >>>>
Sample % Abs (% Abs Ref)
Sample #1 = 3.5340 (-0.0140)
Sample #2 = 3.5280 (-0.0080)
Sample #3 = 3.4760 (0.0300)
Sample #4 = 3.5320 (0.0190)
Avg % Abs = 3.5120 (0.0137)
STD DEV = 0.0312 (0.0196)
REL STD DEV = 0.890 (143.073)

Sol Value = 0.040 g/210L ***
Fit value = 0.1905 mg/l %%%
Samples Taken = 4, Discarded = 1
Sum Io = 12670, Sum Io = 13291

<<<< CHANNEL 1 >>>>
Sample % Abs (% Abs Ref)
Sample #1 = 0.7590 (-0.0100)
Sample #2 = 0.7910 (-0.0180)
Sample #3 = 0.7550 (0.0180)
Sample #4 = 0.7470 (0.0420)
Avg % Abs = 0.7643 (0.0140)
STD DEV = 0.0234 (0.0302)
REL STD DEV = 3.066 (215.710)

Solution Stats Quadratic Fit Chan 2
Act Fit Residual
g/210L g/210L g/210L
0.000 -0.000 -0.0001
0.040 0.040 0.0001
0.100 0.100 0.0001
0.200 0.200 -0.0001
0.300 0.300 0.0001

Optical Calibration Adjustment

By: TDG

Post-Cal Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																																
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083																																																																																																																																																
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>																																																																																																																																																
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>																																																																																																																																																
<p>HENDRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 10/13/2023 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>10:36</td></tr> <tr><td>Control Test</td><td>0.049</td><td>10:37</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:37</td></tr> <tr><td>Control Test</td><td>0.049</td><td>10:38</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:38</td></tr> <tr><td>Control Test</td><td>0.049</td><td>10:39</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:40</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0490</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	10:36	Control Test	0.049	10:37	Air Blank	0.000	10:37	Control Test	0.049	10:38	Air Blank	0.000	10:38	Control Test	0.049	10:39	Air Blank	0.000	10:40	Control Test Stats			Average	0.0490		Std Dev	0.0000		Rel Std Dev(%)	0.0000		<p>HENDRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 10/13/2023 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>10:42</td></tr> <tr><td>Control Test</td><td>0.080</td><td>10:43</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:44</td></tr> <tr><td>Control Test</td><td>0.079</td><td>10:44</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:45</td></tr> <tr><td>Control Test</td><td>0.078</td><td>10:46</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:46</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0790</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.2658</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	10:42	Control Test	0.080	10:43	Air Blank	0.000	10:44	Control Test	0.079	10:44	Air Blank	0.000	10:45	Control Test	0.078	10:46	Air Blank	0.000	10:46	Control Test Stats			Average	0.0790		Std Dev	0.0010		Rel Std Dev(%)	1.2658		<p>HENDRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 10/13/2023 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>10:55</td></tr> <tr><td>Control Test</td><td>0.199</td><td>10:56</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:56</td></tr> <tr><td>Control Test</td><td>0.199</td><td>10:57</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:58</td></tr> <tr><td>Control Test</td><td>0.199</td><td>10:58</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:59</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1990</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	10:55	Control Test	0.199	10:56	Air Blank	0.000	10:56	Control Test	0.199	10:57	Air Blank	0.000	10:58	Control Test	0.199	10:58	Air Blank	0.000	10:59	Control Test Stats			Average	0.1990		Std Dev	0.0000		Rel Std Dev(%)	0.0000		<p>HENDRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 10/13/2023 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>10:27</td></tr> <tr><td>Control Test</td><td>0.078</td><td>10:28</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:28</td></tr> <tr><td>Control Test</td><td>0.078</td><td>10:29</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:29</td></tr> <tr><td>Control Test</td><td>0.078</td><td>10:29</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:30</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0780</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	10:27	Control Test	0.078	10:28	Air Blank	0.000	10:28	Control Test	0.078	10:29	Air Blank	0.000	10:29	Control Test	0.078	10:29	Air Blank	0.000	10:30	Control Test Stats			Average	0.0780		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:36																																																																																																																																																	
Control Test	0.049	10:37																																																																																																																																																	
Air Blank	0.000	10:37																																																																																																																																																	
Control Test	0.049	10:38																																																																																																																																																	
Air Blank	0.000	10:38																																																																																																																																																	
Control Test	0.049	10:39																																																																																																																																																	
Air Blank	0.000	10:40																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0490																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:42																																																																																																																																																	
Control Test	0.080	10:43																																																																																																																																																	
Air Blank	0.000	10:44																																																																																																																																																	
Control Test	0.079	10:44																																																																																																																																																	
Air Blank	0.000	10:45																																																																																																																																																	
Control Test	0.078	10:46																																																																																																																																																	
Air Blank	0.000	10:46																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0790																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	1.2658																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:55																																																																																																																																																	
Control Test	0.199	10:56																																																																																																																																																	
Air Blank	0.000	10:56																																																																																																																																																	
Control Test	0.199	10:57																																																																																																																																																	
Air Blank	0.000	10:58																																																																																																																																																	
Control Test	0.199	10:58																																																																																																																																																	
Air Blank	0.000	10:59																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1990																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:27																																																																																																																																																	
Control Test	0.078	10:28																																																																																																																																																	
Air Blank	0.000	10:28																																																																																																																																																	
Control Test	0.078	10:29																																																																																																																																																	
Air Blank	0.000	10:29																																																																																																																																																	
Control Test	0.078	10:29																																																																																																																																																	
Air Blank	0.000	10:30																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0780																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
<p>Operator's Signature</p>	<p>Operator's Signature</p>	<p>Operator's Signature</p>	<p>Operator's Signature</p>																																																																																																																																																

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: HENDRY COUNTY SO
Time of Inspection: 13:34

Date of Inspection: 10/13/2023

Serial Number: 80-006386
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test* (g/210L) Lot#:AG223802 Exp: 08/26/2024
0.000	0.049	0.078	0.199	0.078
0.000	0.049	0.078	0.199	0.078
0.000	0.049	0.077	0.199	0.078
0.000	0.049	0.078	0.199	0.078
0.000	0.049	0.078	0.199	0.078
0.000	0.049	0.078	0.199	0.078
0.000	0.049	0.078	0.199	0.078
0.000	0.049	0.078	0.200	0.078
0.000	0.049	0.078	0.199	0.077
0.000	0.049	0.078	0.198	0.077

Standard Deviations	0.0000	0.0003	0.0004	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



TAYLOR D GUTSCHOW

Signature and Printed Name

10/13/2023
Date

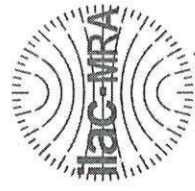

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006386, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-006386</u>	UNCERTAINTY* \pm
Owning Agency:	<u>HENDRY COUNTY SO</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>10/13/2023</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>13:34</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).


The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

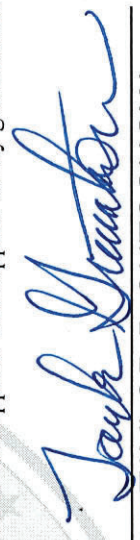

10/13/2023

Date
TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality