

Stability Checks

0.05g/210L 0.047 to 0.053	0.08g/210L 0.077 to 0.083	0.20g/210L 0.194 to 0.206	DGS 0.08g/210L 0.077 to 0.083																																																																																																																																				
COLLIER COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/06/2023 Software: 8100.27 SN 80-007077 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>11:32</td></tr> <tr><td>Control Test</td><td>0.051</td><td>11:32</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:33</td></tr> <tr><td>Control Test</td><td>0.050</td><td>11:34</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:34</td></tr> <tr><td>Control Test</td><td>0.050</td><td>11:35</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:35</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0503</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.1471</td><td></td></tr> </table>	Air Blank	0.000	11:32	Control Test	0.051	11:32	Air Blank	0.000	11:33	Control Test	0.050	11:34	Air Blank	0.000	11:34	Control Test	0.050	11:35	Air Blank	0.000	11:35	Control Test Stats			Average	0.0503		Std Dev	0.0006		Rel Std Dev(%)	1.1471		COLLIER COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/06/2023 Software: 8100.27 SN 80-007077 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>11:40</td></tr> <tr><td>Control Test</td><td>0.079</td><td>11:40</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:41</td></tr> <tr><td>Control Test</td><td>0.078</td><td>11:41</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:42</td></tr> <tr><td>Control Test</td><td>0.078</td><td>11:43</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:43</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0783</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7370</td><td></td></tr> </table>	Air Blank	0.000	11:40	Control Test	0.079	11:40	Air Blank	0.000	11:41	Control Test	0.078	11:41	Air Blank	0.000	11:42	Control Test	0.078	11:43	Air Blank	0.000	11:43	Control Test Stats			Average	0.0783		Std Dev	0.0006		Rel Std Dev(%)	0.7370		COLLIER COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/06/2023 Software: 8100.27 SN 80-007077 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>11:46</td></tr> <tr><td>Control Test</td><td>0.199</td><td>11:47</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:47</td></tr> <tr><td>Control Test</td><td>0.200</td><td>11:48</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:49</td></tr> <tr><td>Control Test</td><td>0.198</td><td>11:49</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:50</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1990</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.5025</td><td></td></tr> </table>	Air Blank	0.000	11:46	Control Test	0.199	11:47	Air Blank	0.000	11:47	Control Test	0.200	11:48	Air Blank	0.000	11:49	Control Test	0.198	11:49	Air Blank	0.000	11:50	Control Test Stats			Average	0.1990		Std Dev	0.0010		Rel Std Dev(%)	0.5025		COLLIER COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/06/2023 Software: 8100.27 SN 80-007077 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>11:22</td></tr> <tr><td>Control Test</td><td>0.078</td><td>11:22</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:23</td></tr> <tr><td>Control Test</td><td>0.078</td><td>11:23</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:24</td></tr> <tr><td>Control Test</td><td>0.078</td><td>11:24</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:25</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0780</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </table>	Air Blank	0.000	11:22	Control Test	0.078	11:22	Air Blank	0.000	11:23	Control Test	0.078	11:23	Air Blank	0.000	11:24	Control Test	0.078	11:24	Air Blank	0.000	11:25	Control Test Stats			Average	0.0780		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Air Blank	0.000	11:32																																																																																																																																					
Control Test	0.051	11:32																																																																																																																																					
Air Blank	0.000	11:33																																																																																																																																					
Control Test	0.050	11:34																																																																																																																																					
Air Blank	0.000	11:34																																																																																																																																					
Control Test	0.050	11:35																																																																																																																																					
Air Blank	0.000	11:35																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0503																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	1.1471																																																																																																																																						
Air Blank	0.000	11:40																																																																																																																																					
Control Test	0.079	11:40																																																																																																																																					
Air Blank	0.000	11:41																																																																																																																																					
Control Test	0.078	11:41																																																																																																																																					
Air Blank	0.000	11:42																																																																																																																																					
Control Test	0.078	11:43																																																																																																																																					
Air Blank	0.000	11:43																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0783																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.7370																																																																																																																																						
Air Blank	0.000	11:46																																																																																																																																					
Control Test	0.199	11:47																																																																																																																																					
Air Blank	0.000	11:47																																																																																																																																					
Control Test	0.200	11:48																																																																																																																																					
Air Blank	0.000	11:49																																																																																																																																					
Control Test	0.198	11:49																																																																																																																																					
Air Blank	0.000	11:50																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.1990																																																																																																																																						
Std Dev	0.0010																																																																																																																																						
Rel Std Dev(%)	0.5025																																																																																																																																						
Air Blank	0.000	11:22																																																																																																																																					
Control Test	0.078	11:22																																																																																																																																					
Air Blank	0.000	11:23																																																																																																																																					
Control Test	0.078	11:23																																																																																																																																					
Air Blank	0.000	11:24																																																																																																																																					
Control Test	0.078	11:24																																																																																																																																					
Air Blank	0.000	11:25																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0780																																																																																																																																						
Std Dev	0.0000																																																																																																																																						
Rel Std Dev(%)	0.0000																																																																																																																																						
Operator's Signature 	Operator's Signature 	Operator's Signature 	Operator's Signature 																																																																																																																																				

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: COLLIER COUNTY SO
Time of Inspection: 14:06

Date of Inspection: 09/06/2023

Serial Number: 80-007077
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202201C Exp: 01/11/2024	0.08g/210L Test (g/210L) Lot#:202201D Exp: 01/18/2024	0.20g/210L Test (g/210L) Lot#:202201E Exp: 01/18/2024	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG223802 Exp: 08/26/2024
0.000	0.050	0.079	0.200	0.078
0.000	0.050	0.079	0.200	0.078
0.000	0.050	0.079	0.199	0.078
0.000	0.051	0.079	0.200	0.078
0.000	0.051	0.079	0.200	0.078
0.000	0.050	0.080	0.199	0.078
0.000	0.050	0.079	0.200	0.078
0.000	0.051	0.079	0.200	0.078
0.000	0.051	0.079	0.200	0.077
0.000	0.051	0.079	0.199	0.078

Standard Deviations	0.0005	0.0003	0.0004	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5


Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

 TAYLOR D GUTSCHOW
Signature and Printed Name

09/06/2023
Date

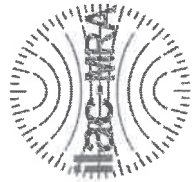

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-007077, manufactured by CML, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-007077</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>COLLIER COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>09/06/2023</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>14:06</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

09/06/2023

Date

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality