



INSTRUMENT PROCESSING SHEET

Agency Miami PDS/N 80-000873Florida Department of
Law EnforcementDate In 07/15/2024 DI Completion Date 07/22/2024☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>TDG</u> Date <u>07/15/2024</u>		Quality Checks By <u>TDG</u> Date <u>07/22/2024</u>		Flow Calibration By <u>TDG</u> Date <u>07/22/2024</u>																																									
☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>Missing back right foot.</u>		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>247</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.156</u> (.156 - .190) 53 mm <u>0.226*</u> (.228 - .278) 103 mm <u>0.472</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks		Flow Column # <u>ATP106</u> ☒ 5L/min - 17mm ☒ 15L/min - 53mm ☒ 30L/min - 103mm ☒ R-Value <u>247</u> ☒ Post Calibration Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.503</u> (.447 - .547)																																									
		<table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP6287</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP6288</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>01923080A3 02/05/2025</td></tr></tbody></table>		Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K 03/29/2025	0.080	MP6287	202303L 03/29/2025	0.200	MP6288	202304C 04/05/2025	0.080 DGS	N/A	01923080A3 02/05/2025	Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____																										
Simulator	Serial #	Lot #/Exp																																											
0.050	MP6286	202303K 03/29/2025																																											
0.080	MP6287	202303L 03/29/2025																																											
0.200	MP6288	202304C 04/05/2025																																											
0.080 DGS	N/A	01923080A3 02/05/2025																																											
Calibration Adjustment By _____				Department Inspection By <u>TDG</u>																																									
Barometric Pressure Gauge _____ ID # _____				Barometric Pressure ID# <u>26932</u> Gauge <u>1017</u> Instrument <u>1015</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2023-B</u>																																									
<table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>				Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			<table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table>		Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288
Simulator	Serial #	Lot #	Expiration																																										
0.000		N/A	N/A																																										
0.040																																													
0.100																																													
0.200																																													
0.300																																													
0.080 DGS	N/A																																												
Simulator	Serial Number																																												
0.000	MP6284																																												
Interferent	MP6285																																												
0.050	MP6286																																												
0.080	MP6287																																												
0.200	MP6288																																												
☐ Post Calibration Adjustment Stability Checks				Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment																																									
<table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>				Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			☐ Post-Stability Checks ☒ Flow Calibration ☐ Form 40 ☐ Other _____																					
Simulator	Serial #	Lot #	Expiration																																										
0.050																																													
0.080																																													
0.200																																													
0.080 DGS	N/A																																												
Notes/Suggested Service: <u>*Outside nominal range. (TDG)</u>				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																									
				Phil Nicodemo Digitally signed by Phil Nicodemo Date: 2024.07.30 12:36:00 -0400																																									
				Shayla Platt Digitally signed by Shayla Platt Date: 2024.07.31 15:26:12 -0400																																									
				Tech Review / Date _____ Admin Review / Date _____																																									

Flow Calibration Adjustment(s)

Performed by TDG

MIAMI PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000873
07/22/2024
Software: 8100.27

Flow Rate Calibration*****

1: Rate (Liters/min) = 5

SQRT(Diff) = 6.926

2: Rate (Liters/min) = 15

SQRT(Diff) = 11.789

3: Rate (Liters/min) = 30

SQRT(Diff) = 21.094

Dependent Data Scale Factor = 100000 L/min

Independent Data Scale Factor = 256

Rounded Slope = 681

Rounded Intercept = -6479.14

Correlation = 0.99795

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
✓	✓	✓	✓
MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 07/22/2024 Software: 8100.27 SN 80-000873 Test g/210L Time Air Blank 0.000 11:04 Control Test 0.052 11:04 Air Blank 0.000 11:05 Control Test 0.051 11:06 Air Blank 0.000 11:06 Control Test 0.050 11:07 Air Blank 0.000 11:08 Control Test Stats Average 0.0510 Std Dev 0.0010 Rel Std Dev(%) 1.9508	MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 07/22/2024 Software: 8100.27 SN 80-000873 Test g/210L Time Air Blank 0.000 11:10 Control Test 0.081 11:11 Air Blank 0.000 11:12 Control Test 0.081 11:12 Air Blank 0.000 11:13 Control Test 0.080 11:14 Air Blank 0.000 11:14 Control Test Stats Average 0.0807 Std Dev 0.0006 Rel Std Dev(%) 0.7157	MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 07/22/2024 Software: 8100.27 SN 80-000873 Test g/210L Time Air Blank 0.000 11:17 Control Test 0.200 11:18 Air Blank 0.000 11:18 Control Test 0.199 11:19 Air Blank 0.000 11:20 Control Test 0.199 11:20 Air Blank 0.000 11:21 Control Test Stats Average 0.1993 Std Dev 0.0006 Rel Std Dev(%) 0.2896	MIAMI PD Intoxilyzer - Alcohol Analyzer Model 8000 07/22/2024 Software: 8100.27 SN 80-000873 Test g/210L Time Air Blank 0.000 10:52 Control Test 0.080 10:52 Air Blank 0.000 10:52 Control Test 0.080 10:53 Air Blank 0.000 10:53 Control Test 0.079 10:54 Air Blank 0.000 10:54 Control Test Stats Average 0.0797 Std Dev 0.0006 Rel Std Dev(%) 0.7247
Operator's Signature <i>W</i>	Operator's Signature <i>W</i>	Operator's Signature <i>W</i>	Operator's Signature <i>W</i>

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI PD

Serial Number: 80-000873

Time of Inspection: 13:24

Date of Inspection: 07/22/2024

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: 202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#: 202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#: 202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: 01923080A3 Exp: 02/05/2025
0.000	0.050	0.081	0.202	0.078
0.000	0.051	0.081	0.202	0.079
0.000	0.050	0.081	0.201	0.078
0.000	0.051	0.081	0.202	0.078
0.000	0.051	0.080	0.201	0.078
0.000	0.051	0.081	0.201	0.078
0.000	0.051	0.081	0.201	0.078
0.000	0.051	0.081	0.201	0.078
0.000	0.051	0.081	0.201	0.078
0.000	0.051	0.082	0.201	0.078


Standard Deviations	0.0004	0.0004	0.0004	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

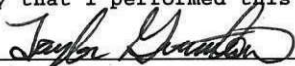

TAYLOR D GUTSCHOW

Signature and Printed Name

07/22/2024
Date

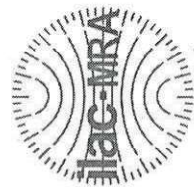

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000873, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000873</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>MIAMI PD</u>	0.050 g/210 L	0.004
Calibration Date:	<u>07/22/2024</u>	0.080 g/210 L	0.004
Calibration Time:	<u>13:24</u>	0.200 g/210 L	0.007
		0.080 g/210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

07/22/2024

Date

TAYLOR D GUTSCHOW,

Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality