



INSTRUMENT PROCESSING SHEET

Agency Charlotte CSOS/N 80-000946Florida Department of
Law EnforcementDate In 08/22/2024 DI Completion Date N/A☐ Ship ☐ P/U ☐ H/D ☒ CMI ☐ EE

Intake By <u>TDG</u> Date <u>08/23/2024</u>	Quality Checks By _____ Date _____	Flow Calibration By _____ Date _____																																																											
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: <u>DGS</u>	☐ Breath Tube Screen ☐ Replace External O-Rings ☐ Instrument Set Up Verified ☐ R-Value _____ ☐ Flow Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) ☐ Barometric Pressure Check Gauge ID # _____ ☐ Stability Checks <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050			0.080			0.200			0.080 DGS	N/A		Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____																																												
Simulator	Serial #	Lot #/Exp																																																											
0.050																																																													
0.080																																																													
0.200																																																													
0.080 DGS	N/A																																																												
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> Notes/Suggested Service: <u>Instrument sent to FDLE for upload prior to repair. Modem port became recessed and cannot be accessed. Uploaded via direct-connect on 8/23. (TDG)</u>	Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection By _____ Barometric Pressure ID# _____ Gauge _____ Instrument _____ Mouth Alcohol Solution Lot # _____ Acetone Stock Solution Lot # _____ <table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td></td></tr><tr><td>Interferent</td><td></td></tr><tr><td>0.050</td><td></td></tr><tr><td>0.080</td><td></td></tr><tr><td>0.200</td><td></td></tr></tbody></table> Attachments ☐ Form 41 ☐ Post-Stability Checks ☐ Stability Checks ☐ Flow Calibration ☐ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☒ Other <u>Form 51</u> ☐ Instrument Complies with Chapter 11D-8, FAC ☒ Instrument Does Not Comply with Chapter 11D-8, FAC ☐ Return to/Place into Evidentiary Use ☒ Remain Out of Evidentiary Use ☐ Conduct an Agency Inspection Before Evidentiary Use Tech Review / Date _____ Admin Review / Date _____	Simulator	Serial Number	0.000		Interferent		0.050		0.080		0.200	
Simulator	Serial #	Lot #	Expiration																																																										
0.000		N/A	N/A																																																										
0.040																																																													
0.100																																																													
0.200																																																													
0.300																																																													
0.080 DGS	N/A																																																												
Simulator	Serial #	Lot #	Expiration																																																										
0.050																																																													
0.080																																																													
0.200																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number																																																												
0.000																																																													
Interferent																																																													
0.050																																																													
0.080																																																													
0.200																																																													

Return Material Authorization

Ship to:

☒ CMI, Inc.

☐ Enforcement Electronics

Shipment to repair facility authorized by: Michael Kern on 08/22/2024

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-000946

Bill To Address:

Charlotte County Sheriff's Office

Attn: Michael Kern

Ship to Address:

Florida Department of Law Enforcement

Fort Myers Regional Operations Center

Attn: Taylor Gutschow

4700 Terminal Drive, Suite 1

Fort Myers, FL 33907

Reason for Return:

Modem port became pushed up inside the instrument and can't be accessed. Records were
uploaded via direct-connect by FDLE.

Please choose one of the following options:

- ☐ 1. I _____, authorize all repairs.
- ☐ 2. I _____, authorize repairs up to \$ _____.
- ☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Michael Kern

Phone #: 941-979-2864

Email: MKern@ccsofl.net

ATP Contact Name: Taylor Gutschow

ATP Email: TaylorGutschow@fdle.state.fl.us



INSTRUMENT PROCESSING SHEET

Agency Charlotte CSOS/N 80-000946Florida Department of
Law EnforcementDate In 06/10/2024 DI Completion Date 07/03/2024☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>TDG</u> Date <u>06/12/2024</u> ☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____	Quality Checks By <u>TDG</u> Date <u>07/03/2024</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>154</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP106</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.500</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP6287</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP6288</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>01923080A3 02/05/2025</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K 03/29/2025	0.080	MP6287	202303L 03/29/2025	0.200	MP6288	202304C 04/05/2025	0.080 DGS	N/A	01923080A3 02/05/2025	Flow Calibration By _____ Date _____ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By <u>TDG</u> Date <u>07/02/2024</u> ☒ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____																																	
Simulator	Serial #	Lot #/Exp																																																
0.050	MP6286	202303K 03/29/2025																																																
0.080	MP6287	202303L 03/29/2025																																																
0.200	MP6288	202304C 04/05/2025																																																
0.080 DGS	N/A	01923080A3 02/05/2025																																																
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> Notes/Suggested Service: <u>Checked breath tube screen, replaced o-rings, and verified setup on 7/2. Finished the rest of the Quality Checks on 7/3. (TDG)</u> _____ _____ _____ _____ _____			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A		
Simulator	Serial #	Lot #	Expiration																																															
0.000		N/A	N/A																																															
0.040																																																		
0.100																																																		
0.200																																																		
0.300																																																		
0.080 DGS	N/A																																																	
Simulator	Serial #	Lot #	Expiration																																															
0.050																																																		
0.080																																																		
0.200																																																		
0.080 DGS	N/A																																																	
Department Inspection By <u>TDG</u> Barometric Pressure ID# <u>26932</u> Gauge <u>1019</u> Instrument <u>1019</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2023-B</u> <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table> Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☒ Form 40 ☐ Calibration Adjustment ☐ Other _____ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use Conduct an Agency Inspection Before Evidentiary Use Phil Nicodemio Digitally signed by Phil Nicodemio Date: 2024.07.15 09:56:37 -0400 Shayla Platt Digitally signed by Shayla Platt Date: 2024.07.15 10:26:25 -0400 Tech Review / Date _____ Admin Review / Date _____			Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288																																				
Simulator	Serial Number																																																	
0.000	MP6284																																																	
Interferent	MP6285																																																	
0.050	MP6286																																																	
0.080	MP6287																																																	
0.200	MP6288																																																	

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT – INTOXILYZER 8000

Agency: CHARLOTTE COUNTY SO
Time of Inspection: 10:27

Date of Inspection: 07/03/2024

Serial Number: 80-000946
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

AI NOT CONDUCTED. BYPASSED TO BRING OUT OF DISABLED MODE.

Not determined ^{mg} 7/3/24

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Taylor D Gutschow

TAYLOR D GUTSCHOW

Signature and Printed Name

07/03/2024
Date

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
CHARLOTTE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 07/03/2024 SN 80-000946 Software: 8100.27 Test g/210L Time Air Blank 0.000 10:32 Control Test 0.050 10:33 Air Blank 0.000 10:33 Control Test 0.050 10:34 Air Blank 0.000 10:34 Control Test 0.050 10:35 Air Blank 0.000 10:36 Control Test Stats Average 0.0500 Std Dev 0.0000 Rel Std Dev(%) 0.0000 Operator's Signature 	CHARLOTTE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 07/03/2024 SN 80-000946 Software: 8100.27 Test g/210L Time Air Blank 0.000 10:43 Control Test 0.078 10:43 Air Blank 0.000 10:44 Control Test 0.079 10:45 Air Blank 0.000 10:45 Control Test 0.079 10:46 Air Blank 0.000 10:46 Control Test Stats Average 0.0787 Std Dev 0.0006 Rel Std Dev(%) 0.7339 Operator's Signature 	CHARLOTTE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 07/03/2024 SN 80-000946 Software: 8100.27 Test g/210L Time Air Blank 0.000 10:57 Control Test 0.198 10:58 Air Blank 0.000 10:59 Control Test 0.197 10:59 Air Blank 0.000 11:00 Control Test 0.196 11:00 Air Blank 0.000 11:01 Control Test Stats Average 0.1970 Std Dev 0.0010 Rel Std Dev(%) 0.5076 Operator's Signature 	CHARLOTTE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 07/03/2024 SN 80-000946 Software: 8100.27 Test g/210L Time Air Blank 0.000 11:13 Control Test 0.079 11:13 Air Blank 0.000 11:13 Control Test 0.079 11:14 Air Blank 0.000 11:14 Control Test 0.078 11:15 Air Blank 0.000 11:15 Control Test Stats Average 0.0787 Std Dev 0.0006 Rel Std Dev(%) 0.7339 Operator's Signature 

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: CHARLOTTE COUNTY SO
Time of Inspection: 14:25

Date of Inspection: 07/03/2024

Serial Number: 80-000946
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:01923080A3 Exp: 02/05/2025
0.000	0.050	0.079	0.198	0.079
0.000	0.050	0.078	0.197	0.079
0.000	0.050	0.080	0.196	0.079
0.000	0.049	0.078	0.197	0.078
0.000	0.050	0.078	0.198	0.078
0.000	0.050	0.079	0.196	0.079
0.000	0.049	0.079	0.197	0.078
0.000	0.050	0.079	0.197	0.078
0.000	0.050	0.078	0.197	0.079
0.000	0.049	0.079	0.197	0.078


Standard Deviations	0.0004	0.0006	0.0006	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0005 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

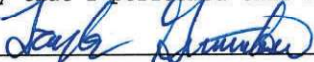

TAYLOR D GUTSCHOW

Signature and Printed Name

07/03/2024
Date

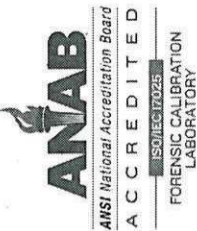

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000946, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000946</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>CHARLOTTE COUNTY SO</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>07/03/2024</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>14:25</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.



07/03/2024

Date

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality