



INSTRUMENT PROCESSING SHEET

Agency New Port Richey PDS/N 80-001064Florida Department of
Law EnforcementDate In 08/22/2024 DI Completion Date 08/28/2024☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>TDG</u> Date <u>08/26/2024</u> ☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____	Quality Checks By <u>TDG</u> Date <u>08/27/2024</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>194</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.230</u> (.228 - .278) 103 mm <u>0.503</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5099</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP5100</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP5101</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>01923080A3 02/05/2025</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP5099	202303K 03/29/2025	0.080	MP5100	202303L 03/29/2025	0.200	MP5101	202304C 04/05/2025	0.080 DGS	N/A	01923080A3 02/05/2025	Flow Calibration By _____ Date _____ Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____ _____																																	
Simulator	Serial #	Lot #/Exp																																																
0.050	MP5099	202303K 03/29/2025																																																
0.080	MP5100	202303L 03/29/2025																																																
0.200	MP5101	202304C 04/05/2025																																																
0.080 DGS	N/A	01923080A3 02/05/2025																																																
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A		
Simulator	Serial #	Lot #	Expiration																																															
0.000		N/A	N/A																																															
0.040																																																		
0.100																																																		
0.200																																																		
0.300																																																		
0.080 DGS	N/A																																																	
Simulator	Serial #	Lot #	Expiration																																															
0.050																																																		
0.080																																																		
0.200																																																		
0.080 DGS	N/A																																																	
Department Inspection By <u>TDG</u> Barometric Pressure ID# <u>26932</u> Gauge <u>1018</u> Instrument <u>1020</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2023-B</u> <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5097</td></tr><tr><td>Interferent</td><td>MP5098</td></tr><tr><td>0.050</td><td>MP5099</td></tr><tr><td>0.080</td><td>MP5100</td></tr><tr><td>0.200</td><td>MP5101</td></tr></tbody></table> Attachments <table border="1" style="width:100%"><tr><td>☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment</td><td>☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 ☐ Other _____</td></tr></table> ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use <table border="1" style="width:100%"><tr><td>Phil Nicodem Digitally signed by Phil Nicodem Date: 2024.09.09 10:53:02 -0400</td><td>Shayla Platt Digitally signed by Shayla Platt Date: 2024.09.09 14:13:56 -0400</td></tr></table> Tech Review / Date _____ Admin Review / Date _____			Simulator	Serial Number	0.000	MP5097	Interferent	MP5098	0.050	MP5099	0.080	MP5100	0.200	MP5101	☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 ☐ Other _____	Phil Nicodem Digitally signed by Phil Nicodem Date: 2024.09.09 10:53:02 -0400	Shayla Platt Digitally signed by Shayla Platt Date: 2024.09.09 14:13:56 -0400																																
Simulator	Serial Number																																																	
0.000	MP5097																																																	
Interferent	MP5098																																																	
0.050	MP5099																																																	
0.080	MP5100																																																	
0.200	MP5101																																																	
☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 ☐ Other _____																																																	
Phil Nicodem Digitally signed by Phil Nicodem Date: 2024.09.09 10:53:02 -0400	Shayla Platt Digitally signed by Shayla Platt Date: 2024.09.09 14:13:56 -0400																																																	

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: NEW PORT RICHEY PD

Time of Inspection: 14:42

Date of Inspection: 08/27/2024

Serial Number: 80-001064

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____


Remarks:

AI NOT CONDUCTED. BYPASSED TO BRING OUT OF DISABLED MODE.

Not determined ^{ML} 08/27/2024

The above instrument complies (☒) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

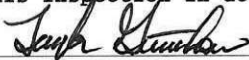

TAYLOR D GUTSCHOW

Signature and Printed Name

08/27/2024
Date

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
✓	✓	✓	✓
NEW PORT RICHEV PD Intoxilyzer - Alcohol Analyzer Model 8000 08/27/2024 Software: 8100.27 SN 80-001064 Test g/210L Time Air Blank 0.000 15:18 Control Test 0.150 15:19 Air Blank 0.000 15:20 Control Test 0.050 15:20 Air Blank 0.000 15:21 Control Test 0.050 15:21 Air Blank 0.000 15:22 Control Test Stats Average 0.0500 Std Dev 0.0000 Rel Std Dev(%) 0.0000	NEW PORT RICHEV PD Intoxilyzer - Alcohol Analyzer Model 8000 08/27/2024 Software: 8100.27 SN 80-001064 Test g/210L Time Air Blank 0.000 14:59 Control Test 0.080 15:00 Air Blank 0.000 15:00 Control Test 0.079 15:01 Air Blank 0.000 15:02 Control Test 0.080 15:02 Air Blank 0.000 15:03 Control Test Stats Average 0.0797 Std Dev 0.0006 Rel Std Dev(%) 0.7247	NEW PORT RICHEV PD Intoxilyzer - Alcohol Analyzer Model 8000 08/27/2024 Software: 8100.27 SN 80-001064 Test g/210L Time Air Blank 0.000 15:05 Control Test 0.197 15:06 Air Blank 0.000 15:06 Control Test 0.198 15:07 Air Blank 0.000 15:07 Control Test 0.197 15:08 Air Blank 0.000 15:09 Control Test Stats Average 0.1973 Std Dev 0.0006 Rel Std Dev(%) 0.2926	NEW PORT RICHEV PD Intoxilyzer - Alcohol Analyzer Model 8000 08/27/2024 Software: 8100.27 SN 80-001064 Test g/210L Time Air Blank 0.000 15:11 Control Test 0.080 15:12 Air Blank 0.000 15:12 Control Test 0.080 15:13 Air Blank 0.000 15:13 Control Test 0.080 15:13 Air Blank 0.000 15:14 Control Test Stats Average 0.0800 Std Dev 0.0000 Rel Std Dev(%) 0.0000
Operator's Signature	Operator's Signature	Operator's Signature	Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: NEW PORT RICHEY PD
Time of Inspection: 13:23

Date of Inspection: 08/28/2024

Serial Number: 80-001064
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:01923080A3 Exp: 02/05/2025
0.000	0.049	0.079	0.197	0.080
0.000	0.050	0.079	0.197	0.080
0.000	0.050	0.079	0.197	0.080
0.000	0.050	0.079	0.197	0.080
0.000	0.050	0.079	0.197	0.080
0.000	0.050	0.080	0.196	0.080
0.000	0.050	0.079	0.197	0.080
0.000	0.050	0.079	0.197	0.080
0.000	0.050	0.080	0.197	0.080
0.000	0.050	0.079	0.197	0.080


Standard Deviations	0.0003	0.0004	0.0003	0.0000
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

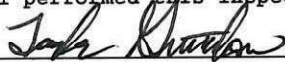

TAYLOR D GUTSCHOW

Signature and Printed Name

08/28/2024
Date

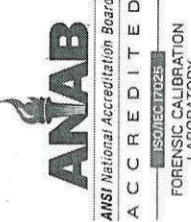

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001064, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001064</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>NEW PORT RICHEY PD</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>08/28/2024</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>13:23</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

08/28/2024

Date

TAYLOR D GUTSCHOW,

Department Inspector

Service • Integrity • Respect • Quality



INSTRUMENT PROCESSING SHEET

Agency New Port Richey PDS/N 80-001064Florida Department of
Law EnforcementDate In 05/28/2024 DI Completion Date 06/05/2024☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>TDG</u> Date <u>06/04/2024</u>		Quality Checks By <u>TDG</u> Date <u>06/05/2024</u>		Flow Calibration By _____ Date _____																													
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>Upload evaluation.</u>		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>209</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.152</u> (.139 - .169) 36 mm <u>0.171</u> (.156 - .190) 53 mm <u>0.246</u> (.228 - .278) 103 mm <u>0.500</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks		☐ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																													
		<table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP6287</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP6288</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG222203 08/10/2024</td></tr></tbody></table>		Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K 03/29/2025	0.080	MP6287	202303L 03/29/2025	0.200	MP6288	202304C 04/05/2025	0.080 DGS	N/A	AG222203 08/10/2024	Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____														
Simulator	Serial #	Lot #/Exp																															
0.050	MP6286	202303K 03/29/2025																															
0.080	MP6287	202303L 03/29/2025																															
0.200	MP6288	202304C 04/05/2025																															
0.080 DGS	N/A	AG222203 08/10/2024																															
Calibration Adjustment By _____				Department Inspection By <u>TDG</u>																													
Barometric Pressure Gauge _____ ID # _____				Barometric Pressure ID# <u>26932</u>																													
<table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>				Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge <u>1012</u> Instrument <u>1013</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2023-B</u>	
Simulator	Serial #	Lot #	Expiration																														
0.000		N/A	N/A																														
0.040																																	
0.100																																	
0.200																																	
0.300																																	
0.080 DGS	N/A																																
☐ Post Calibration Adjustment Stability Checks				<table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table>		Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288																
Simulator	Serial Number																																
0.000	MP6284																																
Interferent	MP6285																																
0.050	MP6286																																
0.080	MP6287																																
0.200	MP6288																																
<table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>				Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☒ Form 40 ☐ Calibration Adjustment ☐ Other _____									
Simulator	Serial #	Lot #	Expiration																														
0.050																																	
0.080																																	
0.200																																	
0.080 DGS	N/A																																
Notes/Suggested Service: <u>Checked breath tube screen and replaced o-rings on 6/4. Finished the rest of the Quality Checks on 6/5. Instrument uploaded on the first attempt. Added plastic caps for shipping. (TDG)</u>				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																													
				Digitally signed by Phil Nicodemo Digitally signed by Phil Nicodemo Date: 2024.06.19 09:52:14 -0400 Shayla Platt Shayla Platt Date: 2024.06.20 12:11:22 -04'00'																													
Tech Review / Date _____				Admin Review / Date _____																													

Original Instrument Processing Sheet lost during digital signature process.
SP 6/13/24.

Shayla
Platt

Digitally signed by
Shayla Platt
Date: 2024.06.07
13:07:47 -04'00'

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: NEW PORT RICHEY PD
Time of Inspection: 09:18

Date of Inspection: 06/05/2024

Serial Number: 80-001064
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

AI NOT CONDUCTED. BYPASSED TO BRING OUT OF DISABLED MODE.

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Taylor D Gutschow

TAYLOR D GUTSCHOW

Signature and Printed Name

06/05/2024
Date

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
Printed to external printer. Results within acceptable range. Will attach. ML 6/15/24	NEW PORT RICHEY PD Intoxilyzer - Alconol Analyzer Model: 8000 SN 80-001064 06/05/2024 Software: 8100.27 Test g/210L Time Air Blank 0.000 09:28 Control Test 0.078 09:29 Air Blank 0.000 09:30 Control Test 0.078 09:31 Air Blank 0.000 09:31 Control Test 0.079 09:32 Air Blank 0.000 09:32 Control Test Stats Average 0.0783 Std Dev 0.0006 Rel Std Dev(%) 0.7370 Operator's Signature 	NEW PORT RICHEY PD Intoxilyzer - Alconol Analyzer Model: 8000 SN 80-001064 06/05/2024 Software: 8100.27 Test g/210L Time Air Blank 0.000 09:35 Control Test 0.196 09:36 Air Blank 0.000 09:37 Control Test 0.196 09:38 Air Blank 0.000 09:39 Control Test 0.196 09:39 Air Blank 0.000 09:39 Control Test Stats Average 0.1960 Std Dev 0.0000 Rel Std Dev(%) 0.0000 Operator's Signature 	NEW PORT RICHEY PD Intoxilyzer - Alconol Analyzer Model: 8000 SN 80-001064 06/05/2024 Software: 8100.27 Test g/210L Time Air Blank 0.000 09:39 Control Test 0.180 09:40 Air Blank 0.000 09:40 Control Test 0.180 09:41 Air Blank 0.000 09:41 Control Test 0.180 09:41 Air Blank 0.000 09:41 Control Test Stats Average 0.1800 Std Dev 0.0000 Rel Std Dev(%) 0.0000 Operator's Signature 

NEW PORT RICHEY PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001064
06/05/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:22
Control Test	0.049	09:22
Air Blank	0.000	09:23
Control Test	0.050	09:23
Air Blank	0.000	09:24
Control Test	0.049	09:25
Air Blank	0.000	09:25
Control Test Stats		
Average	0.0493	
Std Dev	0.0006	
Rel Std Dev(%)	1.1703	



Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT – INTOXILYZER 8000

Agency: NEW PORT RICHEY PD
Time of Inspection: 11:57

Date of Inspection: 06/05/2024

Serial Number: 80-001064
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG222203 Exp: 08/10/2024
0.000	0.050	0.079	0.196	0.080
0.000	0.050	0.079	0.196	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.050	0.079	0.196	0.080
0.000	0.049	0.079	0.197	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.050	0.079	0.196	0.080


Standard Deviations	0.0005	0.0000	0.0003	0.0000
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

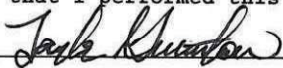

TAYLOR D GUTSCHOW

Signature and Printed Name

06/05/2024
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001064, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-001064	UNCERTAINTY* $\pm$
Owning Agency:	NEW PORT RICHEY PD	0.050 g/ 210 L 0.004
Calibration Date:	06/05/2024	0.080 g/ 210 L 0.004
Calibration Time:	11:57	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

06/05/2024

Date

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Shayla
Platt

Digitally signed by
Shayla Platt
Date: 2024.06.07
13:07:47 -04'00'

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: NEW PORT RICHEY PD
Time of Inspection: 09:18

Date of Inspection: 06/05/2024

Serial Number: 80-001064
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

AI NOT CONDUCTED. BYPASSED TO BRING OUT OF DISABLED MODE.

Not determined *ML 6/5/24*

The above instrument complies (☒) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Taylor D Gutschow

TAYLOR D GUTSCHOW

Signature and Printed Name

06/05/2024
Date

Stability Checks

0.05g/210L		0.08g/210L		0.20g/210L		DGS 0.08g/210L																																																																																																										
0.047 to 0.053		0.077 to 0.083		0.194 to 0.206		0.077 to 0.083																																																																																																										
✓		✓		✓		✓																																																																																																										
Printed to external printer. Results within acceptable range. Will attach. ML 6/15/24		NEW PORT RICHEY PD Intoxilyzer - Alcotest Analyzer Model: 8000 SN 80-001064 06/05/2024 Software: 8100.27 <table><thead><tr><th>Test</th><th>9/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>09:28</td></tr><tr><td>Control Test</td><td>0.078</td><td>09:29</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:30</td></tr><tr><td>Control Test</td><td>0.078</td><td>09:30</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:31</td></tr><tr><td>Control Test</td><td>0.079</td><td>09:31</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:32</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0783</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7370</td><td></td></tr></tbody></table> Operator's Signature ML	Test	9/210L	Time	Air Blank	0.000	09:28	Control Test	0.078	09:29	Air Blank	0.000	09:30	Control Test	0.078	09:30	Air Blank	0.000	09:31	Control Test	0.079	09:31	Air Blank	0.000	09:32	Control Test Stats			Average	0.0783		Std Dev	0.0006		Rel Std Dev(%)	0.7370		NEW PORT RICHEY PD Intoxilyzer - Alcotest Analyzer Model: 8000 SN 80-001064 06/05/2024 Software: 8100.27 <table><thead><tr><th>Test</th><th>9/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>09:35</td></tr><tr><td>Control Test</td><td>0.196</td><td>09:36</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:36</td></tr><tr><td>Control Test</td><td>0.196</td><td>09:37</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:38</td></tr><tr><td>Control Test</td><td>0.196</td><td>09:38</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:39</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.1960</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></tbody></table> Operator's Signature ML	Test	9/210L	Time	Air Blank	0.000	09:35	Control Test	0.196	09:36	Air Blank	0.000	09:36	Control Test	0.196	09:37	Air Blank	0.000	09:38	Control Test	0.196	09:38	Air Blank	0.000	09:39	Control Test Stats			Average	0.1960		Std Dev	0.0000		Rel Std Dev(%)	0.0000		NEW PORT RICHEY PD Intoxilyzer - Alcotest Analyzer Model: 8000 SN 80-001064 06/05/2024 Software: 8100.27 <table><thead><tr><th>Test</th><th>9/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>09:39</td></tr><tr><td>Control Test</td><td>0.080</td><td>09:40</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:40</td></tr><tr><td>Control Test</td><td>0.080</td><td>09:41</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:41</td></tr><tr><td>Control Test</td><td>0.080</td><td>09:41</td></tr><tr><td>Air Blank</td><td>0.000</td><td>09:42</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0800</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></tbody></table> Operator's Signature ML	Test	9/210L	Time	Air Blank	0.000	09:39	Control Test	0.080	09:40	Air Blank	0.000	09:40	Control Test	0.080	09:41	Air Blank	0.000	09:41	Control Test	0.080	09:41	Air Blank	0.000	09:42	Control Test Stats			Average	0.0800		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Test	9/210L	Time																																																																																																														
Air Blank	0.000	09:28																																																																																																														
Control Test	0.078	09:29																																																																																																														
Air Blank	0.000	09:30																																																																																																														
Control Test	0.078	09:30																																																																																																														
Air Blank	0.000	09:31																																																																																																														
Control Test	0.079	09:31																																																																																																														
Air Blank	0.000	09:32																																																																																																														
Control Test Stats																																																																																																																
Average	0.0783																																																																																																															
Std Dev	0.0006																																																																																																															
Rel Std Dev(%)	0.7370																																																																																																															
Test	9/210L	Time																																																																																																														
Air Blank	0.000	09:35																																																																																																														
Control Test	0.196	09:36																																																																																																														
Air Blank	0.000	09:36																																																																																																														
Control Test	0.196	09:37																																																																																																														
Air Blank	0.000	09:38																																																																																																														
Control Test	0.196	09:38																																																																																																														
Air Blank	0.000	09:39																																																																																																														
Control Test Stats																																																																																																																
Average	0.1960																																																																																																															
Std Dev	0.0000																																																																																																															
Rel Std Dev(%)	0.0000																																																																																																															
Test	9/210L	Time																																																																																																														
Air Blank	0.000	09:39																																																																																																														
Control Test	0.080	09:40																																																																																																														
Air Blank	0.000	09:40																																																																																																														
Control Test	0.080	09:41																																																																																																														
Air Blank	0.000	09:41																																																																																																														
Control Test	0.080	09:41																																																																																																														
Air Blank	0.000	09:42																																																																																																														
Control Test Stats																																																																																																																
Average	0.0800																																																																																																															
Std Dev	0.0000																																																																																																															
Rel Std Dev(%)	0.0000																																																																																																															

NEW PORT RICHEY PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001064
06/05/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	09:22
Control Test	0.049	09:22
Air Blank	0.000	09:23
Control Test	0.050	09:23
Air Blank	0.000	09:24
Control Test	0.049	09:25
Air Blank	0.000	09:25
Control Test Stats		
Average	0.0493	
Std Dev	0.0006	
Rel Std Dev(%)	1.1703	



Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: NEW PORT RICHEY PD
Time of Inspection: 11:57

Date of Inspection: 06/05/2024

Serial Number: 80-001064
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG222203 Exp: 08/10/2024
0.000	0.050	0.079	0.196	0.080
0.000	0.050	0.079	0.196	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.050	0.079	0.196	0.080
0.000	0.049	0.079	0.197	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.050	0.079	0.196	0.080


Standard Deviations	0.0005	0.0000	0.0003	0.0000
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

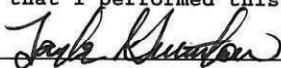

TAYLOR D GUTSCHOW

Signature and Printed Name

06/05/2024
Date

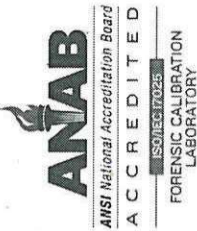

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001064, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-001064	UNCERTAINTY* $\pm$
Owning Agency:	NEW PORT RICHEY PD	0.050 g/ 210 L 0.004
Calibration Date:	06/05/2024	0.080 g/ 210 L 0.004
Calibration Time:	11:57	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

06/05/2024

Date

TAYLOR D GUTSCHOW,

Department Inspector

Service • Integrity • Respect • Quality



Agency New Port Richey PD

S/N 80-001064

Date In 12-22-2023 DI Completion Date 01-23-2024

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

FDLE/ATP Form 48 January 2022
Issuing Authority: Alcohol Testing Program

Page 1 of 1

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: NEW PORT RICHEY PD
Time of Inspection: 10:59

Date of Inspection: 01/23/2024

Serial Number: 80-001064
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG229803 Exp: 10/25/2024
0.000	0.000 / 0.050	0.079	0.197	0.079
0.000	0.000 / 0.050	0.080	0.197	0.079
0.000	0.050 / 0.050	0.080	0.198	0.079
0.000	0.050 / 0.050	0.080	0.198	0.079
0.000	0.050 / 0.050	0.079	0.198	0.079
0.000	0.050 / 0.050	0.080	0.197	0.079
0.000	0.050 / 0.051	0.080	0.198	0.079
0.000	0.050 / 0.050	0.080	0.198	0.079
0.000	0.051 / 0.051	0.080	0.198	0.079
0.000	0.050 / 0.050	0.080	0.198	0.079
Standard Deviations	0.0211 / 0.0004	0.0004	0.0004	0.0000

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

05: Control Outside Tolerance SIM NOT CONNECTED.

connected sim and repeated 0.050g/210L test
successfully ✓✓

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Israel Soto

ISRAEL SOTO

Signature and Printed Name

01/23/2024
Date

Stability checks

01-04-24

80-001064

NEW PORT RICHEY PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001064
01/04/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:00
Control Test	0.051	12:01
Air Blank	0.000	12:01
Control Test	0.051	12:02
Air Blank	0.000	12:03
Control Test	0.050	12:03
Air Blank	0.000	12:04
Control Test Stats		
Average	0.0507	
Std Dev	0.0006	
Rel Std Dev(%)	1.1395	

Operator's Signature

AZ

NEW PORT RICHEY PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001064
01/04/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:22
Control Test	0.081	12:22
Air Blank	0.000	12:23
Control Test	0.081	12:23
Air Blank	0.000	12:23
Control Test	0.081	12:24
Air Blank	0.000	12:24
Control Test Stats		
Average	0.0810	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

AZ

NEW PORT RICHEY PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001064
01/04/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:13
Control Test	0.199	12:14
Air Blank	0.000	12:14
Control Test	0.197	12:15
Air Blank	0.000	12:15
Control Test	0.197	12:16
Air Blank	0.000	12:17
Control Test Stats		
Average	0.1977	
Std Dev	0.0012	
Rel Std Dev(%)	0.5842	

Operator's Signature

AZ

NEW PORT RICHEY PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001064
01/04/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:06
Control Test	0.000	12:07
Air Blank	0.000	12:07
Control Test	0.080	12:08
Air Blank	0.000	12:08
Control Test	0.080	12:09
Air Blank	0.000	12:10
Control Test Stats		
Average	0.0800	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

AZ

065



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001064, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001064</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>NEW PORT RICHEY PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>01/23/2024</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>10:59</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 March 2022

Issuing Authority: Alcohol Testing Program

01/23/2024

Date

Israel Soto

ISRAEL SOTO,

Department Inspector

Service • Integrity • Respect • Quality