



INSTRUMENT PROCESSING SHEET

Agency Clermont PDS/N 80-001174Florida Department of
Law EnforcementDate In 07/26/2024 DI Completion Date 08/06/2024☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By	TDG	Date	08/05/2024	Quality Checks	By	TDG	Date	08/05/2024	Flow Calibration	By	Date																																							
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: <u>Shipped in pelican case.</u>					☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>227</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.152</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.500</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks					Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																									
					<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5099</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP5100</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP5101</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>01923080A3 02/05/2025</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP5099	202303K 03/29/2025	0.080	MP5100	202303L 03/29/2025	0.200	MP5101	202304C 04/05/2025	0.080 DGS	N/A	01923080A3 02/05/2025	<table border="1"><thead><tr><th colspan="2">Maintenance By _____ Date _____</th></tr><tr><td>☐ Battery Replacement</td><td></td></tr><tr><td>☐ Dry Gas Regulator Replacement</td><td></td></tr><tr><td>☐ Breath Tube Replacement</td><td></td></tr><tr><td>☐ Other _____</td><td></td></tr></thead></table>				Maintenance By _____ Date _____		☐ Battery Replacement		☐ Dry Gas Regulator Replacement		☐ Breath Tube Replacement		☐ Other _____																
Simulator	Serial #	Lot #/Exp																																																	
0.050	MP5099	202303K 03/29/2025																																																	
0.080	MP5100	202303L 03/29/2025																																																	
0.200	MP5101	202304C 04/05/2025																																																	
0.080 DGS	N/A	01923080A3 02/05/2025																																																	
Maintenance By _____ Date _____																																																			
☐ Battery Replacement																																																			
☐ Dry Gas Regulator Replacement																																																			
☐ Breath Tube Replacement																																																			
☐ Other _____																																																			
Calibration Adjustment By _____					Department Inspection By <u>TDG</u>																																														
Barometric Pressure Gauge _____ ID # _____					Barometric Pressure ID# <u>26932</u> Gauge <u>1014</u> Instrument <u>1009</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2023-B</u>																																														
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>					Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5097</td></tr><tr><td>Interferent</td><td>MP5098</td></tr><tr><td>0.050</td><td>MP5099</td></tr><tr><td>0.080</td><td>MP5100</td></tr><tr><td>0.200</td><td>MP5101</td></tr></tbody></table>							Simulator	Serial Number	0.000	MP5097	Interferent	MP5098	0.050	MP5099	0.080	MP5100	0.200	MP5101
Simulator	Serial #	Lot #	Expiration																																																
0.000		N/A	N/A																																																
0.040																																																			
0.100																																																			
0.200																																																			
0.300																																																			
0.080 DGS	N/A																																																		
Simulator	Serial Number																																																		
0.000	MP5097																																																		
Interferent	MP5098																																																		
0.050	MP5099																																																		
0.080	MP5100																																																		
0.200	MP5101																																																		
☐ Post Calibration Adjustment Stability Checks					<table border="1"><thead><tr><th colspan="2">Attachments</th></tr></thead><tbody><tr><td>☒ Form 41</td><td>☐ Post-Stability Checks</td></tr><tr><td>☒ Stability Checks</td><td>☐ Flow Calibration</td></tr><tr><td>☒ Calibration Certificate</td><td>☐ Form 40</td></tr><tr><td>☐ Calibration Adjustment</td><td>☐ Other _____</td></tr></tbody></table>							Attachments		☒ Form 41	☐ Post-Stability Checks	☒ Stability Checks	☐ Flow Calibration	☒ Calibration Certificate	☐ Form 40	☐ Calibration Adjustment	☐ Other _____																														
Attachments																																																			
☒ Form 41	☐ Post-Stability Checks																																																		
☒ Stability Checks	☐ Flow Calibration																																																		
☒ Calibration Certificate	☐ Form 40																																																		
☐ Calibration Adjustment	☐ Other _____																																																		
Notes/Suggested Service: <u>Added static bag for shipping back to the agency. (TDG)</u>					☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																														
					Phil Nicodemo Tech Review / Date _____ Digitally signed by Phil Nicodemo Date: 2024.08.08 14:47:24 -04'00'																																														
					Shayla Platt Admin Review / Date _____ Digitally signed by Shayla Platt Date: 2024.08.08 23:18:58 -04'00'																																														

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																				
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083																																																																																																																																				
CLERMONT PD Intoxilyzer - Alcohol Analyzer Model 8000 08/05/2024 Software: 8100.27 SN 80-001174 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>15:29</td></tr> <tr><td>Control Test</td><td>0.050</td><td>15:30</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:30</td></tr> <tr><td>Control Test</td><td>0.049</td><td>15:31</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:31</td></tr> <tr><td>Control Test</td><td>0.049</td><td>15:32</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:33</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0493</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.1703</td><td></td></tr> </table>	Air Blank	0.000	15:29	Control Test	0.050	15:30	Air Blank	0.000	15:30	Control Test	0.049	15:31	Air Blank	0.000	15:31	Control Test	0.049	15:32	Air Blank	0.000	15:33	Control Test Stats			Average	0.0493		Std Dev	0.0006		Rel Std Dev(%)	1.1703		CLERMONT PD Intoxilyzer - Alcohol Analyzer Model 8000 08/05/2024 Software: 8100.27 SN 80-001174 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>15:35</td></tr> <tr><td>Control Test</td><td>0.078</td><td>15:36</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:36</td></tr> <tr><td>Control Test</td><td>0.079</td><td>15:37</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:37</td></tr> <tr><td>Control Test</td><td>0.078</td><td>15:38</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:39</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0783</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7370</td><td></td></tr> </table>	Air Blank	0.000	15:35	Control Test	0.078	15:36	Air Blank	0.000	15:36	Control Test	0.079	15:37	Air Blank	0.000	15:37	Control Test	0.078	15:38	Air Blank	0.000	15:39	Control Test Stats			Average	0.0783		Std Dev	0.0006		Rel Std Dev(%)	0.7370		CLERMONT PD Intoxilyzer - Alcohol Analyzer Model 8000 08/05/2024 Software: 8100.27 SN 80-001174 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>15:42</td></tr> <tr><td>Control Test</td><td>0.199</td><td>15:42</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:43</td></tr> <tr><td>Control Test</td><td>0.199</td><td>15:44</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:44</td></tr> <tr><td>Control Test</td><td>0.199</td><td>15:45</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:45</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1990</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </table>	Air Blank	0.000	15:42	Control Test	0.199	15:42	Air Blank	0.000	15:43	Control Test	0.199	15:44	Air Blank	0.000	15:44	Control Test	0.199	15:45	Air Blank	0.000	15:45	Control Test Stats			Average	0.1990		Std Dev	0.0000		Rel Std Dev(%)	0.0000		CLERMONT PD Intoxilyzer - Alcohol Analyzer Model 8000 08/05/2024 Software: 8100.27 SN 80-001174 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>15:47</td></tr> <tr><td>Control Test</td><td>0.081</td><td>15:47</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:47</td></tr> <tr><td>Control Test</td><td>0.080</td><td>15:48</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:48</td></tr> <tr><td>Control Test</td><td>0.081</td><td>15:49</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:49</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0807</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7157</td><td></td></tr> </table>	Air Blank	0.000	15:47	Control Test	0.081	15:47	Air Blank	0.000	15:47	Control Test	0.080	15:48	Air Blank	0.000	15:48	Control Test	0.081	15:49	Air Blank	0.000	15:49	Control Test Stats			Average	0.0807		Std Dev	0.0006		Rel Std Dev(%)	0.7157	
Air Blank	0.000	15:29																																																																																																																																					
Control Test	0.050	15:30																																																																																																																																					
Air Blank	0.000	15:30																																																																																																																																					
Control Test	0.049	15:31																																																																																																																																					
Air Blank	0.000	15:31																																																																																																																																					
Control Test	0.049	15:32																																																																																																																																					
Air Blank	0.000	15:33																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0493																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	1.1703																																																																																																																																						
Air Blank	0.000	15:35																																																																																																																																					
Control Test	0.078	15:36																																																																																																																																					
Air Blank	0.000	15:36																																																																																																																																					
Control Test	0.079	15:37																																																																																																																																					
Air Blank	0.000	15:37																																																																																																																																					
Control Test	0.078	15:38																																																																																																																																					
Air Blank	0.000	15:39																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0783																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.7370																																																																																																																																						
Air Blank	0.000	15:42																																																																																																																																					
Control Test	0.199	15:42																																																																																																																																					
Air Blank	0.000	15:43																																																																																																																																					
Control Test	0.199	15:44																																																																																																																																					
Air Blank	0.000	15:44																																																																																																																																					
Control Test	0.199	15:45																																																																																																																																					
Air Blank	0.000	15:45																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.1990																																																																																																																																						
Std Dev	0.0000																																																																																																																																						
Rel Std Dev(%)	0.0000																																																																																																																																						
Air Blank	0.000	15:47																																																																																																																																					
Control Test	0.081	15:47																																																																																																																																					
Air Blank	0.000	15:47																																																																																																																																					
Control Test	0.080	15:48																																																																																																																																					
Air Blank	0.000	15:48																																																																																																																																					
Control Test	0.081	15:49																																																																																																																																					
Air Blank	0.000	15:49																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0807																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.7157																																																																																																																																						
Operator's Signature	Operator's Signature	Operator's Signature	Operator's Signature																																																																																																																																				

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: CLERMONT PD
Time of Inspection: 12:46

Date of Inspection: 08/06/2024

Serial Number: 80-001174
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:01923080A3 Exp: 02/05/2025
0.000	0.050	0.080	0.200	0.081
0.000	0.050	0.080	0.200	0.081
0.000	0.050	0.080	0.200	0.081
0.000	0.050	0.080	0.201	0.080
0.000	0.050	0.079	0.200	0.081
0.000	0.050	0.080	0.200	0.080
0.000	0.050	0.080	0.200	0.080
0.000	0.050	0.080	0.200	0.081
0.000	0.050	0.080	0.200	0.080
0.000	0.050	0.079	0.200	0.080

Standard Deviations	0.0000	0.0004	0.0003	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



TAYLOR D GUTSCHOW

Signature and Printed Name

08/06/2024
Date

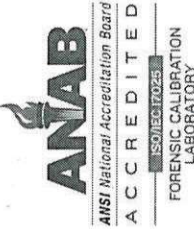

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001174, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001174</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>CLERMONT PD</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>08/06/2024</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>12:46</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

08/06/2024

Date

Taylor D Gutschow

TAYLOR D GUTSCHOW,

Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality