



INSTRUMENT PROCESSING SHEET

Agency Hialeah PDS/N 80-002462Florida Department of
Law EnforcementDate In 07/11/2024 DI Completion Date 07/15/2024☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>TDG</u> Date <u>07/12/2024</u> ☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____	Quality Checks By <u>TDG</u> Date <u>07/12/2024</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>211</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.218*</u> (.228 - .278) 103 mm <u>0.468</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP6287</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP6288</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>01923080A3 02/05/2025</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K 03/29/2025	0.080	MP6287	202303L 03/29/2025	0.200	MP6288	202304C 04/05/2025	0.080 DGS	N/A	01923080A3 02/05/2025	Flow Calibration By <u>TDG</u> Date <u>07/12/2024</u> Flow Column # <u>ATP106</u> ☒ 5L/min - 17mm ☒ 15L/min - 53mm ☒ 30L/min - 103mm ☒ R-Value <u>211</u> ☒ Post Calibration Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.507</u> (.447 - .547) Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____																																	
Simulator	Serial #	Lot #/Exp																																																
0.050	MP6286	202303K 03/29/2025																																																
0.080	MP6287	202303L 03/29/2025																																																
0.200	MP6288	202304C 04/05/2025																																																
0.080 DGS	N/A	01923080A3 02/05/2025																																																
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> Notes/Suggested Service: <u>*Outside nominal range. (TDG)</u> _____ _____ _____ _____ _____ _____ _____ _____			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A		
Simulator	Serial #	Lot #	Expiration																																															
0.000		N/A	N/A																																															
0.040																																																		
0.100																																																		
0.200																																																		
0.300																																																		
0.080 DGS	N/A																																																	
Simulator	Serial #	Lot #	Expiration																																															
0.050																																																		
0.080																																																		
0.200																																																		
0.080 DGS	N/A																																																	
Department Inspection By <u>TDG</u> Barometric Pressure ID# <u>26932</u> Gauge <u>1019</u> Instrument <u>1017</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2023-B</u> <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table> Attachments <table border="1" style="width:100%"><tr><td>☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment</td><td>☐ Post-Stability Checks ☒ Flow Calibration ☐ Form 40 ☐ Other _____</td></tr></table> ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Digitally signed by <u>Shayla Platt</u> Date: 2024.07.26 14:56:13 -04'00' Phil Nicodemus Digitally signed by Phil Nicodemus Date: 2024.07.25 09:01:19 -04'00' Platt Tech Review / Date _____ Admin Review / Date _____			Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288	☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☒ Flow Calibration ☐ Form 40 ☐ Other _____																																		
Simulator	Serial Number																																																	
0.000	MP6284																																																	
Interferent	MP6285																																																	
0.050	MP6286																																																	
0.080	MP6287																																																	
0.200	MP6288																																																	
☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☒ Flow Calibration ☐ Form 40 ☐ Other _____																																																	

Flow Calibration Adjustment(s)

Performed by TDG

HIALEAH PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-002462
07/12/2024
Software: 8100.27

Flow Rate Calibration*****

1: Rate (Liters/min) = 5
SQRT(Diff) = 7.348
2: Rate (Liters/min) = 15
SQRT(Diff) = 11.703
3: Rate (Liters/min) = 30
SQRT(Diff) = 20.832

Dependent Data Scale Factor = 100000 L/min

Independent Data Scale Factor = 256

Rounded Slope = 712

Rounded Intercept = -755094

Correlation = 0.99625

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																				
0.047 to 0.053 ✓	0.077 to 0.083 ✓	0.194 to 0.206 ✓	0.077 to 0.083 ✓ ≤0.003 of Wet ✓																																																																																																																																				
HiLEAH PD Intoxilyzer - Alcohol Analyzer Model 8000 07/12/2024 Software: 8100.27 SN 80-002462 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>13:36</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:37</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:37</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:38</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:39</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:39</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:40</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0490</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </table>	Air Blank	0.000	13:36	Control Test	0.049	13:37	Air Blank	0.000	13:37	Control Test	0.049	13:38	Air Blank	0.000	13:39	Control Test	0.049	13:39	Air Blank	0.000	13:40	Control Test Stats			Average	0.0490		Std Dev	0.0000		Rel Std Dev(%)	0.0000		HiLEAH PD Intoxilyzer - Alcohol Analyzer Model 8000 07/12/2024 Software: 8100.27 SN 80-002462 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>13:43</td></tr> <tr><td>Control Test</td><td>0.078</td><td>13:44</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:45</td></tr> <tr><td>Control Test</td><td>0.078</td><td>13:45</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:46</td></tr> <tr><td>Control Test</td><td>0.078</td><td>13:46</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:47</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0780</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </table>	Air Blank	0.000	13:43	Control Test	0.078	13:44	Air Blank	0.000	13:45	Control Test	0.078	13:45	Air Blank	0.000	13:46	Control Test	0.078	13:46	Air Blank	0.000	13:47	Control Test Stats			Average	0.0780		Std Dev	0.0000		Rel Std Dev(%)	0.0000		HiLEAH PD Intoxilyzer - Alcohol Analyzer Model 8000 07/12/2024 Software: 8100.27 SN 80-002462 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>13:49</td></tr> <tr><td>Control Test</td><td>0.196</td><td>13:50</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:51</td></tr> <tr><td>Control Test</td><td>0.196</td><td>13:51</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:52</td></tr> <tr><td>Control Test</td><td>0.195</td><td>13:52</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:53</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1957</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2951</td><td></td></tr> </table>	Air Blank	0.000	13:49	Control Test	0.196	13:50	Air Blank	0.000	13:51	Control Test	0.196	13:51	Air Blank	0.000	13:52	Control Test	0.195	13:52	Air Blank	0.000	13:53	Control Test Stats			Average	0.1957		Std Dev	0.0006		Rel Std Dev(%)	0.2951		HiLEAH PD Intoxilyzer - Alcohol Analyzer Model 8000 07/12/2024 Software: 8100.27 SN 80-002462 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>13:54</td></tr> <tr><td>Control Test</td><td>0.079</td><td>13:54</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:55</td></tr> <tr><td>Control Test</td><td>0.079</td><td>13:55</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:56</td></tr> <tr><td>Control Test</td><td>0.078</td><td>13:56</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:56</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0787</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7339</td><td></td></tr> </table>	Air Blank	0.000	13:54	Control Test	0.079	13:54	Air Blank	0.000	13:55	Control Test	0.079	13:55	Air Blank	0.000	13:56	Control Test	0.078	13:56	Air Blank	0.000	13:56	Control Test Stats			Average	0.0787		Std Dev	0.0006		Rel Std Dev(%)	0.7339	
Air Blank	0.000	13:36																																																																																																																																					
Control Test	0.049	13:37																																																																																																																																					
Air Blank	0.000	13:37																																																																																																																																					
Control Test	0.049	13:38																																																																																																																																					
Air Blank	0.000	13:39																																																																																																																																					
Control Test	0.049	13:39																																																																																																																																					
Air Blank	0.000	13:40																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0490																																																																																																																																						
Std Dev	0.0000																																																																																																																																						
Rel Std Dev(%)	0.0000																																																																																																																																						
Air Blank	0.000	13:43																																																																																																																																					
Control Test	0.078	13:44																																																																																																																																					
Air Blank	0.000	13:45																																																																																																																																					
Control Test	0.078	13:45																																																																																																																																					
Air Blank	0.000	13:46																																																																																																																																					
Control Test	0.078	13:46																																																																																																																																					
Air Blank	0.000	13:47																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0780																																																																																																																																						
Std Dev	0.0000																																																																																																																																						
Rel Std Dev(%)	0.0000																																																																																																																																						
Air Blank	0.000	13:49																																																																																																																																					
Control Test	0.196	13:50																																																																																																																																					
Air Blank	0.000	13:51																																																																																																																																					
Control Test	0.196	13:51																																																																																																																																					
Air Blank	0.000	13:52																																																																																																																																					
Control Test	0.195	13:52																																																																																																																																					
Air Blank	0.000	13:53																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.1957																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.2951																																																																																																																																						
Air Blank	0.000	13:54																																																																																																																																					
Control Test	0.079	13:54																																																																																																																																					
Air Blank	0.000	13:55																																																																																																																																					
Control Test	0.079	13:55																																																																																																																																					
Air Blank	0.000	13:56																																																																																																																																					
Control Test	0.078	13:56																																																																																																																																					
Air Blank	0.000	13:56																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0787																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.7339																																																																																																																																						
Operator's Signature <i>TK</i>	Operator's Signature <i>TK</i>	Operator's Signature <i>TK</i>	Operator's Signature <i>TK</i>																																																																																																																																				

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: HIALEAH PD
Time of Inspection: 13:09

Date of Inspection: 07/15/2024

Serial Number: 80-002462
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:01923080A3 Exp: 02/05/2025
0.000	0.048	0.077	0.195	0.079
0.000	0.049	0.078	0.195	0.079
0.000	0.050	0.077	0.195	0.080
0.000	0.049	0.078	0.197	0.079
0.000	0.049	0.078	0.195	0.080
0.000	0.049	0.078	0.195	0.080
0.000	0.050	0.078	0.196	0.080
0.000	0.049	0.078	0.196	0.080
0.000	0.049	0.078	0.196	0.079
0.000	0.050	0.078	0.195	0.080

Standard Deviations	0.0006	0.0004	0.0007	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0005 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



TAYLOR D GUTSCHOW

Signature and Printed Name

07/15/2024
Date

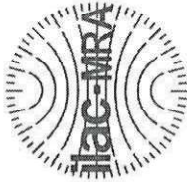

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-002462, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-002462</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>HIALEAH PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>07/15/2024</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>13:09</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

07/15/2024

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality



INSTRUMENT PROCESSING SHEET

Agency Hialeah PDS/N 80-002462Florida Department of
Law EnforcementDate In 05/03/2024 DI Completion Date _____☐ Ship ☐ P/U ☐ H/D ☒ CMI ☐ EE

Intake By <u>TDG</u> Date <u>05/17/2024</u>		Quality Checks By <u>TDG</u> Date <u>05/17/2024</u>		Flow Calibration By _____ Date _____																																																													
☒ Annual ☐ Registration ☒ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>213</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP104</u> 32 mm <u>0.152</u> (.139 - .169) 36 mm <u>0.167</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.503</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP6287</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP6288</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG222203 08/10/2024</td></tr></tbody></table>		Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K 03/29/2025	0.080	MP6287	202303L 03/29/2025	0.200	MP6288	202304C 04/05/2025	0.080 DGS	N/A	AG222203 08/10/2024	Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____ _____																																														
Simulator	Serial #	Lot #/Exp																																																															
0.050	MP6286	202303K 03/29/2025																																																															
0.080	MP6287	202303L 03/29/2025																																																															
0.200	MP6288	202304C 04/05/2025																																																															
0.080 DGS	N/A	AG222203 08/10/2024																																																															
Calibration Adjustment By _____		Department Inspection By _____																																																															
Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Barometric Pressure ID# _____ Gauge _____ Instrument _____ Mouth Alcohol Solution Lot # _____ Acetone Stock Solution Lot # _____ <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td></td></tr><tr><td>Interferent</td><td></td></tr><tr><td>0.050</td><td></td></tr><tr><td>0.080</td><td></td></tr><tr><td>0.200</td><td></td></tr></tbody></table> Attachments ☐ Form 41 ☒ Stability Checks ☐ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other <u>Extra Stabilities / Form 51</u>				Simulator	Serial Number	0.000		Interferent		0.050		0.080		0.200	
Simulator	Serial #	Lot #	Expiration																																																														
0.000		N/A	N/A																																																														
0.040																																																																	
0.100																																																																	
0.200																																																																	
0.300																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial #	Lot #	Expiration																																																														
0.050																																																																	
0.080																																																																	
0.200																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial Number																																																																
0.000																																																																	
Interferent																																																																	
0.050																																																																	
0.080																																																																	
0.200																																																																	
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☐ Instrument Complies with Chapter 11D-8, FAC ☒ Instrument Does Not Comply with Chapter 11D-8, FAC ☐ Return to/Place into Evidentiary Use ☒ Remain Out of Evidentiary Use ☐ Conduct an Agency Inspection Before Evidentiary Use Tech Review / Date _____ Admin Review / Date _____																																																															

Stability Checks

0.05g/210L 0.047 to 0.053 ✓	0.08g/210L 0.077 to 0.083 ✓	0.20g/210L 0.194 to 0.206 ✓	DGS 0.08g/210L 0.077 to 0.083 ✓ ≤0.003 of Wet ✓																																																																																																																																																
HIALEAH PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-002462 05/17/2024 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:17</td></tr> <tr><td>Control Test</td><td>0.050</td><td>12:18</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:18</td></tr> <tr><td>Control Test</td><td>0.049</td><td>12:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:19</td></tr> <tr><td>Control Test</td><td>0.049</td><td>12:20</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:20</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0493</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.1703</td><td></td></tr> </tbody> </table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	12:17	Control Test	0.050	12:18	Air Blank	0.000	12:18	Control Test	0.049	12:19	Air Blank	0.000	12:19	Control Test	0.049	12:20	Air Blank	0.000	12:20	Control Test Stats			Average	0.0493		Std Dev	0.0006		Rel Std Dev(%)	1.1703		HIALEAH PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-002462 05/17/2024 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:24</td></tr> <tr><td>Control Test</td><td>0.077</td><td>12:25</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:25</td></tr> <tr><td>Control Test</td><td>0.078</td><td>12:26</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:26</td></tr> <tr><td>Control Test</td><td>0.079</td><td>12:27</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:27</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0780</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.2821</td><td></td></tr> </tbody> </table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	12:24	Control Test	0.077	12:25	Air Blank	0.000	12:25	Control Test	0.078	12:26	Air Blank	0.000	12:26	Control Test	0.079	12:27	Air Blank	0.000	12:27	Control Test Stats			Average	0.0780		Std Dev	0.0010		Rel Std Dev(%)	1.2821		HIALEAH PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-002462 05/17/2024 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:30</td></tr> <tr><td>Control Test</td><td>0.195</td><td>12:31</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:32</td></tr> <tr><td>Control Test</td><td>0.195</td><td>12:32</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:33</td></tr> <tr><td>Control Test</td><td>0.195</td><td>12:33</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:34</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1950</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	12:30	Control Test	0.195	12:31	Air Blank	0.000	12:32	Control Test	0.195	12:32	Air Blank	0.000	12:33	Control Test	0.195	12:33	Air Blank	0.000	12:34	Control Test Stats			Average	0.1950		Std Dev	0.0000		Rel Std Dev(%)	0.0000		DLS HIALEAH PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-002462 05/17/2024 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:07</td></tr> <tr><td>Control Test</td><td>0.079</td><td>12:07</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:08</td></tr> <tr><td>Control Test</td><td>0.078</td><td>12:08</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:08</td></tr> <tr><td>Control Test</td><td>0.079</td><td>12:09</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:09</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0787</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7339</td><td></td></tr> </tbody> </table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	12:07	Control Test	0.079	12:07	Air Blank	0.000	12:08	Control Test	0.078	12:08	Air Blank	0.000	12:08	Control Test	0.079	12:09	Air Blank	0.000	12:09	Control Test Stats			Average	0.0787		Std Dev	0.0006		Rel Std Dev(%)	0.7339	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:17																																																																																																																																																	
Control Test	0.050	12:18																																																																																																																																																	
Air Blank	0.000	12:18																																																																																																																																																	
Control Test	0.049	12:19																																																																																																																																																	
Air Blank	0.000	12:19																																																																																																																																																	
Control Test	0.049	12:20																																																																																																																																																	
Air Blank	0.000	12:20																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0493																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	1.1703																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:24																																																																																																																																																	
Control Test	0.077	12:25																																																																																																																																																	
Air Blank	0.000	12:25																																																																																																																																																	
Control Test	0.078	12:26																																																																																																																																																	
Air Blank	0.000	12:26																																																																																																																																																	
Control Test	0.079	12:27																																																																																																																																																	
Air Blank	0.000	12:27																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0780																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	1.2821																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:30																																																																																																																																																	
Control Test	0.195	12:31																																																																																																																																																	
Air Blank	0.000	12:32																																																																																																																																																	
Control Test	0.195	12:32																																																																																																																																																	
Air Blank	0.000	12:33																																																																																																																																																	
Control Test	0.195	12:33																																																																																																																																																	
Air Blank	0.000	12:34																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1950																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:07																																																																																																																																																	
Control Test	0.079	12:07																																																																																																																																																	
Air Blank	0.000	12:08																																																																																																																																																	
Control Test	0.078	12:08																																																																																																																																																	
Air Blank	0.000	12:08																																																																																																																																																	
Control Test	0.079	12:09																																																																																																																																																	
Air Blank	0.000	12:09																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0787																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7339																																																																																																																																																		

0.200
ARS

Extra Stabilities

HIALEAH PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-002462
05/17/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:35
Control Test	0.194	12:35
Air Blank	0.000	12:36
Control Test	0.193	12:37
Air Blank	0.000	12:37
Control Test	0.194	12:38
Air Blank	0.000	12:38
Control Test	0.191	12:39
Air Blank	0.000	12:40
Control Test	0.196	12:40
Air Blank	0.000	12:41
Control Test	0.190	12:41
Air Blank	0.000	12:42
Control Test	0.189	12:43
Air Blank	0.000	12:43
Control Test	0.192	12:44
Air Blank	0.000	12:44
Control Test	0.195	12:45
Air Blank	0.000	12:46
Control Test	0.189	12:46
Air Blank	0.000	12:47
Control Test	0.191	12:47
Air Blank	0.000	12:48
Control Test	0.190	12:49
Air Blank	0.000	12:49
Control Test	0.194	12:50
Air Blank	0.000	12:50
Control Test	0.194	12:51
Air Blank	0.000	12:52
Control Test	0.188	12:52
Air Blank	0.000	12:53
Control Test	0.192	12:53
Air Blank	0.000	12:54
Control Test	0.191	12:55
Air Blank	0.000	12:55
Control Test	0.188	12:56
Air Blank	0.000	12:56
Control Test	0.189	12:57
Air Blank	0.000	12:58
Control Test	0.189	12:58
Air Blank	0.000	12:59
Control Test Stats		
Average	0.1914	
Std Dev	0.0025	
Rel Std Dev(%)	1.2848	

During the Quality Check Stabilities, it was noticed The values would continually bounce around without stabilizing. Instrument was originally sent to repair for imprecise 0.20 measurements. Will return to CMI.

ML
05/17/2024

ML
Operator's Signature

Return Material Authorization

Ship to: ☒ CMI, Inc.

☐ Enforcement Electronics

Shipment to repair facility authorized by: Det. Cesar Correa on 5/17/2024

Items Returned: Instrument ☒ Supplies ☐ Other ☐ Describe: _____

Instrument Model: Intoxilyzer 8000 Serial Number: 80-002462

Bill To Address:

Hialeah Police Department

Attn: Cesar Correa

Ship to Address:

Florida Department of Law Enforcement

Fort Myers Regional Operations Center

Attn: Taylor Gutschow

4700 Terminal Drive, Suite 1

Fort Myers, FL 33907

Reason for Return:

Instrument was sent to CMI in 2023 for imprecise measurements at the 0.20 concentration that could not be corrected by optical cal adjust. Instrument returned from CMI (work order 408255) but continues to have bouncing values. Per the above work order, CMI performed a calibration adjustment on the instrument, but it already requires another calibration adjustment. See Extra Stabilities. Records were uploaded by FDLE.

Please choose one of the following options:

☐ 1. I _____, authorize all repairs.

☐ 2. I _____, authorize repairs up to \$_____.

☒ 3. I require an estimate **BEFORE** any repairs will be authorized and/ or conducted.

Please contact: Name: Det. Cesar Correa

Phone #: 954-790-2809

Email: ccorrea@hialeahfl.gov

ATP Contact Name: Taylor Gutschow

ATP Email: TaylorGutschow@fdle.state.fl.us