



INSTRUMENT PROCESSING SHEET

Agency Hillsborough County SO

S/N 80-005113

Florida Department of
Law Enforcement

Date In 4/12/2024

DI Completion Date 5/17/2024

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

5/30/2024 PN 5/30/2024

Intake	By BS	Date 4/12/2024	Quality Checks	By DA	Date 4/15/2024	Flow Calibration	By	Date																																																
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: Breath tube rubber tape is worn against brace. Instrument takes several minutes and several attempts to turn on DA Barometric pressure out of range needs calibration adjustment DA 4/15/24			☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value 200 ☒ Flow Verification (L/s) Flow Column # ATP105 32 mm 0.160 (.139 - .169) 36 mm 0.171 (.156 - .190) 53 mm 0.242 (.228 - .278) 103 mm 0.500 (.447 - .547) ☒ Barometric Pressure Check Gauge ID # 30793 ☒ Stability Checks			Flow Column # ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value ☐ Post Calibration Verification (L/s) Flow Column # 32 mm (.139 - .169) 36 mm (.156 - .190) 53 mm (.228 - .278) 103 mm (.447 - .547)																																																		
			<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP5089</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP5090</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>06723080A5 04/05/2025</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP5088	202303K 03/29/2025	0.080	MP5089	202303L 03/29/2025	0.200	MP5090	202304C 04/05/2025	0.080 DGS	N/A	06723080A5 04/05/2025	<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP5089</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP5090</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>06723080A5 04/05/2025</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP5088	202303K 03/29/2025	0.080	MP5089	202303L 03/29/2025	0.200	MP5090	202304C 04/05/2025	0.080 DGS	N/A	06723080A5 04/05/2025	<table border="1"><thead><tr><th colspan="2">Maintenance By</th><th>Date</th></tr></thead><tbody><tr><td>☐ Battery Replacement</td><td></td><td></td></tr><tr><td>☐ Dry Gas Regulator Replacement</td><td></td><td></td></tr><tr><td>☐ Breath Tube Replacement</td><td></td><td></td></tr><tr><td>☐ Other</td><td></td><td></td></tr></tbody></table>			Maintenance By		Date	☐ Battery Replacement			☐ Dry Gas Regulator Replacement			☐ Breath Tube Replacement			☐ Other		
Simulator	Serial #	Lot #/Exp																																																						
0.050	MP5088	202303K 03/29/2025																																																						
0.080	MP5089	202303L 03/29/2025																																																						
0.200	MP5090	202304C 04/05/2025																																																						
0.080 DGS	N/A	06723080A5 04/05/2025																																																						
Simulator	Serial #	Lot #/Exp																																																						
0.050	MP5088	202303K 03/29/2025																																																						
0.080	MP5089	202303L 03/29/2025																																																						
0.200	MP5090	202304C 04/05/2025																																																						
0.080 DGS	N/A	06723080A5 04/05/2025																																																						
Maintenance By		Date																																																						
☐ Battery Replacement																																																								
☐ Dry Gas Regulator Replacement																																																								
☐ Breath Tube Replacement																																																								
☐ Other																																																								
Calibration Adjustment By PN			Department Inspection By PN																																																					
Barometric Pressure Gauge 1009 ID # 28421			Barometric Pressure ID# 28427																																																					
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td>MP6294</td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td>MP6295</td><td>23400</td><td>10/24/2025</td></tr><tr><td>0.100</td><td>MP6296</td><td>23450</td><td>12/5/2025</td></tr><tr><td>0.200</td><td>MP6297</td><td>24080</td><td>2/13/2026</td></tr><tr><td>0.300</td><td>MP6298</td><td>23410</td><td>11/1/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG310901</td><td>04/19/2025</td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.000	MP6294	N/A	N/A	0.040	MP6295	23400	10/24/2025	0.100	MP6296	23450	12/5/2025	0.200	MP6297	24080	2/13/2026	0.300	MP6298	23410	11/1/2025	0.080 DGS	N/A	AG310901	04/19/2025	Gauge 1007 /1014 Instrument 1009 /1016																									
Simulator	Serial #	Lot #	Expiration																																																					
0.000	MP6294	N/A	N/A																																																					
0.040	MP6295	23400	10/24/2025																																																					
0.100	MP6296	23450	12/5/2025																																																					
0.200	MP6297	24080	2/13/2026																																																					
0.300	MP6298	23410	11/1/2025																																																					
0.080 DGS	N/A	AG310901	04/19/2025																																																					
			Mouth Alcohol Solution Lot # 2024-A																																																					
			Acetone Stock Solution Lot # 2023-B																																																					
			<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5086 /MP6294</td></tr><tr><td>Interferent</td><td>MP5087 /MP6295</td></tr><tr><td>0.050</td><td>MP5088 /MP6296</td></tr><tr><td>0.080</td><td>MP5089 /MP6297</td></tr><tr><td>0.200</td><td>MP5090 /MP6298</td></tr></tbody></table>					Simulator	Serial Number	0.000	MP5086 /MP6294	Interferent	MP5087 /MP6295	0.050	MP5088 /MP6296	0.080	MP5089 /MP6297	0.200	MP5090 /MP6298																																					
Simulator	Serial Number																																																							
0.000	MP5086 /MP6294																																																							
Interferent	MP5087 /MP6295																																																							
0.050	MP5088 /MP6296																																																							
0.080	MP5089 /MP6297																																																							
0.200	MP5090 /MP6298																																																							
☒ Post Calibration Adjustment Stability Checks			Attachments																																																					
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202303K</td><td>03/29/2025</td></tr><tr><td>0.080</td><td>MP5089</td><td>202303L</td><td>03/29/2025</td></tr><tr><td>0.200</td><td>MP5090</td><td>202304C</td><td>04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>06723080A5</td><td>04/05/2025</td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.050	MP5088	202303K	03/29/2025	0.080	MP5089	202303L	03/29/2025	0.200	MP5090	202304C	04/05/2025	0.080 DGS	N/A	06723080A5	04/05/2025	<table border="1"><tbody><tr><td>☒ Form 41 X2</td><td>☒ Post-Stability Checks</td></tr><tr><td>☒ Stability Checks</td><td>☐ Flow Calibration</td></tr><tr><td>☒ Calibration Certificate</td><td>☐ Form 40</td></tr><tr><td>☒ Calibration Adjustment</td><td>☒ Other Dry Control Test</td></tr></tbody></table>					☒ Form 41 X2	☒ Post-Stability Checks	☒ Stability Checks	☐ Flow Calibration	☒ Calibration Certificate	☐ Form 40	☒ Calibration Adjustment	☒ Other Dry Control Test																					
Simulator	Serial #	Lot #	Expiration																																																					
0.050	MP5088	202303K	03/29/2025																																																					
0.080	MP5089	202303L	03/29/2025																																																					
0.200	MP5090	202304C	04/05/2025																																																					
0.080 DGS	N/A	06723080A5	04/05/2025																																																					
☒ Form 41 X2	☒ Post-Stability Checks																																																							
☒ Stability Checks	☐ Flow Calibration																																																							
☒ Calibration Certificate	☐ Form 40																																																							
☒ Calibration Adjustment	☒ Other Dry Control Test																																																							
Notes/Suggested Service: Accidentally performed a Dry Control Test during the post calibration adjustment stability checks. see 'Other' attachment. PN 5/17/2024 Tech Review: The same DGS lot was used during calibration adjustment and inspection on 5/17. Inspection repeated using a different DGS lot on 5/30. Notified AI of instrument power issue described in the intake notes. AI will monitor instrument. PN 5/30/2024			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2024.06.05 10:57:20 -04'00' Shayla Platt Digitally signed by Shayla Platt Date: 2024.06.05 12:39:26 -04'00'																																																					
Tech Review / Date			Admin Review / Date																																																					

Stability Checks 80-005113 DA 4/15/24

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005113
04/15/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:50
Control Test	0.050	10:51
Air Blank	0.000	10:51
Control Test	0.050	10:52
Air Blank	0.000	10:53
Control Test	0.049	10:53
Air Blank	0.000	10:54
Control Test Stats		
Average	0.0497	
Std Dev	0.0006	
Rel Std Dev(%)	1.1625	


Operator's Signature

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005113
04/15/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:55
Control Test	0.079	10:56
Air Blank	0.000	10:57
Control Test	0.079	10:57
Air Blank	0.000	10:58
Control Test	0.079	10:58
Air Blank	0.000	10:59
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator's Signature

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005113
04/15/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:03
Control Test	0.198	11:04
Air Blank	0.000	11:04
Control Test	0.197	11:05
Air Blank	0.000	11:05
Control Test	0.197	11:06
Air Blank	0.000	11:07
Control Test Stats		
Average	0.1973	
Std Dev	0.0006	
Rel Std Dev(%)	0.2926	


Operator's Signature

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005113
04/15/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:08
Control Test	0.081	11:08
Air Blank	0.000	11:09
Control Test	0.084	11:09
Air Blank	0.000	11:09
Control Test	0.080	11:10
Air Blank	0.000	11:10
Control Test Stats		
Average	0.0817	
Std Dev	0.0021	
Rel Std Dev(%)	2.5490	

DGS


Operator's Signature

OPTICAL CALIBRATION ADJUSTMENT



HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005113
05/16/2024 08:26:09

Auto Calibration

pg 1 of 2

<<<< 3um >>>>			<<<< 9um >>>>		

Solution = 0.000 g/210L or 0.0000 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	0.1240	(-0.0150)	0.1130	(-0.0190)	
Sample #2	0.0830	(0.0090)	0.0920	(-0.0020)	
Sample #3	0.0870	(-0.0010)	0.1050	(-0.0080)	
Sample #4	0.0820	(0.0270)	0.0910	(0.0080)	
Avg % Abs	0.0840	(0.0117)	0.0960	(-0.0007)	
STD DEV	0.0026	(0.0142)	0.0078	(0.0081)	
REL STD DEV	3.150	(121.622)	8.136	(1212.436)	

Solution = 0.040 g/210L or 0.1905 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	0.9190	(-0.0090)	1.6090	(-0.0060)	
Sample #2	0.8860	(0.0070)	1.5910	(0.0170)	
Sample #3	0.8740	(0.0190)	1.5720	(0.0190)	
Sample #4	0.8600	(0.0330)	1.5700	(0.0240)	
Avg % Abs	0.8733	(0.0197)	1.5777	(0.0200)	
STD DEV	0.0130	(0.0130)	0.0116	(0.0036)	
REL STD DEV	1.490	(66.167)	0.735	(18.028)	

Solution = 0.100 g/210L or 0.4762 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	2.0740	(-0.0080)	3.7680	(-0.0040)	
Sample #2	2.0250	(0.0130)	3.7260	(0.0240)	
Sample #3	2.0400	(0.0010)	3.7210	(0.0240)	
Sample #4	2.0330	(0.0250)	3.7130	(0.0390)	
Avg % Abs	2.0327	(0.0130)	3.7200	(0.0290)	
STD DEV	0.0075	(0.0120)	0.0066	(0.0087)	
REL STD DEV	0.369	(92.308)	0.176	(29.863)	

Solution = 0.200 g/210L or 0.9524 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	3.9050	(0.0000)	7.1430	(-0.0070)	
Sample #2	3.8710	(0.0420)	7.0720	(0.0570)	
Sample #3	3.8630	(0.0500)	7.0530	(0.0760)	
Sample #4	3.8980	(0.0320)	7.0860	(0.0570)	
Avg % Abs	3.8773	(0.0413)	7.0703	(0.0633)	
STD DEV	0.0183	(0.0090)	0.0166	(0.0110)	
REL STD DEV	0.473	(21.819)	0.234	(17.321)	

Solution = 0.300 g/210L or 1.4286 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	5.7250	(-0.0030)	10.3580	(0.0020)	
Sample #2	5.7210	(0.0240)	10.3130	(0.0680)	
Sample #3	5.7240	(0.0320)	10.2930	(0.0810)	
Sample #4	5.6940	(0.0540)	10.2630	(0.1040)	
Avg % Abs	5.7130	(0.0367)	10.2897	(0.0843)	
STD DEV	0.0165	(0.0155)	0.0252	(0.0182)	
REL STD DEV	0.289	(42.368)	0.245	(21.617)	

OPTICAL CALIBRATION ADJUSTMENT

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005113
05/16/2024 08:26:09

Auto Calibration

pg 2 of 2

<<<< 3um >>>>			<<<< 9um >>>>		
-----			-----		
Zero Order Coef	-219.29			-132.22	
First Order Coef	2423.94			1275.03	
Second Order Coef	20.51			12.34	
-----			-----		
Act	Fit	Residual	Act	Fit	Residual
(g/210L)	(g/210L)	(g/210L)	(g/210L)	(g/210L)	(g/210L)
0.000	-0.000	0.0003	0.000	-0.000	0.0002
0.040	0.040	-0.0002	0.040	0.040	-0.0001
0.100	0.101	-0.0006	0.100	0.100	-0.0004
0.200	0.199	0.0008	0.200	0.199	0.0005
0.300	0.300	-0.0003	0.300	0.300	-0.0002
-----			-----		

<<<< 3um >>>>		<<<< 9um >>>>	
-----		-----	
Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1			
Sample			
Sample #1	3214.00		3410.00
Sample #2	3186.00		3434.00
Sample #3	3189.00		3423.00
Sample #4	3172.00		3427.00
Avg	3182.3333		3428.0000
STD DEV	9.0738		5.5678
REL STD DEV	0.285		0.162
H2O adjust (mg/l*10k)	627		381

Barometric Pressure = 1009

*****CALIBRATION SUCCESSFUL*****

POST STABILITY CHECKS

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000
05/16/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:01
Control Test	0.050	11:02
Air Blank	0.000	11:02
Control Test	0.050	11:03
Air Blank	0.000	11:03
Control Test	0.050	11:04
Air Blank	0.000	11:05
Control Test Stats		
Average	0.0500	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000
05/16/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:06
Control Test	0.080	11:07
Air Blank	0.000	11:07
Control Test	0.080	11:08
Air Blank	0.000	11:09
Control Test	0.080	11:09
Air Blank	0.000	11:10
Control Test Stats		
Average	0.0800	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000
05/16/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:12
Control Test	0.200	11:12
Air Blank	0.000	11:13
Control Test	0.199	11:13
Air Blank	0.000	11:14
Control Test	0.198	11:15
Air Blank	0.000	11:15
Control Test Stats		
Average	0.1990	
Std Dev	0.0010	
Rel Std Dev(%)	0.5025	

Operator's Signature

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000
05/16/2024
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:19
Control Test	0.080	11:19
Air Blank	0.000	11:19
Control Test	0.080	11:20
Air Blank	0.000	11:20
Control Test	0.080	11:20
Air Blank	0.000	11:21
Control Test Stats		
Average	0.0800	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

DGS

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: HILLSBOROUGH CO SO

Serial Number: 80-005113

Time of Inspection: 11:36

Date of Inspection: 05/30/2024

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: 202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#: 202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#: 202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: 06723080A5 Exp: 04/05/2025
0.000	0.050	0.080	0.199	0.079
0.000	0.050	0.080	0.199	0.079
0.000	0.050	0.080	0.199	0.079
0.000	0.050	0.081	0.199	0.079
0.000	0.049	0.080	0.199	0.079
0.000	0.049	0.080	0.198	0.079
0.000	0.050	0.080	0.199	0.078
0.000	0.050	0.080	0.198	0.079
0.000	0.050	0.080	0.198	0.079
0.000	0.050	0.080	0.198	0.079


Standard Deviations	0.0004	0.0003	0.0005	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply () with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

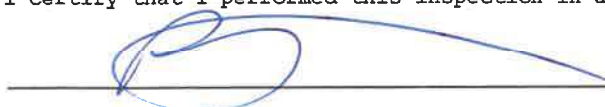

PHIL NICODEMO

Signature and Printed Name

05/30/2024

Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road.
Suite B1032
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-005113, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-005113	UNCERTAINTY* $\pm$	
Owning Agency:	HILLSBOROUGH CO SO	0.050 g/ 210 L	0.004
Calibration Date:	05/30/2024	0.080 g/ 210 L	0.004
Calibration Time:	11:36	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within $\pm$ 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

05/30/2024

Date

PHIL NICODEMO,

Department Inspector

FDLE/ATP Form 69 March 2022

Issuing Authority: Alcohol Testing Program


Service • Integrity • Respect • Quality

DRY CONTROL TEST

HILLSBOROUGH CO SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005113
05/16/2024
Software: 8100.27

: DRY CONTROL TEST :

Test	g/210L	Time
Air Blank	0.000	10:59
Control Test	0.000	11:00
Air Blank	0.000	11:00

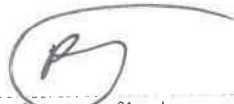

Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: HILLSBOROUGH CO SO
Time of Inspection: 13:58

Date of Inspection: 05/17/2024

Serial Number: 80-005113
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG310901 Exp: 04/19/2025
0.000	0.050	0.080	0.199	0.080
0.000	0.050	0.079	0.199	0.080
0.000	0.050	0.079	0.199	0.080
0.000	0.050	0.080	0.199	0.079
0.000	0.050	0.079	0.199	0.080
0.000	0.050	0.080	0.199	0.080
0.000	0.050	0.080	0.199	0.079
0.000	0.051	0.080	0.199	0.079
0.000	0.051	0.080	0.199	0.079
0.000	0.051	0.080	0.199	0.079


Standard Deviations	0.0004	0.0004	0.0000	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

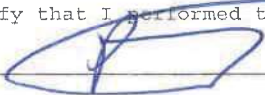

PHIL NICODEMO

Signature and Printed Name

05/17/2024
Date