



INSTRUMENT PROCESSING SHEET

Agency Miami Dade PDS/N 80-005549Florida Department of
Law EnforcementDate In 02/06/2024 DI Completion Date 02/12/2024☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By TDG	Date	Quality Checks By TDG	Date	Flow Calibration By	Date																																																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>Missing one screw on the DGS regulator.</u>	<u>02/08/2024</u>	☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>183</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP106</u> 32 mm <u>0.160</u> (.139 - .169) 36 mm <u>0.175</u> (.156 - .190) 53 mm <u>0.242</u> (.228 - .278) 103 mm <u>0.496</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5094</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP5095</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP5096</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>01923080A3 02/05/2025</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP5094	202303K 03/29/2025	0.080	MP5095	202303L 03/29/2025	0.200	MP5096	202304C 04/05/2025	0.080 DGS	N/A	01923080A3 02/05/2025	<u>02/12/2024</u>	Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By TDG <u>02/08/2024</u> ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☒ Other <u>Added screw to regulator.</u>																																														
Simulator	Serial #	Lot #/Exp																																																															
0.050	MP5094	202303K 03/29/2025																																																															
0.080	MP5095	202303L 03/29/2025																																																															
0.200	MP5096	202304C 04/05/2025																																																															
0.080 DGS	N/A	01923080A3 02/05/2025																																																															
Calibration Adjustment By _____		Department Inspection By TDG _____																																																															
Barometric Pressure Gauge _____ ID # _____ <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Barometric Pressure ID# <u>26932</u> Gauge <u>1012</u> Instrument <u>1007</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2022-B</u> <table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5092</td></tr><tr><td>Interferent</td><td>MP5093</td></tr><tr><td>0.050</td><td>MP5094</td></tr><tr><td>0.080</td><td>MP5095</td></tr><tr><td>0.200</td><td>MP5096</td></tr></tbody></table> Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☐ Other _____ ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Phil Nicodemo Digitally signed by Phil Nicodemo Date: 2024.02.13 11:36:07 +05'00' <u>Shayla Platt</u> Digitally signed by Shayla Platt Date: 2024.02.16 10:26:59 -05'00' Tech Review / Date _____ Admin Review / Date _____				Simulator	Serial Number	0.000	MP5092	Interferent	MP5093	0.050	MP5094	0.080	MP5095	0.200	MP5096
Simulator	Serial #	Lot #	Expiration																																																														
0.000		N/A	N/A																																																														
0.040																																																																	
0.100																																																																	
0.200																																																																	
0.300																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial #	Lot #	Expiration																																																														
0.050																																																																	
0.080																																																																	
0.200																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial Number																																																																
0.000	MP5092																																																																
Interferent	MP5093																																																																
0.050	MP5094																																																																
0.080	MP5095																																																																
0.200	MP5096																																																																
Notes/Suggested Service: _____ _____ _____ _____ _____ _____																																																																	

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
✓	✓	✓	✓
MIAMI DDOE PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-005549 02/12/2024 Software: 8100.27 Test g/210L Time Air Blank 0.000 10:24 Control Test 0.050 10:25 Air Blank 0.000 10:25 Control Test 0.049 10:26 Air Blank 0.000 10:26 Control Test 0.050 10:27 Air Blank 0.000 10:28 Control Test Stats Average 0.0497 Std Dev 0.0006 Rel Std Dev(%) 1.1625 Operator's Signature 	MIAMI DDOE PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-005549 02/12/2024 Software: 8100.27 Test g/210L Time Air Blank 0.000 10:31 Control Test 0.079 10:32 Air Blank 0.000 10:32 Control Test 0.079 10:33 Air Blank 0.000 10:34 Control Test 0.079 10:34 Air Blank 0.000 10:35 Control Test Stats Average 0.0790 Std Dev 0.0000 Rel Std Dev(%) 0.0000 Operator's Signature 	MIAMI DDOE PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-005549 02/12/2024 Software: 8100.27 Test g/210L Time Air Blank 0.000 10:39 Control Test 0.198 10:39 Air Blank 0.000 10:40 Control Test 0.197 10:41 Air Blank 0.000 10:41 Control Test 0.197 10:42 Air Blank 0.000 10:42 Control Test Stats Average 0.1973 Std Dev 0.0006 Rel Std Dev(%) 0.2926 Operator's Signature 	MIAMI DDOE PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-005549 02/12/2024 Software: 8100.27 Test g/210L Time Air Blank 0.000 10:44 Control Test 0.078 10:44 Air Blank 0.000 10:45 Control Test 0.079 10:45 Air Blank 0.000 10:45 Control Test 0.079 10:46 Air Blank 0.000 10:46 Control Test Stats Average 0.0787 Std Dev 0.0006 Rel Std Dev(%) 0.7339 Operator's Signature 

065

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI DADE PD
Time of Inspection: 13:15

Date of Inspection: 02/12/2024

Serial Number: 80-005549
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test* (g/210L) Lot#:01923080A3 Exp: 02/05/2025
0.000	0.050	0.080	0.197	0.078
0.000	0.050	0.080	0.198	0.078
0.000	0.050	0.079	0.197	0.079
0.000	0.050	0.080	0.197	0.078
0.000	0.050	0.080	0.196	0.079
0.000	0.050	0.079	0.197	0.079
0.000	0.050	0.079	0.197	0.079
0.000	0.050	0.079	0.198	0.078
0.000	0.050	0.079	0.196	0.079
0.000	0.050	0.080	0.197	0.079

Standard Deviations	0.0000	0.0005	0.0006	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



Signature and Printed Name

TAYLOR D GUTSCHOW

02/12/2024
Date

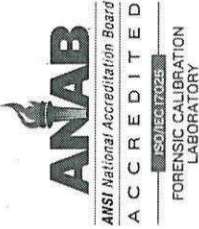

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-005549, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-005549	UNCERTAINTY* $\pm$	
Owning Agency:	MIAMI DADE PD	0.050 g/ 210 L	0.004
Calibration Date:	02/12/2024	0.080 g/ 210 L	0.004
Calibration Time:	13:15	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

02/12/2024

Date


TAYLOR D GUTSCHOW,

Department Inspector

Service • Integrity • Respect • Quality