



INSTRUMENT PROCESSING SHEET

Agency Miami Dade PD

S/N 80-005551

Florida Department of
Law Enforcement

Date In 04/16/2024 DI Completion Date 05/15/2024

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By TDG	Date 05/09/2024	Quality Checks	By TDG	Date 05/15/2024	Flow Calibration	By	Date																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: AI reports the instrument has lost the forms and agency codes. Missing back right foot on DGS shelf.			☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value 123 ☒ Flow Verification (L/s) Flow Column # ATP104 32 mm 0.148 (.139 - .169) 36 mm 0.164 (.156 - .190) 53 mm 0.234 (.228 - .278) 103 mm 0.492 (.447 - .547) ☒ Barometric Pressure Check Gauge ID # 26932 ☒ Stability Checks			Flow Column # ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value ☐ Post Calibration Verification (L/s) Flow Column # 32 mm (.139 - .169) 36 mm (.156 - .190) 53 mm (.228 - .278) 103 mm (.447 - .547)																														
			<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP6287</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP6288</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG222203 08/10/2024</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K 03/29/2025	0.080	MP6287	202303L 03/29/2025	0.200	MP6288	202304C 04/05/2025	0.080 DGS	N/A	AG222203 08/10/2024				Maintenance By TDG Date 05/14/2024 ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☒ Other CMI remotely reloaded forms and agency codes.												
Simulator	Serial #	Lot #/Exp																																		
0.050	MP6286	202303K 03/29/2025																																		
0.080	MP6287	202303L 03/29/2025																																		
0.200	MP6288	202304C 04/05/2025																																		
0.080 DGS	N/A	AG222203 08/10/2024																																		
Calibration Adjustment			By			Department Inspection By TDG																														
Barometric Pressure Gauge ID #			Barometric Pressure ID# 26932			Gauge 1013 Instrument 1012																														
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A						Mouth Alcohol Solution Lot # 2023-A Acetone Stock Solution Lot # 2023-B		
Simulator	Serial #	Lot #	Expiration																																	
0.000		N/A	N/A																																	
0.040																																				
0.100																																				
0.200																																				
0.300																																				
0.080 DGS	N/A																																			
☐ Post Calibration Adjustment Stability Checks						<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table>			Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288																
Simulator	Serial Number																																			
0.000	MP6284																																			
Interferent	MP6285																																			
0.050	MP6286																																			
0.080	MP6287																																			
0.200	MP6288																																			
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A						Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other Form 38 (Test)										
Simulator	Serial #	Lot #	Expiration																																	
0.050																																				
0.080																																				
0.200																																				
0.080 DGS	N/A																																			
Notes/Suggested Service: Checked the breath tube screen and replaced the o-rings on 5/9. Finished the rest of the Quality Checks on 5/15. (TDG)						☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																														
						Digitally signed by Shayla Platt Shayla Platt Date: 2024.05.16 10:36:18 -0400																														
			Phil Nicodemus Digitally signed by Phil Nicodemus Date: 2024.05.16 07:45:35 -0400			Tech Review / Date Admin Review / Date																														

0.05g/210L		0.08g/210L		0.20g/210L		DGS 0.08g/210L																																																																																																									
0.047 to 0.053		0.077 to 0.083		0.194 to 0.206		0.077 to 0.083																																																																																																									
MIAWI DAGE PD Intoxilyzer - Alcotest Analyzer Model 8000 SN 80-005551 05/15/2024 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>12:18</td></tr><tr><td>Control Test</td><td>0.050</td><td>12:18</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:19</td></tr><tr><td>Control Test</td><td>0.050</td><td>12:20</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:20</td></tr><tr><td>Control Test</td><td>0.050</td><td>12:21</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:21</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0500</td><td></td></tr><tr><td>Std Dev</td><td>0.0000</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr></tbody></table> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	12:18	Control Test	0.050	12:18	Air Blank	0.000	12:19	Control Test	0.050	12:20	Air Blank	0.000	12:20	Control Test	0.050	12:21	Air Blank	0.000	12:21	Control Test Stats			Average	0.0500		Std Dev	0.0000		Rel Std Dev(%)	0.0000		MIAWI DAGE PD Intoxilyzer - Alcotest Analyzer Model 8000 SN 80-005551 05/15/2024 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>12:24</td></tr><tr><td>Control Test</td><td>0.081</td><td>12:25</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:26</td></tr><tr><td>Control Test</td><td>0.080</td><td>12:26</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:27</td></tr><tr><td>Control Test</td><td>0.079</td><td>12:28</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:28</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0800</td><td></td></tr><tr><td>Std Dev</td><td>0.0010</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>1.2500</td><td></td></tr></tbody></table> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	12:24	Control Test	0.081	12:25	Air Blank	0.000	12:26	Control Test	0.080	12:26	Air Blank	0.000	12:27	Control Test	0.079	12:28	Air Blank	0.000	12:28	Control Test Stats			Average	0.0800		Std Dev	0.0010		Rel Std Dev(%)	1.2500		MIAWI DAGE PD Intoxilyzer - Alcotest Analyzer Model 8000 SN 80-005551 05/15/2024 Software: 8100.27 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>12:31</td></tr><tr><td>Control Test</td><td>0.203</td><td>12:31</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:32</td></tr><tr><td>Control Test</td><td>0.201</td><td>12:33</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:33</td></tr><tr><td>Control Test</td><td>0.202</td><td>12:34</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:35</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.2020</td><td></td></tr><tr><td>Std Dev</td><td>0.0010</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.4950</td><td></td></tr></tbody></table> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	12:31	Control Test	0.203	12:31	Air Blank	0.000	12:32	Control Test	0.201	12:33	Air Blank	0.000	12:33	Control Test	0.202	12:34	Air Blank	0.000	12:35	Control Test Stats			Average	0.2020		Std Dev	0.0010		Rel Std Dev(%)	0.4950		DGS 0.08g/210L ☒ 0.077 to 0.083 ☒ ≤0.003 of Wet ✓ Inadvertently printed to external printer. Will attach. 05/15/2024
Test	g/210L	Time																																																																																																													
Air Blank	0.000	12:18																																																																																																													
Control Test	0.050	12:18																																																																																																													
Air Blank	0.000	12:19																																																																																																													
Control Test	0.050	12:20																																																																																																													
Air Blank	0.000	12:20																																																																																																													
Control Test	0.050	12:21																																																																																																													
Air Blank	0.000	12:21																																																																																																													
Control Test Stats																																																																																																															
Average	0.0500																																																																																																														
Std Dev	0.0000																																																																																																														
Rel Std Dev(%)	0.0000																																																																																																														
Test	g/210L	Time																																																																																																													
Air Blank	0.000	12:24																																																																																																													
Control Test	0.081	12:25																																																																																																													
Air Blank	0.000	12:26																																																																																																													
Control Test	0.080	12:26																																																																																																													
Air Blank	0.000	12:27																																																																																																													
Control Test	0.079	12:28																																																																																																													
Air Blank	0.000	12:28																																																																																																													
Control Test Stats																																																																																																															
Average	0.0800																																																																																																														
Std Dev	0.0010																																																																																																														
Rel Std Dev(%)	1.2500																																																																																																														
Test	g/210L	Time																																																																																																													
Air Blank	0.000	12:31																																																																																																													
Control Test	0.203	12:31																																																																																																													
Air Blank	0.000	12:32																																																																																																													
Control Test	0.201	12:33																																																																																																													
Air Blank	0.000	12:33																																																																																																													
Control Test	0.202	12:34																																																																																																													
Air Blank	0.000	12:35																																																																																																													
Control Test Stats																																																																																																															
Average	0.2020																																																																																																														
Std Dev	0.0010																																																																																																														
Rel Std Dev(%)	0.4950																																																																																																														

MIAMI DADE PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005551
05/15/2024
Software: 8100.27

065

Test	g/210L	Time
Air Blank	0.000	12:06
Control Test	0.081	12:06
Air Blank	0.000	12:07
Control Test	0.081	12:07
Air Blank	0.000	12:08
Control Test	0.082	12:08
Air Blank	0.000	12:09
Control Test Stats		
Average	0.0813	
Std Dev	0.0006	
Rel Std Dev(%)	0.7099	



Operator's Signature

**FLORIDA DEPARTMENT OF LAW ENFORCEMENT
ALCOHOL TESTING PROGRAM
BREATH ALCOHOL TEST AFFIDAVIT**

Instrument Type: Intoxilyzer 8000
Instrument Registered To: MIAMI DADE PD
Instrument Serial Number: 80-005551 Software: 8100.27
Date of Test: 05/15/2024

Date of Last Agency Inspection: 04/11/2024

Observation Period Began: 11:30

Subject's Name: TEST TEST

DOB: 12/12/2000 Sex: M

The subject was observed for at-least twenty-minutes prior to the administration of the breath test to ensure that the subject did not take anything orally and did not regurgitate.

Results:	Test	g/210L	Time
	Diagnostics Check	OK	11:54
	Air Blank	0.000	11:55
	Control Test	0.081	11:55
	Air Blank	0.000	11:56
	Subject Sample #1	ABT*	11:56
	Air Blank	0.000	11:56

*Sequence Aborted

Test of forms and agency codes. Both successful.

ML

05/15/2024

Cylinder Lot: AG222203
Exp: 08/10/2024

State of Florida, County of _____,

Personally appeared before me the undersigned authority, who () is personally known to me or () produced _____ as identification, and who after being placed under oath, states:

I, TAYLOR D GUTSCHOW, hold a valid Breath Test Operator permit issued by the Florida Department of Law Enforcement, I administered the above breath test to the subject named above in accordance with Chapter 11D-8, Florida Administrative Code, and this form is a true and accurate report of that breath test.

Breath Test Operator: _____ Date: _____
Signature

Sworn to (or affirmed) before me this _____ day of _____, _____

Signature of Notary Public-State of Florida Printed Name of Notary Public-State of Florida

Note: Pursuant to section 117.10, Florida Statutes, law enforcement officers, correctional officers, traffic accident investigation officers and traffic infraction enforcement officers are notaries public when engaged in the performance of official duties. In accordance with section 316.1934(5), F.S., this completed form is admissible without further authentication and is presumptive proof of the results herein. To be used in accordance with Section 316.1934(5), F.S., and in administrative proceedings pursuant to 322.2615, F.S.

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI DADE PD
Time of Inspection: 15:02

Date of Inspection: 05/15/2024

Serial Number: 80-005551
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG222203 Exp: 08/10/2024
0.000	0.050	0.080	0.202	0.080
0.000	0.050	0.080	0.201	0.079
0.000	0.050	0.080	0.201	0.079
0.000	0.049	0.080	0.201	0.079
0.000	0.049	0.080	0.201	0.079
0.000	0.050	0.080	0.202	0.079
0.000	0.050	0.080	0.202	0.079
0.000	0.050	0.080	0.202	0.079
0.000	0.050	0.080	0.202	0.079
0.000	0.050	0.081	0.202	0.079

Standard Deviations	0.0004	0.0003	0.0005	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

Inadvertently repeated the Minimum Sample Volume Check.
Both the original ^{TDL} check and the repeated check were passing.

TDL
05/15/2024

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Taylor D Gutschow

TAYLOR D GUTSCHOW

Signature and Printed Name

05/15/2024
Date

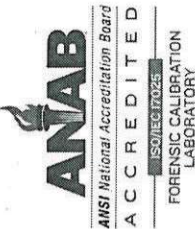

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-005551, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-005551</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>MIAMI DADE PD</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>05/15/2024</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>15:02</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

05/15/2024

Date

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality