



INSTRUMENT PROCESSING SHEET

Agency Monroe CSOS/N 80-007167Florida Department of
Law EnforcementDate In 03/13/2024 DI Completion Date 03/13/2024☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake By TDG _____ Date <u>03/13/2024</u>		Quality Checks By TDG _____ Date <u>03/13/2024</u>		Flow Calibration By _____ Date _____																																																															
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>204</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP106</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.230</u> (.228 - .278) 103 mm <u>0.496</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP6287</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP6288</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>01923080A3 02/05/2025</td></tr></tbody></table>		Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K 03/29/2025	0.080	MP6287	202303L 03/29/2025	0.200	MP6288	202304C 04/05/2025	0.080 DGS	N/A	01923080A3 02/05/2025	Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____ _____																																																
Simulator	Serial #	Lot #/Exp																																																																	
0.050	MP6286	202303K 03/29/2025																																																																	
0.080	MP6287	202303L 03/29/2025																																																																	
0.200	MP6288	202304C 04/05/2025																																																																	
0.080 DGS	N/A	01923080A3 02/05/2025																																																																	
Calibration Adjustment By _____		Department Inspection By TDG _____																																																																	
Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Barometric Pressure ID# <u>26932</u> Gauge <u>1019</u> Instrument <u>1018</u> Mouth Alcohol Solution Lot # <u>2023-A</u> Acetone Stock Solution Lot # <u>2022-B</u> <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table> Attachments <table border="1" style="width:100%"><tr><td>☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment</td><td>☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____</td></tr></table>				Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288	☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____
Simulator	Serial #	Lot #	Expiration																																																																
0.000		N/A	N/A																																																																
0.040																																																																			
0.100																																																																			
0.200																																																																			
0.300																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial #	Lot #	Expiration																																																																
0.050																																																																			
0.080																																																																			
0.200																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial Number																																																																		
0.000	MP6284																																																																		
Interferent	MP6285																																																																		
0.050	MP6286																																																																		
0.080	MP6287																																																																		
0.200	MP6288																																																																		
☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____																																																																		
Notes/Suggested Service: <u>Received a DSP Fail on startup. Cycled power on instrument and did not receive a DSP Fail. AI has not observed DSP Fails in the field and will monitor. (TDG)</u> _____ _____ _____ _____ _____		☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use <table border="1" style="width:100%"><tr><td>Phil Nicodemo Tech Review / Date _____</td><td>Digitally signed by Phil Nicodemo Date: 2024.03.13 14:20:07 -04'00' by Shayla Platt Date: 2024.03.13 15:48:13 -04'00' Admin Review / Date _____</td></tr></table>				Phil Nicodemo Tech Review / Date _____	Digitally signed by Phil Nicodemo Date: 2024.03.13 14:20:07 -04'00' by Shayla Platt Date: 2024.03.13 15:48:13 -04'00' Admin Review / Date _____																																																												
Phil Nicodemo Tech Review / Date _____	Digitally signed by Phil Nicodemo Date: 2024.03.13 14:20:07 -04'00' by Shayla Platt Date: 2024.03.13 15:48:13 -04'00' Admin Review / Date _____																																																																		

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083
✓	✓	✓	✓
✓	✓	✓	✓
MONROE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 03/13/2024 Software: 8100.27 SN 80-007167	MONROE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 03/13/2024 Software: 8100.27 SN 80-007167	MONROE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 03/13/2024 Software: 8100.27 SN 80-007167	MONROE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 03/13/2024 Software: 8100.27 SN 80-007167
Test g/210L Time Air Blank 0.000 10:04 Control Test 0.049 10:05 Air Blank 0.000 10:05 Control Test 0.049 10:06 Air Blank 0.000 10:07 Control Test 0.049 10:07 Air Blank 0.000 10:08 Control Test Stats Average 0.0490 Std Dev 0.0000 Rel Std Dev(%) 0.0000	Test g/210L Time Air Blank 0.000 10:11 Control Test 0.078 10:12 Air Blank 0.000 10:12 Control Test 0.079 10:13 Air Blank 0.000 10:13 Control Test 0.078 10:14 Air Blank 0.000 10:15 Control Test Stats Average 0.0783 Std Dev 0.0006 Rel Std Dev(%) 0.7370	Test g/210L Time Air Blank 0.000 10:18 Control Test 0.197 10:18 Air Blank 0.000 10:19 Control Test 0.196 10:19 Air Blank 0.000 10:20 Control Test 0.196 10:21 Air Blank 0.000 10:21 Control Test Stats Average 0.1963 Std Dev 0.0006 Rel Std Dev(%) 0.2941	Test g/210L Time Air Blank 0.000 09:48 Control Test 0.081 09:49 Air Blank 0.000 09:49 Control Test 0.080 09:50 Air Blank 0.000 09:50 Control Test 0.080 09:50 Air Blank 0.000 09:51 Control Test Stats Average 0.0803 Std Dev 0.0006 Rel Std Dev(%) 0.7187
Operator's Signature <i>ML</i>	Operator's Signature <i>ML</i>	Operator's Signature <i>ML</i>	Operator's Signature <i>ML</i>

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MONROE COUNTY SO
Time of Inspection: 12:32

Date of Inspection: 03/13/2024

Serial Number: 80-007167
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test* (g/210L) Lot#:01923080A3 Exp: 02/05/2025
0.000	0.049	0.078	0.197	0.081
0.000	0.049	0.079	0.197	0.081
0.000	0.050	0.078	0.197	0.080
0.000	0.049	0.079	0.197	0.081
0.000	0.049	0.078	0.197	0.081
0.000	0.049	0.078	0.196	0.081
0.000	0.049	0.078	0.196	0.080
0.000	0.050	0.079	0.196	0.080
0.000	0.049	0.078	0.197	0.081
0.000	0.050	0.079	0.197	0.081

Standard Deviations	0.0004	0.0005	0.0004	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5


Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Taylor D Gutschow TAYLOR D GUTSCHOW
Signature and Printed Name

03/13/2024
Date

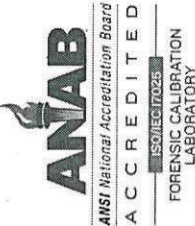

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-007167, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-007167</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>MONROE COUNTY SO</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>03/13/2024</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>12:32</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

03/13/2024

Date

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality