



INSTRUMENT PROCESSING SHEET

Agency FDLES/N 80-000219Florida Department of
Law EnforcementDate In 6/9/2025DI Completion Date 06/13/2025☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>SLH</u> Date <u>6/9/2025</u>		Quality Checks By <u>SLH</u> Date <u>6/10/2025</u>		Flow Calibration By _____ Date _____																																																													
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>226</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 105</u> 32 mm <u>0.152</u> (.139 - .169) 36 mm <u>0.171</u> (.156 - .190) 53 mm <u>0.238</u> (.228 - .278) 103 mm <u>0.507</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28662</u> ☒ Stability Checks		☐ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																													
		<table border="1" style="width:100%"><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr><tr><td>0.050</td><td>MP5088</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP5089</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP5090</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602 10/22/2026</td></tr></table>		Simulator	Serial #	Lot #/Exp	0.050	MP5088	202406K 06/19/2026	0.080	MP5089	202406L 06/19/2026	0.200	MP5090	202406N 06/20/2026	0.080 DGS	N/A	AG429602 10/22/2026	Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____ _____																																														
Simulator	Serial #	Lot #/Exp																																																															
0.050	MP5088	202406K 06/19/2026																																																															
0.080	MP5089	202406L 06/19/2026																																																															
0.200	MP5090	202406N 06/20/2026																																																															
0.080 DGS	N/A	AG429602 10/22/2026																																																															
Calibration Adjustment By _____				Department Inspection By <u>SLH</u>																																																													
Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%"><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%"><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></table>				Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Barometric Pressure ID# <u>28427</u> Gauge <u>1013</u> Instrument <u>1016</u> Mouth Alcohol Solution Lot # <u>2025A</u> Acetone Stock Solution Lot # <u>2024B</u> <table border="1" style="width:100%"><tr><th>Simulator</th><th>Serial Number</th></tr><tr><td>0.000</td><td>MP6294</td></tr><tr><td>Interferent</td><td>MP5087</td></tr><tr><td>0.050</td><td>MP5088</td></tr><tr><td>0.080</td><td>MP5089</td></tr><tr><td>0.200</td><td>MP5090</td></tr></table>		Simulator	Serial Number	0.000	MP6294	Interferent	MP5087	0.050	MP5088	0.080	MP5089	0.200	MP5090
Simulator	Serial #	Lot #	Expiration																																																														
0.000		N/A	N/A																																																														
0.040																																																																	
0.100																																																																	
0.200																																																																	
0.300																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial #	Lot #	Expiration																																																														
0.050																																																																	
0.080																																																																	
0.200																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial Number																																																																
0.000	MP6294																																																																
Interferent	MP5087																																																																
0.050	MP5088																																																																
0.080	MP5089																																																																
0.200	MP5090																																																																
Notes/Suggested Service: _____ Internal printer producing faint receipt, recommend using external printer. SLH 6/17/2025 _____ _____ _____ _____ _____				Attachments <table border="1" style="width:100%"><tr><td>☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment</td><td>☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 ☐ Other _____</td></tr></table> ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use		☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 ☐ Other _____																																																										
☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 ☐ Other _____																																																																
Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2025.06.19 13:14:04 -04'00'				Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2025.06.23 15:04:20 -04'00'																																																													
Tech Review / Date _____				Admin Review / Date _____																																																													

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: FDLE

Time of Inspection: 10:32

Date of Inspection: 06/10/2025

Serial Number: 80-000219

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

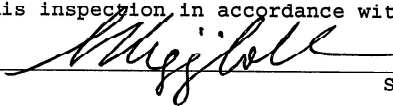
Number of Simulators Used: _____

Remarks:

BYPASS AGENCY INSPECTION COMPLIANCE UNDETERMINED

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.


 LEANDRA HIGGINBOTHAM
 Signature and Printed Name

06/10/2025
 Date

Stability Checks

Sut 6/10/2025
80-000219

0.050 g/210L	0.080 g/210L	0.200 g/210L	DGS 0.080 g/210L																																																																								
0.047 to 0.053 g/210L	0.077 to 0.083 g/210L	0.194 to 0.206 g/210L	0.077 to 0.083 g/210L																																																																								
Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis																																																																								
Reran 0.050 g/210L for improved print receipt. Still light print so move to external printer. Sut #2																																																																											
FDLE 6/10/2025 Printed: 6/10/2025 11:02:00 06/10/2025 ✓ Software: 8100.27 ✓ 80-000219 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr> <td>Air Blank</td> <td>0.000</td> <td>11:02</td> </tr> <tr> <td>Control Test</td> <td>0.049 ✓</td> <td>11:02</td> </tr> <tr> <td>Air Blank</td> <td>0.000</td> <td>11:03</td> </tr> <tr> <td>Control Test</td> <td>0.048 ✓</td> <td>11:04</td> </tr> <tr> <td>Air Blank</td> <td>0.000</td> <td>11:04</td> </tr> <tr> <td>Control Test</td> <td>0.051 ✓</td> <td>11:05</td> </tr> <tr> <td>Air Blank</td> <td>0.000</td> <td>11:05</td> </tr> <tr> <td colspan="3">Control Test Stats</td> </tr> <tr> <td>Average</td> <td>0.0493</td> <td></td> </tr> <tr> <td>Std Dev</td> <td>0.0015</td> <td></td> </tr> <tr> <td>Rel Std Dev(%)</td> <td>3.0963</td> <td></td> </tr> </tbody> </table> Operator's Signature: <i>Miguel</i>	Test	g/210L	Time	Air Blank	0.000	11:02	Control Test	0.049 ✓	11:02	Air Blank	0.000	11:03	Control Test	0.048 ✓	11:04	Air Blank	0.000	11:04	Control Test	0.051 ✓	11:05	Air Blank	0.000	11:05	Control Test Stats			Average	0.0493		Std Dev	0.0015		Rel Std Dev(%)	3.0963		FDLE 6/10/2025 Printed: 6/10/2025 11:02:00 06/10/2025 ✓ Software: 8100.27 ✓ 80-000219 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr> <td>Air Blank</td> <td>0.000</td> <td>11:01</td> </tr> <tr> <td>Control Test</td> <td>0.049 ✓</td> <td>11:02</td> </tr> <tr> <td>Air Blank</td> <td>0.000</td> <td>11:02</td> </tr> <tr> <td>Control Test</td> <td>0.049 ✓</td> <td>11:03</td> </tr> <tr> <td>Air Blank</td> <td>0.000</td> <td>11:03</td> </tr> <tr> <td>Control Test</td> <td>0.050 ✓</td> <td>11:04</td> </tr> <tr> <td>Air Blank</td> <td>0.000</td> <td>11:04</td> </tr> <tr> <td colspan="3">Control Test Stats</td> </tr> <tr> <td>Average</td> <td>0.0493</td> <td></td> </tr> <tr> <td>Std Dev</td> <td>0.0016</td> <td></td> </tr> <tr> <td>Rel Std Dev(%)</td> <td>1.1703</td> <td></td> </tr> </tbody> </table> Operator's Signature: <i>Miguel</i>	Test	g/210L	Time	Air Blank	0.000	11:01	Control Test	0.049 ✓	11:02	Air Blank	0.000	11:02	Control Test	0.049 ✓	11:03	Air Blank	0.000	11:03	Control Test	0.050 ✓	11:04	Air Blank	0.000	11:04	Control Test Stats			Average	0.0493		Std Dev	0.0016		Rel Std Dev(%)	1.1703			
Test	g/210L	Time																																																																									
Air Blank	0.000	11:02																																																																									
Control Test	0.049 ✓	11:02																																																																									
Air Blank	0.000	11:03																																																																									
Control Test	0.048 ✓	11:04																																																																									
Air Blank	0.000	11:04																																																																									
Control Test	0.051 ✓	11:05																																																																									
Air Blank	0.000	11:05																																																																									
Control Test Stats																																																																											
Average	0.0493																																																																										
Std Dev	0.0015																																																																										
Rel Std Dev(%)	3.0963																																																																										
Test	g/210L	Time																																																																									
Air Blank	0.000	11:01																																																																									
Control Test	0.049 ✓	11:02																																																																									
Air Blank	0.000	11:02																																																																									
Control Test	0.049 ✓	11:03																																																																									
Air Blank	0.000	11:03																																																																									
Control Test	0.050 ✓	11:04																																																																									
Air Blank	0.000	11:04																																																																									
Control Test Stats																																																																											
Average	0.0493																																																																										
Std Dev	0.0016																																																																										
Rel Std Dev(%)	1.1703																																																																										

FDLE

Intoxilyzer - Alcohol Analyzer

Model 8000

SN 80-000219

06/10/2025

Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:58
Control Test	0.078	11:58
Air Blank	0.000	11:59
Control Test	0.079	12:00
Air Blank	0.000	12:00
Control Test	0.079	12:01
Air Blank	0.000	12:01
Control Test Stats		
Average	0.0787 ✓	
Std Dev	0.0006	
Rel Std Dev(%)	0.7339	

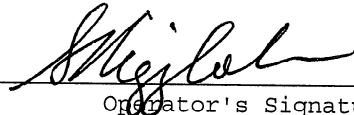


Operator's Signature

FDLE
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000219
06/10/2025
Software: 8100.27

DGS

Test	g/210L	Time
Air Blank	0.000	12:06
Control Test	0.078	12:07
Air Blank	0.000	12:07
Control Test	0.078	12:08
Air Blank	0.000	12:08
Control Test	0.077	12:08
Air Blank	0.000	12:09
Control Test Stats		
Average	0.0777 ✓	
Std Dev	0.0006	
Rel Std Dev(%)	0.7434	

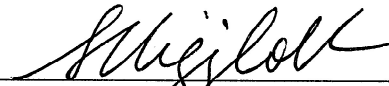


Operator's Signature

FDLE
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000219
06/10/2025
Software: 8100.27

#1

Test	g/210L	Time
Air Blank	0.000	11:41
Control Test	0.163	11:41
Air Blank	0.000	11:42
Control Test	0.163	11:43
Air Blank	0.000	11:43
Control Test	0.162	11:44
Air Blank	0.000	11:44
Control Test Stats		
Average	0.1627	
Std Dev	0.0006	
Rel Std Dev(%)	0.3549	


Operator's Signature

RCA = root cause analysis — retightened
and checked seal sub

FDLE
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000219
06/10/2025
Software: 8100.27

#2

Test	g/210L	Time
Air Blank	0.000	11:47
Control Test	0.197 ✓	11:48
Air Blank	0.000 ✓	11:48
Control Test	0.197 ✓	11:49
Air Blank	0.000 ✓	11:49
Control Test	0.197 ✓	11:50
Air Blank	0.000	11:51
Control Test Stats		
Average	0.1970	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: FDLE

Time of Inspection: 17:28

Date of Inspection: 06/13/2025

Serial Number: 80-000219

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.047 / 0.049	0.079	0.198	0.077
0.000	0.047 / 0.049	0.078	0.197	0.077
0.000	0.046 / 0.048	0.079	0.197	0.078
0.000	0.046 / 0.048	0.079	0.197	0.077
0.000	0.045 / 0.048	0.078	0.197	0.077
0.000	0.045 / 0.047	0.078	0.197	0.077
0.000	0.045 / 0.047	0.078	0.196	0.077
0.000	0.045 / 0.047	0.077	0.196	0.077
0.000	0.044 / 0.046	0.077	0.196	0.076
0.000	0.044 / 0.046	0.077	0.196	0.077

Standard Deviations	0.0010 / 0.0010	0.0008	0.0006	0.0004
---------------------	-----------------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0007 Number of Simulators Used: 5

Remarks:

05: Control Outside Tolerance. - tightened lid of simulator and re-ran. Sat 6/13/2025

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

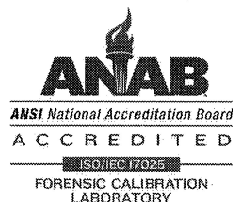
I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



LEANDRA HIGGINBOTHAM

Signature and Printed Name

06/13/2025
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000219, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000219</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>FDLE</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>06/13/2025</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>17:28</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

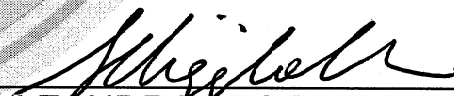
Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

06/13/2025

Date


LEANDRA HIGGINBOTHAM,
Department Inspector

FDLE/ATP Form 69 October 2024

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1