



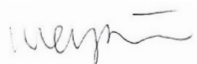
INSTRUMENT PROCESSING SHEET

Agency UMATILLA PD(LAKE SO)S/N 80-000836Florida Department of
Law EnforcementDate In 7/8/25DI Completion Date 07/17/2025☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By <u>KTS</u>	Date <u>7/8/25</u>	Quality Checks	By <u>WKP</u>	Date <u>7/10/2025</u>	Flow Calibration	By <u>SLH</u>	Date <u>07/16/2025</u>																												
☐ Annual ☒ Registration <u>SP 7/21/25</u> ☐ Return from CMI / EE Visual Inspection: ☐ Case ☐ Handle ☐ Keyboard ☐ Dry Gas Shelf ☒ Feet ☐ Breath Tube ☐ Ports ☐ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: <u>Missing right foot.</u>			☐ Breath Tube Screen ☐ Replace External O-Rings ☐ Instrument Set Up Verified ☐ R-Value <u>136</u> ☐ Flow Verification (L/s) Flow Column # <u>ATP 102</u> 32 mm <u>0.140</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.222*</u> (.228 - .278) 103 mm <u>0.488</u> (.447 - .547) ☐ Barometric Pressure Check Gauge ID # <u>28662</u> ☐ Stability Checks			☐ Flow Column # <u>ATP103</u> ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value <u>137</u> ☐ Post Calibration Verification (L/s) Flow Column # <u>ATP102</u> 32 mm <u>0.152</u> (.139 - .169) 36 mm <u>0.171</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.523</u> (.447 - .547)																														
			<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP5089</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP5090</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602 10/22/2026</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP5088	202406K 06/19/2026	0.080	MP5089	202406L 06/19/2026	0.200	MP5090	202406N 06/20/2026	0.080 DGS	N/A	AG429602 10/22/2026	<table border="1"><thead><tr><th colspan="2">Maintenance By _____ Date _____</th></tr></thead><tbody><tr><td>☐ Battery Replacement</td><td></td></tr><tr><td>☐ Dry Gas Regulator Replacement</td><td></td></tr><tr><td>☐ Breath Tube Replacement</td><td></td></tr><tr><td>☐ Other _____</td><td></td></tr></tbody></table>			Maintenance By _____ Date _____		☐ Battery Replacement		☐ Dry Gas Regulator Replacement		☐ Breath Tube Replacement		☐ Other _____				
Simulator	Serial #	Lot #/Exp																																		
0.050	MP5088	202406K 06/19/2026																																		
0.080	MP5089	202406L 06/19/2026																																		
0.200	MP5090	202406N 06/20/2026																																		
0.080 DGS	N/A	AG429602 10/22/2026																																		
Maintenance By _____ Date _____																																				
☐ Battery Replacement																																				
☐ Dry Gas Regulator Replacement																																				
☐ Breath Tube Replacement																																				
☐ Other _____																																				
Calibration Adjustment By _____			Department Inspection By <u>SLH</u>																																	
Barometric Pressure Gauge _____ ID # _____			Barometric Pressure ID# <u>28421</u>																																	
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge <u>1016</u> Instrument <u>1012</u> Mouth Alcohol Solution Lot # <u>2025-A</u> Acetone Stock Solution Lot # <u>2024-B</u>					
Simulator	Serial #	Lot #	Expiration																																	
0.000		N/A	N/A																																	
0.040																																				
0.100																																				
0.200																																				
0.300																																				
0.080 DGS	N/A																																			
☐ Post Calibration Adjustment Stability Checks			<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6294</td></tr><tr><td>Interferent</td><td>MP5087</td></tr><tr><td>0.050</td><td>MP5088</td></tr><tr><td>0.080</td><td>MP5089</td></tr><tr><td>0.200</td><td>MP5090</td></tr></tbody></table>						Simulator	Serial Number	0.000	MP6294	Interferent	MP5087	0.050	MP5088	0.080	MP5089	0.200	MP5090																
Simulator	Serial Number																																			
0.000	MP6294																																			
Interferent	MP5087																																			
0.050	MP5088																																			
0.080	MP5089																																			
0.200	MP5090																																			
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Attachments ☐ Form 41 ☐ Post-Stability Checks ☐ Stability Checks ☐ Flow Calibration ☐ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☒ Other <u>Form 47 and 0.08 stability check SP 7/21/25</u>													
Simulator	Serial #	Lot #	Expiration																																	
0.050																																				
0.080																																				
0.200																																				
0.080 DGS	N/A																																			
Notes/Suggested Service: <u>Performed root cause analysis on failed flow verification, repeated and flow verification was outside of range. WKP 7/10/2025*</u>			☐ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☐ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☐ Conduct an Agency Inspection Before Evidentiary Use																																	
Tech Review: <u>Added Form 47 and stability check performed on 7/21/25. SP 7/21/25 Checked box for registration updated DI date on Form 47. SP 7/21/25</u>			Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2025.07.21 13:05:29 -0400 Phil Nicodemo Digitally signed by Phil Nicodemo Date: 2025.07.21 13:17:47 -0400																																	
Tech Review / Date			Admin Review / Date																																	

*Tech Review: Flow verification repeated measurement remained the same value as the initial measurement. WKP 7/21/2025

Stability Checks

0.05g/210L			0.08g/210L			0.20g/210L			DGS 0.08g/210L		
0.047 to 0.053			0.077 to 0.083			0.194 to 0.206			0.077 to 0.083		
✓			✓			✓			✓		
Performed Root Case Analysis			Performed Root Case Analysis			Performed Root Case Analysis			Performed Root Case Analysis		
UMATILLA POLYMER S01 Intoxilyzer - Alcohol Analyzer Model 8000 SN 60-000836 07/10/2025 Software: 8100.27			UMATILLA POLYMER S01 Intoxilyzer - Alcohol Analyzer Model 8000 SN 60-000836 07/10/2025 Software: 8100.27			UMATILLA POLYMER S01 Intoxilyzer - Alcohol Analyzer Model 8000 SN 60-000836 07/10/2025 Software: 8100.27			UMATILLA POLYMER S01 Intoxilyzer - Alcohol Analyzer Model 8000 SN 60-000836 07/10/2025 Software: 8100.27		
Test	g/210L	Time	Test	g/210L	Time	Test	g/210L	Time	Test	g/210L	Time
Air Blank	0.000	13:25	Air Blank	0.000	13:36	Air Blank	0.000	13:43	Air Blank	0.000	13:48
Control Test	0.050	13:26	Control Test	0.079	13:37	Control Test	0.198	13:43	Control Test	0.080	13:49
Air Blank	0.000	13:27	Air Blank	0.000	13:38	Air Blank	0.000	13:44	Air Blank	0.000	13:49
Control Test	0.049	13:27	Control Test	0.079	13:38	Control Test	0.198	13:45	Control Test	0.080	13:50
Air Blank	0.000	13:28	Air Blank	0.000	13:39	Air Blank	0.000	13:45	Air Blank	0.000	13:50
Control Test	0.049	13:28	Control Test	0.079	13:39	Control Test	0.197	13:46	Control Test	0.079	13:50
Air Blank	0.000	13:29	Air Blank	0.000	13:40	Air Blank	0.000	13:47	Air Blank	0.000	13:51
Control Test Stats			Control Test Stats			Control Test Stats			Control Test Stats		
Average	0.0493		Average	0.0790		Average	0.1977		Average	0.0797	
Std Dev	0.0006		Std Dev	0.0000		Std Dev	0.0006		Std Dev	0.0006	
Rel Std Dev(%)	1.1703		Rel Std Dev(%)	0.0000		Rel Std Dev(%)	0.2921		Rel Std Dev(%)	0.7247	
 Operator's Signature			 Operator's Signature			 Operator's Signature			 Operator's Signature		

UMATILLA PD(LAKE SO)
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000836
07/16/2025
Software: 8100.27

Flow Rate Calibration*****

1: Rate (Liters/min) = 5
SQRT(Diff)) = 6.781
2: Rate (Liters/min) = 15
SQRT(Diff)) = 11.000
3: Rate (Liters/min) = 30
SQRT(Diff)) = 20.047

Dependent Data Scale Factor = 100000 L/min

Independent Data Scale Factor = 256

Rounded Slope = 722

Rounded Intercept = -664388

Correlation = 0.99576

SLS
7/8/16/25

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: UMATILLA PD (LAKE SO)

Serial Number: 80-000836

Time of Inspection: 15:53

Date of Inspection: 07/17/2025

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.049	0.079	0.197	0.079
0.000	0.049	0.079	0.198	0.078
0.000	0.049	0.079	0.198	0.078
0.000	0.049	0.079	0.198	0.078
0.000	0.049	0.079	0.198	0.078
0.000	0.049	0.079	0.198	0.078
0.000	0.050	0.079	0.198	0.078
0.000	0.049	0.079	0.198	0.078
0.000	0.050	0.079	0.198	0.078
0.000	0.050	0.079	0.198	0.078

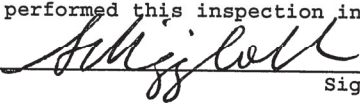
Standard Deviations	0.0004	0.0000	0.0003	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

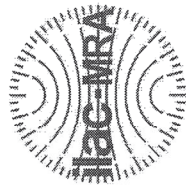
Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.


 LEANDRA HIGGINBOTHAM
 Signature and Printed Name

07/17/2025
 Date



ANAB

ANSI National Accreditation Board
ACCREDITED
ISO/IEC 17025
FORENSIC CALIBRATION
LABORATORY

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000836, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000836</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>UMATILLA PD(LAKE SO)</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>07/17/2025</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>15:53</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 October 2024

Issuing Authority: Alcohol Testing Program

07/17/2025

Date

Leandra Higginbotham

LEANDRA HIGGINBOTHAM,

Department Inspector

Service • Integrity • Respect • Quality

0.08 Stability Check

80-000836

LAKE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 6100 SN 80-000836
07/21/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:23
Control Test	0.079	10:23
Air Blank	0.000	10:24
Control Test	0.079	10:24
Air Blank	0.000	10:25
Control Test	0.079	10:25
Air Blank	0.000	10:25
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

SP
Operator's Signature

Agency Name Updated
Post Inspection.
Ran stability check to
demonstrate and
confirm name
change.

SP 7/21/25



Florida Department of
Law Enforcement

REQUEST FOR REGISTRATION

MAKE AND MODEL OF INSTRUMENT: INTOXILYZER 8000

SERIAL NUMBER: 80-000836

OWNING AGENCY: LAKE COUNTY SHERIFF'S OFFICE

DATE OF DEPARTMENT INSPECTION: JULY 17, 2025
SP 7/21/25

AGENCY INSPECTOR: Sgt. JERRY S. DeSANTIS 818
CPL MATTHEW BOWDEN 1100

ADDRESS: 360 W. RUBY ST. TAVANES FL 32778

CITY, STATE, ZIP: TAVANES FL 32778

TELEPHONE NUMBER: 352-602-9722

FAX NUMBER: NA

EMAIL ADDRESS (if available): JEFF.DeSANTIS@LCSO.ORG

For Program Office Use Only:

- ☒ Registration Issued
- ☒ Instrument Added to Evidentiary Instrument Database
- ☐ Instrument Added to Monthly Statistics Database
- ☒ Contact Information Added to Instrument Database