



INSTRUMENT PROCESSING SHEET

Agency Indian River County Sheriff's OfficeS/N 80-001329Florida Department of
Law EnforcementDate In 07/15/2025 DI Completion Date 08/13/2025☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By	TDG	Date	08/08/2025	Quality Checks	By	TDG	Date	08/12/2025	Flow Calibration	By	Date																											
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>Agency Inspector reports the instrument failed the monthly inspection at the 0.20 ARS Test.</u>					☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>178</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP106</u> 32 mm <u>0.140</u> (.139 - .169) 36 mm <u>0.156</u> (.156 - .190) 53 mm <u>0.230</u> (.228 - .278) 103 mm <u>0.496</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28199</u> ☒ Stability Checks					Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																													
					<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP6287</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP6288</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602 10/22/2026</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP6286	202406K 06/19/2026	0.080	MP6287	202406L 06/19/2026	0.200	MP6288	202406N 06/20/2026	0.080 DGS	N/A	AG429602 10/22/2026				<table border="1"><thead><tr><th colspan="2">Maintenance By _____ Date _____</th></tr><tr><td>☐ Battery Replacement</td><td></td></tr><tr><td>☐ Dry Gas Regulator Replacement</td><td></td></tr><tr><td>☐ Breath Tube Replacement</td><td></td></tr><tr><td>☐ Other _____</td><td></td></tr></thead></table>			Maintenance By _____ Date _____		☐ Battery Replacement		☐ Dry Gas Regulator Replacement		☐ Breath Tube Replacement		☐ Other _____		
Simulator	Serial #	Lot #/Exp																																					
0.050	MP6286	202406K 06/19/2026																																					
0.080	MP6287	202406L 06/19/2026																																					
0.200	MP6288	202406N 06/20/2026																																					
0.080 DGS	N/A	AG429602 10/22/2026																																					
Maintenance By _____ Date _____																																							
☐ Battery Replacement																																							
☐ Dry Gas Regulator Replacement																																							
☐ Breath Tube Replacement																																							
☐ Other _____																																							
Calibration Adjustment By _____					Department Inspection By <u>TDG</u>																																		
Barometric Pressure Gauge _____ ID # _____					Barometric Pressure ID# <u>28199</u>																																		
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>					Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge <u>1018</u> Instrument <u>1009</u> Mouth Alcohol Solution Lot # <u>2024-A</u> Acetone Stock Solution Lot # <u>2024-B</u>						
Simulator	Serial #	Lot #	Expiration																																				
0.000		N/A	N/A																																				
0.040																																							
0.100																																							
0.200																																							
0.300																																							
0.080 DGS	N/A																																						
☐ Post Calibration Adjustment Stability Checks					<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table>							Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288																
Simulator	Serial Number																																						
0.000	MP6284																																						
Interferent	MP6285																																						
0.050	MP6286																																						
0.080	MP6287																																						
0.200	MP6288																																						
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>					Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☒ Form 40 ☐ Calibration Adjustment ☐ Other _____														
Simulator	Serial #	Lot #	Expiration																																				
0.050																																							
0.080																																							
0.200																																							
0.080 DGS	N/A																																						
Notes/Suggested Service: <u>Added plastic caps for return to agency. Could not duplicate the reported low readings at the 0.20 concentration. (TDG 8/13/25)</u>					☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																		
					Digitally signed by Phil Nicodemo Date: 2025.08.13 15:18:12 -0400 Tech Review / Date _____																																		
					Digitally signed by Shayla Platt Date: 2025.08.14 12:28:47 -0400 Admin Review / Date _____																																		

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: INDIAN RIVER CSO
Time of Inspection: 13:36

Date of Inspection: 08/12/2025

Serial Number: 80-001329
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

AI NOT CONDUCTED. COMPLIANCE NOT DETERMINED.

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



TAYLOR D GUTSCHOW

Signature and Printed Name

08/12/2025
Date

0.05g/210L 0.047 to 0.053	0.08g/210L 0.077 to 0.083	0.20g/210L 0.194 to 0.206	DGS 0.08g/210L 0.077 to 0.083 ☒ ≤0.003 of Wet ☒																																																																																																												
INDIAN RIVER CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001329 08/12/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:01</td></tr> <tr><td>Control Test</td><td>0.049</td><td>14:02</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:03</td></tr> <tr><td>Control Test</td><td>0.048</td><td>14:03</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:04</td></tr> <tr><td>Control Test</td><td>0.048</td><td>14:05</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:05</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0483</td><td></td></tr> <tr><td>Std Dev</td><td>0.0016</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.1945</td><td></td></tr> </tbody> </table> Operator's Signature <i>ML</i>	Test	g/210L	Time	Air Blank	0.000	14:01	Control Test	0.049	14:02	Air Blank	0.000	14:03	Control Test	0.048	14:03	Air Blank	0.000	14:04	Control Test	0.048	14:05	Air Blank	0.000	14:05	Control Test Stats			Average	0.0483		Std Dev	0.0016		Rel Std Dev(%)	1.1945		INDIAN RIVER CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001329 08/12/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank*</td><td>0.000</td><td>14:08</td></tr> <tr><td>Control Test</td><td>0.080</td><td>14:09</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:09</td></tr> <tr><td>Control Test</td><td>0.080</td><td>14:10</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:10</td></tr> <tr><td>Control Test</td><td>0.079</td><td>14:11</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:12</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0797</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7247</td><td></td></tr> </tbody> </table> Operator's Signature <i>ML</i>	Test	g/210L	Time	Air Blank*	0.000	14:08	Control Test	0.080	14:09	Air Blank	0.000	14:09	Control Test	0.080	14:10	Air Blank	0.000	14:10	Control Test	0.079	14:11	Air Blank	0.000	14:12	Control Test Stats			Average	0.0797		Std Dev	0.0006		Rel Std Dev(%)	0.7247		INDIAN RIVER CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001329 08/12/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:54</td></tr> <tr><td>Control Test</td><td>0.199</td><td>13:54</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:55</td></tr> <tr><td>Control Test</td><td>0.198</td><td>13:56</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:56</td></tr> <tr><td>Control Test</td><td>0.198</td><td>13:57</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:57</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1983</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2911</td><td></td></tr> </tbody> </table> Operator's Signature <i>ML</i>	Test	g/210L	Time	Air Blank	0.000	13:54	Control Test	0.199	13:54	Air Blank	0.000	13:55	Control Test	0.198	13:56	Air Blank	0.000	13:56	Control Test	0.198	13:57	Air Blank	0.000	13:57	Control Test Stats			Average	0.1983		Std Dev	0.0006		Rel Std Dev(%)	0.2911		Printed to external printer. Will attach. ML 8/12/25
Test	g/210L	Time																																																																																																													
Air Blank	0.000	14:01																																																																																																													
Control Test	0.049	14:02																																																																																																													
Air Blank	0.000	14:03																																																																																																													
Control Test	0.048	14:03																																																																																																													
Air Blank	0.000	14:04																																																																																																													
Control Test	0.048	14:05																																																																																																													
Air Blank	0.000	14:05																																																																																																													
Control Test Stats																																																																																																															
Average	0.0483																																																																																																														
Std Dev	0.0016																																																																																																														
Rel Std Dev(%)	1.1945																																																																																																														
Test	g/210L	Time																																																																																																													
Air Blank*	0.000	14:08																																																																																																													
Control Test	0.080	14:09																																																																																																													
Air Blank	0.000	14:09																																																																																																													
Control Test	0.080	14:10																																																																																																													
Air Blank	0.000	14:10																																																																																																													
Control Test	0.079	14:11																																																																																																													
Air Blank	0.000	14:12																																																																																																													
Control Test Stats																																																																																																															
Average	0.0797																																																																																																														
Std Dev	0.0006																																																																																																														
Rel Std Dev(%)	0.7247																																																																																																														
Test	g/210L	Time																																																																																																													
Air Blank	0.000	13:54																																																																																																													
Control Test	0.199	13:54																																																																																																													
Air Blank	0.000	13:55																																																																																																													
Control Test	0.198	13:56																																																																																																													
Air Blank	0.000	13:56																																																																																																													
Control Test	0.198	13:57																																																																																																													
Air Blank	0.000	13:57																																																																																																													
Control Test Stats																																																																																																															
Average	0.1983																																																																																																														
Std Dev	0.0006																																																																																																														
Rel Std Dev(%)	0.2911																																																																																																														

INDIAN RIVER CSO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001329
08/12/2025
Software: 8100.27

D65

Test	g/210L	Time
Air Blank	0.000	13:44
Control Test	0.080	13:45
Air Blank	0.000	13:45
Control Test	0.078	13:46
Air Blank	0.000	13:46
Control Test	0.078	13:46
Air Blank	0.000	13:47
Control Test Stats		
Average	0.0787	
Std Dev	0.0012	
Rel Std Dev(%)	1.4678	

ML

Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: INDIAN RIVER CSO
Time of Inspection: 12:19

Date of Inspection: 08/13/2025

Serial Number: 80-001329
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.048	0.078	0.197	0.081
0.000	0.048	0.078	0.197	0.082
0.000	0.048	0.078	0.198	0.081
0.000	0.048	0.079	0.198	0.081
0.000	0.048	0.079	0.197	0.081
0.000	0.048	0.078	0.198	0.081
0.000	0.048	0.079	0.197	0.080
0.000	0.048	0.079	0.198	0.081
0.000	0.048	0.079	0.198	0.081
0.000	0.048	0.079	0.198	0.081


Standard Deviations	0.0000	0.0005	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

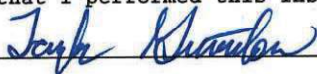

TAYLOR D GUTSCHOW

Signature and Printed Name

08/13/2025
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001329, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001329</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>INDIAN RIVER CSO</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>08/13/2025</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>12:19</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Taylor
Gutschow
Digitally signed by Taylor
Gutschow
Date: 2025.08.13 14:42:53
+04'00'

08/13/2025

Date

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality



INSTRUMENT PROCESSING SHEET

Agency Indian River CSOS/N 80-001329Florida Department of
Law EnforcementDate In 02/06/2025 DI Completion Date 02/10/2025☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>TDG</u> Date <u>02/06/2025</u> ☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____	Quality Checks By <u>TDG</u> Date <u>02/07/2025</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>184</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP106</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.503</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP6287</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP6288</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602 10/22/2026</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K 03/29/2025	0.080	MP6287	202303L 03/29/2025	0.200	MP6288	202304C 04/05/2025	0.080 DGS	N/A	AG429602 10/22/2026	Flow Calibration By _____ Date _____ Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____																																	
Simulator	Serial #	Lot #/Exp																																																
0.050	MP6286	202303K 03/29/2025																																																
0.080	MP6287	202303L 03/29/2025																																																
0.200	MP6288	202304C 04/05/2025																																																
0.080 DGS	N/A	AG429602 10/22/2026																																																
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A		
Simulator	Serial #	Lot #	Expiration																																															
0.000		N/A	N/A																																															
0.040																																																		
0.100																																																		
0.200																																																		
0.300																																																		
0.080 DGS	N/A																																																	
Simulator	Serial #	Lot #	Expiration																																															
0.050																																																		
0.080																																																		
0.200																																																		
0.080 DGS	N/A																																																	
Department Inspection By <u>TDG</u> Barometric Pressure ID# <u>26932</u> Gauge <u>1023</u> Instrument <u>1018</u> Mouth Alcohol Solution Lot # <u>2024-A</u> Acetone Stock Solution Lot # <u>2023-B</u> <table border="1" style="width:100%"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table> Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☐ Other _____ ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Phil _____ Digitally signed by Shayla Platt Nicodemo _____ Date: 2025.02.12 11:24:46 -05'00' Tech Review / Date _____ Admin: Review / Date _____			Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288																																				
Simulator	Serial Number																																																	
0.000	MP6284																																																	
Interferent	MP6285																																																	
0.050	MP6286																																																	
0.080	MP6287																																																	
0.200	MP6288																																																	

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																				
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083																																																																																																																																				
INDIAN RIVER CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001329 02/07/2025 Software: 8100.27 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>13:54</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:55</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:56</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:56</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:57</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:57</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:58</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0490</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </table> Operator's Signature 	Air Blank	0.000	13:54	Control Test	0.049	13:55	Air Blank	0.000	13:56	Control Test	0.049	13:56	Air Blank	0.000	13:57	Control Test	0.049	13:57	Air Blank	0.000	13:58	Control Test Stats			Average	0.0490		Std Dev	0.0000		Rel Std Dev(%)	0.0000		INDIAN RIVER CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001329 02/07/2025 Software: 8100.27 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>14:01</td></tr> <tr><td>Control Test</td><td>0.078</td><td>14:02</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:02</td></tr> <tr><td>Control Test</td><td>0.078</td><td>14:03</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:04</td></tr> <tr><td>Control Test</td><td>0.079</td><td>14:04</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:05</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0783</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7370</td><td></td></tr> </table> Operator's Signature 	Air Blank	0.000	14:01	Control Test	0.078	14:02	Air Blank	0.000	14:02	Control Test	0.078	14:03	Air Blank	0.000	14:04	Control Test	0.079	14:04	Air Blank	0.000	14:05	Control Test Stats			Average	0.0783		Std Dev	0.0006		Rel Std Dev(%)	0.7370		INDIAN RIVER CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001329 02/07/2025 Software: 8100.27 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>13:47</td></tr> <tr><td>Control Test</td><td>0.194</td><td>13:48</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:49</td></tr> <tr><td>Control Test</td><td>0.195</td><td>13:49</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:50</td></tr> <tr><td>Control Test</td><td>0.195</td><td>13:51</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:51</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1947</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2965</td><td></td></tr> </table> Operator's Signature 	Air Blank	0.000	13:47	Control Test	0.194	13:48	Air Blank	0.000	13:49	Control Test	0.195	13:49	Air Blank	0.000	13:50	Control Test	0.195	13:51	Air Blank	0.000	13:51	Control Test Stats			Average	0.1947		Std Dev	0.0006		Rel Std Dev(%)	0.2965		INDIAN RIVER CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001329 02/07/2025 Software: 8100.27 Test g/210L Time <table border="1"> <tr><td>Air Blank</td><td>0.000</td><td>13:32</td></tr> <tr><td>Control Test</td><td>0.081</td><td>13:33</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:33</td></tr> <tr><td>Control Test</td><td>0.079</td><td>13:34</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:34</td></tr> <tr><td>Control Test</td><td>0.080</td><td>13:34</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:35</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0800</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.2500</td><td></td></tr> </table> Operator's Signature 	Air Blank	0.000	13:32	Control Test	0.081	13:33	Air Blank	0.000	13:33	Control Test	0.079	13:34	Air Blank	0.000	13:34	Control Test	0.080	13:34	Air Blank	0.000	13:35	Control Test Stats			Average	0.0800		Std Dev	0.0010		Rel Std Dev(%)	1.2500	
Air Blank	0.000	13:54																																																																																																																																					
Control Test	0.049	13:55																																																																																																																																					
Air Blank	0.000	13:56																																																																																																																																					
Control Test	0.049	13:56																																																																																																																																					
Air Blank	0.000	13:57																																																																																																																																					
Control Test	0.049	13:57																																																																																																																																					
Air Blank	0.000	13:58																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0490																																																																																																																																						
Std Dev	0.0000																																																																																																																																						
Rel Std Dev(%)	0.0000																																																																																																																																						
Air Blank	0.000	14:01																																																																																																																																					
Control Test	0.078	14:02																																																																																																																																					
Air Blank	0.000	14:02																																																																																																																																					
Control Test	0.078	14:03																																																																																																																																					
Air Blank	0.000	14:04																																																																																																																																					
Control Test	0.079	14:04																																																																																																																																					
Air Blank	0.000	14:05																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0783																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.7370																																																																																																																																						
Air Blank	0.000	13:47																																																																																																																																					
Control Test	0.194	13:48																																																																																																																																					
Air Blank	0.000	13:49																																																																																																																																					
Control Test	0.195	13:49																																																																																																																																					
Air Blank	0.000	13:50																																																																																																																																					
Control Test	0.195	13:51																																																																																																																																					
Air Blank	0.000	13:51																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.1947																																																																																																																																						
Std Dev	0.0006																																																																																																																																						
Rel Std Dev(%)	0.2965																																																																																																																																						
Air Blank	0.000	13:32																																																																																																																																					
Control Test	0.081	13:33																																																																																																																																					
Air Blank	0.000	13:33																																																																																																																																					
Control Test	0.079	13:34																																																																																																																																					
Air Blank	0.000	13:34																																																																																																																																					
Control Test	0.080	13:34																																																																																																																																					
Air Blank	0.000	13:35																																																																																																																																					
Control Test Stats																																																																																																																																							
Average	0.0800																																																																																																																																						
Std Dev	0.0010																																																																																																																																						
Rel Std Dev(%)	1.2500																																																																																																																																						

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: INDIAN RIVER CSO
Time of Inspection: 11:31

Date of Inspection: 02/10/2025

Serial Number: 80-001329
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.048	0.079	0.197	0.081
0.000	0.048	0.079	0.197	0.082
0.000	0.048	0.079	0.197	0.082
0.000	0.048	0.079	0.197	0.081
0.000	0.049	0.079	0.197	0.082
0.000	0.048	0.079	0.197	0.081
0.000	0.048	0.078	0.197	0.081
0.000	0.048	0.079	0.197	0.082
0.000	0.049	0.079	0.197	0.081
0.000	0.049	0.079	0.197	0.081

Standard Deviations	0.0004	0.0003	0.0000	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



Signature and Printed Name

TAYLOR D GUTSCHOW

02/10/2025
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001329, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-001329	UNCERTAINTY* $\pm$
Owning Agency:	INDIAN RIVER CSO	0.050 g/ 210 L 0.004
Calibration Date:	02/10/2025	0.080 g/ 210 L 0.004
Calibration Time:	11:31	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

02/10/2025

Date

Taylor D Gutschow,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality