



INSTRUMENT PROCESSING SHEET


Agency ORANGE COUNTY SOS/N 80-001420Florida Department of
Law EnforcementDate In 5/29/2025DI Completion Date 6/3/2025☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By WKP	Date	Quality Checks	By WKP	Date	Flow Calibration	By	Date																																																											
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>DGS port plastic connector was</u> <u>broken upon arrival.WKP 5/29/25</u> 		<u>5/29/2025</u>	☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>212</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP-105</u> 32 mm <u>0.167</u> (.139 - .169) 36 mm <u>0.183</u> (.156 - .190) 53 mm <u>0.250</u> (.228 - .278) 103 mm <u>0.500</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28427</u> ☒ Stability Checks		<u>5/29/2025</u>	☐ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																																													
			<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP5089</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP5090</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602 10/22/2026</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP5088	202406K 06/19/2026	0.080	MP5089	202406L 06/19/2026	0.200	MP5090	202406N 06/20/2026	0.080 DGS	N/A	AG429602 10/22/2026	<table border="1"><thead><tr><th colspan="2">Maintenance By SLH</th><th>Date</th></tr></thead><tbody><tr><td colspan="2">☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☒ Other <u>Replaced DGS plastic elbow</u> <u>connector.SLH</u></td><td><u>05/29/2025</u></td></tr></tbody></table>			Maintenance By SLH		Date	☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☒ Other <u>Replaced DGS plastic elbow</u> <u>connector.SLH</u>		<u>05/29/2025</u>																																						
Simulator	Serial #	Lot #/Exp																																																																	
0.050	MP5088	202406K 06/19/2026																																																																	
0.080	MP5089	202406L 06/19/2026																																																																	
0.200	MP5090	202406N 06/20/2026																																																																	
0.080 DGS	N/A	AG429602 10/22/2026																																																																	
Maintenance By SLH		Date																																																																	
☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☒ Other <u>Replaced DGS plastic elbow</u> <u>connector.SLH</u>		<u>05/29/2025</u>																																																																	
Calibration Adjustment By _____			Department Inspection By <u>SLH</u>																																																																
Barometric Pressure Gauge _____ ID # _____ <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Barometric Pressure ID# <u>28662</u> Gauge <u>1013</u> Instrument <u>1013</u> Mouth Alcohol Solution Lot # <u>2025A</u> Acetone Stock Solution Lot # <u>2024B</u> <table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5086</td></tr><tr><td>Interferent</td><td>MP5087</td></tr><tr><td>0.050</td><td>MP5088</td></tr><tr><td>0.080</td><td>MP5089</td></tr><tr><td>0.200</td><td>MP5090</td></tr></tbody></table> Attachments <table border="1"><tbody><tr><td>☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment</td><td>☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____</td></tr></tbody></table>			Simulator	Serial Number	0.000	MP5086	Interferent	MP5087	0.050	MP5088	0.080	MP5089	0.200	MP5090	☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____
Simulator	Serial #	Lot #	Expiration																																																																
0.000		N/A	N/A																																																																
0.040																																																																			
0.100																																																																			
0.200																																																																			
0.300																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial #	Lot #	Expiration																																																																
0.050																																																																			
0.080																																																																			
0.200																																																																			
0.080 DGS	N/A																																																																		
Simulator	Serial Number																																																																		
0.000	MP5086																																																																		
Interferent	MP5087																																																																		
0.050	MP5088																																																																		
0.080	MP5089																																																																		
0.200	MP5090																																																																		
☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment	☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____																																																																		
Notes/Suggested Service: <u>Inadvertently left DGS unhooked</u> <u>when performing stability check, identified user error and</u> <u>repeated analysis. WKP 5/29/2025</u> 			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2025.06.09 13:34:23 +04'00' Shayla Platt Digitally signed by Shayla Platt Date: 2025.06.09 13:59:07 -04'00' Tech Review / Date _____ Admin Review / Date _____																																																																

0.05g/210L

ORANGE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001420
05/29/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	15:13
Control Test	0.049	15:13
Air Blank	0.000	15:14
Control Test	0.049	15:14
Air Blank	0.000	15:15
Control Test	0.048	15:16
Air Blank	0.000	15:16
Control Test Stats		
Average	0.0487	
Std Dev	0.0006	
Rel Std Dev(%)	1.1863	

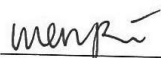


Operator's Signature

0.08g/210L
WET

ORANGE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001420
05/29/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	15:17
Control Test	0.078	15:18
Air Blank	0.000	15:19
Control Test	0.078	15:19
Air Blank	0.000	15:20
Control Test	0.079	15:21
Air Blank	0.000	15:21
Control Test Stats		
Average	0.0783	
Std Dev	0.0006	
Rel Std Dev(%)	0.7370	

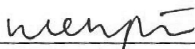


Operator's Signature

0.200g/210L

ORANGE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001420
05/29/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	15:22
Control Test	0.197	15:23
Air Blank	0.000	15:24
Control Test	0.196	15:24
Air Blank	0.000	15:25
Control Test	0.197	15:26
Air Blank	0.000	15:26
Control Test Stats		
Average	0.1967	
Std Dev	0.0006	
Rel Std Dev(%)	0.2936	

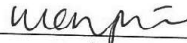

Operator's Signature

0.08g/210L
DGS #1

ORANGE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001420
05/29/2025
Software: 8100.27

* Inadvertently left
DGS unhooked
when performing
stability check,
identified user
error and
repeated analysis

Test	g/210L	Time
Air Blank	0.000	15:28
Control Test	0.000	15:28
Air Blank	0.000	15:28
Control Test	0.080	15:29
Air Blank	0.000	15:29
Control Test	0.079	15:30
Air Blank	0.000	15:30
Control Test Stats		
Average	0.0530	
Std Dev	0.0459	
Rel Std Dev(%)	86.6077	

menpi

Operator's Signature

0.08g/210L
DAS #2

ORANGE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001420
05/29/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	15:33
Control Test	0.079	15:33
Air Blank	0.000	15:34
Control Test	0.080	15:34
Air Blank	0.000	15:35
Control Test	0.080	15:35
Air Blank	0.000	15:35
Control Test Stats		
Average	0.0797	
Std Dev	0.0006	
Rel Std Dev(%)	0.7247	

Wenph

Operator's Signature

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: ORANGE COUNTY SO
Time of Inspection: 14:41

Date of Inspection: 06/03/2025

Serial Number: 80-001420
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.049	0.078	0.196	0.081
0.000	0.049	0.079	0.196	0.080
0.000	0.050	0.079	0.197	0.080
0.000	0.049	0.079	0.197	0.080
0.000	0.049	0.079	0.196	0.080
0.000	0.049	0.079	0.198	0.080
0.000	0.049	0.079	0.197	0.080
0.000	0.049	0.079	0.197	0.080
0.000	0.049	0.079	0.197	0.080
0.000	0.049	0.079	0.197	0.080
0.000	0.049	0.079	0.197	0.080
Standard Deviations	0.0003	0.0003	0.0006	0.0003

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.



I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

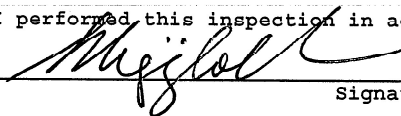

LEANDRA HIGGINBOTHAM

Signature and Printed Name

06/03/2025
Date

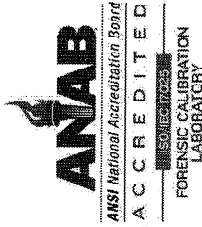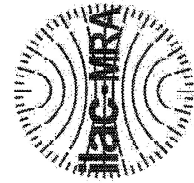

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001420, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-001420	UNCERTAINTY* $\pm$	
Owning Agency:	ORANGE COUNTY SO	0.050 g/210 L	0.004
Calibration Date:	06/03/2025	0.080 g/210 L	0.004
Calibration Time:	14:41	0.200 g/210 L	0.007
		0.080 g/210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.


TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

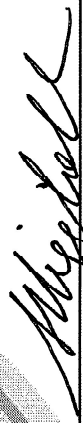
LEANDRA HIGGINBOTHAM,
Department Inspector

06/03/2025

Date

FDLE/ATP Form 69 October 2024
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality