

Stability Checks

0.05g/210L 0.047 to 0.053	0.08g/210L 0.077 to 0.083	0.20g/210L 0.194 to 0.206	DGS 0.08g/210L 0.077 to 0.083	≤0.003 of Wet																																																																																																																																																
HOVESTEAD PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006199 06/04/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>16:03</td></tr> <tr><td>Control Test</td><td>0.047</td><td>16:04</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:04</td></tr> <tr><td>Control Test</td><td>0.046</td><td>16:05</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:06</td></tr> <tr><td>Control Test</td><td>0.046</td><td>16:06</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:07</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0463</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.2461</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	16:03	Control Test	0.047	16:04	Air Blank	0.000	16:04	Control Test	0.046	16:05	Air Blank	0.000	16:06	Control Test	0.046	16:06	Air Blank	0.000	16:07	Control Test Stats			Average	0.0463		Std Dev	0.0006		Rel Std Dev(%)	1.2461		HOVESTEAD PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006199 06/04/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>16:25</td></tr> <tr><td>Control Test</td><td>0.077</td><td>16:26</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:26</td></tr> <tr><td>Control Test</td><td>0.077</td><td>16:27</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:27</td></tr> <tr><td>Control Test</td><td>0.077</td><td>16:28</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:29</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0770</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	16:25	Control Test	0.077	16:26	Air Blank	0.000	16:26	Control Test	0.077	16:27	Air Blank	0.000	16:27	Control Test	0.077	16:28	Air Blank	0.000	16:29	Control Test Stats			Average	0.0770		Std Dev	0.0000		Rel Std Dev(%)	0.0000		HOVESTEAD PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006199 06/04/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>16:12</td></tr> <tr><td>Control Test</td><td>0.195</td><td>16:13</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:14</td></tr> <tr><td>Control Test</td><td>0.194</td><td>16:14</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:15</td></tr> <tr><td>Control Test</td><td>0.195</td><td>16:16</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:16</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1947</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2966</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	16:12	Control Test	0.195	16:13	Air Blank	0.000	16:14	Control Test	0.194	16:14	Air Blank	0.000	16:15	Control Test	0.195	16:16	Air Blank	0.000	16:16	Control Test Stats			Average	0.1947		Std Dev	0.0006		Rel Std Dev(%)	0.2966		HOVESTEAD PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006199 06/04/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>16:19</td></tr> <tr><td>Control Test</td><td>0.078</td><td>16:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:19</td></tr> <tr><td>Control Test</td><td>0.077</td><td>16:20</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:20</td></tr> <tr><td>Control Test</td><td>0.077</td><td>16:21</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>16:21</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0773</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7466</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	16:19	Control Test	0.078	16:19	Air Blank	0.000	16:19	Control Test	0.077	16:20	Air Blank	0.000	16:20	Control Test	0.077	16:21	Air Blank	0.000	16:21	Control Test Stats			Average	0.0773		Std Dev	0.0006		Rel Std Dev(%)	0.7466		 Operator's Signature  Operator's Signature  Operator's Signature  Operator's Signature
Test	g/210L	Time																																																																																																																																																		
Air Blank	0.000	16:03																																																																																																																																																		
Control Test	0.047	16:04																																																																																																																																																		
Air Blank	0.000	16:04																																																																																																																																																		
Control Test	0.046	16:05																																																																																																																																																		
Air Blank	0.000	16:06																																																																																																																																																		
Control Test	0.046	16:06																																																																																																																																																		
Air Blank	0.000	16:07																																																																																																																																																		
Control Test Stats																																																																																																																																																				
Average	0.0463																																																																																																																																																			
Std Dev	0.0006																																																																																																																																																			
Rel Std Dev(%)	1.2461																																																																																																																																																			
Test	g/210L	Time																																																																																																																																																		
Air Blank	0.000	16:25																																																																																																																																																		
Control Test	0.077	16:26																																																																																																																																																		
Air Blank	0.000	16:26																																																																																																																																																		
Control Test	0.077	16:27																																																																																																																																																		
Air Blank	0.000	16:27																																																																																																																																																		
Control Test	0.077	16:28																																																																																																																																																		
Air Blank	0.000	16:29																																																																																																																																																		
Control Test Stats																																																																																																																																																				
Average	0.0770																																																																																																																																																			
Std Dev	0.0000																																																																																																																																																			
Rel Std Dev(%)	0.0000																																																																																																																																																			
Test	g/210L	Time																																																																																																																																																		
Air Blank	0.000	16:12																																																																																																																																																		
Control Test	0.195	16:13																																																																																																																																																		
Air Blank	0.000	16:14																																																																																																																																																		
Control Test	0.194	16:14																																																																																																																																																		
Air Blank	0.000	16:15																																																																																																																																																		
Control Test	0.195	16:16																																																																																																																																																		
Air Blank	0.000	16:16																																																																																																																																																		
Control Test Stats																																																																																																																																																				
Average	0.1947																																																																																																																																																			
Std Dev	0.0006																																																																																																																																																			
Rel Std Dev(%)	0.2966																																																																																																																																																			
Test	g/210L	Time																																																																																																																																																		
Air Blank	0.000	16:19																																																																																																																																																		
Control Test	0.078	16:19																																																																																																																																																		
Air Blank	0.000	16:19																																																																																																																																																		
Control Test	0.077	16:20																																																																																																																																																		
Air Blank	0.000	16:20																																																																																																																																																		
Control Test	0.077	16:21																																																																																																																																																		
Air Blank	0.000	16:21																																																																																																																																																		
Control Test Stats																																																																																																																																																				
Average	0.0773																																																																																																																																																			
Std Dev	0.0006																																																																																																																																																			
Rel Std Dev(%)	0.7466																																																																																																																																																			

HOMESTEAD PD
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-006199
 06/20/2025 12:57:49

Auto Calibration

pg 1 of 2

<<<<< 3um >>>>>			<<<<< 9um >>>>>		

Solution = 0.000 g/210L or 0.0000 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	0.1500	(-0.0150)	0.1650	(0.0520)	
Sample #2	0.1220	(0.0050)	0.2260	(0.0240)	
Sample #3	0.1470	(0.0120)	0.2320	(0.0130)	
Sample #4	0.1210	(0.0240)	0.2300	(0.0230)	
Avg % Abs	0.1300	(0.0137)	0.2293	(0.0200)	
STD DEV	0.0147	(0.0096)	0.0031	(0.0061)	
REL STD DEV	11.331	(70.310)	1.332	(30.414)	

Solution = 0.040 g/210L or 0.1905 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	0.8300	(0.0000)	1.5600	(-0.0050)	
Sample #2	0.8540	(-0.0080)	1.6080	(-0.0090)	
Sample #3	0.8620	(-0.0080)	1.5640	(-0.0220)	
Sample #4	0.8140	(0.0320)	1.5390	(0.0000)	
Avg % Abs	0.8433	(0.0053)	1.5703	(-0.0103)	
STD DEV	0.0257	(0.0231)	0.0349	(0.0111)	
REL STD DEV	3.049	(433.013)	2.225	(107.037)	

Solution = 0.100 g/210L or 0.4762 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	1.8790	(-0.0130)	3.5080	(-0.0030)	
Sample #2	1.8300	(0.0150)	3.4810	(0.0160)	
Sample #3	1.8800	(0.0070)	3.5240	(-0.0270)	
Sample #4	1.8400	(0.0370)	3.4540	(-0.0110)	
Avg % Abs	1.8500	(0.0197)	3.4863	(-0.0073)	
STD DEV	0.0265	(0.0155)	0.0353	(0.0217)	
REL STD DEV	1.430	(78.991)	1.013	(296.362)	

Solution = 0.200 g/210L or 0.9524 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	3.5430	(0.0000)	6.6750	(-0.0250)	
Sample #2	3.5810	(0.0000)	6.6480	(0.0180)	
Sample #3	3.5190	(0.0300)	6.6620	(-0.0020)	
Sample #4	3.5770	(0.0120)	6.6750	(0.0110)	
Avg % Abs	3.5590	(0.0140)	6.6617	(0.0090)	
STD DEV	0.0347	(0.0151)	0.0135	(0.0101)	
REL STD DEV	0.975	(107.855)	0.203	(112.765)	

Solution = 0.300 g/210L or 1.4286 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	5.2360	(-0.0070)	9.7050	(-0.0100)	
Sample #2	5.2060	(0.0390)	9.6960	(0.0140)	
Sample #3	5.1740	(0.0470)	9.6430	(0.0070)	
Sample #4	5.2310	(0.0300)	9.6470	(0.0270)	
Avg % Abs	5.2037	(0.0387)	9.6620	(0.0160)	
STD DEV	0.0286	(0.0085)	0.0295	(0.0101)	
REL STD DEV	0.549	(21.995)	0.305	(63.431)	

HOMESTEAD PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006199
06/20/2025 12:57:49

Auto Calibration

pg 2 of 2

<<<<< 3um >>>>>

Zero Order Coef -365.10
First Order Coef 2716.88
Second Order Coef 18.76

<<<<< 9um >>>>>

Zero Order Coef -326.94
First Order Coef 1416.50
Second Order Coef 9.84

Act Fit Residual
(g/210L) (g/210L) (g/210L)
0.000 -0.000 0.0002
0.040 0.041 -0.0007
0.100 0.099 0.0008
0.200 0.200 -0.0004
0.300 0.300 0.0001

Act Fit Residual
(g/210L) (g/210L) (g/210L)
0.000 -0.000 0.0000
0.040 0.040 -0.0004
0.100 0.099 0.0006
0.200 0.200 -0.0005
0.300 0.300 0.0002

<<<<< 3um >>>>> <<<<< 9um >>>>>

Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1
Sample
Sample #1 3083.00 3255.00
Sample #2 3091.00 3246.00
Sample #3 3157.00 3202.00
Sample #4 3193.00 3285.00
Avg 3147.0000 3244.3333
STD DEV 51.7301 41.5251
REL STD DEV 1.644 1.280
H2O adjust (mg/l*10k) 662 565

Barometric Pressure = 1017

*****CALIBRATION SUCCESSFUL*****

0.05g/210L 0.047 to 0.053	0.08g/210L 0.077 to 0.083	0.20g/210L 0.194 to 0.206	DGS 0.08g/210L 0.077 to 0.083																																																																																																																																																
HOMESTAD PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006199 06/20/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:30</td></tr> <tr><td>Control Test</td><td>0.049</td><td>15:31</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:31</td></tr> <tr><td>Control Test</td><td>0.050</td><td>15:32</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:32</td></tr> <tr><td>Control Test</td><td>0.048</td><td>15:33</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:34</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0490</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>2.0408</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	15:30	Control Test	0.049	15:31	Air Blank	0.000	15:31	Control Test	0.050	15:32	Air Blank	0.000	15:32	Control Test	0.048	15:33	Air Blank	0.000	15:34	Control Test Stats			Average	0.0490		Std Dev	0.0010		Rel Std Dev(%)	2.0408		HOMESTAD PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006199 06/20/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:22</td></tr> <tr><td>Control Test</td><td>0.080</td><td>15:23</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:23</td></tr> <tr><td>Control Test</td><td>0.081</td><td>15:24</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:24</td></tr> <tr><td>Control Test</td><td>0.080</td><td>15:25</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:26</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0803</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7187</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	15:22	Control Test	0.080	15:23	Air Blank	0.000	15:23	Control Test	0.081	15:24	Air Blank	0.000	15:24	Control Test	0.080	15:25	Air Blank	0.000	15:26	Control Test Stats			Average	0.0803		Std Dev	0.0006		Rel Std Dev(%)	0.7187		HOMESTAD PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006199 06/20/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:14</td></tr> <tr><td>Control Test</td><td>0.200</td><td>15:15</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:16</td></tr> <tr><td>Control Test</td><td>0.199</td><td>15:16</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:17</td></tr> <tr><td>Control Test</td><td>0.198</td><td>15:18</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:18</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1990</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.5025</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	15:14	Control Test	0.200	15:15	Air Blank	0.000	15:16	Control Test	0.199	15:16	Air Blank	0.000	15:17	Control Test	0.198	15:18	Air Blank	0.000	15:18	Control Test Stats			Average	0.1990		Std Dev	0.0010		Rel Std Dev(%)	0.5025		HOMESTAD PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006199 06/20/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:09</td></tr> <tr><td>Control Test</td><td>0.079</td><td>15:09</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:10</td></tr> <tr><td>Control Test</td><td>0.080</td><td>15:10</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:10</td></tr> <tr><td>Control Test</td><td>0.079</td><td>15:11</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:11</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0793</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7277</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	15:09	Control Test	0.079	15:09	Air Blank	0.000	15:10	Control Test	0.080	15:10	Air Blank	0.000	15:10	Control Test	0.079	15:11	Air Blank	0.000	15:11	Control Test Stats			Average	0.0793		Std Dev	0.0006		Rel Std Dev(%)	0.7277	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:30																																																																																																																																																	
Control Test	0.049	15:31																																																																																																																																																	
Air Blank	0.000	15:31																																																																																																																																																	
Control Test	0.050	15:32																																																																																																																																																	
Air Blank	0.000	15:32																																																																																																																																																	
Control Test	0.048	15:33																																																																																																																																																	
Air Blank	0.000	15:34																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0490																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	2.0408																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:22																																																																																																																																																	
Control Test	0.080	15:23																																																																																																																																																	
Air Blank	0.000	15:23																																																																																																																																																	
Control Test	0.081	15:24																																																																																																																																																	
Air Blank	0.000	15:24																																																																																																																																																	
Control Test	0.080	15:25																																																																																																																																																	
Air Blank	0.000	15:26																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0803																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7187																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:14																																																																																																																																																	
Control Test	0.200	15:15																																																																																																																																																	
Air Blank	0.000	15:16																																																																																																																																																	
Control Test	0.199	15:16																																																																																																																																																	
Air Blank	0.000	15:17																																																																																																																																																	
Control Test	0.198	15:18																																																																																																																																																	
Air Blank	0.000	15:18																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1990																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	0.5025																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:09																																																																																																																																																	
Control Test	0.079	15:09																																																																																																																																																	
Air Blank	0.000	15:10																																																																																																																																																	
Control Test	0.080	15:10																																																																																																																																																	
Air Blank	0.000	15:10																																																																																																																																																	
Control Test	0.079	15:11																																																																																																																																																	
Air Blank	0.000	15:11																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0793																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7277																																																																																																																																																		
Operator's Signature	Operator's Signature	Operator's Signature	Operator's Signature																																																																																																																																																

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: HOMESTEAD PD
Time of Inspection: 13:42

Date of Inspection: 06/25/2025

Serial Number: 80-006199
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.049	0.079	0.199	0.079
0.000	0.048	0.079	0.199	0.078
0.000	0.048	0.079	0.198	0.078
0.000	0.048	0.080	0.199	0.078
0.000	0.048	0.079	0.199	0.078
0.000	0.048	0.080	0.199	0.078
0.000	0.048	0.079	0.198	0.079
0.000	0.048	0.079	0.198	0.078
0.000	0.048	0.079	0.198	0.078
0.000	0.048	0.079	0.199	0.078

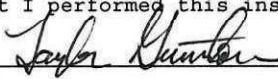
Standard Deviations	0.0003	0.0004	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.


 TAYLOR D GUTSCHOW
 Signature and Printed Name

06/25/2025
 Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006199 , manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-006199</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>HOMESTEAD PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>06/25/2025</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>13:42</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Taylor
Gutschow
Digitally signed by Taylor
Gutschow
Date: 2025.06.27 12:22:54
-04'00'

06/25/2025

Date

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality