



INSTRUMENT PROCESSING SHEET

Agency Manatee County Sheriff's OfficeS/N 80-006631Florida Department of
Law EnforcementDate In 05/30/2025 DI Completion Date 06/25/2025☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>TDG</u> Date <u>06/06/2025</u>		Quality Checks By <u>TDG</u> Date <u>06/13/2025</u>		Flow Calibration By _____ Date _____																																																											
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>212</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP106</u> 32 mm <u>0.160</u> (.139 - .169) 36 mm <u>0.171</u> (.156 - .190) 53 mm <u>0.246</u> (.228 - .278) 103 mm <u>0.507</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28199</u> ☒ Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP6287</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP6288</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602 10/22/2026</td></tr></tbody></table>		Simulator	Serial #	Lot #/Exp	0.050	MP6286	202406K 06/19/2026	0.080	MP6287	202406L 06/19/2026	0.200	MP6288	202406N 06/20/2026	0.080 DGS	N/A	AG429602 10/22/2026	Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____ _____																																												
				Simulator	Serial #	Lot #/Exp																																																									
				0.050	MP6286	202406K 06/19/2026																																																									
				0.080	MP6287	202406L 06/19/2026																																																									
0.200	MP6288	202406N 06/20/2026																																																													
0.080 DGS	N/A	AG429602 10/22/2026																																																													
Calibration Adjustment By <u>TDG</u>		Department Inspection By <u>TDG</u>																																																													
Barometric Pressure Gauge <u>1019</u> ID # <u>26932</u>		Barometric Pressure ID# <u>28199</u> Gauge <u>1020</u> Instrument <u>1020</u> Mouth Alcohol Solution Lot # <u>2024-A</u> Acetone Stock Solution Lot # <u>2023-B</u>																																																													
<table border="1" style="width:100%; border-collapse: collapse;"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td>MP6285</td><td>23400</td><td>10/24/2025</td></tr><tr><td>0.100</td><td>MP6286</td><td>24110</td><td>03/05/2026</td></tr><tr><td>0.200</td><td>MP6287</td><td>24080</td><td>02/13/2026</td></tr><tr><td>0.300</td><td>MP6288</td><td>23410</td><td>11/01/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>28424080A3</td><td>11/05/2026</td></tr></tbody></table> ☒ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td>MP4863</td><td>202406K</td><td>06/19/2026</td></tr><tr><td>0.080</td><td>MP4864</td><td>202406L</td><td>06/19/2026</td></tr><tr><td>0.200</td><td>MP6284</td><td>202406N</td><td>06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602</td><td>10/22/2026</td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000	MP6284	N/A	N/A	0.040	MP6285	23400	10/24/2025	0.100	MP6286	24110	03/05/2026	0.200	MP6287	24080	02/13/2026	0.300	MP6288	23410	11/01/2025	0.080 DGS	N/A	28424080A3	11/05/2026	Simulator	Serial #	Lot #	Expiration	0.050	MP4863	202406K	06/19/2026	0.080	MP4864	202406L	06/19/2026	0.200	MP6284	202406N	06/20/2026	0.080 DGS	N/A	AG429602	10/22/2026	<table border="1" style="width:100%; border-collapse: collapse;"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP5097</td></tr><tr><td>Interferent</td><td>MP5098</td></tr><tr><td>0.050</td><td>MP5099</td></tr><tr><td>0.080</td><td>MP5100</td></tr><tr><td>0.200</td><td>MP5101</td></tr></tbody></table> Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☒ Calibration Adjustment ☒ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____		Simulator	Serial Number	0.000	MP5097	Interferent	MP5098	0.050	MP5099	0.080	MP5100	0.200	MP5101
Simulator	Serial #	Lot #	Expiration																																																												
0.000	MP6284	N/A	N/A																																																												
0.040	MP6285	23400	10/24/2025																																																												
0.100	MP6286	24110	03/05/2026																																																												
0.200	MP6287	24080	02/13/2026																																																												
0.300	MP6288	23410	11/01/2025																																																												
0.080 DGS	N/A	28424080A3	11/05/2026																																																												
Simulator	Serial #	Lot #	Expiration																																																												
0.050	MP4863	202406K	06/19/2026																																																												
0.080	MP4864	202406L	06/19/2026																																																												
0.200	MP6284	202406N	06/20/2026																																																												
0.080 DGS	N/A	AG429602	10/22/2026																																																												
Simulator	Serial Number																																																														
0.000	MP5097																																																														
Interferent	MP5098																																																														
0.050	MP5099																																																														
0.080	MP5100																																																														
0.200	MP5101																																																														
Notes/Suggested Service: <u>Re-used MP6284 for post-cal stability. Rinsed and dried prior to adding new solution. (TDG 6/20/25)</u> _____ _____ _____ _____ _____		☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																																													
				Phil Nicodemo Digitally signed by Phil Nicodemo Date: 2025.07.01 10:19:29 +04'00' LeAndra Higginbotham Digitally signed by LeAndra Higginbotham Date: 2025.07.08 08:44:51 -04'00' Tech Review / Date _____ Admin Review / Date _____																																																											

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																																
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083																																																																																																																																																
MANATEE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006631 05/13/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:36</td></tr> <tr><td>Control Test</td><td>0.048</td><td>13:37</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:38</td></tr> <tr><td>Control Test</td><td>0.048</td><td>13:38</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:39</td></tr> <tr><td>Control Test</td><td>0.048</td><td>13:40</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:40</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0480</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel. Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> Operator's Signature <i>ML</i>	Test	g/210L	Time	Air Blank	0.000	13:36	Control Test	0.048	13:37	Air Blank	0.000	13:38	Control Test	0.048	13:38	Air Blank	0.000	13:39	Control Test	0.048	13:40	Air Blank	0.000	13:40	Control Test Stats			Average	0.0480		Std Dev	0.0000		Rel. Std Dev(%)	0.0000		MANATEE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006631 05/13/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:30</td></tr> <tr><td>Control Test</td><td>0.077</td><td>13:31</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:32</td></tr> <tr><td>Control Test</td><td>0.078</td><td>13:32</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:33</td></tr> <tr><td>Control Test</td><td>0.078</td><td>13:34</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:34</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0777</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel. Std Dev(%)</td><td>0.7434</td><td></td></tr> </tbody> </table> Operator's Signature <i>ML</i>	Test	g/210L	Time	Air Blank	0.000	13:30	Control Test	0.077	13:31	Air Blank	0.000	13:32	Control Test	0.078	13:32	Air Blank	0.000	13:33	Control Test	0.078	13:34	Air Blank	0.000	13:34	Control Test Stats			Average	0.0777		Std Dev	0.0006		Rel. Std Dev(%)	0.7434		MANATEE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006631 06/13/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:04</td></tr> <tr><td>Control Test</td><td>0.196</td><td>14:04</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:05</td></tr> <tr><td>Control Test</td><td>0.196</td><td>14:06</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:06</td></tr> <tr><td>Control Test</td><td>0.196</td><td>14:07</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:07</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1960</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel. Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> Operator's Signature <i>ML</i>	Test	g/210L	Time	Air Blank	0.000	14:04	Control Test	0.196	14:04	Air Blank	0.000	14:05	Control Test	0.196	14:06	Air Blank	0.000	14:06	Control Test	0.196	14:07	Air Blank	0.000	14:07	Control Test Stats			Average	0.1960		Std Dev	0.0000		Rel. Std Dev(%)	0.0000		MANATEE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006631 06/13/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:15</td></tr> <tr><td>Control Test</td><td>0.083</td><td>13:16</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:16</td></tr> <tr><td>Control Test</td><td>0.083</td><td>13:16</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:17</td></tr> <tr><td>Control Test</td><td>0.082</td><td>13:17</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:18</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0827</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel. Std Dev(%)</td><td>0.6984</td><td></td></tr> </tbody> </table> Operator's Signature <i>ML</i>	Test	g/210L	Time	Air Blank	0.000	13:15	Control Test	0.083	13:16	Air Blank	0.000	13:16	Control Test	0.083	13:16	Air Blank	0.000	13:17	Control Test	0.082	13:17	Air Blank	0.000	13:18	Control Test Stats			Average	0.0827		Std Dev	0.0006		Rel. Std Dev(%)	0.6984	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	13:36																																																																																																																																																	
Control Test	0.048	13:37																																																																																																																																																	
Air Blank	0.000	13:38																																																																																																																																																	
Control Test	0.048	13:38																																																																																																																																																	
Air Blank	0.000	13:39																																																																																																																																																	
Control Test	0.048	13:40																																																																																																																																																	
Air Blank	0.000	13:40																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0480																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel. Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	13:30																																																																																																																																																	
Control Test	0.077	13:31																																																																																																																																																	
Air Blank	0.000	13:32																																																																																																																																																	
Control Test	0.078	13:32																																																																																																																																																	
Air Blank	0.000	13:33																																																																																																																																																	
Control Test	0.078	13:34																																																																																																																																																	
Air Blank	0.000	13:34																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0777																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel. Std Dev(%)	0.7434																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	14:04																																																																																																																																																	
Control Test	0.196	14:04																																																																																																																																																	
Air Blank	0.000	14:05																																																																																																																																																	
Control Test	0.196	14:06																																																																																																																																																	
Air Blank	0.000	14:06																																																																																																																																																	
Control Test	0.196	14:07																																																																																																																																																	
Air Blank	0.000	14:07																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1960																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel. Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	13:15																																																																																																																																																	
Control Test	0.083	13:16																																																																																																																																																	
Air Blank	0.000	13:16																																																																																																																																																	
Control Test	0.083	13:16																																																																																																																																																	
Air Blank	0.000	13:17																																																																																																																																																	
Control Test	0.082	13:17																																																																																																																																																	
Air Blank	0.000	13:18																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0827																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel. Std Dev(%)	0.6984																																																																																																																																																		

MANATEE COUNTY SO
 Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-006631
 06/20/2025 12:51:38

Auto Calibration

pg 1 of 2

<<<<< 3um >>>>>			<<<<< 9um >>>>>		

Solution = 0.000 g/210L or 0.0000 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	0.0880	(-0.0070)	0.2090	(-0.0020)	
Sample #2	0.1050	(-0.0020)	0.1980	(0.0140)	
Sample #3	0.1020	(0.0090)	0.1900	(0.0070)	
Sample #4	0.0530	(0.0360)	0.1790	(0.0200)	
Avg % Abs	0.0867	(0.0143)	0.1890	(0.0137)	
STD DEV	0.0292	(0.0196)	0.0095	(0.0065)	
REL STD DEV	33.686	(136.419)	5.047	(47.608)	

Solution = 0.040 g/210L or 0.1905 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	0.7920	(-0.0190)	1.5780	(-0.0120)	
Sample #2	0.7730	(0.0000)	1.6020	(-0.0150)	
Sample #3	0.7810	(0.0080)	1.5810	(-0.0020)	
Sample #4	0.7490	(0.0320)	1.5510	(0.0130)	
Avg % Abs	0.7677	(0.0133)	1.5780	(-0.0013)	
STD DEV	0.0167	(0.0167)	0.0256	(0.0140)	
REL STD DEV	2.169	(124.900)	1.624	(1050.893)	

Solution = 0.100 g/210L or 0.4762 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	1.7560	(-0.0200)	3.5600	(-0.0030)	
Sample #2	1.7790	(-0.0240)	3.5680	(-0.0080)	
Sample #3	1.7680	(0.0000)	3.5650	(0.0050)	
Sample #4	1.7680	(0.0110)	3.5700	(0.0100)	
Avg % Abs	1.7717	(-0.0043)	3.5677	(0.0023)	
STD DEV	0.0064	(0.0179)	0.0025	(0.0093)	
REL STD DEV	0.358	(413.027)	0.071	(398.210)	

Solution = 0.200 g/210L or 0.9524 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	3.4470	(-0.0080)	6.8400	(-0.0170)	
Sample #2	3.4470	(0.0000)	6.8100	(0.0080)	
Sample #3	3.4010	(0.0300)	6.7960	(0.0100)	
Sample #4	3.4060	(0.0330)	6.7790	(0.0270)	
Avg % Abs	3.4180	(0.0210)	6.7950	(0.0150)	
STD DEV	0.0252	(0.0182)	0.0155	(0.0104)	
REL STD DEV	0.738	(86.897)	0.228	(69.602)	

Solution = 0.300 g/210L or 1.4286 mg/l, Samples = 4, Discarded = 1					
Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)	
Sample #1	5.0560	(-0.0190)	9.8720	(-0.0020)	
Sample #2	5.0370	(0.0000)	9.8640	(0.0040)	
Sample #3	5.0350	(0.0020)	9.8460	(0.0170)	
Sample #4	4.9890	(0.0320)	9.8380	(0.0150)	
Avg % Abs	5.0203	(0.0113)	9.8493	(0.0120)	
STD DEV	0.0272	(0.0179)	0.0133	(0.0070)	
REL STD DEV	0.541	(158.169)	0.135	(58.333)	

MANATEE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006631
06/20/2025 12:51:38

Auto Calibration

pg 2 of 2

<<<< 3um >>>>

Zero Order Coef -241.59
First Order Coef 2782.35
Second Order Coef 22.13

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	-0.000	0.0000
0.040	0.040	-0.0001
0.100	0.100	0.0001
0.200	0.200	-0.0001
0.300	0.300	0.0000

<<<< 9um >>>>

-261.71
1362.58
11.59

Act (g/210L)	Fit (g/210L)	Residual (g/210L)
0.000	-0.000	0.0001
0.040	0.040	-0.0003
0.100	0.100	0.0003
0.200	0.200	-0.0002
0.300	0.300	0.0001

<<<< 3um >>>>

Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1
Sample

Sample #1	3195.00	3193.00
Sample #2	3195.00	3200.00
Sample #3	3178.00	3214.00
Sample #4	3170.00	3166.00
Avg	3181.0000	3193.3333
STD DEV	12.7671	24.6847
REL STD DEV	0.401	0.773
H2O adjust (mg/l*10k)	628	616

Barometric Pressure = 1019

*****CALIBRATION SUCCESSFUL*****

Post-Cal Stability Checks

0.05g/210L 0.047 to 0.053	0.08g/210L 0.077 to 0.083	0.20g/210L 0.194 to 0.206	DGS 0.08g/210L 0.077 to 0.083	≤0.003 of Wet																																																																																																																																																
MANATEE COUNTY SO Intoxilyzer - Alcotest Analyzer Model 8000 06/20/2025 Software: 8100.27 SN 80-006631 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:28</td></tr> <tr><td>Control Test</td><td>0.049</td><td>14:29</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:29</td></tr> <tr><td>Control Test</td><td>0.048</td><td>14:30</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:31</td></tr> <tr><td>Control Test</td><td>0.049</td><td>14:31</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:32</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0487</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.1863</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	14:28	Control Test	0.049	14:29	Air Blank	0.000	14:29	Control Test	0.048	14:30	Air Blank	0.000	14:31	Control Test	0.049	14:31	Air Blank	0.000	14:32	Control Test Stats			Average	0.0487		Std Dev	0.0006		Rel Std Dev(%)	1.1863		MANATEE COUNTY SO Intoxilyzer - Alcotest Analyzer Model 8000 06/20/2025 Software: 8100.27 SN 80-006631 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:36</td></tr> <tr><td>Control Test</td><td>0.080</td><td>14:37</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:37</td></tr> <tr><td>Control Test</td><td>0.079</td><td>14:38</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:39</td></tr> <tr><td>Control Test</td><td>0.080</td><td>14:39</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:40</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0797</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7247</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	14:36	Control Test	0.080	14:37	Air Blank	0.000	14:37	Control Test	0.079	14:38	Air Blank	0.000	14:39	Control Test	0.080	14:39	Air Blank	0.000	14:40	Control Test Stats			Average	0.0797		Std Dev	0.0006		Rel Std Dev(%)	0.7247		MANATEE COUNTY SO Intoxilyzer - Alcotest Analyzer Model 8000 06/20/2025 Software: 8100.27 SN 80-006631 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:51</td></tr> <tr><td>Control Test</td><td>0.199</td><td>14:52</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:53</td></tr> <tr><td>Control Test</td><td>0.198</td><td>14:53</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:54</td></tr> <tr><td>Control Test</td><td>0.198</td><td>14:55</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:55</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1983</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2911</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	14:51	Control Test	0.199	14:52	Air Blank	0.000	14:53	Control Test	0.198	14:53	Air Blank	0.000	14:54	Control Test	0.198	14:55	Air Blank	0.000	14:55	Control Test Stats			Average	0.1983		Std Dev	0.0006		Rel Std Dev(%)	0.2911		MANATEE COUNTY SO Intoxilyzer - Alcotest Analyzer Model 8000 06/20/2025 Software: 8100.27 SN 80-006631 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:18</td></tr> <tr><td>Control Test</td><td>0.079</td><td>14:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:19</td></tr> <tr><td>Control Test</td><td>0.079</td><td>14:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:20</td></tr> <tr><td>Control Test</td><td>0.079</td><td>14:20</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:21</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0790</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	14:18	Control Test	0.079	14:19	Air Blank	0.000	14:19	Control Test	0.079	14:19	Air Blank	0.000	14:20	Control Test	0.079	14:20	Air Blank	0.000	14:21	Control Test Stats			Average	0.0790		Std Dev	0.0000		Rel Std Dev(%)	0.0000		Operator's Signature Operator's Signature Operator's Signature Operator's Signature
Test	g/210L	Time																																																																																																																																																		
Air Blank	0.000	14:28																																																																																																																																																		
Control Test	0.049	14:29																																																																																																																																																		
Air Blank	0.000	14:29																																																																																																																																																		
Control Test	0.048	14:30																																																																																																																																																		
Air Blank	0.000	14:31																																																																																																																																																		
Control Test	0.049	14:31																																																																																																																																																		
Air Blank	0.000	14:32																																																																																																																																																		
Control Test Stats																																																																																																																																																				
Average	0.0487																																																																																																																																																			
Std Dev	0.0006																																																																																																																																																			
Rel Std Dev(%)	1.1863																																																																																																																																																			
Test	g/210L	Time																																																																																																																																																		
Air Blank	0.000	14:36																																																																																																																																																		
Control Test	0.080	14:37																																																																																																																																																		
Air Blank	0.000	14:37																																																																																																																																																		
Control Test	0.079	14:38																																																																																																																																																		
Air Blank	0.000	14:39																																																																																																																																																		
Control Test	0.080	14:39																																																																																																																																																		
Air Blank	0.000	14:40																																																																																																																																																		
Control Test Stats																																																																																																																																																				
Average	0.0797																																																																																																																																																			
Std Dev	0.0006																																																																																																																																																			
Rel Std Dev(%)	0.7247																																																																																																																																																			
Test	g/210L	Time																																																																																																																																																		
Air Blank	0.000	14:51																																																																																																																																																		
Control Test	0.199	14:52																																																																																																																																																		
Air Blank	0.000	14:53																																																																																																																																																		
Control Test	0.198	14:53																																																																																																																																																		
Air Blank	0.000	14:54																																																																																																																																																		
Control Test	0.198	14:55																																																																																																																																																		
Air Blank	0.000	14:55																																																																																																																																																		
Control Test Stats																																																																																																																																																				
Average	0.1983																																																																																																																																																			
Std Dev	0.0006																																																																																																																																																			
Rel Std Dev(%)	0.2911																																																																																																																																																			
Test	g/210L	Time																																																																																																																																																		
Air Blank	0.000	14:18																																																																																																																																																		
Control Test	0.079	14:19																																																																																																																																																		
Air Blank	0.000	14:19																																																																																																																																																		
Control Test	0.079	14:19																																																																																																																																																		
Air Blank	0.000	14:20																																																																																																																																																		
Control Test	0.079	14:20																																																																																																																																																		
Air Blank	0.000	14:21																																																																																																																																																		
Control Test Stats																																																																																																																																																				
Average	0.0790																																																																																																																																																			
Std Dev	0.0000																																																																																																																																																			
Rel Std Dev(%)	0.0000																																																																																																																																																			

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MANATEE COUNTY SO
Time of Inspection: 12:11

Date of Inspection: 06/25/2025

Serial Number: 80-006631
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.048	0.079	0.199	0.079
0.000	0.048	0.079	0.200	0.079
0.000	0.048	0.079	0.199	0.079
0.000	0.048	0.079	0.200	0.079
0.000	0.048	0.080	0.200	0.079
0.000	0.049	0.079	0.200	0.079
0.000	0.048	0.079	0.199	0.078
0.000	0.048	0.079	0.199	0.079
0.000	0.048	0.079	0.199	0.079
0.000	0.048	0.079	0.199	0.078


Standard Deviations	0.0003	0.0003	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

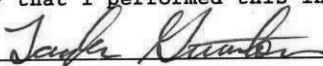

TAYLOR D GUTSCHOW

Signature and Printed Name

06/25/2025
Date

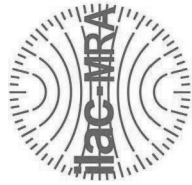

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

Calibration Certificate

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006631, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-006631</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>MANATEE COUNTY SO</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>06/25/2025</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>12:11</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Taylor
Gutschow
Date: 2025.06.27 12:25:34
-04'00'

Date
06/25/2025
TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality